



## Laboratory Procedure Manual

*Analytes:* **Antimony, Arsenic, Barium, Beryllium, Cadmium, Cesium, Cobalt, Lead, Manganese, Molybdenum, Platinum, Strontium, Thallium, Tin, Tungsten, and Uranium**

*Matrix:* **Urine**

*Method:* **Urine Multi-Element ICP-DRC-MS**  
Renamed from "Inductively Coupled Plasma-Mass Spectrometry (ICP-DRC-MS)"

*Method No:* **3018.3 (15 element panel) and 3018A.2 (total arsenic)**

*Revised:* **March 19, 2012**

*As performed by:* Inorganic Radionuclides and Toxicology  
Division of Laboratory Sciences  
National Center for Environmental Health

*Contact:* Dr. Kathleen L. Caldwell  
Phone: 770-488-7990  
Fax: 770-488-4097  
Email: [KCaldwell@cdc.gov](mailto:KCaldwell@cdc.gov)

James L. Pirkle, M.D., Ph.D.  
Director, Division of Laboratory Sciences

### **Important Information for Users**

The Centers for Disease Control and Prevention (CDC) periodically refines these laboratory methods. It is the responsibility of the user to contact the person listed on the title page of each write-up before using the analytical method to find out whether any changes have been made and what revisions, if any, have been incorporated.

**Multi-Element in urine  
NHANES 2011-2012**

**Public Release Data Set Information**

This document details the Lab Protocol for testing the items listed in the following table:

<b>Data File Name</b>	<b>Variable Name</b>	<b>SAS Label</b>
UHM_G	URXUAS	Urinary Total Arsenic (µg/L)
	URXUBA	Barium, urine (ng/mL)
	URDUCD	Cadmium, urine (ng/mL)
	URXUCO	Cobalt, urine (ng/mL)
	URXUCS	Cesium, urine (ng/mL)
	URXUMN	Manganese, urine (ng/mL)
	URXUMO	Molybdenum, urine (ng/mL)
	URXUPB	Lead, urine (ng/mL)
	URXUSN	Tin, urine (ng/mL)
	URXUSR	Strontium, urine (ng/mL)h
	URXUSB	Antimony, urine (ng/mL)
	URXUTL	Thallium, urine (ng/mL)
	URXUTU	Tungsten, urine (ng/mL)
	URXUUR	Uranium, urine (ng/mL)

## **1) Clinical Relevance & Summary of Test Principle**

### **a. Clinical Relevance:**

This method is used to achieve rapid and accurate quantification of fifteen elements of toxicological and nutritional interest including Antimony (Sb), Arsenic (As), Barium (Ba), Beryllium (Be), Cadmium (Cd), Cesium (Cs), Cobalt (Co), Lead (Pb), Manganese (Mn), Molybdenum (Mo), Platinum (Pt), Strontium (Sr), Thallium (Tl), Tin (Sn), Tungsten (W), and Uranium (U). The method may be used to screen urine when people are suspected to be acutely exposed to these elements or to evaluate chronic environmental or other non-occupational exposure. [1-4].

### **b. Test Principle:**

Inductively coupled plasma dynamic reaction cell mass spectrometry (ICP-DRC-MS) is a multi-element analytical technique capable of trace level elemental analysis [1-4]. This ICP-DRC-MS method is used to measure either arsenic, a 15 element panel (Antimony, Barium, Beryllium, Cadmium, Cesium, Cobalt, Lead, Manganese, Molybdenum, Platinum, Strontium, Thallium, Tin, Tungsten, and Uranium), or any subgroup of the 15 element panel.

Liquid samples are introduced into the ICP through a nebulizer and spray chamber carried by a flowing argon stream. By coupling radio-frequency power into flowing argon, plasma is created in which the predominant species are positive argon ions and electrons and has a temperature of 6000-8000 K. The sample passes through a region of the plasma and the thermal energy atomizes the sample and then ionizes the atoms. The ions, along with the argon, enter the mass spectrometer through an interface that separates the ICP (at atmospheric pressure, ~760 torr) from the mass spectrometer (operating at a pressure of  $10^{-5}$  torr). The ions pass through a focusing region, the dynamic reaction cell, the quadrupole mass filter, and finally are counted in rapid sequence at the detector allowing individual isotopes of an element to be determined. The dynamic reaction cell operates in one of two modes. In 'standard' mode the cell is not pressurized and ions pass through the cell to the quadrupole mass filter unaffected. In 'DRC' mode the cell is pressurized with a gas which will collide or react with the incoming ions to either eliminate an interfering ion or change the ion of interest to a new mass which is free from interference. In this method the instrument is operated in DRC mode when analyzing for cadmium, manganese and arsenic, but in standard mode when analyzing for all of the other analytes. For arsenic, the reaction cell is pressurized with a mixture of hydrogen (10%) and argon (90%) which causes the breakup of the  $^{40}\text{Ar}^{35}\text{Cl}^+$  ion which would otherwise interfere with detection of  $^{75}\text{As}$  at  $m/z$  75. When analyzing for cadmium, the reaction cell is pressurized with oxygen. The  $^{98}\text{Mo}^{16}\text{O}^+$  ions which would normally interfere with detection of  $^{114}\text{Cd}$  at  $m/z$  114 react with the oxygen in the cell creating  $^{98}\text{Mo}^{16}\text{O}_2^+$  and  $^{98}\text{Mo}^{16}\text{O}_3^+$  at masses which no longer

represent an interference to  $^{114}\text{Cd}$  analysis. The DRC is also pressurized with oxygen gas when analyzing for  $^{55}\text{Mn}$ . The  $^{39}\text{K}^{16}\text{O}^+$  ions which would normally interfere with the detection of  $^{55}\text{Mn}$  at  $m/z$  55 react with the oxygen in the cell and no longer represent interference to  $^{55}\text{Mn}$  analysis. Electrical signals resulting from the detection of ions are processed into digital information that is used to indicate first the intensity of the ions and then the concentration of the element. This method was originally based on the method by Mulligan et al. [5]. The DRC portions of the method are based on work published by Tanner et al. [2, 3]. Urine samples are diluted 1+ 9 with 2% (v/v) concentrated nitric acid (and 1.5% ethanol in the case of arsenic). The diluent for the 15 element panel contains iridium (Ir), rhodium (Rh) for multi-internal standardization. The diluent for arsenic contains gallium (Ga) for internal standardization. Nitric acid is used for the purpose of solubilizing and stabilizing metals in solution. Internal standards are a constant concentration in all blanks, calibrators and samples. Monitoring the instrument signal ratio of a metal to its internal standard allows correction for instrument noise and drift, and sample-to-sample matrix differences. Ethanol is used in the case of arsenic for the purpose of providing a constant amount of signal enhancement (carbon effect) across all blanks, calibrators, and samples.

## **2) Limitations of Method; Interfering Substances and Conditions**

### **a. Interferences Addressed by This Method**

i. Breakup of Argon Chloride ( $^{40}\text{Ar}^{35}\text{Cl}$ ) Interference on Arsenic ( $^{75}\text{As}$ ) Using DRC: The dynamic reaction cell of the ELAN ICP-DRC-MS is used in this method to break apart the argon chloride ( $^{40}\text{Ar}^{35}\text{Cl}$ ) interference on arsenic at  $m/z$  75 [6] which is common to urine analysis by ICP-MS (see Section 1.b for an explanation of this process).

ii. Correction & Elimination of Interferences ( $^{114}\text{Sn}$ ,  $^{98}\text{Mo}^{16}\text{O}$ ) on Cadmium ( $^{114}\text{Cd}$ ).

1. Mathematical Correction for Tin ( $^{114}\text{Sn}$ ) Interference:

The correction equation ( $-0.026826*\text{Sn}118$ ) is used in the "Equations" tab of the method to correct the counts observed as  $m/z$  114 to exclude counts due to  $^{114}\text{Sn}$ .

2. Elimination of Molybdenum Oxide ( $^{98}\text{Mo}^{16}\text{O}$ ) Interference Using DRC:

The dynamic reaction cell of the ELAN ICP-DRC-MS is used in this method to eliminate interference from molybdenum oxide ( $^{98}\text{Mo}^{16}\text{O}$ ) onto cadmium at  $m/z$  114 [7]. See Section 1.b for an explanation of this process.

iii. Elimination of interference ( $^{39}\text{K}^{16}\text{O}$ ) on manganese  $^{55}\text{Mn}$  using DRC:

The dynamic reaction cell of the ELAN ICP-DRC-MS is used in this method to reduce the potassium oxide ( $^{39}\text{K}^{16}\text{O}$ ) interference on manganese at  $m/z$  55. See section 1.b for an explanation of this process.

iv. Matrix Enhancement of Arsenic Signal:

Matrix induced signal enhancement in ICP-MS analysis from carbon on arsenic has been previously reported in the literature [8, 9]. When arsenic is being determined by this method, ethanol (1.5% v/v) is added in the diluent and rinse solutions to “normalize” the arsenic signal enhancement in all blanks, calibrators, and samples.

b. Limitations of Method (Interferences Remaining in Method)

i. Calcium Chloride ( $^{40}\text{Ca}^{35}\text{Cl}$ ) Interference on Arsenic ( $^{75}\text{As}$ ):

It has been determined that a small interference remains at m/z 75 when the urine matrix contains **both** high chloride **and** high calcium levels [6]. Even at extreme calcium and chloride levels, this interference is has not been found to be significant (approximately 0.4  $\mu\text{g/L}$ ).

ii. Gallium Oxide ( $^{71}\text{Ga}^{17}\text{O}$ ) Interference on Strontium ( $^{88}\text{Sr}$ ):

Arsenic only analysis requires Ga to be used as the internal standard and is added to the diluent at a concentration of 10  $\mu\text{g/L}$  (see section 6.c for the preparation of diluent). Gallium should not be added to the diluent (for the purpose of being used as an internal standard for As) when strontium is measured (15 element panel) due to the formation of  $^{71}\text{Ga}^{17}\text{O}^+$ , which occurs at the same m/z as  $^{88}\text{Sr}^+$ . A 5  $\mu\text{g/L}$  solution of Ga in 2%  $\text{HNO}_3$  resulted in a background equivalent concentration (BEC) of 50  $\text{ng/L}$  for  $^{88}\text{Sr}$ . Based on these results, the expected increase in  $^{88}\text{Sr}$  concentration in a diluted urine sample is 1  $\mu\text{g/L}$ .

**3) Procedures for Collecting, Storing, and Handling Specimens; Criteria for Specimen Rejection; Specimen Accountability and Tracking**

- a. Procedures for Collecting, Storing, and Handling Specimens: Specimen handling conditions, special requirements, and procedures for collection and transport are discussed in the division (DLS) Policies and Procedures Manual [10]. Copies are available in branch, laboratory, and special activities specimen-handling offices. An electronic copy is available at:

[http://intranet.nceh.cdc.gov/dls/pdf/policiesprocedures/Policy\\_and\\_Procedures\\_Manual.DLS.2002mod.pdf](http://intranet.nceh.cdc.gov/dls/pdf/policiesprocedures/Policy_and_Procedures_Manual.DLS.2002mod.pdf). In general,

- i. No fasting or special diets are required before collection of urine.
- ii. Use sterile, lot screened collectors for specimen acquisition.
- iii. Urine specimens should be transported frozen (packed in dry ice during shipment is preferred when possible).

**Multi-Element in urine  
NHANES 2011-2012**

- iv. Once received, store long term at  $\leq -20^{\circ}\text{C}$  until time for analysis. Short-term storage at  $2-4^{\circ}\text{C}$  is acceptable. Refreeze at  $\leq -20^{\circ}\text{C}$  portions of the sample that remain after analytical aliquots are withdrawn. Thawing and refreezing samples has not been found to compromise sample results.
  - v. Acceptable containers for analytical aliquots include lot screened polypropylene (PP) cryovials or tubes (i.e. 5 mL cryogenic vial or 15mL centrifuge tube).
  - b. Criteria for Specimen Rejection: Specimen characteristics that may compromise test results are indicated above. Reasons for rejection of a sample for analysis include
    - i. Low volume: Optimal amount of urine is 1.8+ mL. The volume of urine used for one analysis is 0.5 mL.
    - ii. Contamination: Improper collection procedures or collection devices can contaminate the urine by contact with dust, dirt, etc.
- In all cases, request a second urine specimen.
- c. Transfer or Referral of Specimens; Procedures for Specimen Accountability and Tracking: Location, status, and final disposition of the specimens will be tracked at least by paper document in the "Study Folder" (created before analysts receive the samples). Apart from this specimen tracking form, this folder will also contain the paper print outs of results from analysis of the specimens. Maintain records for a minimum of 3 years. Use only numerical identifiers for samples within the laboratory (e.g., case ID numbers) in order to safeguard confidentiality. Only the medical supervisor (MS) or project coordinator (PC) i.e. non CDC personnel should have access to the personal identifiers.

#### **4) Safety Precautions**

- a. General Safety
  - i. Observe all safety regulations as detailed in the Division (DLS) Safety Manual. Additional information can be found in your lab's chemical hygiene plan.
  - ii. Observe Universal Precautions when working with urine.
  - iii. Wear appropriate gloves, lab coat, and safety glasses while handling all solutions. Consult the laboratory chemical hygiene plan.
  - iv. Exercise special care when handling and dispensing concentrated nitric acid. Add acid to water. Nitric acid is a caustic chemical that is capable of causing severe eye and skin damage. ***If nitric acid comes in contact with any part***

***of the body, quickly wash the affected area with copious quantities of water for at least 15 minutes.***

- v. Use secondary containment for containers of biological or corrosive liquids.
  - vi. The use of the foot pedal on the Micromedic Digiflex™ is recommended because it reduces analyst contact with work surfaces that have been in contact urine and also keeps the analyst's hands free to hold the specimen cups and autosampler tubes and to wipe off the tip of Micromedic Digiflex™.
  - vii. Training will be given before operating the ICP-DRC-MS, as there are many possible hazards including ultraviolet radiation, high voltages, radio-frequency radiation, and high temperatures. This information is also detailed in the PerkinElmer ELAN® ICP-DRC-MS System Safety Manual.
  - viii. Use flash arrestors on oxygen and argon / hydrogen gas cylinders and properly secure gas cylinders with safety harnesses.
  - ix. Wipe down all work surfaces at the end of the day with bleach-rite spray or freshly prepared 10% (v/v) sodium-hypochlorite solution.
- b. Radiation Safety: All personnel performing this method must successfully meet requirements of a CDC-OHS radiation worker (RW) due to the use of natural uranium in this method and observe all necessary radiation safety considerations indicated in the CDC Radiation Safety Manual [11].
- c. Waste Disposal: Operators of this method should take the CDC-OHS Hazardous Chemical Waste Management Course (initial and yearly refreshers).
- i. Waste to be Placed Into Biohazard Autoclave Bags & Pans:
    - 1. All biological samples and diluted specimens (after analysis run).
    - 2. All disposable plastic and paper which contact urine (autosampler tubes, gloves, etc.).
    - 3. Used non-glass/quartz ICP-MS consumables (i.e. probes, tubing, cones, ion lenses).
  - ii. Waste to be Placed Into Sharps Containers: Pipette Tips, broken glass or quartz instrument consumables (broken spray chambers, torches, nebulizers, etc. . .). Large broken glass which will not fit in the sharps container should be placed in a separate autoclave pan from other waste and labeled as "broken glass" (see the "Autoclaving" section of the CDC safety policies and practices manual located in the laboratory).

iii. Liquid Waste

1. Waste discarded down sink: Only liquid waste from the ICP-DRC-MS instrument can be discarded at the sink. Flush the sink with copious amounts of water.
2. Waste to be Picked up by the Radiation Safety Office: Contact the laboratory radiation inventory person and the CDC Radiation Safety Office for disposal of any single element uranium standard, intermediate stock standard, or intermediate working standard solutions.
3. Waste to be Picked up by Hazardous Waste Program: Submit request for hazardous waste removal of all other liquid waste.

**5) Instrument & Material Sources**

a. Sources for ICP-MS Instrumentation

i. ICP-MS: Inductively Coupled Plasma Dynamic Reaction Cell Mass Spectrometer (ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> or ELAN<sup>®</sup> DRC II) (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).

1. DXi-FAST upgrade: Standard peristaltic pump replaced by DXi\_FAST micro-peristaltic pump / FAST actuator and valve combination unit. For ELAN DRCII, part # DXI-54-P4-F6.

ii. Recirculating chiller / heat exchanger for ICP-MS: Refrigerated chiller (PolyScience 6105PE for ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> instruments) or heat exchanger (PolyScience 3370 for ELAN<sup>®</sup> DRC II instruments) (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).

iii. Autosampler: ESI SC-4 autosampler (Elemental Scientific Inc., Omaha, NE) or equivalent.

iv. FAST Sample Introduction System (Elemental Scientific Inc., Omaha, NE).

1. FAST controller
2. FAST actuator: CTFE high-flow valve head like part number SC-0599-1210 (part number includes lines and probes).

b. Sources for ICP-MS Parts & Consumables

NOTE: The minimum number of spares recommended before reordering (if owning one instrument) are listed as “# Spares =” in the descriptions below.



**Multi-Element in urine  
NHANES 2011-2012**

- i. Adapter, PEEK: Securely connects 1.6mm O.D. PFA tubing to 0.03" I.D. peristaltic tubing. Composed of three PEEK parts.
  1. Female nut for 1.6mm O.D. (1/16") tubing. Like part P-420 (Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)).
  2. PEEK ferrule. Like part P-260x (10pk SuperFlangeless ferrule, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)).
  3. Conical Adapter Body. Like part P-692 (Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)).
- ii. Bottles (for rinse solution): Four liter screw-cap polypropylene container with 2 luer connections (like catalog# SC-0305-1, Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
- iii. Carboy and cap assembly for waste collection: 10-15 L, polypropylene widemouth carboy (100 mm neck size) with handles and no spigot (Like part #7BE-25126, Lab Safety Supply, Janesville, WI, [www.lss.com](http://www.lss.com)) with cap assembly like part # N0690271 (PerkinElmer, Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).
- iv. Coolant, for Polyscience chiller or heat exchanger: Only PerkinElmer part # WE01-6558 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) is approved for use by PerkinElmer. # Spares = 6.
- v. Cone, sampler (nickel/platinum): PerkinElmer part # WE021140/WE027802 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Part # SC2011-Ni (Testing has also found Spectron, Ventura, CA, [www.spectronus.com](http://www.spectronus.com) cones to be comparable). # Spares = 4.
- vi. Cone, skimmer (nickel/platinum): PerkinElmer part # WE021137/WE027803 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Part # SC2012-Ni (Testing has also found Spectron, Ventura, CA, [www.spectronus.com](http://www.spectronus.com) cones to be comparable) # Spares = 4.
- vii. Connector (for tubing): Use to connect 1/8" I.D. PVC tubing to 0.125" I.D. peristaltic pump tubing. Use part # 3140715 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 4.
- viii. Detector, electron multiplier: Like part # N8125001 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer (part # 14210, SGE Incorporated, Austin, Texas, <http://www.etpsci.com>) or various distributors. # Spares = 1.
- ix. FAST accessories
  1. Valve: CTFE High-flow valve head for SC-FAST (uses ¼-28 fittings). Like part # SC-0599-1010 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
  2. Stator: CTFE Stator for 6 port SC-FAST high flow valve (¼-28 fittings). Like part # SC-0599-1010-01 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).

**Multi-Element in urine  
NHANES 2011-2012**

3. Rotor: Composite rotor for 6 port SC-FAST high flow valve (1/4-28 fittings). Like part # SC-0599-1010-05 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
4. Sample Loop:
  - a. DLS 3018: 3 mL Teflon, white connector-nuts for high flow valve head. Like part # SC-0315-30 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
  - b. DLS 3018 subset: A subset of 3 elements (mixed mode of DRC and Standard groups) can be analyzed using a 1.25 mL Teflon sample loop with white nut connectors for high flow valve head of FAST sample introduction system.
  - c. DLS 3018A or DLS 3018 single element: Default volume is 0.5mL Teflon sample loop with white nut connectors for high flow valve head of FAST sample introduction system. This volume loop can be created by cutting 25% off the length of a 2mL Teflon sample loop. Like part # SC-0315-20 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). Volumes larger than 0.5mL can be used, but will require longer loop fill (ESI software) and sample flush (ELAN software) times, and proportionally larger volumes used in sample preparation.
5. Probe, Autosampler: Teflon, carbon fiber support, 0.8mm i.d., blue marker, 1/4-28 fittings. Like part number SC-5037-3751 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 2.
6. Probe, Carrier Solution: Teflon, carbon fiber support, 0.5mm i.d., orange marker, 1/4-28 fittings. Like part number SC-5037-3501 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 2.
7. Tubing, FAST vacuum: Vacuum line for SC-FAST high flow valve, connects to port #6, black nut for connection to valve head, natural brown color nut on other end for connection to SC autosampler vacuum port. Like part # SC-0321 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
8. Tubing, connects nebulizer to valve: See "Nebulizer, PolyPro-ST micro flow"
- x. Hose, for connection to chiller: Push on hose. I.D. = 1/2", O.D. = 3/4". Use part # PB-8 (per inch, Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Do not normally need spare hose (unless moving instrument into a new location).
- xi. Hose, for exhaust of ELAN: Available as part of ELAN installation kit from Perkin Elmer (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer as part # S-LP-10 air connector (Thermafex, Abbeville, SC, [www.thermafex.net](http://www.thermafex.net)). Equivalent part may be substituted. # Spares = 10 feet of 4" diameter and 10 feet of 6" diameter hose.

**Multi-Element in urine  
NHANES 2011-2012**

- xii. Injector, quartz with ball joint: I.D. = 2.0 mm. PerkinElmer part # WE023948 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer as part # 400-30 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or from various distributors. # Spares = 2.
- xiii. Injector support (for pass-through injector): PerkinElmer part # WE023951 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Available direct from manufacturer as part # 400-37 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or from various distributors. # Spares = 2.
- xiv. Ion Lens: PerkinElmer part # WE018034 (PerkinElmer Norwalk, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). # Spares = 3.
- xv. Nebulizer: PolyPro-ST micro flow polypropylene nebulizer with external 1/4-28 threaded connector for liquid delivery, low pressure version or equivalent. Like part # ES-4040-7010 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 1. Different nebulizers may be used, however, the nebulizer gas flow rate, sample flush time, read delay time, loop fill time, loop size, urine sample dilution preparation volume, and sample-to-sample carry-over must be evaluated and optimized.
1. Gas connection:
    - a. Teflon tubing: 4mm o.d., 2.4mm i.d. Teflon tubing (like part # ES-2502, Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 1.
    - b. Adapter kit: Plastic adapters to connect Teflon tubing (2.4mm i.d) to ¼" male Swagelok (compression) port on ICP-DRC-MS. Parts can be obtained as components in a "gas fittings kit for microflow nebulizer", kit part # ES-2501-1000, Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 1.
  2. Liquid connection: Connects nebulizer to port #3 of high flow FAST valve head with green, 1/4- 28 fitting. Like part # SC-0317-0250 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 2.
- xvi. Nut: (for flanged connections of 1.59mm (1/16") o.d. PFA tubing) Flanged, for 1/16" o.d. tubing, 1/4-28 threads. Use part # P-406x (pkg. of 10, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)) or equivalent. Use a Teflon-coated Viton o-ring with this nut instead of the stainless steel washer that comes with part # P-406x). # Spares = 10.
- xvii. Nut and Ferrule set, 1/8" Swagelok: Such as part # SS-200-NFSET (stainless steel) or part # B-200-NFSET (brass) (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. For part numbers listed here a quantity of 1 means 1 nut, 1 front ferrule, and 1 back ferrule. Spares = 20.
- xviii. Nut and Ferrule set, 1/4" Swagelok: Such as part # SS-400-NFSET (stainless steel) or part # B-400-NFSET (brass) (Georgia Valve and Fitting, Atlanta,

**Multi-Element in urine  
NHANES 2011-2012**

- GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. For part numbers listed here a quantity of 1 means 1 nut, 1 front ferrule, and 1 back ferrule. *Spares = 20.*
- xix. Oil, Welch Director Gold: For roughing pumps. Available direct from manufacturer as part # 8995G-15 (1 gallon, Welch Rietschle Thomas, Skokie, IL, [www.welchvacuum.com](http://www.welchvacuum.com)) or from various distributors. Equivalent oil may be substituted. *# Spares = 4.*
- xx. O-ring: (for sampler cone) PerkinElmer part # N8120511 (pkg. of 5, PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. *# Spares = 20 o-rings.*
- xxi. O-ring: (for skimmer cone) PerkinElmer part # N8120512 (pkg. of 5, PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. *# Spares = 20 o-rings.*
- xxii. O-ring: (for flanged connections of 1.59mm (1/16") o.d. PFA tubing) Teflon-coated Viton o-ring, i.d. = 1/16", thickness = 1/16", o.d. = 3/16". Such as part # V75-003 (O-rings West, Seattle, WA, [www.oringswest.com](http://www.oringswest.com)) or equivalent. *# Spares = 20.*
- xxiii. O-ring: (for injector support).
1. Internal o-rings: ID = 1/4", OD = 3/8", thickness = 1/16". Need 2 o-rings per injector support setup. PerkinElmer part # N8122008 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent (such as part # V75-010, O-rings West, Seattle, WA, [www.oringswest.com](http://www.oringswest.com)). *# Spares = 20.*
  2. External o-rings: ID = 3/8", OD = 1/2", thickness = 1/16". Need 2 o-rings for each injector support setup. PerkinElmer part # N8122009 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent (such as part # V75-012, O-rings West, Seattle, WA, [www.oringswest.com](http://www.oringswest.com)). *# Spares = 20.*
- xxiv. O-ring: (for inside spray chamber at nebulizer port) Such as part # 120-56 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). Additional o-rings can sometimes be obtained free of charge or at reduced price when acquired while purchasing spray chambers. *# Spares = 20.*
- xxv. O-ring: (for inside of torch mount): Part # WE017284 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Do not substitute. The PerkinElmer o-ring is specially metal impregnated to minimize RF leakage though the torch mount. *# Spares = 2.*
- xxvi. Photon Stop: PerkinElmer part # WE018278 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). *# Spares = 1.*
- xxvii. Plugs, Quick Change for Roughing Pump Oil: These plugs will only work on the Varian roughing pumps which come standard on ELAN DRC II ICPMS instruments. These plugs will not fit the Leybold pumps which come standard on the ELAN DRC Plus instruments. Part # W1011013 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). No spares typically needed.

**Multi-Element in urine  
NHANES 2011-2012**

xxviii. Probes

1. for ESI autosampler: Teflon, carbon fiber support, 0.8mm i.d., blue marker, 1/4-28 fittings. Like part number SC-5037-3751 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 2.
  2. for carrier solution of FAST sample introduction system: Teflon, carbon fiber support, 0.5mm i.d., orange marker, 1/4-28 fittings. Like part number SC-5037-3501 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)). # Spares = 2.
- xxix. RF coil. PerkinElmer part # WE02-1816 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 2.
- xxx. Screw, for Torch Mount: PerkinElmer part # WE011870. (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 3.
- xxxi. Spray chamber, quartz concentric: PerkinElmer part # WE025221 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. Available direct from manufacturer as part # 400-20 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or from various distributors. # Spares = 2.
- xxxii. Torch, quartz: PerkinElmer part # N812-2006 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. Available direct from manufacturer as part # 400-10 (Precision Glass Blowing, Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)) or various distributors. Damaged torches can often be repaired for substantially lower cost than purchasing a new one by companies such as Wilmad LabGlass (Buena, NJ, [www.wilmad-labglass.com](http://www.wilmad-labglass.com)) or Precision Glass Blowing (Centennial, CO, [www.precisionglassblowing.com](http://www.precisionglassblowing.com)). # New Spares = 2.
- xxxiii. Tubing and adapter, for SC autosampler rinse station drain: Tygon tubing and adapter to attach to back of SC autosampler for draining rinse station waste (like part # SC-0303-002, Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
- xxxiv. Tubing and adapters, for SC autosampler rinse station filling: Teflon tubing and adapters (to attach to back of SC autosampler for filling rinse stations and to attach to rinse containers). Like part # SC-0302-0500, Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
- xxxv. Tubing and nut, for FAST carrier solution: 0.5mm i.d. Teflon tubing (orange marker) with red 1/4-28 male nut. Connects to high flow FAST valve head, port #2. Like part # SC-0316-0500 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).
- xxxvi. Tubing, FAST vacuum: Vacuum line for SC-FAST high flow valve, connects to port #6, black nut for connection to valve head, natural brown color nut on other end for connection to SC autosampler vacuum port. Like part # SC-0321 (Elemental Scientific Inc., Omaha, NE., [www.elementalscientific.com](http://www.elementalscientific.com)).

**Multi-Element in urine  
NHANES 2011-2012**

xxxvii. Tubing, main argon delivery to instrument: I.D. = 1/8", O.D. = 1/4". Such as part # C-06500-02 (pkg. of 100ft, polypropylene, Fisher Scientific International, Hampton, NH, [www.fishersci.com](http://www.fishersci.com)) or equivalent. # Spares = 50ft.

xxxviii. Tubing, PFA: I.D. = 0.5mm, O.D. = 1.59mm (1/16"). Used to transfer liquid between rinse solution jug and peristaltic pump tubing

The Perfluoroalkoxy (PFA) copolymer is a form of Teflon®. Such as part # 1548 (20ft length, Upchurch Scientific, Oak Harbor, WA, [www.upchurch.com](http://www.upchurch.com)) or equivalent. # Spares = 20ft.

xxxix. Tubing, peristaltic, 0.045" i.d. (rinse station feed): Standard PVC, 2-stop (red / red) peristaltic pump tubing, i.d. = 0.045". PerkinElmer part # N0680375, (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 6 packs of 12 tubes.

xl. Tubing, peristaltic, 0.03" i.d. (carrier solution for ESI autosampler): use either

1. Standard PVC, 2-stop (black / black) peristaltic pump tubing, i.d. = 0.03". PerkinElmer part # 09908587 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 6 packs of 12 tubes.
2. Standard PVC, 3-stop (black/ black/black) peristaltic pump tubing, i.d. 0.76 mm. Spectron part # SC0056 (Spectron, Ventura, CA, [www.spectronus.com](http://www.spectronus.com)) or equivalent. #Spares = 6 packs of 12 tubes. Use this type of tubing with ESI DXi micro-peristaltic pump.

xli. Tubing, peristaltic, 0.125" i.d. (spray chamber drain): use either

1. Standard PVC, 2-stop (black / white) peristaltic pump tubing, i.d. = 0.125" or equivalent. PerkinElmer part # N812-2012 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)) or equivalent. # Spares = 6 packs of 12 tubes.
2. Standard Santoprene, 3-stop (grey/ grey/ grey) peristaltic pump tubing, i.d. 1.30 mm. Spectron part # SC0311 (Spectron, Ventura, CA, [www.spectronus.com](http://www.spectronus.com)) or equivalent. #Spares = 6 packs of 12 tubes. Use this type of tubing with ESI DXi micro-peristaltic pump.

xlii. Tubing, PVC, i.d. = 1/8", o.d. = 3/16". May be used to transfer liquid

1. between spray chamber waste port and peristaltic pump
2. between peristaltic pump and liquid waste jug

Like part # 14-169-7A (pkg. of 50ft, Fisher Scientific International, Hampton, NH, [www.fishersci.com](http://www.fishersci.com)) or equivalent. # Spares = 20ft.

xliii. Tubing, Stainless Steel, o.d. = 1/8", wall thickness = 0.028": Used to connect DRC gas cylinders to ELAN DRC gas ports. Also used to replace plastic tubing in the DRC gas path within the ELAN. Like part # SS-T2-S-028-20 (20ft, Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Spares = 20ft.

**Multi-Element in urine  
NHANES 2011-2012**

- xliv. Tubing, Teflon, corrugated, ¼" o.d.: Connects to the auxiliary and plasma gas side-arms of the torch. Part # WE015903 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). # Spares = 2.
  - xlv. Tubing, vinyl (argon delivery to nebulizer): Vinyl Tubing, 1/8" ID x 1/4" OD. Like part # EW-06405-02 (Cole Parmer, Vernon Hills, Illinois, [www.coleparmer.com](http://www.coleparmer.com)) or equivalent. Equivalent tubing material may be substituted. # Spares = 10ft.
  - xlvi. Union Elbow, PTFE ¼" Swagelok: Connects argon tubing to torch auxiliary gas sidearm on bayonet mount ELAN ICPMS instruments. Like part # T-400-9 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Spares = 2.
  - xlvii. Union Tee, PTFE, ¼" Swagelok: Connects argon tubing to torch plasma gas sidearm and holds igniter inside torch sidearm on bayonet mount ELAN ICPMS instruments. Like part # T-400-3 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)) or equivalent. Spares = 2.
- c. Sources for ICP-MS Maintenance Equipment & Supplies
- i. Anemometer: Like digital wind-vane anemometer (Model 840032, SPER Scientific LTD., Scottsdale, AZ, [www.sperscientific.com](http://www.sperscientific.com)) or equivalent. Use to verify adequate exhaust ventilation for ICP-MS (check with hoses fully disconnected).
  - ii. Pan, for changing roughing pump oil: Like part # 53216 (United States Plastics Corporation, Lima, OH, [www.usplastic.com](http://www.usplastic.com)) or equivalent. # On hand = 1.
  - iii. Container, to hold acid baths for glassware: Polypropylene or polyethylene containers with lids (must be large enough for torch, injector, or spray chamber submersion). May be purchased from laboratory or home kitchen supply companies. # On hand = 4.
  - iv. Cotton swabs: Any vendor. For cleaning of cones and glassware.
  - v. Cutter (for 1/8" o.d. metal tubing): Terry tool with 3 replacement wheels. Like part # TT-1008 (Chrom Tech, Inc., Saint Paul, MN, [www.chromtech.com](http://www.chromtech.com)) or equivalent.
  - vi. Getter Regeneration Kit: Part # WE023257 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)). Use this as needed (at least annually) to clean the getter in the pathway of channel A DRC gas.
  - vii. Magnifying glass: Any 10x + pocket loupe for inspection of cones and other ICP-MS parts. Plastic body is preferred for non-corrosion characteristics. Like part # 5BC-42813 (Lab Safety Supply, Janesville, WI, [www.labsafety.com](http://www.labsafety.com)).
  - viii. Screw Driver, for Ion Lens Removal: Screw driver with long, flexible shaft, and 2mm ball-Allen end for removal of ion lens screws (if lens is not in quick-release mount), part # W1010620. Extra 2mm bits, part # W1010598 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).

**Multi-Element in urine  
NHANES 2011-2012**

- ix. Toothbrush: Any vendor. For cleaning ion lens and glassware.
  - x. Ultrasonic bath: Like ULTRASONIK™ Benchtop Cleaners (NEYTECH, Bloomfield, CT, [www.neytech.com](http://www.neytech.com)) or equivalent.
- d. Sources for General Laboratory Consumable Supplies
- i. Bar Code Scanner: Like Code Reader 2.0 (Code Corporation, Draper, UT, [www.codecorp.com](http://www.codecorp.com)) or equivalent. For scanning sample IDs during analysis setup. Any bar code scanner capable of reading Code 128 encoding at a 3 mil label density can be substituted.
  - ii. Carboy (for preparation of urine quality control pool and waste jug for ICPMS sample introduction system): Polypropylene 10-L carboy (like catalog # 02-960-20C, Fisher Scientific, Pittsburgh, PA, [www.fischersci.com](http://www.fischersci.com)) or equivalent. Carboys with spouts are not advised due to potential for leaking.
  - iii. Containers for diluent and Rinse Solution: Two liter Teflon™ containers (like catalog# 02-923-30E, Fisher Scientific, Pittsburgh, PA., [www.fishersci.com](http://www.fishersci.com)) and 4L polypropylene jugs (like catalog# 02-960-10A, Fisher Scientific, Pittsburgh, PA, [www.fishersci.com](http://www.fishersci.com)) have both been used. Acid rinse before use. Equivalent containers may be substituted.
  - iv. Cups for urine collection: Like polypropylene 4.5 oz cup, catalog # 354013 (Becton Dickinson Labware, Franklin Lakes, NJ, [www.bd.com](http://www.bd.com)) or equivalent. Each lot of cups used must be lot screened (tested to be free of trace metal contamination). Clear plastics tend to have lowest trace metal contamination.
  - v. Gloves: Powder-free, low particulate nitrile (like Best Clean-DEX™ 100% nitrile gloves, any vendor). Equivalent nitrile or latex gloves may be substituted.
  - vi. Paper towels: For general lab use, any low-lint paper wipes such as KIMWIPES®EX-L Delicate Task Wipers or KAYDRY®EX-L Delicate Task Wipers (Kimberly-Clark Professional, Atlanta, GA, [www.kcprofessional.com](http://www.kcprofessional.com)). For sensitive applications in cleanrooms, a wipe designed for cleanroom use may be desired such as the Econowipe or Wetwipe (Liberty, East Berlin, CT, [www.liberty-ind.com](http://www.liberty-ind.com)).
  - vii. Pipette (for preparation of urine dilutions to be analyzed): Micromedic Digiflex-CX Automatic™ pipette equipped with 10.0-mL dispensing syringe, 2 mL sampling syringe, 0.75-mm tip, and foot pedal (Titertek, Huntsville, AL, <http://www.titertek.com/>).
  - viii. Pipettes (for preparation of intermediate stock working standards & other reagents): Like Brinkmann Research Pro Electronic pipettes (Brinkmann Instruments, Inc., Westbury, NY, <http://www.brinkmann.com/home/>). 5-100 µL (catalog #4860 000.070), 20-300 µL (catalog #4860 000.089), 50-1000 µL (catalog #4860 000.097), 100-5000 µL (catalog #4860 000.100). Note: pipette catalog numbers are without individual chargers. Can purchase individual chargers (pipette catalog numbers will differ) or a charging stand that will hold four pipettes (catalog #4860 000.860). When purchasing pipette tips (epTips),



purchase one or more boxes, then “reloads” for those boxes after that: 5-100  $\mu\text{L}$  (box catalog # 22 49 133-4, reload catalog # 22 49 153-9), 20-300  $\mu\text{L}$  (box catalog # 22 49 134-2, reload catalog # 22 49 154-7), 50-1000  $\mu\text{L}$  (box catalog # 22 49 135-1, reload catalog # 22 49 155-5), 100-5000  $\mu\text{L}$  (box catalog # 22 49 138-5, reload catalog # 22 49 198-9, bulk bag catalog # 22 49 208-0). Equivalent pipettes and tips can be substituted.

- ix.  Tubes for sample analysis (for autosampler): Like polypropylene 15-mL conical tubes, BD Falcon model #352097 (Becton Dickinson Labware, Franklin Lakes, NJ, [www.bd.com](http://www.bd.com)). Equivalent tubes may be substituted which are shown by lot screening to be free of trace metal contamination. Clear plastics tend to have lowest trace metal contamination. Blue colored caps have also been used successfully for this method.
  - x.  Tubes for storage of intermediate working stock standards: Like polypropylene 50-mL conical tubes, BD Falcon model #352098 (Becton Dickinson Labware, Franklin Lakes, NJ, [www.bd.com](http://www.bd.com)). For use in storage of intermediate working stock standards. Equivalent tubes may be substituted which are shown by lot screening to be free of trace metal contamination. Clear plastics tend to have lowest trace metal contamination. Blue colored caps have also been used successfully for this method.
  - xi.  Vortexer: Like MV-1 Mini Vortexer (VWR, West Chester, PA, [www.vwr.com](http://www.vwr.com)). Used for vortexing urine specimens before removing an aliquot for analysis. Equivalent item can be substituted.
  - xii.  Water purification system: Like NANOpure Diamond Ultrapure Water System (Barnstead International, Dubuque, Iowa, [www.barnstead.com](http://www.barnstead.com)). For ultra-pure water used in reagent and dilution preparations. An equivalent water purification unit capable of producing  $\geq 18$  Mega-ohm-cm water may be substituted.
- e.  Sources of Chemicals, Gases, and Regulators
- i.  Acid, Hydrochloric acid: Veritas™ double-distilled grade, 30-35% (GFS Chemicals Inc. Columbus, OH, [www.gfschemicals.com](http://www.gfschemicals.com)). This is referred to as “concentrated” hydrochloric acid in this method write-up. For use in preparation of intermediate working stock standards. An equivalent hydrochloric acid product may be substituted, but it must meet or exceed the purity specifications of this product for trace metals content.
  - ii.  Acid, Nitric acid: Veritas™ double-distilled grade, 68-70% (GFS Chemicals Inc. Columbus, OH, [www.gfschemicals.com](http://www.gfschemicals.com)). For use in diluent, rinse solution, intermediate working stock standards, and QC pool preparations. This is referred to as “concentrated” nitric acid in this method write-up. An equivalent nitric acid product may be substituted, but it must meet or exceed the purity specifications of this product for trace metals content.

**Multi-Element in urine  
NHANES 2011-2012**

- iii. Ethanol (EtOH): USP dehydrated 200 proof (Pharmco Products, Inc.) or equivalent.
- iv. Argon Gas (for plasma & nebulizer) and Regulator: High purity argon (>99.999% purity, Specialty Gases Southeast, Atlanta, GA, [www.sgsgas.com](http://www.sgsgas.com)) for torch and nebulizer. Minimum tank source is a dewar of liquid argon (180-250L). Bulk tank (1500<sup>+</sup>L is preferred).
  1. Regulator for argon (at dewar): Stainless steel, single stage, specially cleaned regulator with 3000 psig max inlet, 0-100 outlet pressure range, CGA 580 cylinder connector, and needle valve shutoff on delivery side terminating in a ¼" Swagelok connector. Part number "KPRCFRF415A2/AG10-AR1" (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)). An equivalent regulator from an alternate vendor may be substituted. # Spares = 1.
  2. Regulator for argon (between bulk tank and PerkinElmer filter regulator): Single Stage 316SS Regulator, with 0-300 psi Inlet Gauge, 0-200 psi Outlet Gauge, Outlet Spring Range, 0-250 psi, ¼" Swagelok Inlet Connection, ¼ turn Shut off Valve on Outlet with ¼" Swagelok Connection and Teflon Seals. Part number KPR1GRF412A20000-AR1 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)). An equivalent regulator from an alternate vendor may be substituted. # Spares = 1.
  3. Regulator for argon (PerkinElmer filter regulator on back of ELAN): Argon regulator filter kit. Catalog number N812-0508 (PerkinElmer, Shelton, CT, [www.perkinelmer.com](http://www.perkinelmer.com)).
- v. Argon / hydrogen: Argon (90%) / hydrogen (10%) for DRC channel A. Initial purity of argon = 99.9997<sup>+</sup>% ("Research grade 5.7"). Initial purity of hydrogen = 99.9999<sup>+</sup>% ("Research Grade 6.0"). Mixture is typically purchased in cylinder size 35 (6"x24") (Airgas South, Atlanta, GA, [www.airgas.com](http://www.airgas.com)).
  1. Regulator for argon / hydrogen: Stainless steel, two stage, specially cleaned regulator with 3000 psig max inlet, 0-25 outlet pressure range, CGA 350 cylinder connector, and needle valve shutoff on delivery side terminating in a ¼" Swagelok connector. Like part number KCYADPF412A2AD10 (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)). An equivalent regulator from an alternate vendor may be substituted. # Spares = 1.
  2. Flash Arrestor (Stainless steel): Like part # 6104 (Matheson Tri Gas, Montgomeryville, PA, [www.mathesontrigas.com](http://www.mathesontrigas.com)) or equivalent.
- vi. Disinfectant, for work surfaces: Bleach-rite spray (any distributor). On-site dilutions of bleach (1part bleach + 9 parts water) may be substituted, but must be re-made daily.
- vii. Oxygen: Oxygen ("Research Grade Research Grade 5.0", 99.9999% purity) for DRC channel B. Typically purchased in cylinder size 300 (9.5" x 54") (Airgas South, Atlanta, GA, [www.airgas.com](http://www.airgas.com)).

**Multi-Element in urine  
NHANES 2011-2012**

1. Regulator for oxygen: High purity brass body with monel trim, two stage regulator. Stainless steel is not used for this application due to safety concerns of working with oxygen at high pressure [12]. For one regulator, order the following parts, and ask that they be tested and assembled (Engineered Specialty Products, Kennesaw, GA, [www.espgauges.com](http://www.espgauges.com)).
  - a. Tescom part # 44-3410S24-555  
*Regulator body*: Brass bar stock, two stage, Monel trim, TFE seats, Elgiloy diaphragms, Cv=0.05, 3000 psig max inlet, 1-25 psig outlet range, 1/4 FNPT inlet / outlet / gauge ports, O<sub>2</sub> cleaned to ASTM G93 and CGA4.1.
  - b. Tescom part # 60500-3000N  
*Inlet pressure gauge*: 2" diameter, 0-3000 psig range , O<sub>2</sub> cleaned, 1/4" MNPT bottom, brass.
  - c. Tescom part # 60500-0015N  
*Delivery pressure gauge*: 2" diameter, 0-15 psig range , O<sub>2</sub> cleaned, 1/4" MNPT bottom, brass.
  - d. Tescom part # 63842-540-B  
*NPT to CGA Adaptor*: 1/4" NPT to CGA 540 adapter, brass.
  - e. Swagelok part # B-200-1-4:  
*Adaptor*: Brass male connector, 1/4" MNPT to 1/8" Swagelok (Georgia Valve and Fitting, Atlanta, GA, [www.swagelok.com](http://www.swagelok.com)).

An equivalent regulator from an alternate vendor may be substituted.

# Spares = 1.

2. Flash Arrestor (brass): Like part # 6103 (Matheson Tri Gas, Montgomeryville, PA, [www.mathesontrigas.com](http://www.mathesontrigas.com)) or equivalent.
- viii. Standard, Gallium: Like 1,000 mg/L, item # PLGA2-2Y. (SPEX Industries, Inc., Edison, NJ, [www.spexcsp.com](http://www.spexcsp.com)). Used as an internal standard in diluent. Any vendor whose standards are traceable to the National Institute for Standards and Technology may be substituted. The standard must have low trace metal contamination. Gallium is only used in the diluent for the measurement of arsenic (As).
  - ix. Standard, Iridium: Like 1,000 mg/L iridium, item # PLIR3-2Y (SPEX Industries, Inc., Edison, NJ, [www.spexcsp.com](http://www.spexcsp.com)). Used as an internal standard in diluent. Any vendor whose standards are traceable to the National Institute for Standards and Technology may be substituted. The standard must have low trace metal contamination.
  - x. Standard, Multi-element intermediate stock standard: Item number SM-2107-003 (High Purity Standards, Charleston, SC, <http://www.hps.net/>). This is a custom mix solution. This solution is diluted to prepare the intermediate stock working standards, which are in turn diluted to prepare the working calibrators.

This solution can be prepared in-house from NIST traceable single element stock solutions if necessary.

- xi. Standard, Rhodium: Like 1,000 mg/L, item # PLRH3-2Y. (SPEX Industries, Inc., Edison, NJ, [www.spexcsp.com](http://www.spexcsp.com)). Used as an internal standard in diluent. Any vendor whose standards are traceable to the National Institute for Standards and Technology may be substituted. The standard must have low trace metal contamination.
- xii. Standard, single element stock standards for preparation of urine quality control pools: National Institute of Standards and Technology (NIST) Standard Reference Materials (SRMs) 3103a (As), 3105a (Be), 3113 (Co), 3134 (Mo), 3108 (Cd), 3102a (Sb), 3111a (Cs), 3104a (Ba), 3163 (W), 3128 (Pb), 3140 (Pt), 3158 (Tl), and 3164 (U) (National Institute of Standards and Technology (NIST), Office of Standard Reference Materials, Gaithersburg, MD, [www.nist.gov](http://www.nist.gov)). Other sources of standards can be used if they are NIST traceable.
- xiii. Triton X-100™ surfactant: Like “Baker Analyzed” TritonX-100™ (J.T. Baker Chemical Co., [www.jtbaker.com](http://www.jtbaker.com)). Another source may be substituted, but it must be free of trace-metal contamination.

## **6) Preparation of Reagent and Materials.**

### **a. Internal Standard Intermediate Mixture**

- i. For DLS 3018 (Ir and Rh): Preparation of single intermediate solution containing all internal standards will simplify the addition of the internal standards into the final diluent solution. This solution can be purchased rather than prepared. To prepare 200 mL of the Intermediate internal standard solution:
  1. Partially fill a 200 mL acid-washed volumetric flask (PP, PMP, or Teflon™) with  $\geq 18$  Mega-ohm·cm water (approximately 100-150 mL).
  2. Carefully add 4 mL of double-distilled, concentrated nitric acid. Mix into solution.
  3. Add 0.8 mL of 10,000 ug/mL Rh standard. If initial Rh standard concentration is different, adjust volume proportionally.
  4. Add 0.8 mL of 10,000 ug/mL Ir standard. If initial Ir standard concentration is different, adjust volume proportionally.
  5. Fill to mark (200mL) and mix thoroughly.
  6. Label should include “Internal Standard Intermediate Mixture. 40 ug/mL Rh, Ir, and. 2% (v/v) HNO<sub>3</sub>”, “Store at room temperature”, preparation date, expiration date 1 year from preparation date, and preparer’s initials.
- ii. For DLS 3018A (Ga): Dilution of a Ga standard to be used as an internal standard in the final diluents solution. This solution can be purchased rather

than prepared. To prepare 200 mL of the Intermediate internal standard solution:

1. Partially fill a 200 mL acid-washed volumetric flask (PP, PMP, or Teflon™) with  $\geq 18$  Mega-ohm·cm water (approximately 100-150 mL).
2. Carefully add 4 mL of double-distilled, concentrated nitric acid. Mix into solution.
3. Add 0.8 mL of 10,000 ug/mL Ga standard. If initial Ga standard concentration is different, adjust volume proportionally.
4. Fill to mark (200mL) and mix thoroughly.
5. Label should include "Internal Standard Intermediate Mixture. 40 ug/mL Ga, 2% (v/v) HNO<sub>3</sub>", "Store at room temperature", preparation date, expiration date 1 year from preparation date, and preparer's initials.

b. Diluent and Carrier

i. Purpose: All samples (blanks, calibrators, QC, or patient samples) are combined with the diluent during the sample preparation step before analysis. This is where the internal standards are added which during the analysis will compensate for instrumental variations on the analyte signal. If using the FAST sample introduction system, the diluent is also used as the carrier solution.

ii. Preparation:

1. Diluent Preparation for DLS 3018 method (*not including arsenic*) – NO ETHANOL

- a. Contents: An aqueous solution of 10 microgram/L Rh, Ir, and in 2% (v/v) double-distilled nitric acid.
- b. Preparation (4L) & storage: This solution does not have to be made up in a volumetric flask. The important thing about the concentration of the internal standards is that they be consistent within all samples in one run. To prepare different volumes of diluent, add proportionally larger or smaller volumes of the solution constituents.
  - i. Acid-rinse a 4 L container (material may be polypropylene (PP), polymethylpentene (PMP), or Teflon™).
  - ii. Partially fill the 4 L container with  $\geq 18$  megaohm·cm water.
  - iii. Carefully add 80 mL double-distilled, concentrated nitric acid and mix.
  - iv. Add 1 mL of the 40 ug/mL Rh, Ir, internal standard solution. If other concentrations are used, the volume added should be adjusted proportionally.
  - v. Make up to volume (4 L) with  $\geq 18$  megaohm·cm water.

- vi. Store at room temperature and prepare as needed.
- vii. Label should include “10 µg/L Rh, and Ir,” “2% (v/v) HNO<sub>3</sub>”, “Store at room temperature”, preparation date, expiration date (1 year from prep), and preparer’s initials.

2. Diluent Preparation DLS 3018A (*urine arsenic*) method – CONTAINS ETHANOL

- a. Contents: An aqueous solution of 10 microgram/L Ga in 2% (v/v) double-distilled nitric acid and 1.5% (v/v) ethanol.
- b. Preparation (4L) & storage: This solution does not have to be made up in a volumetric flask. The important thing about the concentration of the internal standards is that they be consistent within all samples in one run. To prepare different volumes of diluent, add proportionally larger or smaller volumes of the solution constituents.
  - i. Acid-rinse a 4 L container (material may be polypropylene (PP), polymethylpentene (PMP), or Teflon™).
  - ii. Partially fill the 4 L container with  $\geq 18$  megaohm·cm water.
  - iii. Carefully add 80 mL double-distilled, concentrated nitric acid and mix.
  - iv. Carefully add 60 mL dehydrated 200 proof ethanol and mix.
  - v. Add 1 mL of the 40 ug/mL Ga internal standard solution. If other concentrations are used, the volume added should be adjusted proportionally.
  - vi. Make up to volume (4 L) with  $\geq 18$  megaohm·cm water.
  - vii. Store at room temperature and prepare as needed.
  - viii. Label should include “10 µg/L Ga”, “2% (v/v) HNO<sub>3</sub>”, “1.5% (v/v) Ethanol”, “Store at room temperature”, preparation date, expiration date (1 year from prep), and preparer’s initials.

c. Base Urine

- i. Purpose: This urine pool material will be mixed with the intermediate working calibrators just prior to analysis to matrix-match the calibration curve to the urine matrix of the unknown samples.
- ii. Contents: A mixture of multiple urine sources collected from anonymous donors are used to approximate an average urine matrix.
- iii. Preparation & Storage:
  - 1. Collect urine anonymously by placing screened containers and collection cups in the restrooms with a sign stating the reason the specimens are being collected, the name of the investigator to contact for additional information, and requesting that people provide a urine specimen

**Multi-Element in urine  
NHANES 2011-2012**

(complete details can be found in CDC protocol #3994, ProTrack # DLSITN0313).

2. Once the urine is collected from donors, it should be analyzed to ensure that concentrations of the analytes in this method are relatively low, so as to not interfere with the proper measurement of calibrators.
3. Once screened, mix the urine collections together in a larger container (i.e. acid washed polypropylene (PP), polymethylpentene (PMP), or Teflon™) and stir for 30+ minutes on a large stir plate (acid wash large Teflon™ stir bar before use).
4. For short term storage, store at 2-4°C. For long-term storage, dispense into smaller-volume tubes (i.e., 50-mL acid-washed or lot screened polypropylene tubes) and store at  $\leq -20^{\circ}\text{C}$ .
5. Labels on 50mL tubes should include “Base Urine for Multi-element Method”, “Store Long Term at  $\leq 20^{\circ}\text{C}$ ”, “Store Short Term at 2-4° C”, preparation date, expiration date 3 years from prep date, and preparer’s initials.

d. ICP-DRC-MS Rinse Solution

i. Purpose: Pump this solution into the sample introduction system between samples to prevent carry-over of the analytes of interest from one sample measurement to the next.

ii. Preparation:

1. Intermediate Triton X-100 Solution: To avoid the process of dissolving pure Triton X-100 on a daily basis, prepare an intermediate 2% Triton X-100™ / 5% (v/v) double-distilled, nitric-acid solution for daily use.

a. To prepare 2L of Intermediate Triton X-100 Solution:

- i. Partially fill a 2 L acid-washed bottle (PP, PMP, or Teflon™) with  $\geq 18$  Mega-ohm·cm water (approximately 1-1.5 L). Use of volumetric flask is not required.
  - ii. Add 20 mL of Triton X-100™ and stir until completely dissolved. Use a Teflon™ stir bar and stir plate if necessary (acid wash stir bar before use).
  - iii. Carefully add 100 mL of double-distilled, concentrated nitric acid.
  - iv. Fill to 2 L and stir thoroughly.
  - v. Label should include “2% Triton X-100™ / 5% (v/v) HNO<sub>3</sub>”, “Store at room temperature”, preparation date, expiration date 1 year from preparation date, and preparer’s initials.
2. Rinse Solution Preparation for DLS 3018 method (*not including arsenic*) – NO ETHANOL

**Multi-Element in urine  
NHANES 2011-2012**

- a. Contents: A 0.002% Triton X-100™, 5% (v/v) double-distilled nitric acid solution.
  - b. Preparation & Storage: To Prepare 4 L of the Final Rinse Solution,
    - i. Partially fill a 4 L acid-washed bottle (PP, PMP, or Teflon™ ) with  $\geq 18$  Mega-ohm·cm water (approximately 2-3 L). Use of volumetric flask is not required.
    - ii. Add 4 mL of the 2% Triton X-100™ / 5% (v/v) double-distilled, nitric-acid intermediate stock solution and mix well.
    - iii. Carefully add 200 mL of double distilled concentrated nitric acid and mix well.
    - iv. Fill to 4 L using  $\geq 18$  Megaohm·cm water.
    - v. Store at room temperature and prepare as needed. To prepare volumes other than specified here, add proportionally larger or smaller volumes of the solution constituents.
    - vi. Label should include “0.002% Triton X-100™ / 5% (v/v) HNO<sub>3</sub>”, “Store at room temperature”, preparation date, expiration date one year from preparation date, and preparer’s initials.
3. Rinse Solution Preparation for DLS 3018A (arsenic) method – INCLUDES ETHANOL
- a. Contents: A 0.002% Triton X-100™, 5% (v/v) double-distilled nitric acid solution and 1.5% (v/v) ethanol.
  - b. Preparation & Storage: To Prepare 4 L of the Final Rinse Solution,
    - i. Partially fill a 4 L acid-washed bottle (PP, PMP, or Teflon™ ) with  $\geq 18$  Mega-ohm·cm water (approximately 2-3 L). Use of volumetric flask is not required.
    - ii. Add 4 mL of the 2% Triton X-100™ / 5% (v/v) double-distilled, nitric-acid intermediate stock solution and mix well.
    - iii. Carefully add 200 mL of double distilled concentrated nitric acid and mix well.
    - iv. Carefully add 60 mL dehydrated 200 proof ethanol and mix well.
    - v. Fill to 4 L using  $\geq 18$  Megaohm·cm water.
    - vi. Store at room temperature and prepare as needed. To prepare volumes other than specified here, add proportionally larger or smaller volumes of the solution constituents.
    - vii. Label should include “0.002% Triton X-100™ / 5% (v/v) HNO<sub>3</sub>, 1.5% (v/v) ethanol”, “Store at room temperature”, preparation date, expiration date one year from preparation date, and preparer’s initials.



e. Standards and Calibrators

i. Multi-element Intermediate Stock Calibration Standard

1. Purpose: This master solution will be diluted to prepare five intermediate working calibrators.
2. Contents: An aqueous solution containing all 16 elements of interest (15 element panel analytes, arsenic, and elements for future R&D (see certificate of analysis), but does not include the internal standards). Matrix is 2% (v/v) HNO<sub>3</sub> and 1% (v/v) HCl with traces of HF in  $\geq 18$  Megaohm·cm water.
3. Preparation (Purchase) & Storage:
  - a. Purchasing from vendors: Either purchased as a NIST-traceable custom mixture, or prepared in-house.

Current vendor & preparation process: Currently purchased from High Purity Standards (Charleston, SC, part number SM-2107-003).
  - b. In-house Preparation: Standard may be made in the lab from NIST-traceable single element standards.
  - c. Storage: Store at room temperature. Label with additional information such as “store at room temperature”, date received, date opened, and initials of person to first open.

ii. Multi-element Intermediate Working Calibration Standards

1. Purpose: Use each day of analysis to prepare the final five working calibrators that will be placed on the autosampler.
2. Content: Five aqueous dilutions of the multi-element intermediate stock calibration standard solution in 2% (v/v) double-distilled nitric acid and 1% (v/v) hydrochloric acid.
3. Preparation & Storage: Different volumes may be prepared by adding proportionally larger or smaller volumes of solution constituents.
  - a. Cleaning flasks: Acid-rinse three 100-mL, one 200-mL, one 500-mL PP, and one 2 L PP (or PMP) volumetric flasks. Check their cleanliness by comparing the counts observed on the ICP-DRC-MS for 1% (v/v) HNO<sub>3</sub> before and after contact with the flasks. Mark each flask according to intended use. Dedicate to purpose.
  - b. HNO<sub>3</sub> & HCl Diluent Preparation: In the cleaned 2L flask, add 1-1.5L  $\geq 18$  Megaohm·cm water, 40 mL high purity concentrated HNO<sub>3</sub>, and 20 mL high purity concentrated HCl. Fill to the mark and mix

thoroughly. Use this diluent to fill the remaining flasks during preparation of the intermediate working calibration standards.

c. Dilutions & Storage:

- i. Partially fill the 100 mL, 200 mL, and 500 mL flasks with the HNO<sub>3</sub> & HCl diluent (50-75% full).
- ii. Using the volumes listed pipette the appropriate volume of the multi-element intermediate stock calibration standard solution into each of the five volumetric flasks. Dilute each to the volumetric mark with the HNO<sub>3</sub> & HCl diluent using a pipette for the final drops. Mix each solution thoroughly.
- iii. Once mixed, transfer to acid-cleaned, labeled, 50-mL containers (PP, PMP, or Teflon™) for storage. Labels should include information such as “Multi-element Urine Working Calibrators”, “2% (v/v) HNO<sub>3</sub>, 1% (v/v) HCl”, date of preparation, expiration date (1 year from date of preparation), “store at room temperature”, initials of preparer, and concentrations for each element.

iii. Working Multi-element Calibrators

1. Purpose: The working multi-element calibrators will be analyzed in each run to provide a signal-to-concentration response curve for each analyte in the method. The concentration of an analyte in a patient urine sample dilution is determined by comparing the observed signal from the dilution of the patient urine sample to the response curve from the working multi-element calibrators.
2. Content: Dilutions (1:100) of the corresponding five intermediate working calibration standards. .
3. Preparation & Use: Made immediately prior to analysis when the intermediate working calibration standards are mixed with base urine (Section 7.b) and diluent (Section 7.a) using a Digiflex automatic pipetter.

iv. Multi-element Intermediate Stock Calibration Verification Standard

1. Purpose: This is the master solution from which all working calibration verification standards will be prepared. It will be diluted to prepare intermediate working calibration verification standards which are in turn diluted and used to verify the accuracy of instrument response to analyte concentrations greater than the calibration range. This stock solution contains all elements needed for both the arsenic and the 12 element panel.
2. Contents: For long shelf life, these four aqueous solutions have different matrices which are optimized to the elements in each (this was recommended for the calibration verification stock standard solutions

because the elemental concentrations were very high compared to the concentrations in the calibration stock standard solution.

- a. Solution A: HNO<sub>3</sub> (10%), HF (0.5%)
- b. Solution B: HCl (10%), trace HNO<sub>3</sub>
- c. Solution C: HCl (1%)
- d. Solution D: HCl (2%)

3. Preparation (Purchase) & Storage:

- a. Purchasing from vendors: The intermediate stock calibration verification standard solutions may be purchased as custom mixtures from any vendor which prepares multi-element solutions that are traceable to the National Institute for Standards and Technology (NIST) for their accuracy.

Current vendor & preparation process: Currently it is purchased from High Purity Standards (Charleston, SC, part number SM-2107-012, solutions A, B, C, and D).

- b. In-house Preparation: If outside laboratories were not available to prepare the intermediate stock calibration standard solution, it is also possible to make it in the laboratory from single element standards which are NIST traceable.
- c. Storage: Due to the uranium content, and in keeping with the guidance of the CDC radiation safety manual [11], the intermediate stock standards must be kept in a lockbox. Store the solutions at room temperature. Label these bottles from HPS with additional information such as "store at room temperature", date received, date opened, and initials of person to first open.

v. Multi-element Intermediate Working Calibration Verification Standards

1. Purpose: Verification of accuracy of instrument response to analyte concentrations greater than the calibration range
2. Content: The intermediate working calibration verification standard solutions used in this method are aqueous dilutions of the multi-element intermediate stock calibration verification standard solution in 2% (v/v) double-distilled nitric acid and 1% (v/v) hydrochloric acid containing all 13 elements of interest (does not include the internal standards).
  - a. Preparation & Storage: Prepare the Intermediate Calibration Verification Standards for analysis just as a Intermediate Working calibrators are prepared.

vi. Internal Quality Control Materials ("Bench" QC)

1. Purpose: Internal (or "bench") quality control (QC) materials are used to evaluate the accuracy and precision of the analysis process, and to determine if the analytical system is "in control" (is producing results that

are acceptably accurate and precise). They are included in the beginning and at the end of each analytical run.

2. Content: The internal (or “bench”) quality control (QC) materials used in this method are pooled human urine, acidified to 1-2% (v/v) HNO<sub>3</sub>, and may have been spiked to reach a desired concentration. The analyte concentrations in the “low QC” are in the low-normal concentration range. The analyte concentrations in the “high QC” are in the high-normal concentration range.
3. Preparation & Storage: Quality control materials can be either prepared by and purchased from an external laboratory or prepared within the CDC laboratories. Quality control must always be traceable to the National Institute for Standards and Technology (NIST). The CDC laboratory currently prepares its own bench QC materials using the following procedures:
  - a. Collection of urine: Collect urine anonymously by placing screened containers and / or collection cups in the restrooms with a sign stating the reason the specimens are being collected, the name of the investigator to contact for additional information, and requesting that people provide a urine specimen (complete details can be found in CDC protocol #3994, ProTrack # DLSITN0313). Volume of urine to collect is dependent on the desired pool size. This write-up will assume a 10-L pool size for both the low and high bench QC.
  - b. Screening Urine: Screen collected samples for metal content before mixing together to make 2 separate base urine pools (for preparing the low and high bench QC materials). Samples can be screened individually or after combining several together (reduces number of analyses).
    - i. Keep urine refrigerated whenever possible to minimize microbial growth.
    - ii. Because this is only a quick screen of the metal content, the number of replicates in the urine method can be reduced to one in order to reduce analysis time.
    - iii. Analyte concentrations in the the final urine pool to be spiked for the low bench QC pool should be in the low-normal population range. Analyte concentrations in the final urine pool to be spiked for the high bench QC pool should be less than some pre-selected target concentration values in the high normal population range. See the Second National Report on Human Exposure to Environmental Chemicals for estimations of the normal population ranges for metals (<http://www.cdc.gov/exposurereport/>).
  - c. Combining Collected Urine: Be attentive not to combine only diluted matrix urine samples into the low pool and only concentrated matrix

urine samples into the high pool. The goal is for combining samples is to approach an 'average' matrix for each pool.

- i. Graduate four acid-washed 10-L carboys (PP or PMP) in 0.5 L increments (two will be used for decanting into).
  - ii. Combine collected urine samples into two separate acid-washed 10-L carboys (PP or PMP), according to their concentrations, for the low bench and high bench QC pools.
  - iii. Mix each urine pool using large acid washed, Teflon™ coated stir bars and large stir plates. Keep urine refrigerated whenever possible.
  - iv. Acidify each urine pool to 1% (v/v) HNO<sub>3</sub> by adding the appropriate volume of double distilled HNO<sub>3</sub>. Stir for 30+ min on large stir plates.
- d. Settling out of solids:
- i. Refrigerate the urine (no stirring) for 1-3 days to allow for settling out of solids.
  - ii. For each urine pool, decant the urine into another of the acid-washed 10-L carboys to remove the urine from the solids settled out on the bottom of the carboy.
  - iii. Repeat steps (i) and (ii) until minimal solids are left at the bottom of the carboy after sitting overnight.
- e. Spiking of urine
- i. Analyze a sample of each urine pool. Record these results for future recovery calculations.
  - ii. Use these results to determine target analyte concentrations possible for the pools
  - iii. Calculate the volume of single element standards needed to spike each pool to the desired concentrations.
  - iv. While stirring the pools on large stir plates, spike each pool with calculated volumes of single element standards (all spiking standards used must be traceable to NIST).
  - v. Continue to stir pools for 30+ minutes after spiking, then reanalyze.
  - vi. Repeat steps 4 and 5 until all analytes reach target concentrations keeping track of the total volume of spiking solution added to each urine pool.
- f. Dispensing and Storage of urine
- i. Container Types: Dispense urine into lot screened containers (i.e. – 5 or 15 mL polypropylene tubes). If possible, prepare tubes of QC which have only enough volume for one typical run + 1 repeat

analysis. This allows for one vial of QC to be used per day of analysis, reducing chances of contamination of QC materials due to multi-day use.

- ii. Labels: Place labels on vials after dispensing and capping if the vials are originally bagged separately from the caps. This minimizes the chance for contamination during the process. Include at least the name of QC pool (text and bar code), date of preparation, and a vial number on the labels.
- iii. Dispensing: Dispensing can be accomplished most easily using a Digiflex automatic pipetter in continuous cycling dispense mode. This process should be done in a clean environment (i.e. a class 100 cleanroom area or hood).
  1. Allow urine pool to reach room temperature before dispensing (to prevent temperature gradients possibly causing concentration gradients across the large number of vials being dispensed and to prevent condensation problems during labeling of vials). This may require leaving the carboy of urine at room temperature overnight before dispensing.
  2. Replace the tubing attached to the dispensing syringe (left when looking at front of Digiflex) with a length of clean Teflon™ tubing long enough to reach into the bottom of the 10L carboy while it is sitting on the stir plate.
  3. Check cleanliness of Digiflex before use by analyzing 1-2% (v/v) HNO<sub>3</sub> which has been flushed through the Digiflex with a portion of the same solution which has not been through the Digiflex.
  4. Approximately one hour before dispensing begins,
    - a. With the large stir plate close to the left side of the Digiflex, begin stirring the urine pool to be dispensed.
    - b. Also during this time, flush the Digiflex with urine from the pool to be dispensed. Place the ends of the tubing attached to both the sample and dispensing syringes into the carboy of urine so that urine won't be used up during this process. Be sure to secure both ends of tubing in the carboy with Parafilm so they will not come out during the flushing process.
  5. After dispensing the urine into the vials, cap the vials and label them. Placing labels on vials after capping minimizes the chance for contamination during the process.
- iv. Homogeneity Testing: After dispensing, check homogeneity of analyte concentrations in pool aliquots by analysis of every Nth sample dispensed (where N ~ 20 - 50 depending on the pool size).

Sample more heavily from the beginning and the ending portions of the tubes dispensed (these are the regions where most homogeneity problems occur). Keep samples pulled for homogeneity analysis in the sequence that they were dispensed for the purpose of looking for trends in concentrations. Once dispensed and homogeneity has been shown to be good throughout the tubes of a pool, store tubes at  $\leq -20^{\circ}\text{C}$  and pull tubes out as needed for analysis.

- v. Storage: Urine pools should be stored long term at  $\leq -20^{\circ}\text{C}$ . Short term storage (several days) at refrigerator temperature ( $\sim 2-4^{\circ}\text{C}$ ).

## **7) Analytical Instrumentation & Parameters**

(see Section 6 for details on hardware used, including sources)

### a. Instrumentation & Equipment Setup:

- i. ICP-DRC-MS: Inductively Coupled Plasma Dynamic Reaction Cell Mass Spectrometer ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> or ELAN<sup>®</sup> DRC II.

#### 1. Modifications made to ICP-DRC-MS

- a. Plastic tubing for between mass flow controllers and dynamic reaction cell have been replaced with stainless steel. Stainless steel tubing is preferred between the reaction gas cylinder / regulator and the back of the ICP-DRC-MS instrument.
- b. A second mass flow controller has been added (channel B) for use with oxygen.
- c. Standard built-in peristaltic pump replaced by DXi-FAST micro-peristaltic pump / FAST actuator unit.

#### 2. Configuration of tubing for liquid handling:

##### a. Sample introduction system:

- i. SC-FAST valve setup: Valve connections must match this description for urine total arsenic analysis.

- 1. Port 1: sample loop (white nut), 3.0 mL loop for DLS 3018 and 0.5 mL loop for DLS 3018A. See "Loop, for FAST valve" in section 5.b. for details.
- 2. Port 2: 0.5 mm ID probe (red nut) for carrier solution. See "Probes" in section 5.b. for details.
- 3. Port 3: nebulizer line (green nut) for transfer of liquid to nebulizer. See "Nebulizers" in section 5.b. for details.

4. Port 4: sample loop (white nut). See “Loop, for FAST valve” in section 5.b. for details.
  5. Port 5: 0.8 mm ID probe (blue nut) for diluted samples. See “Probes, for ESI autosampler” in section 5.b. for details.
  6. Port 6: vacuum line (black nut). See “Tubing, FAST vacuum” in section 5.b. for details.
- b. Tubing connection between autosampler rinse station and rinse solution reservoir: Tubing of different inner diameters can be obtained from Elemental Scientific, their distributors, or custom built in the lab to optimize the rinse station fill rate between samples. Rinse station should not go empty at any point.
  - c. Tubing for autosampler rinse station waste removal: Use minimum drain tubing to make this connection. If this tube is too long, the rinse station will not drain properly.
  - d. Rinse solution jug: Leave one of the caps on the top of the rinse jug loose to allow air venting into the jug as liquid is removed. Otherwise the jug will collapse on itself as the liquid is removed and a vacuum is created inside. Use secondary containment tray and label appropriately (see solution preparation instructions).
  - e. Waste solution jug: Use secondary containment tray and label appropriately (see solution preparation instructions).
  - f. Configuration of tubing and probe for carrier solution: Use a ‘peristaltic to Teflon tubing adapter’ (see consumables descriptions in section 5.b) for connections.
3. Cones used  
Nickel and platinum cones have been used.
  4. Gases & Regulators setup:
    - a. Argon: Argon stored as liquid in a dewar (180-250L) or bulk tank. Gaseous argon used for plasma and nebulizer.
      - i. Regulator for argon source (if a dewar): Keep the inlet pressure (headspace pressure of liquid argon dewar) above 100 psi. Set delivery pressure to 60-100psi to allow for pressure drop across tubing that stretches to the instrument. See Section 6.f. for part numbers and details.
      - ii. Step down regulator (if source of argon is a bulk tank): Place this single stage regulator in the lab so that incoming argon pressure



- can be monitored and adjusted. Set delivery pressure to 60-100 psig. See Section 6.f. for part numbers and details.
- iii. Regulator at ICP-DRC-MS: Single stage “argon regulator filter kit” supplied with the ICP-DRC-MS. Set the delivery pressure depending on the specifications for that model of ELAN ICP-DRC-MS instrument (see the PE Hardware Manual). This will be  $52 \pm 1$  psi for instruments having a 0-60psi gauge and  $60 \pm 1$  for instruments having a 0-100psi gauge. See Section 6.f. for regulator part numbers.
- b. Argon (90%) / hydrogen (10%) mixture for DRC channel A. NOTE: Only for arsenic analysis.
    - i. Regulator for Ar / H<sub>2</sub> gas mixture: Set delivery pressure to 5-7 psig. See Section 6.f. for part numbers and details.
    - ii. Flash arrestor: Stainless steel flash arrestor is used on outlet side of regulator. See Section 6.f. for part numbers and details.
  - c. Oxygen (99.999<sup>±</sup>%) gas for DRC channel B. NOTE: Only for cadmium and manganese analysis
    - i. Regulator for O<sub>2</sub> gas: Set the delivery pressure = 5-7 psig. See Section 6.f. for part numbers and details.
    - ii. Flash arrestor: Brass flash arrestor is used on outlet side of regulator. See Section 6.f. for part numbers and details.
5. Chiller / Heat Exchanger: Refrigerated chiller (for ELAN<sup>®</sup> 6100 DRC<sup>Plus</sup> instruments) or heat exchanger (for ELAN<sup>®</sup> DRC II instruments). For refrigerated chiller, set temperature control to 18°C.
- ii. Computer: Dell Optiplex GX150, GX270, or GX280 have all been used. Processors used have included Pentium III (1 GHz) through Pentium IV (2.8 GHz). Recommend 512Mb - 1Gb RAM. External hard disk drive for nightly backups of data connects via USB port. Software used includes Windows XP Professional, service pack 2 and ELAN v3.3.
  - iii. Autosampler: ESI SC4 autosampler with (DLS 3018 or DLS 3018A) FAST sample introduction. Rack calibration, tubing ID for rinse supply, additional rinse time, probe movement speeds, and probe depth is optimized per autosampler.

## **8) Method Procedures**

a. Quality Control: Quality control procedures implemented in this method are defined by the Division Procedures and Practices Guidelines and is subjected to the complete analytical process. The data from these materials are then used to estimate methodological imprecision and to assess the magnitude of any time-associated trends. The concentrations of these materials should cover the expected concentration range of the analytes for the method. Before QC materials can be used to judge patient analytical runs, acceptable QC concentration limits must be calculated from the concentration results observed in at least 20 characterization runs. During the 20 characterization runs, previously characterized QCs or pools with target values assigned by outside laboratories should be included to evaluate the analysis. The process of limits calculation is performed using the laboratory database and the SAS division QC characterization program.

i. Types of Quality Control:

1. “Bench QC”: The bench QC pools used in this method comprise two levels of concentration spanning the “low-normal” and “high-normal” ranges of the analyte of interest. The intent of bench QC is for the analyst to evaluate the performance of the analytical system on the day of analysis. The analyst inserts both the “low” and the “high” bench QC specimens two times in each analytical run (a set of consecutive assays performed without interruption) so that judgments may be made on the day of analysis. The first analysis of the two bench QC pools is done after the calibration standards are analyzed but before any patient samples are analyzed (so that judgments on the calibration curves may be made before analysis of patient samples). The second analysis of the two bench QC pools is done at the end of the run (approximately 20 patient samples total). If more patient samples are analyzed on the same calibration curve after the second run of the bench QC, both the low-normal and high-normal bench QC must be reanalyzed before and after the additional samples.
2. External Reference Materials: Materials produced by laboratories outside of the CDC which have assigned target concentrations can be helpful in verifying method performance. Some examples include Standard Reference Materials (SRM) from the National Institute of Standards and Technology (NIST) (i.e. SRM 2668 Level 1 & 2668 Level 2) and samples from previous challenges of proficiency testing programs (e.g. Centre de Toxicologie du Quebec (CTQ)). However, only the results for the bench QC materials are used to determine if the run results can be used.

ii. Calibration Verification:

- a. Bi-annual tests as defined in the DLS Policy and Procedures manual: CLIA requires the verification of accuracy of instrument response to analyte concentration be completed at least every 6 months. NIST traceable calibrators are analyzed in each run to define this response up to the concentration of the highest calibrator in the run. To verify accuracy of

instrument response at concentrations higher than the highest calibrator in each run, analyze a NIST traceable standard with very high concentrations at least every 6 months. Prepare the Calibration Verification Standard for analysis just as a working calibrator is prepared. Use the "Urine Blank" as the blank when it is analyzed. If the observed concentrations for the Calibration Verification Standard are not within 10% of the target value the lab supervisor should be notified and the issue should be investigated. Do not substitute external reference materials (i.e. biological samples from a PT program) for the Calibration Verification Standard when performing this. Solutions needed for the Calibration Verification checks can be purchased from standards vendors (i.e. SPEX, High Purity Standards, etc . . .) or prepared in-house from NIST traceable single element standards. Always verify that normal background levels have been re-achieved through adequate rinse time following analysis of elevated standards for calibration verification.

- b. As-needed confirmations (per supervisor discretion): When a sample result is greater than the highest calibrator in the run, the supervisor may request that the result be confirmed in an analysis run which includes a standard or external reference material with equivalent (within 10%) or greater concentration than the sample. In order to avoid needless contamination of the instrument with high concentrations of analytes, the analyst should use the lowest appropriate calibration verification solution concentrations to meet the need.

For *infrequent* verification needs, the calibration verification stock solutions can be used to prepare verification standards to appropriate concentrations. This will, however, introduce elevated concentrations of all elements in the method to the sample introduction system. Frequent measurement of these very high concentrations can result in high background levels in the instrument which are difficult to rinse out and which may limit the ability to measure low concentrations.

For frequent verification needs (i.e. when certain studies have many elevated results on particular elements), use NIST-traceable single element stock standards to prepare single element verification standards. This will limit the exposure of the instrument to elevated concentrations of only the elements needing verification.

An external reference material (i.e. historical proficiency testing sample) can be substituted in place of the Calibration Verification Standard sample in these situations IF

- i. The target value has been assigned by an external source (i.e. NIST, or the proficiency testing program).
- ii. The concentration of the external reference material is within 10% or is higher than the concentration of the material you need it to *confirm*.

- iii. *There is confidence that there is no contamination of previously used external reference material.*
- iv. A note to file is made that this was done.
- v. If the observed concentrations are not within 10% of the target value the lab supervisor should be notified and the issue should be investigated.

Always verify that normal background levels have been re-achieved through adequate rinse time following analysis of calibration verification standard 3 or higher (including samples with those concentrations).

b. Daily Analysis of Samples

i. Preparation of the Analytical Equipment

*For further details on any part of this description, see the ITN Daily Startup SOP for ELAN ICPMS instruments.*

1. Power on the computer, printer, and autosampler, and log into the operating system.
2. Peristaltic pump: Set up the peristaltic pump tubing with proper tension for the sample rinse station.
3. Software: Starting the ESI software before starting the ELAN software may improve stability of software.
4. Daily Pre-Ignition Maintenance Checks: Perform daily maintenance checks as described in the IRAT Daily Startup SOP for ELAN instruments (i.e., Ar supply pressure, interface components cleanliness and positioning, interface pump oil condition, vacuum pressure, etc.). Make appropriate notes in the Daily Maintenance Checklist and Instrument Log Book.
5. Start the Plasma: In the INSTRUMENT window of the software (or on the front of the ELAN), press the "Start" button to ignite the plasma.
6. Aspirate rinse (multi-element method) or carrier (arsenic method) solution: Send Probe to Rinse Station (multi-element method) or manually place carrier probe into carrier solution (arsenic method).
7. Start the peristaltic pump: Start the peristaltic pump by pressing the appropriate arrow in the DEVICES window (make sure that the rotational direction is correct for the way the tubing is set up in the peristaltic pump). Set the pump speed to a slow speed in the DEVICES window during warm-up.
8. Warm-up time: Allow approximately 30 to 45 minutes warm-up time for the ICP-DRC-MS after igniting the plasma. This warm-up time is for the RF generator. There will be another "Stability time" for the DRC later in this procedure.
9. Optimizations and Daily Performance Check: After this warm-up time, perform a daily performance check and any optimizations necessary (as

described in the ITN Daily Startup SOP for ELANs). Include Be (m/z 9) in the daily performance check. Fill in the Daily Maintenance Checklist according to the optimization procedures performed.

- a. Magnesium ( $^{24}\text{Mg}$ ) may have high RSDs due to the use of Triton-X100 in the rinse solution. Avoid this problem by either temporarily using non-Triton-containing rinse solution during the daily check, or repeating the daily check multiple times in succession with no rinse time between.
- i. Saving the Files: Save new tuning (mass calibration) parameters to the file "default.tun." Save new optimization parameters (i.e., detector voltages, autolens values, nebulizer gas flow rate) to the file "default.dac."

10. Software setup for Analysis:

- a. Workspace (files & folders): Click on "Open Workspace" from the "File" menu. Select the workspace file "CDC\_urine multi-element.wrk" (or one customized for user preferences). Select "Review Files" from the "File" menu. Verify & set up the correct files and data directories for your analysis .
- b. Samples / Batch Window: Update the window to reflect the current sample set. The only fields which need to be filled in include the autosampler location, sample identification (id), measurement action, method, sample flush time, sample flush speed, read delay time, read delay & analysis speed, wash time, wash speed. Use a bar code scanner to input data whenever possible. Save the Sample window file and re-use it on other days by simply replacing the sample IDs for the patient samples.

1. DRC Stability Time: Best analyte-to-internal standard ratio stability is obtained after 1-1.5 hrs of analysis of urine samples using the DRC method. Analyze enough "dummy" urine sample dilutions prior to any DRC analysis run to fill 1-1.5 hours of analysis time (not necessary if analyzing only a subgroup of the method containing no DRC analytes). If analyzing the full set of method analytes, 10 samples will be sufficient.

2. Urine vs. Aqueous Method Files:

- a. The difference: There are two method files for this one method. It is necessary to use both to accomplish each run because the current PerkinElmer software will not allow for more than one blank per method file. The ONLY DIFFERENCE between these two files is on the Sampling tab where one lists the autosampler positions of the urine blank and urine calibrators (the "urblk" method file) and the

other lists the autosampler position of the aqueous blank (the “aqblk” method file).

b. Use: The ONLY TIME when it matters which of these files is used is when the measurement action *includes* “Run blank” or “Run standards”. When the measurement action is only ‘run sample’, it does not matter whether the “urblk” or “aqblk” method file is used. Analysts typically follow the pattern below, however, for the sake of consistency and as a reminder of which blank must be used for which type of sample.

i. The “urblk” method file: Use to analyze the initial urine blank (blank for the calibration curve), the urine calibrators, and the urine blank checks (urblkchk1 & urblkchk2) at the very beginning of the run. The urine blank method defines the autosampler location of the urine blank and the urine calibration standards.

ii. The “aqblk” method file must be used to analyze all QC materials and patient samples. The aqueous blank method defines the aqueous blank in autosampler location.

3. Notation of Dilutions: To designate an extra dilution of a sample, edit the sample ID to reflect the level of dilution being performed (i.e., A 1:2 dilution of sample 1 would be reflected in the sample ID “sample 1 (2x dilution)”. This sample ID will be edited during the data-import process to the database so that it is recognized as the appropriate sample.

ii. Preparation of Samples for Analysis

1. Thaw the frozen urine specimens; allow them to reach ambient temperature.
2. DRC stability “dummy urine matrix”. Prepare 50+mL of standard 2 or standard 3 to be analyzed for 1-1.5 hr before the beginning of the run. This can be prepared using 50mL polypropylene tubes or a wide-mouth bottle (which can be put on the autosampler in place of one of the tube trays).
3. Set up a series of 15-mL polypropylene tubes corresponding to the number of blanks, standards, QCs, and patient samples to be analyzed.
4. Prepare the following solutions in the 15-mL falcon tubes using the Micromedic Digiflex™.

a. *Aqueous Blank:* Prepare two aqueous blanks consisting of 1,000  $\mu\text{L}$  of  $\geq 18$  Mega-ohm $\cdot\text{cm}$  water and 9,000  $\mu\text{L}$  of diluent. One will be the actual aqueous blank and the other will be a backup (“Aqueous Blank Check”) in case the original aqueous blank gets contaminated.

- b. *Urine Blank*: Prepare two urine blank dilutions consisting of 900  $\mu\text{L}$  of base urine (same material used to prepare the urine calibration standards), 100  $\mu\text{L}$  of  $\geq 18$  Mega-ohm-cm water, and 9,000  $\mu\text{L}$  of diluent. One of these urine blanks will be the blank for the calibration standards; the other will be analyzed twice after standard 5 as UrBlkChk1 and UrBlkChk2, respectively. Results from the UrBlkChks will be used to determine the method limit of detection.
  - c. *Calibrators*: Prepare the working calibration standards as 100  $\mu\text{L}$  of the appropriate aqueous intermediate working calibration standard, 900  $\mu\text{L}$  of base urine, and 9,000  $\mu\text{L}$  of diluent.
  - d. *Patient & QC Samples*: Before taking an aliquot for analysis, mix the sample so that no particulates remain on the bottom of the tube. Prepare urine sample dilutions for DLS 3018 as 4,500  $\mu\text{L}$  of diluent and 500  $\mu\text{L}$  of the urine sample; for DLS 3018A 2,250  $\mu\text{L}$  of diluent and 250  $\mu\text{L}$  of the urine sample.
  - e. Cap all of the blanks, standards, and samples and mix them well. Uncap them and place them in the autosampler of the ELAN<sup>®</sup> ICPMS in the order that was entered in the Samples / Batch window of the ELAN software.
- iii. Specimen Storage and Handling During Testing: Specimens may be left at room temperature during analysis in case confirmation analyses must be made. Take stringent precautions to avoid external contamination by the metals to be determined. Specimens may be stored short term at refrigerated temperatures, but should be stored long term (>4 weeks) at  $\leq -20$  °C.
- iv. Starting the Analysis: To begin analysis, highlight (click and drag with the mouse) the table rows of the samples that should be included in the run, and then click on “Analyze Batch.”
- v. Monitoring the Analysis: Initiate work in a timely manner so that the run may be monitored. Make every effort to complete analysis within the work day so that the entire run can be monitored. If it is not possible to complete the analysis by the end of the work day, the run may be left to complete itself unattended as long as appropriate planning is made for either overnight operation or Auto Stop (see below).

Monitor the analysis for the following:

1. *DRC stability (analyte / internal standard ratio stability)*

After the analysis of the DRC stability “dummy” samples, these results can be reviewed to determine if sufficient stability of the analyte-to-internal standard ratio was reached before beginning analysis. Importing data into an MS Excel template file is useful to visualizing magnitude of drift.

2. *Proper operation of the instrument.*

3. *Contaminated blanks.*

4. *Linear calibration curves.*

- a. Typical correlation coefficients will be 0.999 to 1.000.
- b. The ELAN software generates a “simple linear” calibration curve (using a least squares calculation) for each of the 13 elements in this method. The curves are generated using the results from analysis of the urine blank and the 5 external urine calibrators whose concentrations are defined in the Calibration tab of the Method file. Specifically, the software plots the “net intensity” (y-axis) versus the analyte concentration (x-axis). The “net intensity” is the blank subtracted **ratio** of the measured intensity for the analyte to the measured intensity of the associated internal standard and is calculated as follows:

$$\text{net intensity} = \frac{\text{Analyte Meas Intensity}_{\text{sample}}}{\text{Internal Std Meas. Intensity}_{\text{sample}}} - \frac{\text{Analyte Meas Intensity}_{\text{Blank}}}{\text{Internal Std Meas Intensity}_{\text{Blank}}}$$

5. *Bench QC results within the acceptable limits.*

If an analyte result for the beginning QC material(s) falls outside of the 99% limits, then the following steps are recommended:

- a. If a particular calibration standard is obviously in error, remake a new dilution at the Digiflex of that working calibrator, reanalyze it, and reprocess the sample analyses using this new result as part of the calibration curve.
- b. Prepare a fresh dilution of the failing QC material and reanalyze it.
- c. Prepare fresh dilutions at the Digiflex of all of the calibration standards (working urine multi-element standards) and reanalyze the entire calibration curve using the freshly prepared standards.

If these three steps do not result in correction of the out-of-control values for QC materials, consult the supervisor for other appropriate corrective actions. Do not report analytical results for runs that are not in statistical control.

6. *Good precision among replicates.*

7. *Consistent measured intensities of the internal standards.*

Some sample-to-sample variations are to be expected. However the intensities should be within a few percent of one another, and should fluctuate around an average value (not drift continuously in one direction).

8. *Elevated patient results.* After any sample having a concentration greater than the third calibration verification standard, verify the instrument background levels have returned to normal before proceeding in analysis to the next sample, adding additional rinsing time if necessary. If this cannot be done in real-time, the sample analyzed immediately after the



elevated sample should be repeated for confirmation (within 10% of the original result) prior to reporting. Report the first analytically verified result.

vi. Records of Results: Run results will be documented daily in both electronic and paper form.

1. Electronic Records:

a. Transfer of Results to the Laboratory Information System / Database: Transfer data electronically between computers or software to reduce errors. When keyboard entry must be used, proofread transcribed data after entry.

b. Long-Term Storage of ELAN software files: Files used and produced by the ELAN software in analyzing samples will be backed up long term on compact disk and kept a minimum of three years.

2. Paper Records: The paper copy of the results from the run should be put into the study folder(s) and should include

a. A summary of the calibration curve statistics.

b. A printout of analysis of each measurement made during the run.

c. Optional, but helpful, is a printout of the DRC stability check measurements in graphical form.

d. On the front sheet of the printed records, write the following

i. Analyst initials

ii. Instrument ID

iii. Date of Analysis

iv. Run # for the day on this instrument

v. Study ID and Group Number

vi. Database batch ID (Not known until the run is imported into the database)

vii. Transfer of Results to the Laboratory Database: Every analytical run performed for the analysis of patient samples should be entered into the laboratory results database unless the run is not useable for obvious reasons (i.e. the run is stopped for some reason before ending QC is analyzed, no internal standard spiked into the diluent, etc. . . ).

1. Data Export Process (from ELAN® software to .TXT file): If the data file was not created during the initial analysis, reprocess the data of interest either with “original conditions” option, or by loading the files and folders used during the analysis. In the ELAN® ICP-DRC-MS software, select “Review Files” from the “File” menu. From this window, you must open the files and directories that were used when collecting the data of the run that

you wish to export. (If the analysis has just ended, all of these files and directories will still be open.) NOTE: A second copy of the ELAN® software can be run as an Edit/Reprocess copy without affecting an ongoing analysis by the first copy of the software running in Windows. After you open the relevant files, go to the “Report” page in the METHOD window. Deselect the box that prints a paper copy of data and select the box that sends data to a file. Select the “Report Options Template” named “CDC\_Database Output.rop” and type in a report filename using a format such as “2005-0714a\_group55.txt” to designate data from analysis of group 55 from July 14, 2005, run #1. Under “Report Format”, choose the “Use Separator” option, and under the “File Write” section choose “Append.” Finally, reprocess the data of interest. (See PerkinElmer ELAN® ICPMS Software Manual.) Make sure you apply the aqueous blank to all sample and quality control material analyses.

2. Data Import Process (from .TXT file to Microsoft Access™ database):

- a. Move the .TXT file to the appropriate subdirectory on the network drive where exported data are stored. Directories for data storage are named according to instrument \ year \ month\, such as I:\Instruments\ELANDRC2A\2005\07\.
- b. Using the ITN Database Frontends, import the instrument file into the database. On the GoTo window, click on “Add Sample Results to Database”, then “Import Instrument Data File”.
- c. Enter the appropriate information to identify the instrument, assay, analysis date & time, run number, analyst, calibrator lot number and prep date used (use the “IS Lot Number” field) and study. If other than default values for Method LOD, High Calibrator, Rep Delta Limit, and units were used in the run, document what was used by clicking on the “View/Set Batch Parameters” button, changing the appropriate values, and then clicking “Back”.
- d. Press the “Import” button, then browse to the correct network folder to select the file which contains the results from the run. Select the file and click “OK”.
- e. In the “Import Instrument Results” table, pressing the “Find X’s” button will show only those samples whose sample ID is not recognized as a valid QC pool ID or sample ID for this study. (Sample IDs are set up when the study is logged into the database.) Corrections to sample IDs and dilution factors can be made in this table (e.g., correction of transcription errors and adjustment for level of dilution). If samples were diluted for analysis, both the sample ID and the dilution factor need to be edited in this table before the values are transferred to the database (the Replace command under the Edit window is helpful in this case). When corrections to sample IDs are made, press the “Check IDs” button to re-evaluate the sample IDs. Any sample or analyte row marked “Not Recognized” will not be transferred to the

database when the "Transfer" button is pressed. Once transferred into the database, the data should be evaluated for QC pass / fail, then set with the with the appropriate settings for QC accept / reject, final value status, and comment(s). See the database programmers for more detail on working in the database.

viii. Analyst Evaluation of Run Results:

1. Bench Quality Control: After completing a run, and importing the results into the database, export the QC results to the SAS program where the run will be judged to be in or out of control. The QC limits are based on the average and standard deviation of the beginning and ending analyses of each of the bench QC pools, so it will not be possible to know if the run is *officially* accepted or rejected until it is completed.
  - a. Quality Control Rules: The SAS program applies the division QC rules to the data as follows:
    - i. If both QC run means (low & high bench QC) are within 2Sm limits and individual results are within 2Si limits, then accept the run.
    - ii. If 1 of the 2 QC run means is outside a 2Sm limit - reject run if:
      1. Extreme Outlier – Run mean is beyond the characterization mean +/- 4Sm
      2. 1 3S Rule - Run mean is outside a 3Sm limit
      3. 2 2S Rule - Both run means are outside the same 2Sm limit
      4. 10 X-bar Rule – Current and previous 9 run means are on same side of the characterization mean
    - iii. If one of the 4 QC individual results is outside a 2Si limit - reject run if:
      1. R 4S Rule – Within-run ranges for all pools in the same run exceed 4Sw (i.e., 95% range limit)

Note: Since runs have multiple results per pool for 2 pools, the R 4S rule is applied within runs only.

Abbreviations:

Si = Standard deviation of individual results (the limits are not shown on the chart unless run results are actually single measurements).

Sm = Standard deviation of the run means (the limits are shown on the chart).

Sw = Within-run standard deviation (the limits are not shown on the chart).

- b. Implications of QC Failures: If the division SAS program declares the run out of control” for any analyte, use the following to determine the implications on usability of the data from the run.
  - i. For 1 or 2 analytes: ONLY the analytes which were “out of control” are invalid for reporting from the run. Set all run results for those 1 or 2 analytes as “QC Rejected” in the database.
  - ii. For 3 or more analytes: All results, regardless of analyte, are invalid for reporting from the run. Set all run results for all analytes as “QC Rejected” in the database. Note in the batch comment field why all results were marked QC rejected.

2. Patient Results:

a. Elevated Results:

i. Boundaries Requiring Confirmatory Measurement:

1. Results Greater than the First Upper Boundary (1UB): Concentrations observed greater than the “first upper boundary” (defined in the laboratory database as the “1UB”) should be confirmed by repeat analysis of a new sample preparation. The concentration assigned to the 1UB for an element is determined by study protocol. Report the original result, as long as the confirmation is within 10% of the original. Continue repeat analysis until a concentration can be confirmed.
  2. Results Greater Than Highest Calibrator: When a sample result is greater than the highest calibrator in the run, the supervisor may request that the result be confirmed in an analysis run which includes a standard or external reference material with equivalent (within 10%) or greater concentration than the sample.
  3. Results Greater Than Calibration Verification Tested: Perform an extra dilution up to the maximum extra dilution levels.
  4. Uranium Isotope Ratio Measurement for Elevated Uranium Concentrations: A uranium 235/238 isotope ratio analysis is performed for all urine uranium samples where the urine total uranium concentration is greater than the 2UB boundary.
- ii. Inadequate Precision in Confirmation of a Measurement: If a sample is reanalyzed to obtain a confirmation of an initially elevated result, the confirmation should be within 10% of the original result.
  - iii. Analyst Reporting of Elevated Results: Concentrations observed greater than the “second upper boundary” (defined in the laboratory database as the “2UB”) should be reported to the QC reviewer as an “elevated result”. The concentration assigned to

the 2UB for an element is determined by study protocol. The analyst should report any patient results confirmed to be greater than the second upper boundary to the QC reviewer as an “elevated result”. There is no routine notification for elevated levels for the metals determined in this method. The protocol for supervisors reporting elevated results to medical personnel is defined according to the study protocol.

- b. Inadequate Precision Within One Measurement: If the range of the three replicate readings (maximum replicate concentration value - minimum replicate concentration value) for a single sample analysis is greater than the criteria **and** the range of the three replicate readings is greater than 10% of the observed concentration, do not use the measurement for reporting. Repeat the analysis of the sample.
- ix. Submitting Final Work for Review: Once results have been imported, reviewed, and set as final in the database by the analyst,
  1. Submit an email to the QC reviewer informing them of the readiness of the data for final review. The email should include
    - a. Instrument ID, run Date, run number, study ID, group ID.
    - b. Any bench QC failures (include reasons if known).
    - c. Any patient sample result greater than the 2UB boundaries.
    - d. Anything out of the ordinary about this analytical work which could have a bearing on the availability (i.e. insufficient sample to analyze), accuracy, or precision of the results.
  2. Include all items called for by the study folder cover sheet in the study folder (i.e. printouts from the ICP-MS, bench QC evaluation) together in the study folder before submitting the folder for review when analysis is complete.
- x. Overnight Operation or Using Auto Stop: Make every effort to complete analysis within the work day so that the entire run can be monitored. If it is not possible to complete the analysis by the end of the work day, the run may be left to complete itself unattended as long as appropriate planning is made for either overnight operation or Auto Stop.
  1. *24 hrs / day operation in DRC mode*:
    - a. To reduce startup time in the mornings, the analyst is encouraged to operate the ELAN in DRC mode 24hrs/day during the work week. This eliminates the need for daily 45 minute RF generator warm-up, and possibly the need for DRC stability time (if the DRC gas is not off for extended periods of time before analysis). To maintain the instrument in DRC mode when not analyzing patient samples, setup multiple sample rows in the Samples / Batch window with autosampler position n zero (rinse station of autosampler) and wash time of 1800s (30 minutes). Repeat this sample row enough times to keep the

instrument in analysis mode overnight (1 sample with 15 minute wash will take ~ 25 minutes).

2. *AutoStop*: If 24 hrs / day ELAN operation is not desired, the instrument can shut the plasma off unattended after analysis. Setup this as follows:
  - a. On the “Auto Start / Stop” tab of the Instrument window, enable the Auto Stop feature.
  - b. Press the “Change” button within the Auto Stop box and set the Delayed shutdown time to 5 minutes. This will rinse the sample introduction system of urine matrix before turning off the plasma.
  - c. It will be necessary to replace the sample peristaltic pump tubing the next day since it will have been clamped shut overnight.
- c. Equipment Maintenance: Analysts are expected to follow a 4-day analysis / 1-day maintenance schedule in the laboratory.
  - i. ICPMS Maintenance: On the maintenance day, perform all maintenance per the Inorganic Toxicology and Nutrition Branch ELAN ICP-MS Weekly Maintenance SOP. All equipment maintenance should be documented in the instrument logbook.
  - ii. Data Backup: Data on the ELAN computer will be backed up via two backup routines.
    1. Daily Backups to External Hard Drive: Automatic backups of the “elandata” directory and all subdirectories should be programmed to occur each night onto an external hard disk.
    2. Weekly Backup to CD: Backup all files in the active “elandata” directory and all subdirectories onto one recordable compact disc during the weekly maintenance SOP. When the active “elandata” directory on the ICP-DRC-MS computer hard drive becomes too large to fit onto a single recordable compact disc, the oldest data can be removed from the computer to make it easier to backup the entire directory weekly. This can usually be done annually.
      - a. Backup the oldest data on the hard drive to two duplicate compact disks and verify that the files on the CD are readable
      - b. Label them with the name of the instrument, the date range of the data, the current date, your name, and “Copy 1 of 2” or “Copy 2 of 2”
      - c. After verifying that the CDs are readable, the oldest, backed up data can be deleted from the ICP-MS computer hard drive.
      - d. It is best to not store duplicate copies in the same location.

## **9) Interpretation of the Results**

- a. Reportable Range: Urine multi-element values are reportable in the range between the method LOD and the highest concentration verified accurate by bi-annual calibration verification tests. For example, if a urine cadmium value is

less than the method LOD of 0.042, report it as  $< 0.042 \mu\text{g/L}$ ). Above the highest concentration verified, extra dilutions are made of the urine sample to bring the concentration within the verified range. If extra dilution has been necessary, the reported value will exceed the upper end of the reportable range.

- b. Reference Ranges (Normal Values): In this method the 95% reference ranges for these elements in urine fall within the range of the calibrators.
- c. Action Levels: Due to the uncertainty of the health implications of elevated concentrations of many of the elements determined with this method, there is no routine notification for elevated levels of every analyte determined with this method. The present NRC standard for workplace removal is  $15 \mu\text{g/L}$  of U in urine [13]. Other action levels for reporting to supervising physicians are determined on a study-by-study basis.

## **10) Method Calculations**

- a. Method Limit of Detection (LODs): The method detection limits for elements in blood specimens are defined as 3 times  $s_0$ , where  $s_0$  is the estimate of the standard deviation at zero analyte concentration.  $S_0$  is taken as the y-intercept of a linear or 2<sup>nd</sup> order polynomial regression of standard deviation versus concentration (4 concentration levels of the analytes in blood each measured 60 times across at least a 2-month timeframe). Method LODs are re-evaluated periodically.
- b. Method Limit of Quantitation (LOQ): The Division of Laboratory Sciences does not currently utilize limits of quantitation in regards to reporting limits [10].
- c. QC Limits: Quality control limits are calculated based on concentration results obtained in at least 20 separate runs. It is preferable to perform separate analyses on separate days and using multiple calibrator lot numbers, instruments, and analysts to best mimic real-life variability. The statistical calculations are performed using the SAS program developed for the Division of Laboratory Sciences (DLS\_QC\_compute\_char\_stats.sas).

## **11) Alternate Methods for Performing Test and Storing Specimens If Test System Fails:**

If the analytical system fails, the analysis may be setup on other ELAN DRC instruments in the laboratory. If no other instrument is available, store the specimens at  $\sim 4^\circ\text{C}$  until the analytical system can be restored to functionality. If interruption longer than 4 weeks is anticipated, then store urine specimens at  $\leq -20^\circ\text{C}$ .


## **12) Summary Statistics and Qc Graphs**

See following pages.

**Summary Statistics for Antimony, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUSH7	88	26AUG11	04JUN12	0.5961	0.0204	3.4
LU10709_UMP3_a_URCUSL7	88	26AUG11	04JUN12	0.1594	0.0140	8.8
HU10710_UMP3_b_URCUSH8	79	06FEB12	11FEB13	0.5873	0.0185	3.2
LU10709_UMP3_b_URCUSL8	79	06FEB12	11FEB13	0.1522	0.0177	11.7

**2011-2012 Antimony, urine (ug/L) Quality Control**

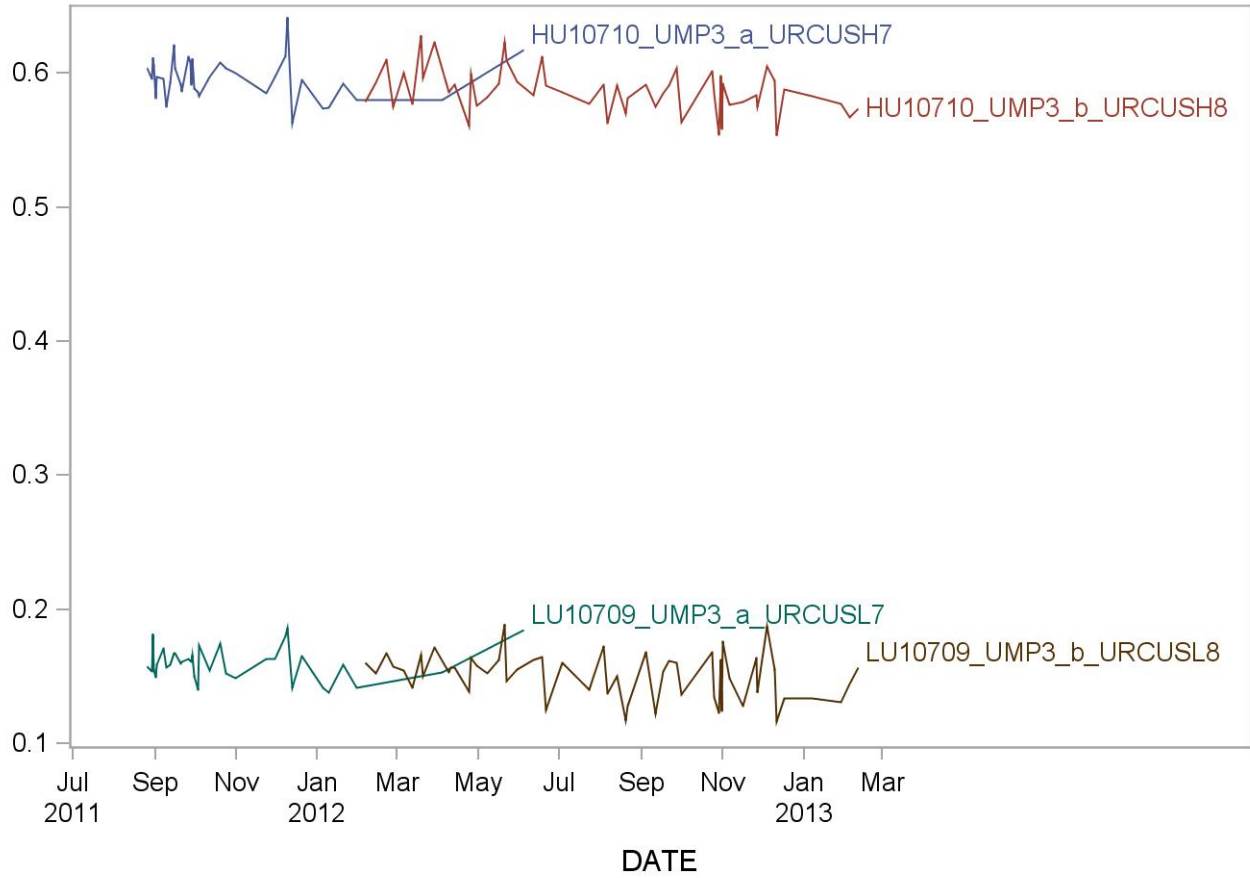




## Summary Statistics for Arsenic, urine (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU-04311_UMP_e_URCTAH6	24	15FEB11	24AUG11	55.6702	1.4081	2.5
LU-04310_UMP_e_URCTAL6	24	15FEB11	24AUG11	3.8173	0.2459	6.4
HU-04311_UMP3_e_URCTAH7	13	31AUG11	12DEC11	55.0215	0.8981	1.6
LU-04310_UMP3_e_URCTAL7	13	31AUG11	12DEC11	3.7238	0.1706	4.6
HU10710_UMP3_a_URCTAH8	6	18OCT11	30JAN12	47.2933	0.6325	1.3
LU10709_UMP3_a_URCTAL8	6	18OCT11	30JAN12	10.8800	0.2009	1.8
HU10710_UMP3_b_URCTAH9	47	06FEB12	29JAN13	47.3248	0.6566	1.4
LU10709_UMP3_b_URCTAL9	47	06FEB12	29JAN13	10.8502	0.2464	2.3


### 2011-2012 Arsenic, urine (ug/L) Quality Control

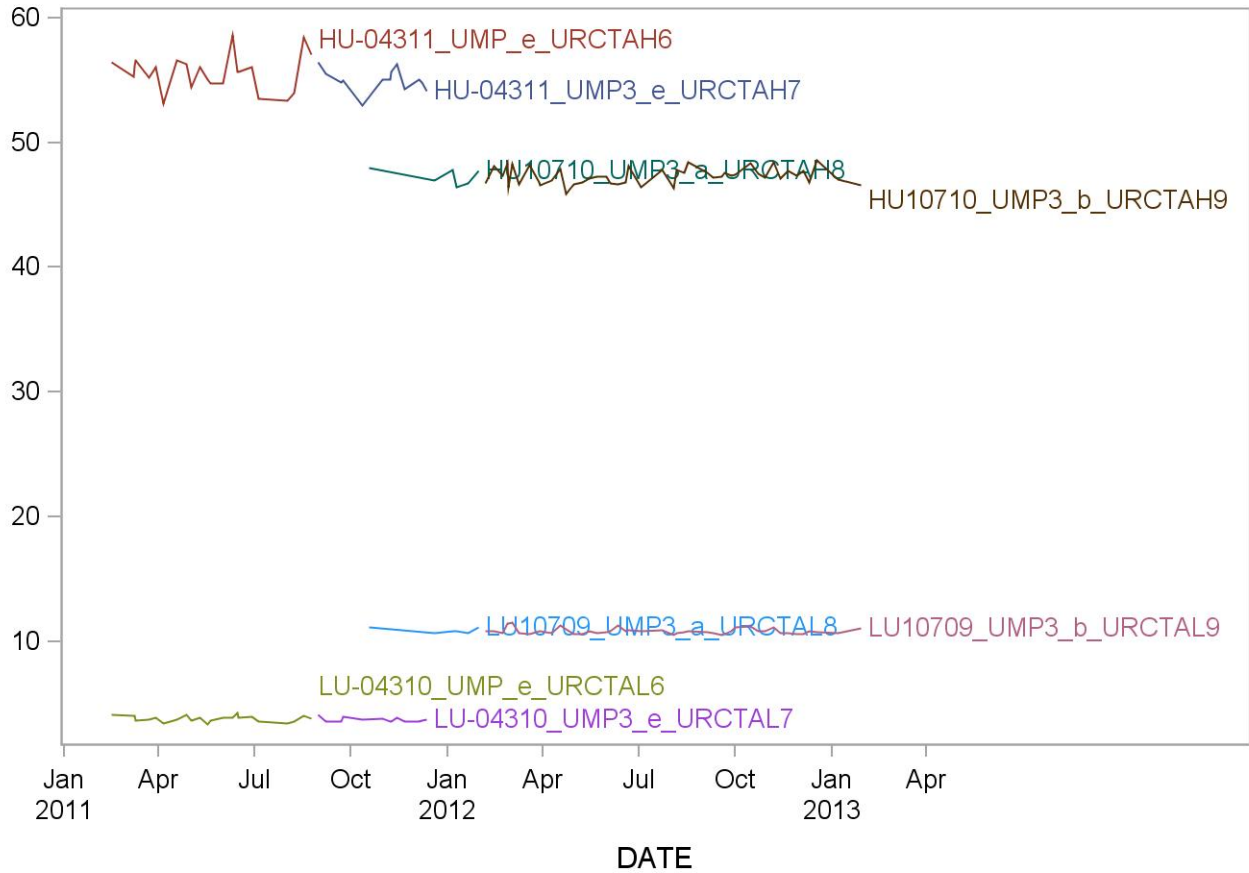

**Summary Statistics for Barium, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUAH7	91	26AUG11	04JUN12	6.362	0.172	2.7
LU10709_UMP3_a_URCUAL7	91	26AUG11	04JUN12	1.222	0.049	4.0
HU10710_UMP3_b_URCUAH8	77	06FEB12	11FEB13	6.425	0.131	2.0
LU10709_UMP3_b_URCUAL8	77	06FEB12	11FEB13	1.229	0.036	2.9


**2011-2012 Barium, urine (ug/L) Quality Control**

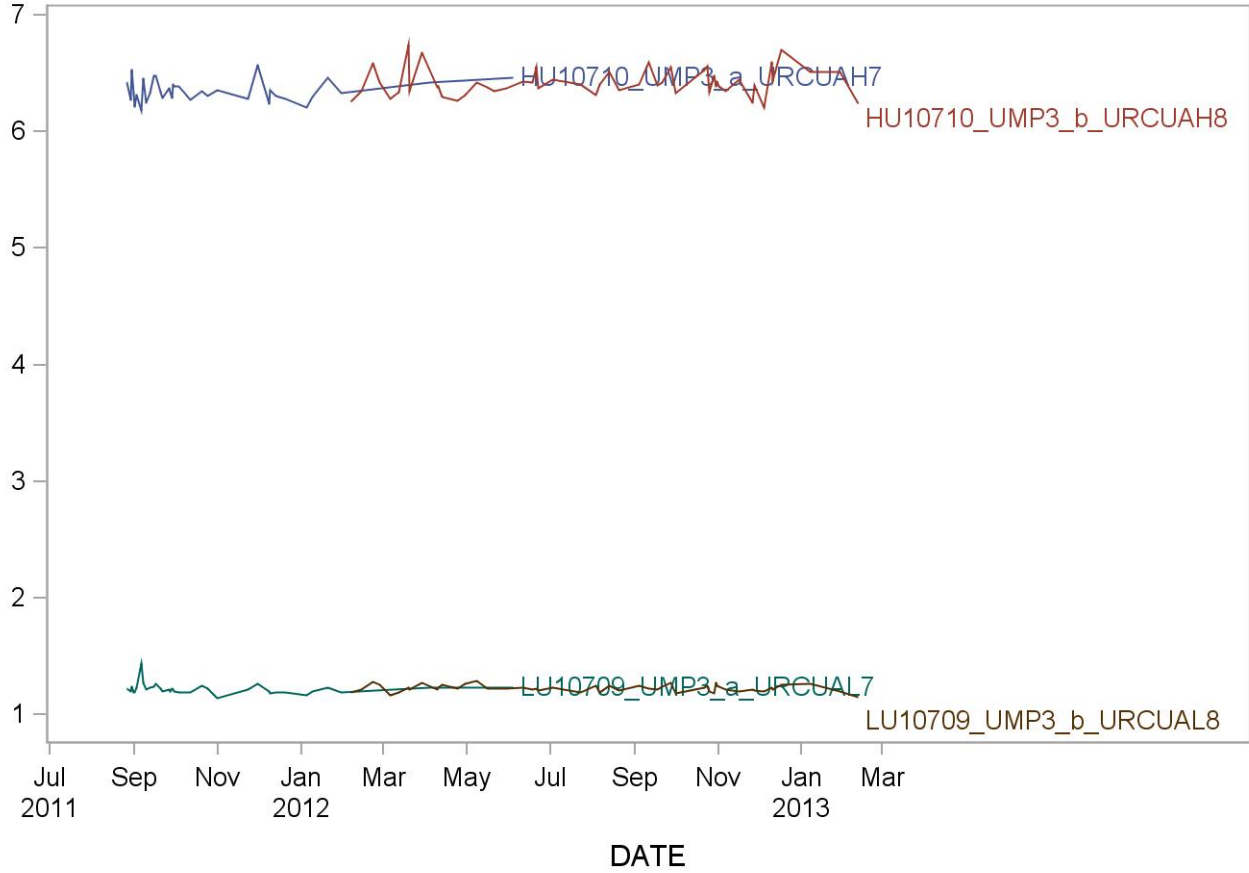

### Summary Statistics for Cadmium, urine (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUDH7	95	26AUG11	04JUN12	1.5247	0.0456	3.0
LU10709_UMP3_a_URCUDL7	95	26AUG11	04JUN12	0.1456	0.0162	11.1
HU10710_UMP3_b_URCUDH8	85	06FEB12	11FEB13	1.5359	0.0287	1.9
LU10709_UMP3_b_URCUDL8	85	06FEB12	11FEB13	0.1467	0.0119	8.1


### 2011-2012 Cadmium, urine (ug/L) Quality Control

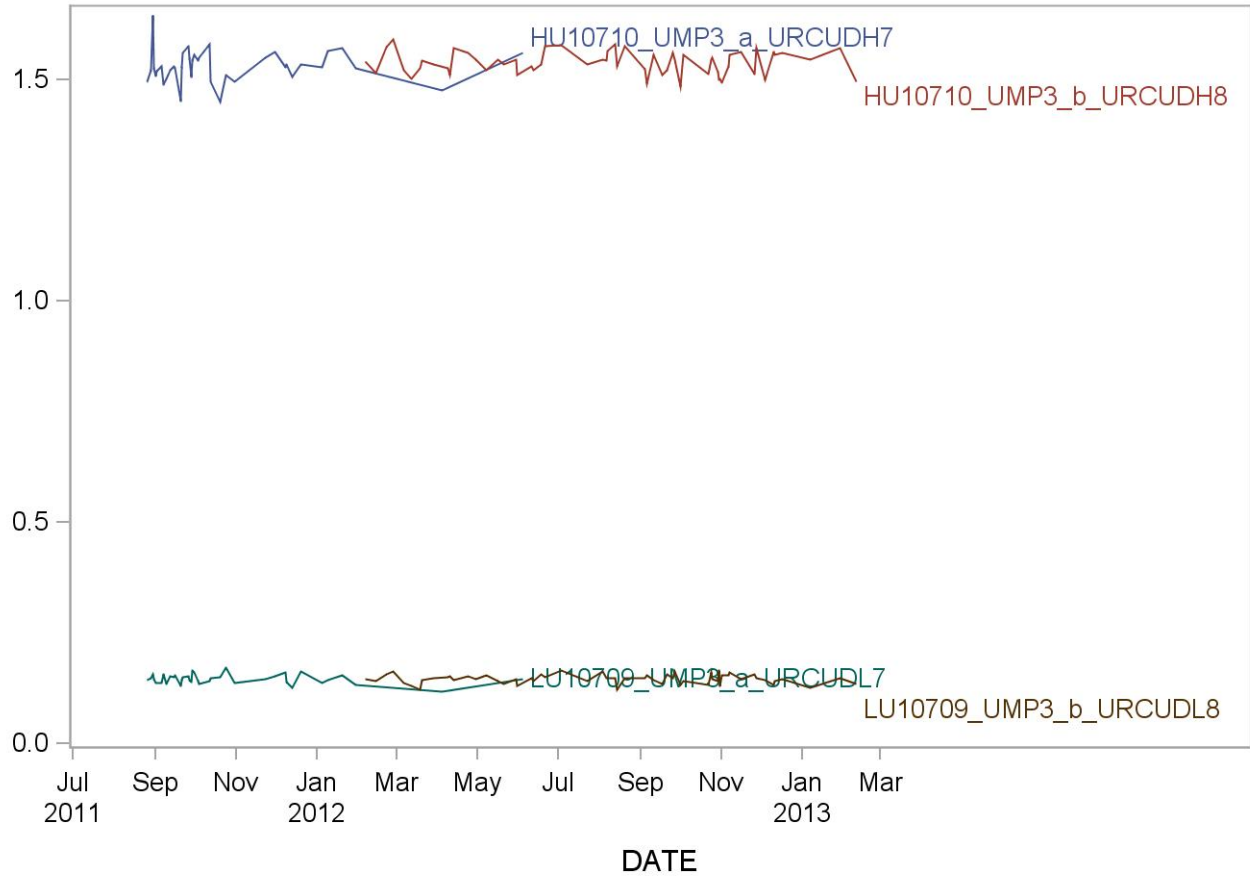

### Summary Statistics for Cesium, urine (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUCH7	88	26AUG11	04JUN12	12.1176	0.2895	2.4
LU10709_UMP3_a_URCUCL7	88	26AUG11	04JUN12	4.3911	0.1154	2.6
HU10710_UMP3_b_URCUCH8	83	06FEB12	12FEB13	12.1892	0.2645	2.2
LU10709_UMP3_b_URCUCL8	83	06FEB12	12FEB13	4.4287	0.0921	2.1


### 2011-2012 Cesium, urine (ug/L) Quality Control

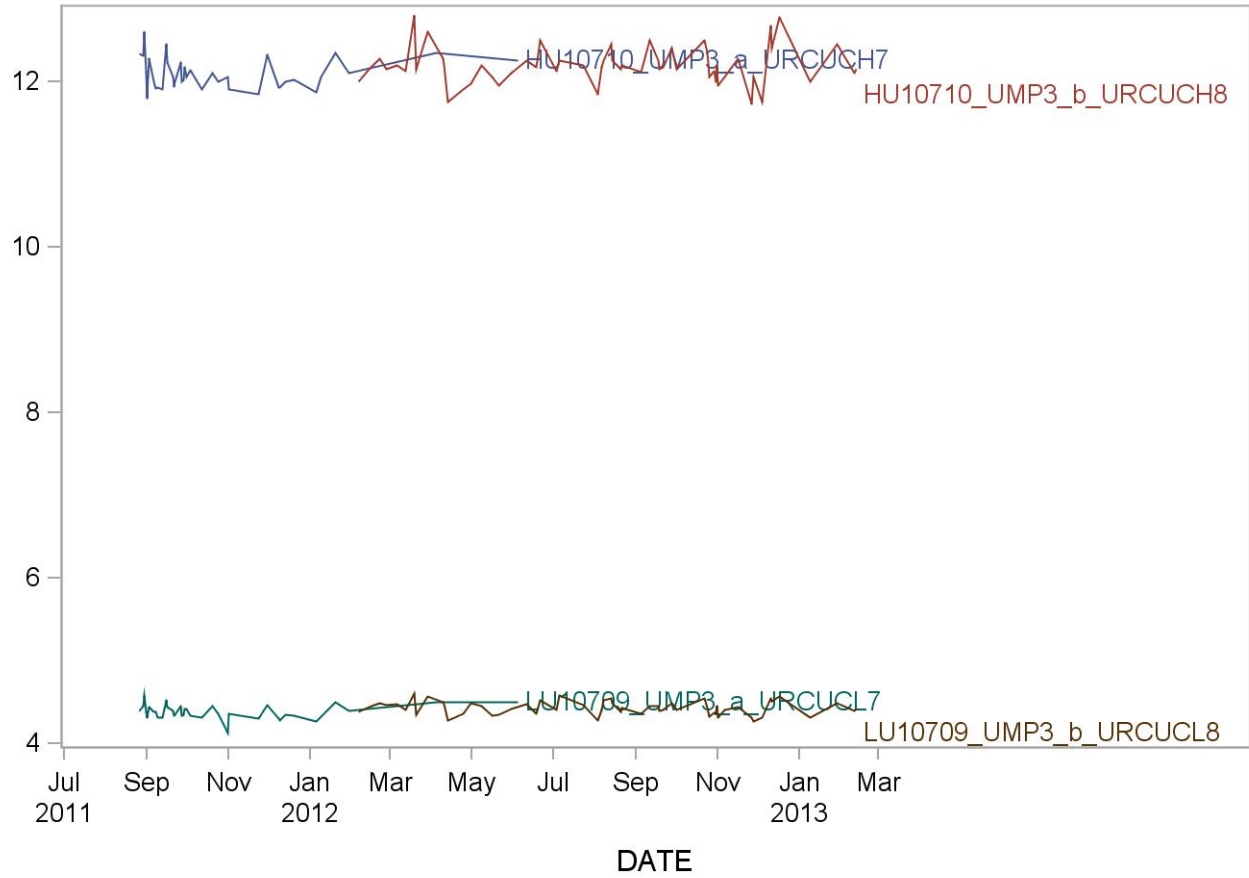

**Summary Statistics for Cobalt, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUOH7	85	26AUG11	04JUN12	2.2669	0.0586	2.6
LU10709_UMP3_a_URCUOL7	85	26AUG11	04JUN12	0.4001	0.0230	5.7
HU10710_UMP3_b_URCUOH8	86	06FEB12	11FEB13	2.2743	0.0662	2.9
LU10709_UMP3_b_URCUOL8	86	06FEB12	11FEB13	0.3915	0.0235	6.0


**2011-2012 Cobalt, urine (ug/L) Quality Control**

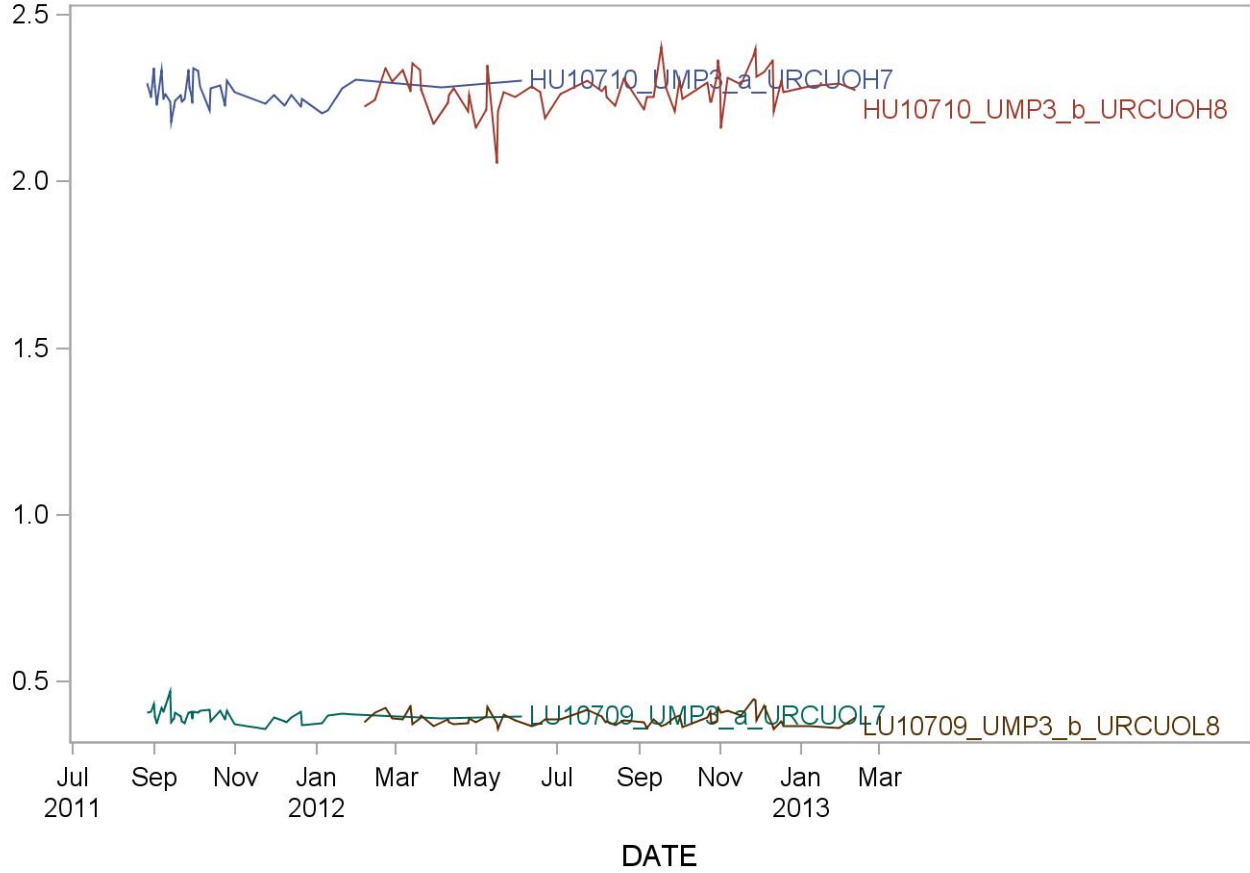

**Summary Statistics for Lead, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUPH7	94	26AUG11	04JUN12	2.776	0.061	2.2
LU10709_UMP3_a_URCUPL7	94	26AUG11	04JUN12	0.382	0.018	4.8
HU10710_UMP3_b_URCUPH8	77	06FEB12	12FEB13	2.772	0.066	2.4
LU10709_UMP3_b_URCUPL8	77	06FEB12	12FEB13	0.371	0.035	9.5


**2011-2012 Lead, urine (ug/L) Quality Control**

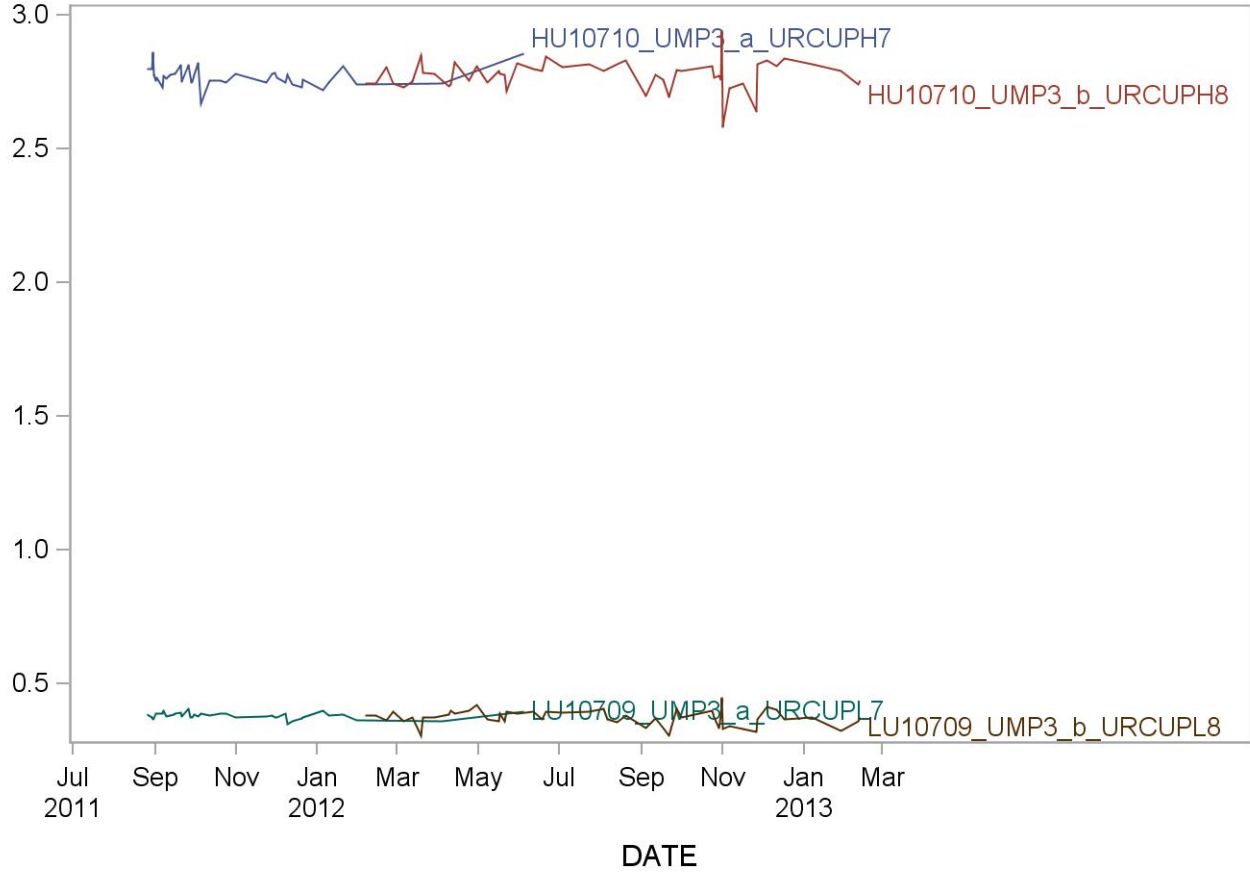

**Summary Statistics for Manganese, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUNH1	83	26AUG11	04JUN12	2.716	0.175	6.4
LU10709_UMP3_a_URCUNL1	83	26AUG11	04JUN12	0.708	0.133	18.8
HU10710_UMP3_b_URCUNH2	78	06FEB12	12FEB13	2.894	0.136	4.7
LU10709_UMP3_b_URCUNL2	78	06FEB12	12FEB13	0.776	0.068	8.7


**2011-2012 Manganese, urine (ug/L) Quality Control**

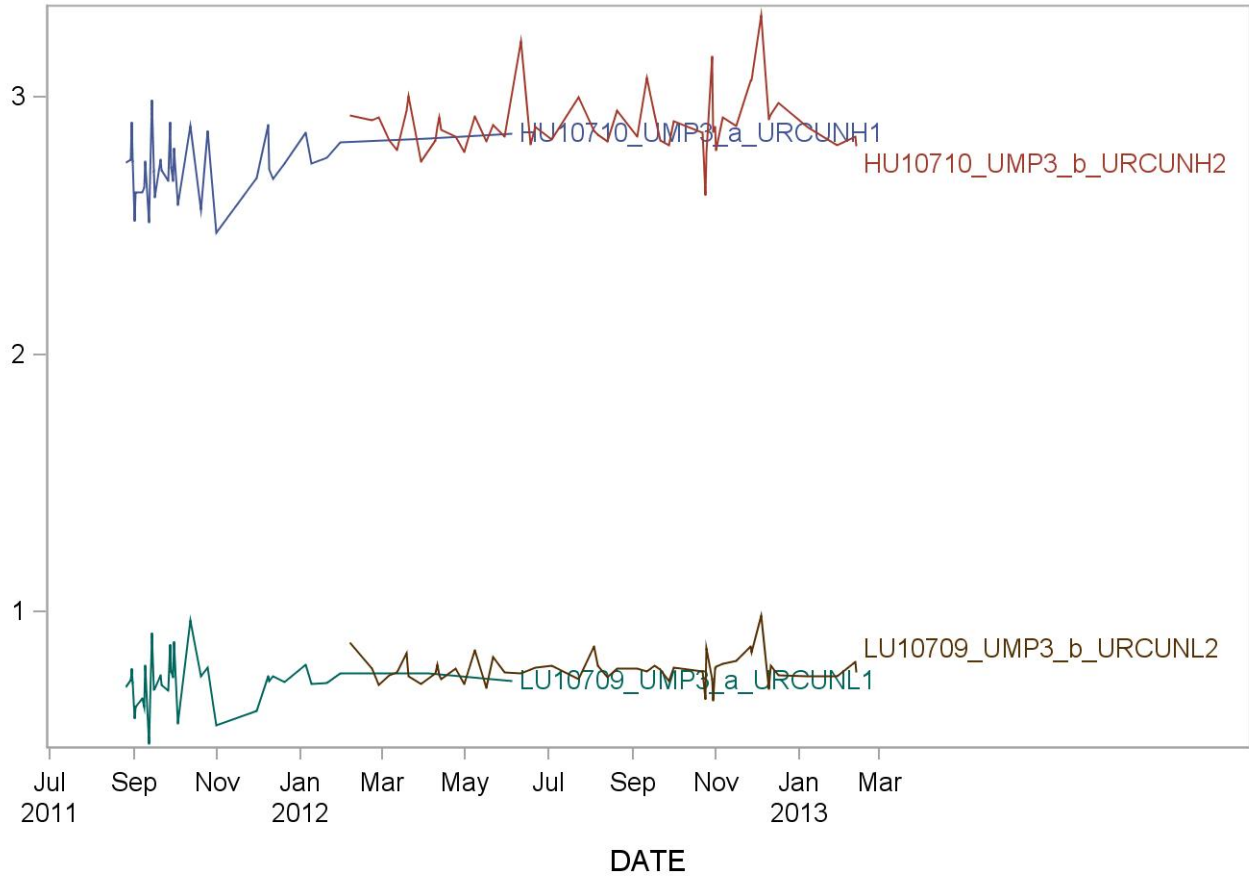

### Summary Statistics for Molybdenum, urine (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUMH7	96	26AUG11	04JUN12	128.755	2.514	2.0
LU10709_UMP3_a_URCUML7	96	26AUG11	04JUN12	38.327	0.743	1.9
HU10710_UMP3_b_URCUMH8	86	06FEB12	11FEB13	129.616	1.738	1.3
LU10709_UMP3_b_URCUML8	86	06FEB12	11FEB13	38.689	0.489	1.3

### 2011-2012 Molybdenum, urine (ug/L) Quality Control

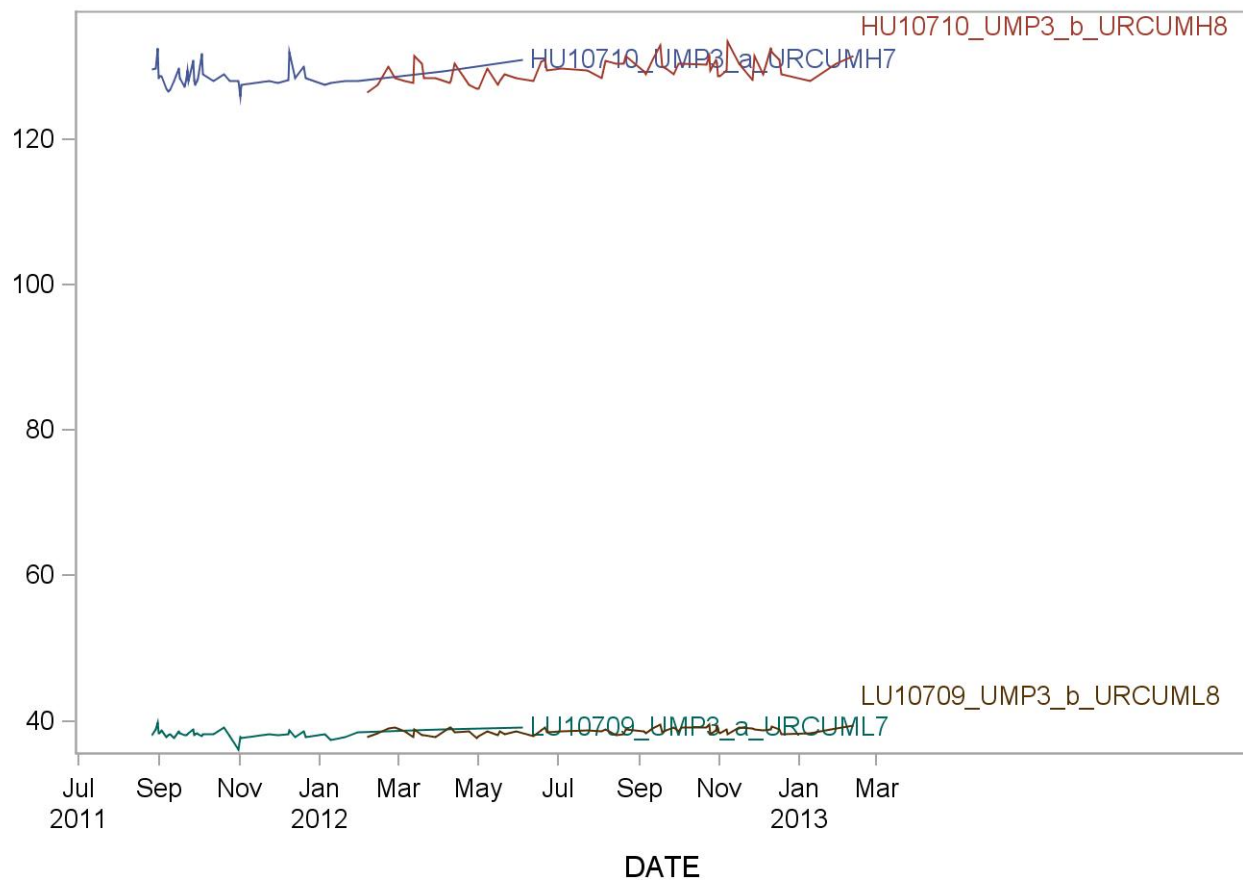




### Summary Statistics for Strontium, urine (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCSRH1	101	26AUG11	04JUN12	164.322	3.794	2.3
LU10709_UMP3_a_URCSRL1	101	26AUG11	04JUN12	79.656	1.898	2.4
HU10710_UMP3_b_URCSRH2	92	06FEB12	12FEB13	165.288	2.771	1.7
LU10709_UMP3_b_URCSRL2	92	06FEB12	12FEB13	80.261	1.449	1.8


### 2011-2012 Strontium, urine (ug/L) Quality Control

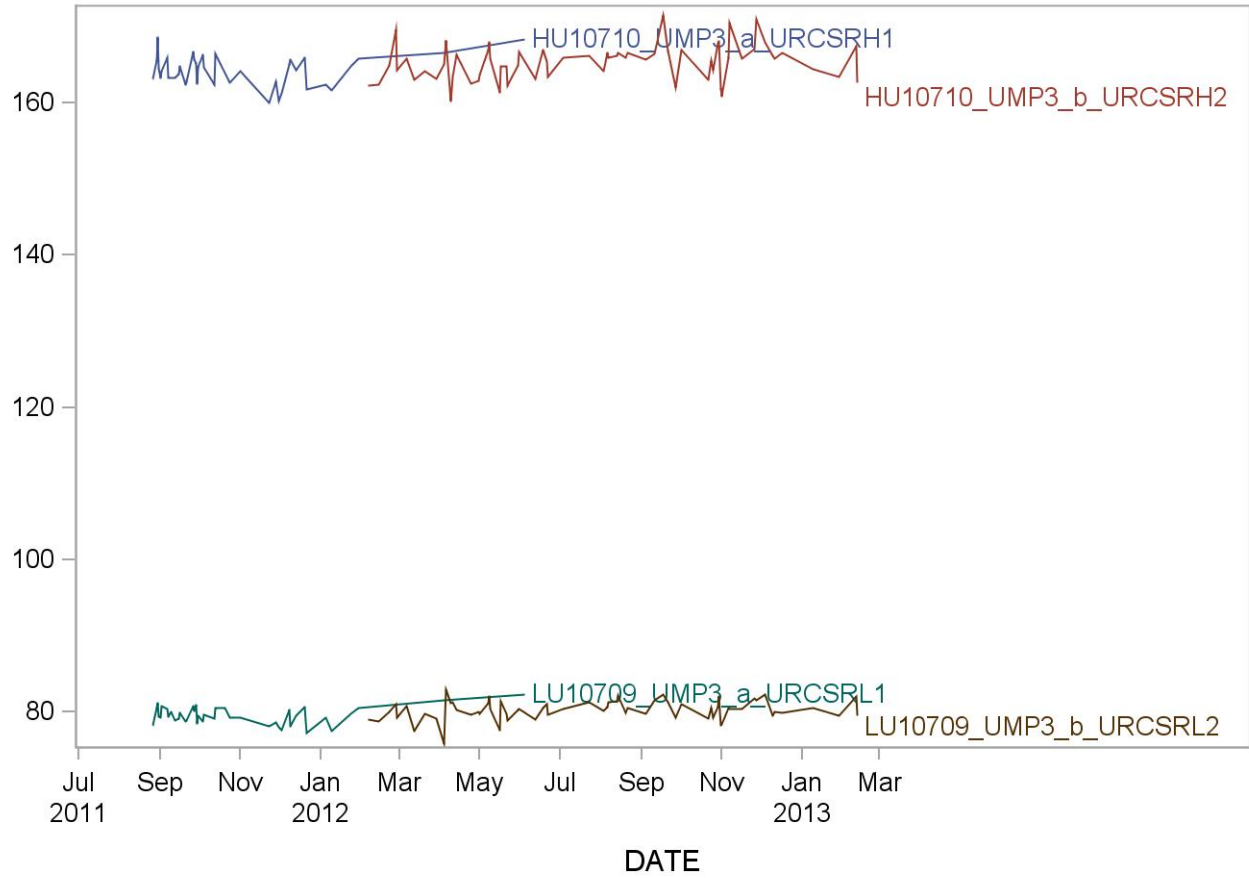

### Summary Statistics for Thallium, urine (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCULH7	86	26AUG11	06JUN12	0.4552	0.0124	2.7
LU10709_UMP3_a_URCULL7	86	26AUG11	06JUN12	0.1686	0.0080	4.8
HU10710_UMP3_b_URCULH8	88	06FEB12	11FEB13	0.4574	0.0091	2.0
LU10709_UMP3_b_URCULL8	88	06FEB12	11FEB13	0.1715	0.0071	4.1


### 2011-2012 Thallium, urine (ug/L) Quality Control

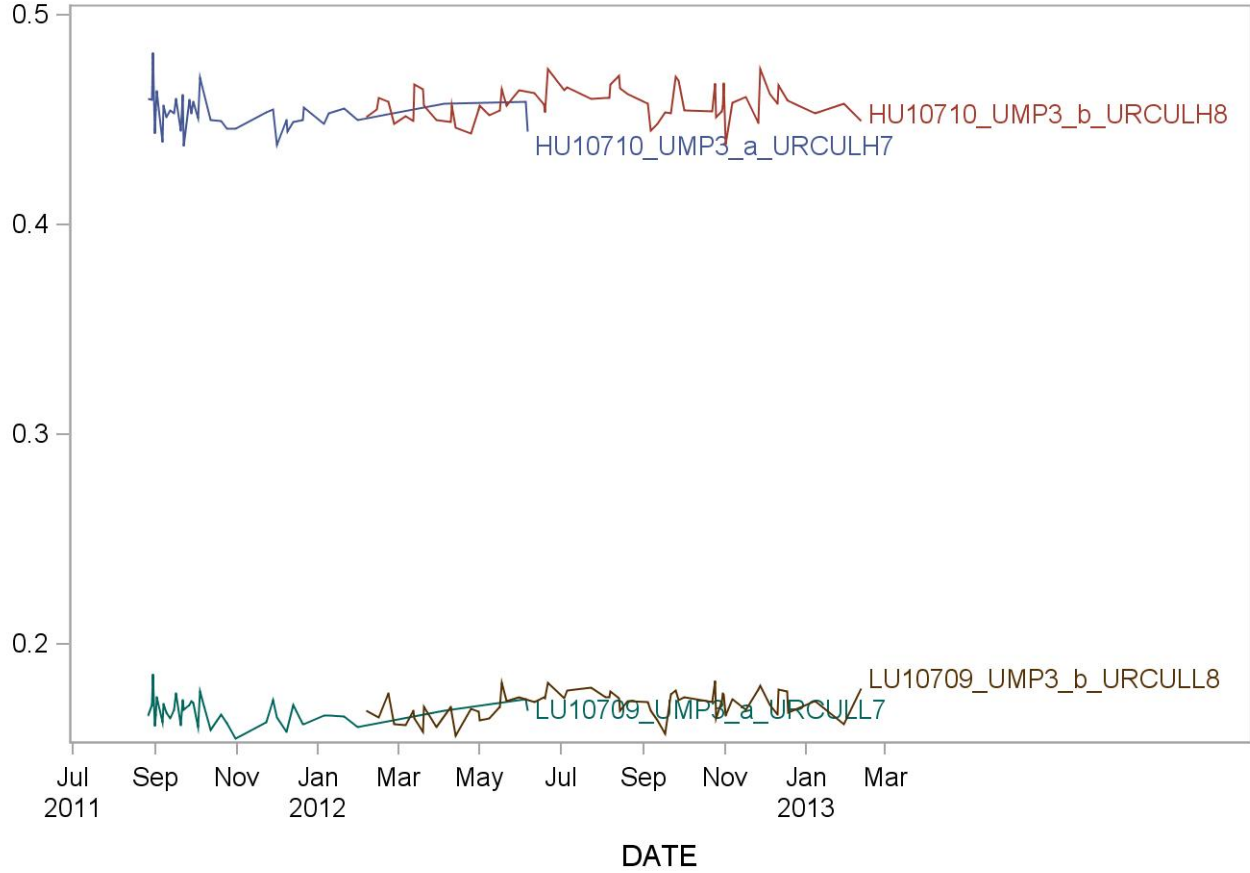

**Summary Statistics for Tin, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCSNH1	93	26AUG11	04JUN12	15.458	0.677	4.4
LU10709_UMP3_a_URCSNL1	93	26AUG11	04JUN12	2.677	0.145	5.4
HU10710_UMP3_b_URCSNH2	81	06FEB12	11FEB13	15.637	0.541	3.5
LU10709_UMP3_b_URCSNL2	81	06FEB12	11FEB13	2.688	0.124	4.6


**2011-2012 Tin, urine (ug/L) Quality Control**

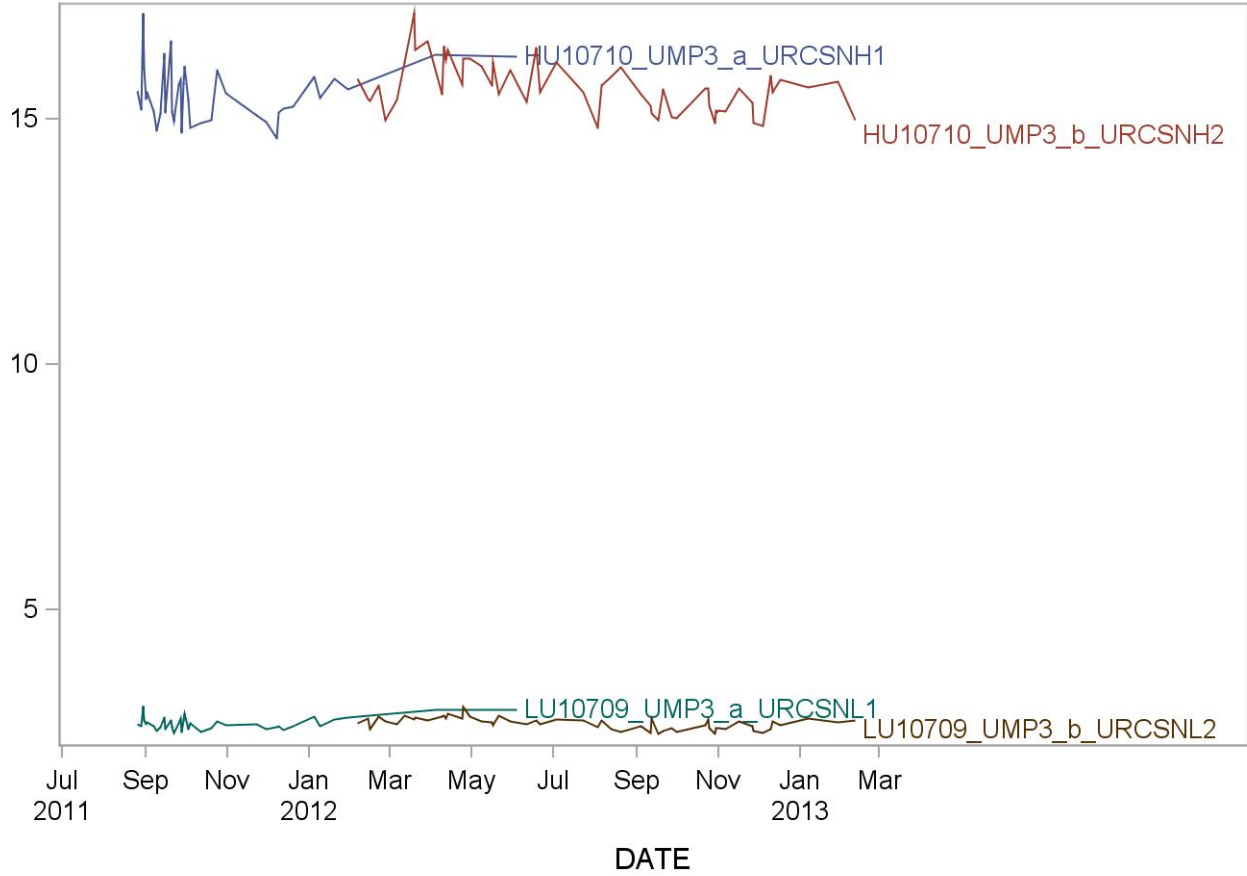

**Summary Statistics for Tungsten, urine (ug/L)**

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCUUH7	92	26AUG11	04JUN12	0.9005	0.0247	2.7
LU10709_UMP3_a_URCUUL7	92	26AUG11	04JUN12	0.1502	0.0100	6.7
HU10710_UMP3_b_URCUUH8	84	06FEB12	11FEB13	0.8921	0.0180	2.0
LU10709_UMP3_b_URCUUL8	84	06FEB12	11FEB13	0.1513	0.0100	6.6


**2011-2012 Tungsten, urine (ug/L) Quality Control**

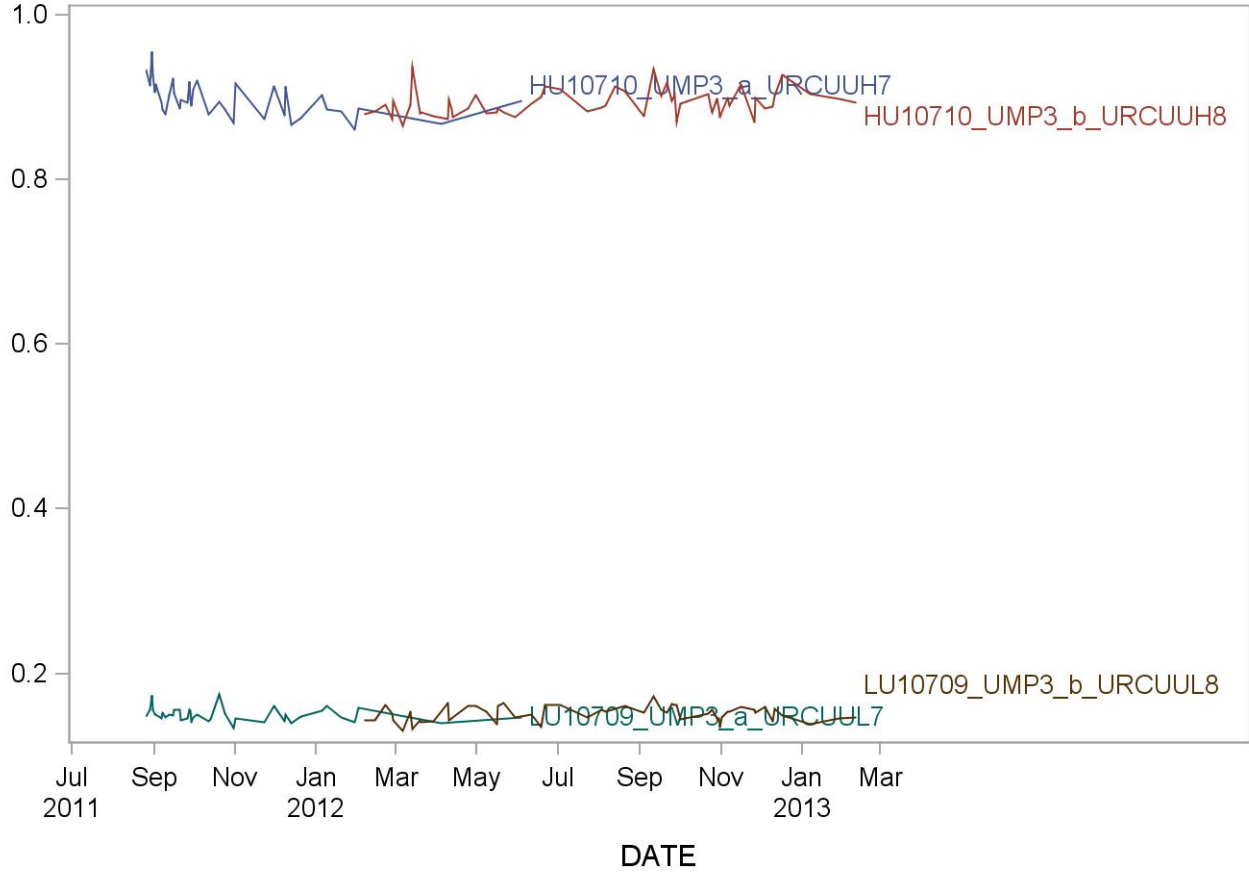

### Summary Statistics for Uranium, urinary (ug/L)

Lot	N	Start Date	End Date	Mean	Standard Deviation	Coefficient of Variation
HU10710_UMP3_a_URCURH7	94	26AUG11	04JUN12	0.11183	0.00632	5.7
LU10709_UMP3_a_URCURL7	94	26AUG11	04JUN12	0.01711	0.00157	9.2
HU10710_UMP3_b_URCURH8	76	06FEB12	11FEB13	0.11691	0.00753	6.4
LU10709_UMP3_b_URCURL8	76	06FEB12	11FEB13	0.01793	0.00128	7.1

### 2011-2012 Uranium, urinary (ug/L) Quality Control

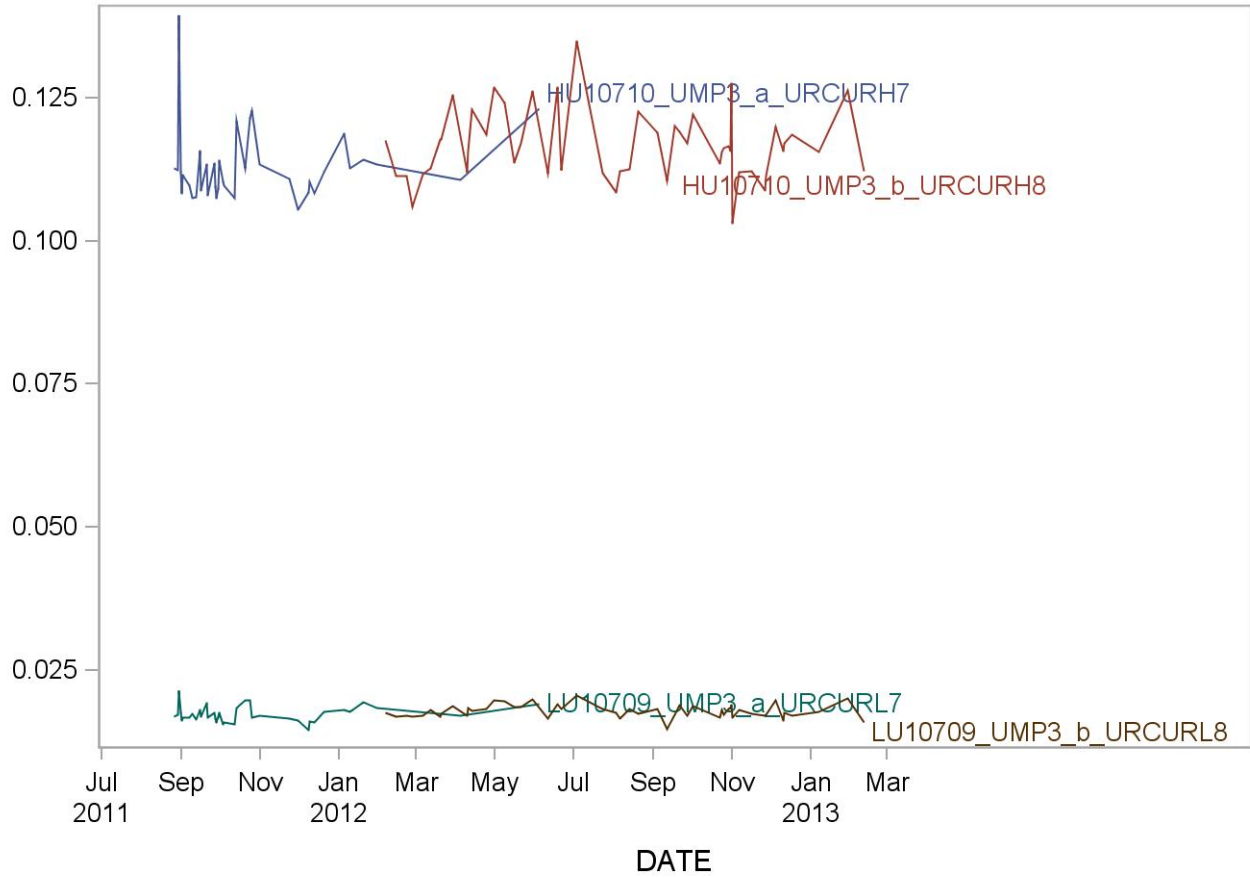

## References

1. Thomas, R., *Practical Guide to ICP-MS (Practical Spectroscopy)*. 2003, New York, NY: Marcel Dekker 336.
2. Tanner, S.D., Baranov, Vladimir I, *Theory, Design, and Operation of a Dynamic Reaction Cell for ICP-MS*. Atomic Spectroscopy, 1999. **20**(2): p. 45-52.
3. Tanner, S.D., V.I. Baranov, and D.R. Bandura, *Reaction cells and collision cells for ICP-MS: a tutorial review*. Spectrochimica Acta Part B-Atomic Spectroscopy, 2002. **57**(9): p. 1361-1452.
4. PerkinElmer SCIEX Instruments, *ELAN DRC II Hardware Guide*. 2001, Canada.
5. Mulligan, K.J., T.M. Davidson, and J.A. Caruso, *Feasibility Of The Direct Analysis Of Urine By Inductively Coupled Argon Plasma Mass-Spectrometry For Biological Monitoring Of Exposure To Metals*. Journal Of Analytical Atomic Spectrometry, 1990. **5**(4): p. 301-306.
6. Jarrett, J.M., *Total Urine Arsenic Biomonitoring Using Inductively Coupled Plasma Mass Spectrometry with a Dynamic Reaction Cell*. 2005, Centers for Disease Control and Prevention.
7. Jarrett, J.M., *Elimination of Molybdenum Oxide Interference In Urine Cadmium Analysis Using Inductively Coupled Plasma Reaction Cell Mass Spectrometry*. 2004, Centers for Disease Control and Prevention.
8. Larsen, E.H. and S. Sturup, *Carbon-enhanced Inductively Coupled Plasma Mass Spectrometric Detection of Arsenic and Selenium and Its Application to Arsenic Speciation*. Journal Of Analytical Atomic Spectrometry, 1994. **9**: p. 1101-1105.
9. Amarasiriwardena, C.J., et al., *Determination of the total arsenic concentration in human urine by inductively coupled plasma mass spectrometry: a comparison of the accuracy of three analytical methods*. Analyst, 1998. **123**(3): p. 441-445.
10. Office of Health and Safety in the Division of Laboratory Sciences, *Policies and Procedures Manual*. 2002, Division of Laboratory Sciences (DLS), National Center for Environmental Health, Centers for Disease Control and Prevention, Public Health Service, Department of Health and Human ServicesCenters for Disease Control and Prevention, .
11. Centers for Disease Control and Prevention (CDC) Radiation Safety Committee, *CDC/ATSDR Occupational Health and Safety Manual (Radiation Safety chapter)*. Centers for Disease Control and Prevention, Public Health Service, Department of Health and Human ServicesCenters for Disease Control and Prevention.
12. Heitland, P. and H.D. Koster, *Biomonitoring of 37 trace elements in blood samples from inhabitants of northern Germany by ICP-MS*. Journal of Trace Elements in Medicine and Biology, 2006. **20**(4): p. 253-262.
13. U.S. Nuclear Regulatory Commission, *Regulatory guide 8.22 (revision 1). Bioassay at uranium mills*. 1988: Atlanta, GA.
14. Centers for Disease Control and Prevention, *Third National Report on Human Exposure to Environmental Chemicals*, <http://www.cdc.gov/exposurereport>. 2005.
15. Stokinger, H.E., *The metals*, in *Patty's industrial hygiene and toxicology*

**Multi-Element in urine  
NHANES 2011-2012**

- G. Clayton and F. Clayton, Editors. 1981, John Wiley and Sons: New York. p. 1493-2060.
16. Fowler, B.A., in *Toxicology of trace elements*
  - R. Goyer and M. Mehlman, Editors. 1977, John Wiley and Sons: New York. p. p. 79.
  17. Iffland, R., *Arsenic*, in *Handbook on metals in clinical and analytical chemistry*, H. Seiler, A. Sigel, and H. Sigel, Editors. 1994, Marcel Dekker, Inc.: New York. p. 238-250.
  18. Gerhardsson, L. and S. Skerfving, *Concepts on biological markers and biomonitoring for metal toxicity*, in *Toxicology of metals*, L. Chang, Editor. 1996, CRC Press: Boca Raton, Florida. p. 98.
  19. *Report on Human Biomonitoring of Environmental Chemicals in Canada*. 2010, Health Canada: Ottawa.
  20. Heitland, P. and H.D. Koster, *Biomonitoring of 30 trace elements in urine of children and adults by ICP-MS*. *Clinica Chimica Acta*, 2006. **365**(1-2): p. 310-318.
  21. Paschal, D.C., et al., *Trace metals in urine of United States residents: Reference range concentrations*. *Environmental Research*, 1998. **76**(1): p. 53-59.
  22. Agency for Toxic Substances and Disease Registry (ATSDR). 2000. Toxicological profile for Manganese. Atlanta, G.U.S.D.o.H.a.H.S., Public Health Service. , *Toxicological Profile for Manganese*, ATSDR, Editor. 2000. p. 15.
  23. Moreno, M.E., et al., *Biomonitoring of metal in children living in a mine tailings zone in Southern Mexico: A pilot study*. *International Journal of Hygiene and Environmental Health*, 2010. **213**(4): p. 252-258.
  24. Gil, F., et al., *Biomonitorization of cadmium, chromium, manganese, nickel and lead in whole blood, urine, axillary hair and saliva in an occupationally exposed population*. *Science of the Total Environment*, 2011. **409**(6): p. 1172-1180.
  25. Wang, H.M., et al., *Urinary heavy metal levels and relevant factors among people exposed to e-waste dismantling*. *Environment International*, 2011. **37**(1): p. 80-85.
  26. Juliao, L., et al., *Exposure of workers in a mineral processing industry in Brazil*. *Radiation Protection Dosimetry*, 2007. **125**(1-4): p. 513-515.