

# National Health and Nutrition Examination Survey 2003-2004

---

## Documentation, Codebook, and Frequencies

### MEC Laboratory Component: Total Cholesterol and HDL-Cholesterol

**Survey Years:**  
**2003 to 2004**

**SAS Export File:**  
**L13\_C.XPT**



**First Publish:** June 2006  
**Last Revised:** N/A

# NHANES 2003–2004 Data Documentation

## Laboratory Assessment: Lab 13 – Total Cholesterol and HDL-Cholesterol

Years of Coverage: 2003–2004

First Published: June 2006

Last Revised: N/A

---

### Component Description

The data will be used to monitor the status of hyperlipidemia and the success of the National Cholesterol Education Program.


The main element of the cardiovascular disease laboratory component in NHANES is blood lipid levels. Cardiovascular disease is the leading cause of death in the United States. The data will be used to monitor the status of hyperlipidemia and the success of the National Cholesterol Education Program.

### Eligible Sample

Participants aged 3 years and older were tested.

### Description of Laboratory Methodology

Total cholesterol is measured enzymatically in serum in a series of coupled reactions that hydrolyze cholesteryl esters and oxidize the 3-OH group of cholesterol. One of the reaction byproducts,  $\text{H}_2\text{O}_2$  is measured quantitatively in a peroxidase-catalyzed reaction that produces a color. Absorbance is measured at 500 nm. The color intensity is proportional to cholesterol concentration. The reaction sequence is as follows:

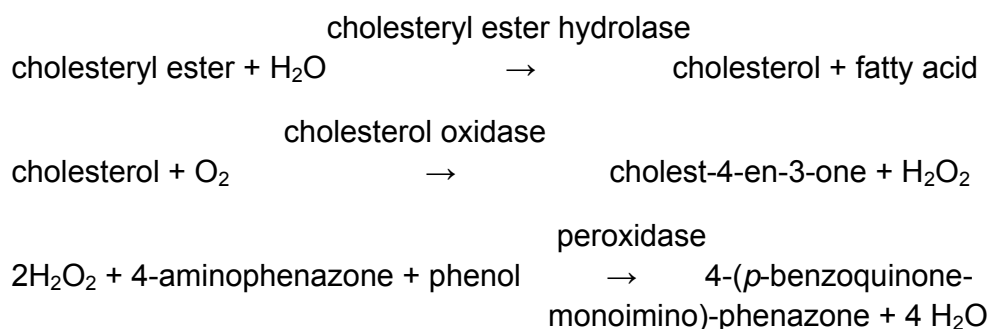

HDL-Cholesterol is measured directly in serum. The apolipoprotein B containing lipoproteins in the specimen are reacted with a blocking reagent that renders them non-reactive with the enzymatic cholesterol reagent under conditions of the assay.

The procedures use the Roche/Boehringer-Mannheim Diagnostics direct HDL method. The method uses sulfated alpha-cyclodextrin in the presence of  $Mg^{+2}$ , which forms complexes with apoB containing lipoproteins, and polyethylene glycol-coupled cholesteryl esterase and cholesterol oxidase for the HDL-cholesterol measurement. The reactions are as follows:

ApoB containing lipoproteins +  $\alpha$ -cyclodextrin +  $Mg^{+2}$  + dextran  $SO_4$  ---  
 > soluble non-reactive complexes with apoB-containing lipoproteins

HDL-cholesteryl esters  $\xrightarrow{\text{PEG-cholesteryl esterase}}$  HDL-unesterified  
 cholesterol + fatty acid

Unesterified cholesterol +  $O_2 \xrightarrow{\text{PEG-cholesterol oxidase}}$  cholestenone +  $H_2O_2$

$H_2O_2$  + 5-aminophenazone + N-ethyl-N-(3-methylphenyl)-N'-succinyl  
 ethylene diamine

+  $H_2O$  +  $H^+$  peroxidase ---> quinoneimine dye +  $H_2O$

Absorbance is measured at 600 nm.

There were no changes to the equipment, lab method, or lab site from the previous 2 years.

A detailed description of the laboratory method used can be found on the NHANES website.

## Laboratory Quality Control and Monitoring

The NHANES quality assurance and quality control (QA/QC) protocols meet the 1988 Clinical Laboratory Improvement Act mandates. Detailed quality control and quality assurance instructions are discussed in the NHANES Laboratory/Medical Technologists Procedures Manual (LPM). Read the LABDOC file for detailed QA/QC protocols. A detailed description of the QA/QC procedures can be found on the NHANES web site.

## Data Processing and Editing

Blood specimens were processed, stored, and shipped to Johns Hopkins Hospital, Baltimore, MD for analysis. Detailed specimen collection and processing instructions are discussed in the NHANES Laboratory/Medical Technologists Procedures Manual (LPM). Read the LABDOC file for detailed data processing and editing protocols. The analytical methods are described in the **Description of the Laboratory Methodology** section.

There was no top coding on this file.

Two derived variables were created in this data file. The formula for their derivation is as follows:

LBDTCSI:

The total cholesterol in mg/dL (LBXTC) was converted to mmol/L (LBDTCSI) by multiplying by 0.02586.

LBDHDDSI:

The HDL-cholesterol in mg/dL (LBXHDD) was converted to mmol/L (LBDHDDSI) by multiplying by 0.02586.

The HDL-cholesterol data was not corrected for 2003-2004 data. This was unlike the previous 1999-2000 and 2001-2002 data where HDL was corrected. The direct immunoassay method for 2003-2004 showed acceptable bias (CV < 4%) and precision (CV < 5%) when compared to HDL-cholesterol quality controls (Solomon Park Research Laboratories, Kirkland, WA) with assigned values established by the Centers for Disease Control and Prevention. Refer to the documentation for 1999-2000 and 2001-2002 for more details on how the data was corrected. (NHANES 1999-2000 Data Files and NHANES 2001-2002 Data Files web site)

Detailed instructions on specimen collection and processing can be found on the NHANES web site.

## **Analytic Notes**

The analysis of NHANES 2003–2004 laboratory data must be conducted with the key survey design and basic demographic variables. The NHANES 2003–2004 Household Questionnaire Data Files contain demographic data, health indicators, and other related information collected during household interviews. They also contain all survey design variables and sample weights for these age groups. The phlebotomy file includes auxiliary information such as the conditions precluding venipuncture. The household questionnaire and phlebotomy files may be linked to the laboratory data file using the unique survey participant identifier SEQN.

#### LBXTC:

The Lab 13 Total Cholesterol data file contains laboratory test results for total cholesterol (LBXTC), which uses the reference analytic method. However, the NHANES Laboratory 40 biochemistry profiles also include measurements of total cholesterol. The Laboratory 40 variable name is LBXSCH. The appropriate variable to use is LBXTC from Laboratory 13.

**References**      N/A

## Locator Fields

**Title:** Total Cholesterol and HDL-Cholesterol

**Contact Number:** 1-866-441-NCHS

**Years of Content:** 2003–2004

**First Published:** June 2006

**Revised:** N/A

**Access Constraints:** None

**Use Constraints:** None

**Geographic Coverage:** National

**Subject:** Cholesterol and HDL-Cholesterol

**Record Source:** NHANES 2003–2004

**Survey Methodology:** NHANES 2003–2004 is a stratified multistage probability sample of the civilian non-institutionalized population of the U.S.

**Medium:** NHANES Web site; SAS transport files

National Health and Nutrition Examination Survey  
Codebook for Data Production (2003-2004)

Total Cholesterol and HDL-Cholesterol (L13\_C)

Person Level Data

June 2006



SEQN	Target
	B(3 Yrs. to 150 Yrs.)
Hard Edits	SAS Label
	Respondent sequence number
English Text: Respondent sequence number.	
English Instructions:	



LBXTC	Target			
	B(3 Yrs. to 150 Yrs.)			
Hard Edits	SAS Label			
	Total Cholesterol (mg/dL)			
English Text: Total Cholesterol (mg/dL)				
English Instructions:				
Code or Value	Description	Count	Cumulative	Skip to Item
68 to 704	Range of Values	7774	7774	
.	Missing	782	8556	

LBXHDD	Target			
	B(3 Yrs. to 150 Yrs.)			
Hard Edits	SAS Label			
	Direct HDL-Cholesterol (mg/dL)			
English Text: Direct HDL-Cholesterol (mg/dL)				
English Instructions:				
Code or Value	Description	Count	Cumulative	Skip to Item
19 to 154	Range of Values	7773	7773	
.	Missing	783	8556	

LBDTCSI	Target			
	B(3 Yrs. to 150 Yrs.)			
Hard Edits	SAS Label			
	Total Cholesterol (mmol/L)			
English Text: Total Cholesterol (mmol/L)				
English Instructions:				
Code or Value	Description	Count	Cumulative	Skip to Item
1.76 to 18.21	Range of Values	7774	7774	
.	Missing	782	8556	

LBDHDDSI	Target			
	B(3 Yrs. to 150 Yrs.)			
Hard Edits	SAS Label			
	Direct HDL-Cholesterol (mmol/L)			
English Text: Direct HDL-Cholesterol (mmol/L)				
English Instructions:				
Code or Value	Description	Count	Cumulative	Skip to Item
0.49 to 3.98	Range of Values	7773	7773	
.	Missing	783	8556	