
National Health
and Nutrition
Examination Survey

CARDIOVASCULAR
FITNESS
PROCEDURES
MANUAL



(Revised November 2000)

TABLE OF CONTENTS

<u>Chapter</u>		<u>Page</u>
1	OVERVIEW OF CARDIOVASCULAR FITNESS	1-1
1.1	Cardiovascular Fitness in NHANES	1-1
1.2	Personnel.....	1-2
1.3	Integrated Survey Information System (ISIS).....	1-2
2	EQUIPMENT/SUPPLIES/MATERIALS.....	2-1
2.1	Quinton MedTrack ST65 Treadmill with Hyperdrive	2-1
2.1.1	Maintenance of the Treadmill	2-2
2.1.2	Calibration of the Treadmill	2-5
2.2	Quinton MedTrack ST Programmable Controller	2-6
2.2.1	Maintenance of the Programmable Controller	2-8
2.2.2	Calibration of the Programmable Controller	2-8
2.3	Colin STBP-780 Automatic Blood Pressure Monitor.....	2-10
2.3.1	Maintenance of the Automatic Blood Pressure Monitor	2-11
2.3.2	Calibration of the Automatic Blood Pressure Monitor.....	2-12
2.3.3	Mercury Manometer (Used for Calibration)	2-14
2.3.4	Mercury Spills and Leaks	2-15
2.4	SECA Floor Scale Model 815.....	2-16
2.4.1	Maintenance and Calibration of the Floor Scales.....	2-17
2.5	SECA Portable Stadiometer Model 222.....	2-18
2.5.1	Maintenance and Calibration of the Stadiometer	2-18
2.6	Description of Supplies	2-19
2.6.1	ECG Electrodes/Alcohol Prep Pads/Abrasive Pads	2-19
2.6.2	Running Shoes.....	2-19
2.6.3	Mesh Vests	2-19
2.7	Contact Information for Equipment Repair or Parts Replacement	2-20
2.8	Equipment Manuals	2-20

TABLE OF CONTENTS (continued)

<u>Chapter</u>		<u>Page</u>
	2.9 Troubleshooting	2-20
	2.9.1 Quinton MedTrack ST 65 Treadmill and Controller	2-20
	2.9.2 Colin STBP-780 Blood Pressure Monitor	2-22
3	PROTOCOL.....	3-1
	3.1 Exclusions to the CV Fitness Component.....	3-1
	3.1.1 Medical Conditions	3-1
	3.1.2 Medications	3-1
	3.1.3 Safety/Exclusion Questions.....	3-3
	3.1.4 Exclusions Based on Heart Rate and Blood Pressure.....	3-4
	3.2 Introduction to Component	3-4
	3.3 Measurement of Height and Weight to Determine Body Mass Index.....	3-5
	3.3.1 Height	3-5
	3.3.2 Weight	3-6
	3.4 Checking SP's Footwear	3-8
	3.5 Determination of Physical Activity Readiness Code	3-9
	3.5.1 Questions to Determine PAR Code: English Version	3-9
	3.5.2 Questions to Determine PAR Code: Spanish Version.....	3-12
	3.6 Determination of Predicted VO ₂ Max.....	3-13
	3.7 Determination of Exercise Protocol	3-14
	3.8 Preparation for Treadmill Test	3-14
	3.8.1 Preparation of Skin and Placement of Electrodes	3-14
	3.8.2 Selection of Arm for Blood Pressure Monitoring.....	3-15
	3.8.3 Selection of Cuff Size and Placement of Cuff.....	3-16
	3.8.4 Room Temperature.....	3-17
	3.8.5 Mesh Vest.....	3-17
	3.9 Treadmill Test	3-18
	3.9.1 Setup of Colin STBP Automatic BP Monitor	3-18
	3.9.2 Setup of Programmable Treadmill Controller.....	3-19
	3.9.3 Obtaining Pretest Blood Pressure and Heart Rate	3-20
	3.9.4 Script for Explaining Treadmill Test.....	3-21

TABLE OF CONTENTS (continued)

<u>Chapter</u>		<u>Page</u>
	3.9.5 Demonstrating How to Walk on a Treadmill	3-21
	3.9.6 Warm-up (00:00 - 02:00 Minutes)	3-23
	3.9.7 RPE Explanation and Assessment.....	3-23
	3.9.8 Stage 1 (00:00 - 03:00 Minutes).....	3-25
	3.9.9 Stage 2 (00:00 - 03:00 Minutes).....	3-26
	3.9.10 Recovery Stages 1, 2, & 3 (If Necessary).....	3-27
	3.9.11 End of Test	3-28
3.10	Procedure if Values Exceed Criteria	3-28
	3.10.1 Exception for RPE greater than 17 at end of Stage 2	3-29
3.11	Procedure if SP is Unable to Adhere to the Protocol	3-30
3.12	Procedure to Change Predicted Exercise Protocol.....	3-31
	3.12.1 Increasing the Exercise Protocol	3-32
	3.12.2 Decreasing the Exercise Protocol	3-33
3.13	Procedure to Stop Treadmill Test Early	3-34
	3.13.1 Priority 1 Stopping Criteria	3-34
	3.13.2 Priority 2 Stopping Criteria	3-35
3.14	Procedure to Recover/Restart an Exam.....	3-37
3.15	Recording Problems with Data Acquisition During the Treadmill Test	3-38
3.16	Component Status	3-38
4	INTEGRATED SURVEY INFORMATION SYSTEM (ISIS) DATA ENTRY SCREENS	4-1
4.1	Height and Weight Measurements to Determine BMI.....	4-1
	4.1.1 Edits on Data Entry for Weight and Height	4-3
4.2	Questions to Determine Exercise Protocol (PAR Code Questions)....	4-8
	4.2.1 First Level Questions.....	4-10
	4.2.2 Second Level Questions (Little or No Regular Physical Activity)	4-11
	4.2.3 Second Level Questions (Regular Modest Physical Activity)	4-12

TABLE OF CONTENTS (continued)

<u>Chapter</u>		<u>Page</u>
	4.2.4 Second Level Questions (Regular Heavy Physical Activity).....	4-14
	4.2.5 Required Data Entry for PAR Questions.....	4-16
4.3	Pretest Measurements.....	4-17
	4.3.1 Digital Thermometer for Measuring Room Temperature ...	4-18
	4.3.2 Exam Preparation –Arm Used, Cuff Size, Mesh Vest Used	4-19
	4.3.3 Capturing Blood Pressure Measurements.....	4-25
	4.3.4 Demonstrating and Explaining How to Walk on a Treadmill	4-26
	4.3.5 Required Data Entry for Arm and Cuff Size	4-27
	4.3.6 Required Data Entry for BP Measurements	4-28
4.4	Treadmill Test	4-29
	4.4.1 Function of Treadmill Test Screen	4-29
	4.4.2 Preparing the Treadmill and Controller.....	4-32
	4.4.3 Getting the SP on the Treadmill	4-33
	4.4.4 Warm-up.....	4-34
	4.4.5 Show Data	4-37
	4.4.6 Stage 1	4-38
	4.4.7 Stage 2	4-44
	4.4.8 Recovery Stages	4-48
	4.4.9 Checking Estimated VO ₂ Max	4-49
	4.4.10 Component Status for Routine Test	4-50
4.5	Changing the Exercise Protocol	4-51
	4.5.1 Changing the Treadmill Controller	4-51
	4.5.2 Decreasing the Exercise Protocol	4-52
	4.5.3 Increasing the Exercise Protocol	4-54
	4.5.4 Lowest and Highest Protocol.....	4-56
4.6	Emergency Stop	4-58
4.7	Priority 1 Stopping Criteria.....	4-59
4.8	Early Termination of Test – Go to Recovery.....	4-62
	4.8.1 Data Not Within Acceptable Limits	4-63
4.9	Priority 2 Stopping Criteria.....	4-68
4.10	Problems With Data Acquisition During Treadmill Test.....	4-71

TABLE OF CONTENTS (continued)

<u>Chapter</u>		<u>Page</u>
4.11	Component Status Categories	4-74
	4.11.1 Complete	4-74
	4.11.2 Partial.....	4-75
	4.11.3 Not Done	4-81
5	REFERRALS AND REPORT OF FINDINGS.....	5-1
5.1	Observation Referrals.....	5-1
5.2	Report of Findings for CV Fitness	5-3
6	QUALITY CONTROL	6-1
6.1	Overview	6-1
6.2	Training.....	6-1
	6.2.1 Initial Training.....	6-1
	6.2.2 Annual Retraining	6-2
6.3	Equipment and Room Set-up Checks.....	6-3
	6.3.1 Daily	6-3
	6.3.2 Weekly.....	6-4
	6.3.3 Stand.....	6-5
6.4	Data Entry Screens for QC on Equipment	6-6
	6.4.1 Quality Control Reminder Message Box.....	6-6
	6.4.2 Utilities Menu to Select Quality Control.....	6-7
	6.4.3 Quality Control Log-on Box.....	6-7
	6.4.4 Data Entry Screen for QC Daily Checks	6-8
	6.4.5 Data Entry Screen for QC Weekly Checks	6-12
	6.4.6 Data Entry Screen for QC Start of Stand Checks	6-16
	6.4.7 Incomplete Entry of Quality Control Data	6-18
6.5	Observational Visits	6-18
	6.5.1 Technician Competency Evaluation Form	6-19
	6.5.2 Quality Control Review of CV Fitness Exclusions in Physician's Exam.....	6-24
	6.5.3 Site Visit Report Form.....	6-26

TABLE OF CONTENTS (continued)

List of Appendixes

Appendix

A	SCRIPTS	A-1
B	EXCLUSIONS BASED ON MEDICAL CONDITIONS.....	B-1
C	EXCLUSIONS BASED ON MEDICATIONS.....	C-1
D	CARDIOVASCULAR SAFETY AND EXCLUSION QUESTIONS.....	D-1
E	DETERMINATION OF PHYSICAL ACTIVITY READINESS CODE.....	E-1
F	CALCULATION OF PREDICTED V_{O_2} MAX.....	F-1
G	CALCULATION OF ESTIMATED V_{O_2} MAX	G-1
H	DETERMINATION OF EXERCISE PROTOCOL.....	H-1
I	REFERENCE RANGES FOR V_{O_2} MAX.....	I-1
J	SPEED AND GRADE TABLE BY PROTOCOL AND STAGE	J-1
K	REVIEW OF EDIT CRITERIA.....	K-1
L	SAMPLE DATA ENTRY SCREENS	L-1

List of Tables

Table

2-1	Troubleshooting problems with the Quinton Med Track ST65 Treadmill and Controller.....	2-21
2-2	Troubleshooting problems with the STBP monitor	2-22
3-1	Determination of exercise protocol based on predicted maximal oxygen consumption	3-11

TABLE OF CONTENTS (continued)

List of Figures

Figure

2-1	Quinton MedTrack ST65 Treadmill.....	2-1
2-2	Treadmill calibration.....	2-4
2-3	Calibration of elevation.....	2-6
2-4	Quinton MedTrack ST Programmable Controller.....	2-7
2-5	Colin STBP-780 Automatic Blood Pressure Monitor.....	2-11
3-1	Correct body placement for height measurement.....	3-7
3-2	Site for placement of electrodes (1)	3-14
3-3	Site for placement of electrodes (2)	3-16
3-4	Treadmill controller - changing the exercise protocol	3-32
5-1	Menu to select observation referral.....	5-1
5-2	Pick list of SP's in current session	5-2
5-3	Observation referral in CV Fitness.....	5-2
5-4	Observation Referral from CV Fitness in physician's referral review box.....	5-3

List of Exhibits

Exhibit

3-1	CV exclusions based on medications in the physician's exam	3-2
4-1	Weight and height measurements from body composition	4-1
4-2	Weight and height measurements in CV Fitness room	4-2
4-3	Edit message on incorrect data entry for weight	4-3
4-4	Edit message on incorrect data entry for weight (below edit range).....	4-4

TABLE OF CONTENTS (continued)

List of Exhibits (continued)

<u>Exhibit</u>		<u>Page</u>
4-5	Edit message on incorrect data entry for height	4-5
4-6	Required data entry for height and weight	4-6
4-7	Cannot calculate VO ₂ max	4-6
4-8	Partial status with comment 'Other, specify'	4-7
4-9	Questions to determine protocol	4-8
4-10	First level questions to determine exercise protocol	4-10
4-11	Little or no regular physical activity (second level questions).....	4-11
4-12	Regular modest physical activity (second level questions).....	4-12
4-13	Second level questions (regular heavy physical activity)	4-14
4-14	Required data entry for first level activity questions	4-16
4-15	Pretest measurements	4-17
4-16	Unable to obtain room temperature.....	4-18
4-17	Exam preparation – arm selected	4-19
4-18	SP excluded if neither arm can be used for BP measurement.....	4-20
4-19	SP excluded due to physical limitation	4-20
4-20	Partial status with physical limitation if CNO selected for arm.....	4-21
4-21	Selection of cuff size.....	4-22
4-22	Problem with cuff size	4-23
4-23	Mesh vest used	4-24
4-24	Blood pressure data capture	4-25
4-25	Demonstrate and explain how to walk on a treadmill.....	4-26

TABLE OF CONTENTS (continued)

List of Exhibits (continued)

<u>Exhibit</u>		<u>Page</u>
4-26	Required data entry for arm and cuff used	4-27
4-27	Required data entry for blood pressure	4-28
4-28	Treadmill test data entry screen	4-29
4-29	Preparing the treadmill and controller.....	4-32
4-30	Starting the timer	4-34
4-31	Explaining the rating of perceived exertion (RPE)	4-35
4-32	Getting RPE for warm-up	4-36
4-33	Show data window	4-37
4-34	Show data message	4-37
4-35	Treadmill test stage 1 – routine test	4-38
4-36	Get rating of perceived exertion (RPE).....	4-39
4-37	Compare heart rate on monitor and screen.....	4-40
4-38	Example 1 – manually entered heart rate - not necessary	4-41
4-39	Show data screen for Example 1	4-41
4-40	Example 2 – manually entered heart rate - necessary	4-42
4-41	Show data window for Example 2.	4-42
4-42	Example 3 - manually entered heart rate – incorrect data entry.....	4-43
4-43	Show data window for Example 3	4-43
4-44	Manual entry fields disabled	4-44
4-45	Treadmill test stage 2 – (19 through 49 years).....	4-45

TABLE OF CONTENTS (continued)

List of Exhibits (continued)

<u>Exhibit</u>		<u>Page</u>
4-46	Treadmill test stage 2 – (12 through 18 years).....	4-46
4-47	Treadmill test – last minute of stage 2	4-47
4-48	Treadmill test, recovery – routine test.....	4-48
4-49	Checking estimated VO ₂ max	4-49
4-50	Checking estimated VO ₂ max	4-49
4-51	Component status - complete	4-50
4-52	Treadmill controller	4-51
4-53	Decreasing the exercise protocol.....	4-52
4-54	Decreased protocol.....	4-53
4-55	Increasing the exercise protocol	4-54
4-56	Changing the protocol in stage 1 or 2 - message.....	4-55
4-57	Protocol 0	4-56
4-58	Protocol 0 message.....	4-56
4-59	Protocol 9	4-57
4-60	Protocol 9 message.....	4-58
4-61	Emergency stop	4-58
4-62	Priority 1 stopping criteria - other	4-59
4-63	Priority 1 stopping criteria - required entry	4-61
4-64	Go to recovery	4-62
4-65	Data not within acceptable limits	4-63

TABLE OF CONTENTS (continued)

List of Exhibits (continued)

<u>Exhibit</u>		<u>Page</u>
4-66	SBP not within acceptable limits	4-64
4-67	DBP not within acceptable limits.....	4-65
4-68	Rating of perceived exertion not within acceptable limits	4-66
4-69	Heart rate greater than 85 percent PMHR.....	4-67
4-70	Priority 2 stopping criteria	4-68
4-71	Component status for priority 2	4-69
4-72	SP unable to stop gripping handrails.....	4-70
4-73	Problems with data acquisition - No problems	4-71
4-74	Problems with data acquisition - BP, data capture, and HR.....	4-72
4-75	Problems with data acquisition - electrodes, BP cuff, and treadmill controller	4-73
4-76	Component status - complete	4-74
4-77	Component status - partial (1).....	4-75
4-78	Component status - partial (2).....	4-76
4-79	Component status - partial (3).....	4-77
4-80	Component status – partial (4)	4-79
4-81	Component status - partial - required comments.....	4-80
4-82	Component status - not done.....	4-81
5-1	Menu to select observation referral.....	5-1
5-2	Pick list of SPs in current session.....	5-2
5-3	Observation referral in CV Fitness.....	5-2

TABLE OF CONTENTS (continued)

List of Exhibits (continued)

<u>Exhibit</u>		<u>Page</u>
5-4	Observation referral from CV Fitness in physician’s referral review box	5-3
6-1	Quality control reminder message box.....	6-6
6-2	Utilities menu to select quality control	6-7
6-3	Quality control log-on	6-7
6-4	Quality control daily checks (1)	6-8
6-5	Quality control daily checks (2)	6-9
6-6	Quality control daily checks (3)	6-10
6-7	Quality control daily checks (4)	6-11
6-8	Quality control weekly checks (1)	6-12
6-9	Quality control weekly checks (2)	6-13
6-10	Quality control weekly checks (3)	6-14
6-11	Quality control weekly checks (4)	6-15
6-12	Quality control ‘start of stand’ checks (1).....	6-16
6-13	Quality control ‘start of stand’ checks (2).....	6-16
6-14	Quality control ‘start of stand’ checks (3).....	6-17
6-15	Quality control incomplete entry.....	6-18
6-16	Technical Competency Evaluation Form	6-20
6-17	Checklist for medication review	6-24
6-18	Checklist for CV Fitness exclusion questions	6-25
6-19	Site Visit Report form	6-26

1. OVERVIEW OF CARDIOVASCULAR FITNESS

1.1 Cardiovascular Fitness in the National Health and Nutrition Examination Survey (NHANES)

Cardiovascular disease is currently the leading cause of death in the United States for both men and women. A sedentary lifestyle (along with elevated blood cholesterol levels, hypertension, and smoking) is a major modifiable risk factor for cardiovascular disease. According to current data however, only about 24 percent of American adults currently meet the physical activity recommendations outlined in the 1996 Surgeon General's Report.

In previous versions of NHANES, no direct cardiovascular fitness assessment components have been included in the data collection. A cardiovascular fitness component will be included with the current NHANES.

Cardiovascular fitness is defined as the body's ability to uptake, transport, and utilize oxygen. Other terminology that is sometimes used to describe cardiovascular fitness includes VO_2 max, maximal oxygen uptake, maximal oxygen consumption, and aerobic power.

A maximal treadmill test is considered to be the most valid method of measuring cardiovascular fitness. By collecting and analyzing expired air during the test, one can directly measure VO_2 max. This type of testing is done in a clinical setting and is extremely expensive and time consuming. A 12-lead electrocardiogram (ECG) is typically utilized during a maximal test. Thus, maximal treadmill testing can be used to diagnose stress induced cardiac arrhythmias and coronary heart disease. Due to time constraints, test setting, expense, and the large number of SPs in NHANES, it is neither feasible nor possible to utilize this method to evaluate cardiovascular fitness.

Submaximal treadmill testing will be utilized as a means to estimate cardiovascular fitness levels during NHANES. Based on variables including gender, age, BMI, and self-reported level of physical activity, SPs will be assigned to one of eight treadmill test protocols, each varying in difficulty. Each of the eight protocols was designed so that the SP could walk at all times. The goal of each protocol is to elicit a heart rate that is approximately 75 percent of the predicted maximum ($220 - \text{age}$) by the end of the test. Each protocol includes a 2-minute warm-up, followed by two 3-minute stages and a 2-minute cool-down period. By monitoring heart rate response to each of the 3-minute stages, VO_2 max can be predicted with a reasonable degree of accuracy. It is important to note that while the heart rate will be monitored during the treadmill test, there will be no electrocardiogram (ECG) on the screen. Therefore,

the NHANES submaximal treadmill test is intended only to estimate cardiovascular fitness levels. It cannot be used to diagnose or identify any cardiovascular disease.

SPs will be carefully screened prior to the treadmill test using questionnaires and measurements of resting heart rate and resting blood pressure. High risk (symptomatic, with disease, etc.) individuals will be excluded, so that only apparently healthy individuals will undergo treadmill testing. The combination of screening, monitoring, and the submaximal nature of the test will all contribute to the safety of this component.

1.2 Personnel

Maximizing safety and obtaining meaningful data are the two primary goals of the technician with regard to the submaximal treadmill test. The demeanor of the treadmill technician should be one of relaxed confidence so as to immediately put the SP at ease. Although there are time constraints to be aware of, treadmill testing procedures should not be rushed, and all procedures must be clearly explained prior to initiating the process. Technicians are expected to pay close attention to the data and the SP at all times, and also be sensitive to the confidentiality of the data. Knowledge of normal vs. abnormal responses to exercise is essential, as well as knowing when a test should be terminated early. Since the submaximal treadmill test is a fitness test (not a diagnostic test), and since technicians are not licensed physicians, no diagnoses will be attempted or made by technicians at any time. Although submaximal treadmill testing is very safe, all technicians must be CPR certified.

1.3 Integrated Survey Information System (ISIS)

The Integrated Survey Information System (ISIS) is a computer-based infrastructure designed to support all survey operations including sample management, data collection, data editing, quality control, analysis, and delivery of NHANES data. Each component in NHANES such as Cardiovascular (CV) Fitness has a computer application for direct data entry. Data collected in the CV Fitness room of the mobile examination center is directly entered in the ISIS system computers. In addition, data from biomedical equipment such as the blood pressure monitor in CV Fitness is directly downloaded to the ISIS system where it is displayed on the computer screen and stored in the system database. For more information about the NHANES and the ISIS system, see the ISIS Manuals and Presentations/ISIS Overview Presentation on the IntraWeb.

2. EQUIPMENT/SUPPLIES/MATERIALS

Below you will find descriptions of the equipment necessary to complete the Cardiovascular Fitness component of the National Health and Nutrition Examination Survey (NHANES).

2.1 Quinton MedTrack ST65 Treadmill with Hyperdrive

The Quinton MedTrack treadmill is a heavy-duty treadmill designed specifically for stress test applications. The Quinton MedTrack Programmable controller is required to operate the treadmill. See Figure 2-1.



Figure 2-1. Quinton MedTrack ST65 Treadmill

2.1.1 Maintenance of the Treadmill

Routine maintenance for the Quinton MedTrack Treadmill involves:

- Confirming that the cable from the treadmill is securely connected to the back of the programmable treadmill controller.
- Daily visual inspection of walking belt for wear and to confirm that the belt is centered on rollers. Always check for potential hazards on treadmill.
- Cleaning surfaces after each test with disinfectant and towel. Make sure the deck does not get wet.
- A weekly check of speed and grade to ensure accuracy of treadmill calibration (see section 2.1.2 on calibration).
- Weekly inspections of treadmill belt tracking and tension (see instructions 2.1.1.1 and 2.1.1.2 below).

2.1.1.1 Belt Tension - Weekly Check

The belt tension should be adjusted whenever the belt slips or moves unsteadily during operation. Notify the Chief Technician if the belt tension needs adjustment. Belt tension is adjusted as follows.

All technicians:

- Use the controls on the controller to start the treadmill at minimum speed and grade, then increase speed to 2.5 mph.
- Hold onto handrail for balance and walk heavily on the treadmill by marching flat-footed. If the belt hesitates or lags noticeably, the Chief Technician will tighten the belt as described in the following steps.

Chief Technician:

- Locate the two tension bolts at the rear of the treadmill.
- Turn both tension bolts clockwise $\frac{1}{4}$ turn. Test the belt tension. Repeat if necessary until the belt runs smoothly without slipping. If more than three adjustments are necessary, a service call may be necessary.
- Stop the treadmill.
- Check the position of the rear roller guards and adjust if necessary (see section on rear roller guards below).

2.1.1.2 Belt Tracking – Weekly Check

All technicians:

- Notify the Chief Technician whenever the belt moves to one side or the other (off center). The Chief Technician will adjust the belt tracking as follows.

Chief Technician:

- **Stay off the belt when adjusting tracking.**
- Use the controls on the controller to start the treadmill at minimum speed and grade.
- Increase speed to 3.7 mph to make adjustments.
- If the belt moves to the right, turn the right tension bolt a quarter turn clockwise.
- If the belt moves to the left, turn the left tension bolt a quarter turn clockwise.
- After making an initial adjustment, run the treadmill for several minutes and observe how the belt tracks. Adjustments to belt tracking take several minutes to become apparent. If the belt continues to move off center, adjust accordingly until it is properly centered. **Do not tighten the belt too much as this will cause accelerated wear on the roller bearings.**
- Stop the treadmill.
- Check the position of the roller guards (see section 2.1.1.3 on rear roller guards).

2.1.1.3 Rear Roller Guards - Weekly Check

All technicians:

- See Figure 2-2. Check the position of the rear roller guards each time you readjust the belt tracking or belt tension. The roller guards are located at the rear of the treadmill between the rear roller and the deck.
- If the gap between the roller and the guard exceeds 3/8 of an inch, notify the Chief Technician.

Chief Technician:

- Adjust the guards when the gap between the roller and the guard exceeds 3/8 of an inch.
- Disconnect the treadmill power cord from the receptacle (outlet).

- Loosen, but do not remove, the two screws attaching the rear roller guard to the deck (hold nuts underside with your finger).
- Slide the rear roller guard towards the rear until the gap between the two is approximately 0.1 inch (the thickness of two quarters).
- Tighten the two mounting screws on the rear roller guards.
- Connect the power cord.

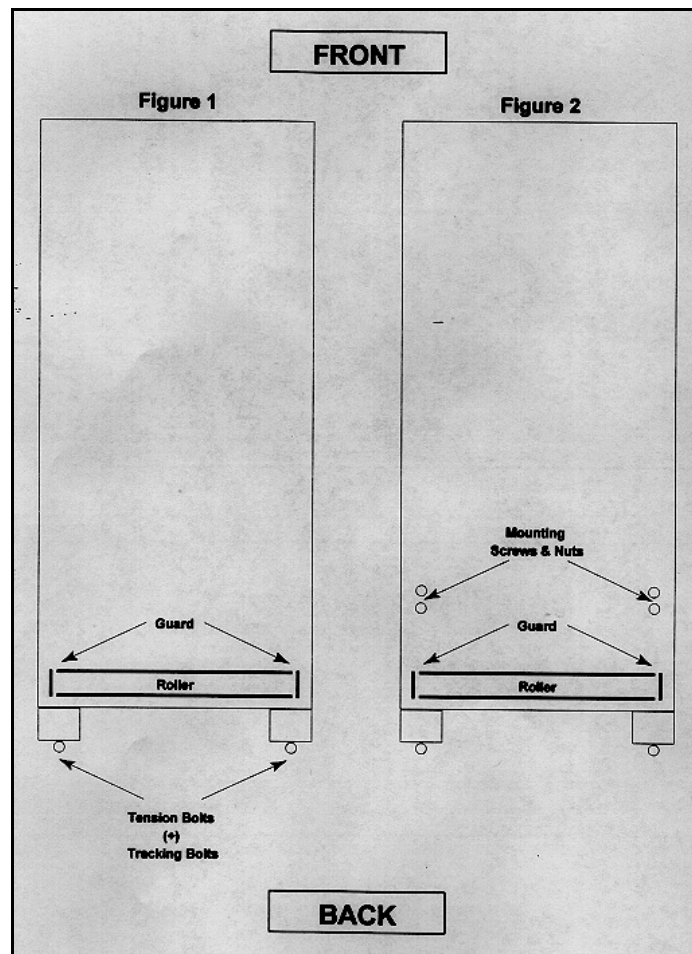


Figure 2-2. Treadmill calibration

2.1.2 Calibration of the Treadmill

Treadmill calibration should be checked weekly to ensure accuracy of testing results. However, if the treadmill is not within calibration specifications, a service technician should be consulted.

2.1.2.1 Speed – Weekly Check

NOTE: This calibration is to confirm the validity of the belt speed. Section 2.2.2 deals with confirming that speed and grade settings on the treadmill controller are programmed correctly.

- Determination of belt length is essential to calculate belt speed. Use two pieces of masking tape to mark one spot on the belt and the second on the edge of the treadmill next to the belt.
- Measure the length of the belt in feet. (Belt length should be determined following any adjustment to the belt.)
- Make sure the belt is free of objects.
- Use the controller to start the treadmill and increase the speed to 2.5 mph.
- Set a timer for 3 minutes.
- Start the timer and begin counting the number of times the two pieces of tape meet (starting with zero).
- Divide the number of revolutions by the time (3 minutes) to obtain revolutions per minute.
- Multiply the revolutions per min. by the belt length and divide the product by 88 ft/min to obtain belt speed.
- **$$[\# \text{ revs/min.} * \text{ belt length}) / 88 \text{ ft/min} = \text{belt speed (mph)}]$$**
- If speed readout on controller exceeds ± 0.2 mph of actual speed, repeat the previous steps.
- See Section 6.4.5 for data entry of these results in the equipment Quality Control screens.
- If treadmill still exceeds range, report this to the MEC manager. The MEC manager will call a service representative.

2.1.2.2 Elevation - Weekly Check

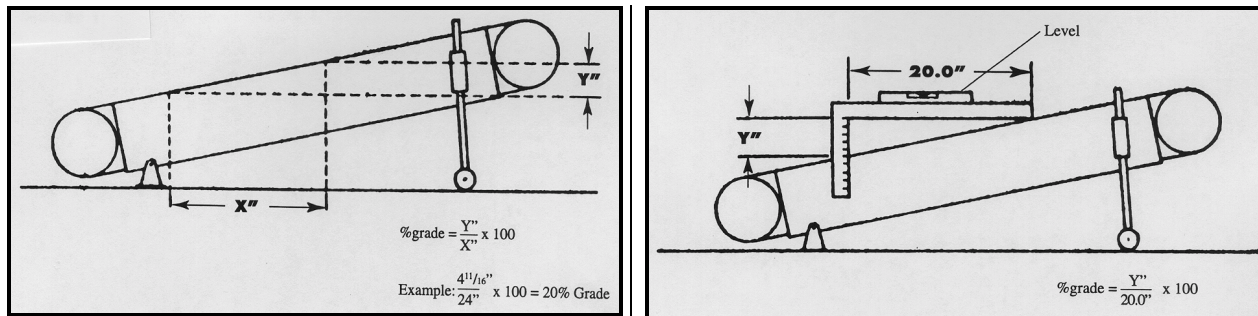


Figure 2-3. Calibration of elevation

All technicians:

- Check to make sure the treadmill is level prior to calibration of elevation using a level. If the treadmill is not level, make adjustments to level treadmill before continuing.
- Raise the treadmill elevation through its entire range and return to zero. This will allow the treadmill to perform a self-check.
- Raise treadmill to 20 percent grade.
- Using the Carpenters Square and a level, measure the change in grade from low to high (rise over run). The change in grade should equal 4 12/16 inches for every 24 inches of length (or 2 2/8 inches for every 11 inches of length. See Figure 2-3.
- If the elevation error exceeds $\pm 1\%$ grade, repeat the previous steps.
- If the elevation still exceeds $\pm 1\%$ grade in error, report this to the MEC manager. The MEC manager will call a service representative.
- See Section 6.4.5 for data entry of these results in the equipment Quality Control screens.

2.2 Quinton MedTrack ST Programmable Controller

The MedTrack ST Programmable Controller is a computerized treadmill controller designed to regulate speed and grade for Quinton treadmills. See Figure 2-4. Eight different protocols have been programmed into the controller for NHANES. Each protocol is programmed such that the actual protocol can be adjusted up or down one level depending on the sample person's (SP's) heart rate (HR) response and set criteria. Noted exceptions are Protocol 1 which can only increase to Protocol 2, and Protocol 8 which can only decrease to Protocol 7. Maximal treadmill speed is set to 3.7 mph to ensure SP safety.

The dimensions of the controller are 22.75" x 8.50" x 2.5" (57.79 x 21.59 x 6.35 cm) and the weight is 11 lbs. (5 kg).



Figure 2-4. Quinton MedTrack ST Programmable Controller

2.2.1 Maintenance of the Programmable Controller

The MedTrack ST Controller requires little maintenance. Wipe surfaces down weekly with a damp cloth to remove dirt and dust from the screen and inspect cables for unusual wear. Routine maintenance of the Colin STBP-780 BP monitor involves:

- Checking that the cable from the treadmill is connected **securely** to the back of the treadmill controller.

2.2.2 Calibration of the Programmable Controller

The Quinton MedTrack Controller does not require any specific calibration. The speed and grade for each of the eight exercise protocols on the treadmill controller will be confirmed on a weekly basis. The **Chief Technician** will be responsible for this procedure.

2.2.2.1 Confirming the Speed and Grade Settings on the Treadmill Controller

NOTE: This calibration check is to determine that the speed and grade settings on the treadmill controller are programmed correctly. Section 2.1.2.1 deals with confirming the validity of the belt speed.

Chief Technician:

- Get Speed and Grade Table by Protocol and Stage (See Appendix J). Use this table to check the settings on the treadmill controller.
- Press the Power button on the treadmill controller to turn the power on.
- To enter programming mode, simultaneously press Stop Belt, Cool Down, and Start Course, then release them.
- The message 'Program' scrolls across the display.
- Press Enter.
- Press 'Input +' to go to Protocol 1.
- Press 'Next Stage'.
- The yellow dot below Warm-up (current) will be flashing to indicate the controller is now in Protocol 1, current warm-up (the assigned exercise protocol).
- Check the speed and grade on the treadmill controller with the speed and grade on the Speed & Grade Table by Protocol and Stage to confirm that the treadmill settings are the same as listed in the table.
- *NOTE: If the numbers are not the same, do not attempt to change the settings. Report this immediately to the component specialist at the Home Office. The component specialist will confirm the settings and give instructions on how to change the settings if necessary.*
- Press 'Next Stage' to advance to Stage 1, Current. Check the speed and grade on the treadmill controller with the speed and grade on the table. (As noted above, the yellow dot will be flashing under the stage that the treadmill is currently set. Each time 'Next Stage' is pressed, the treadmill controller moves to the next stage and the yellow dot flashes under the new stage.)

- Press 'Next Stage' to advance to Stage 2, Current. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Recovery, Current. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Stage 1, Lower. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Stage 2, Lower. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Recovery, Lower. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Stage 1, Higher. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Stage 2, Higher. Confirm the settings for speed and grade.
- Press 'Next Stage' to advance to Recovery, Higher. Confirm the settings for speed and grade.
- Repeat the last 6 steps to check the second set of settings for the next lower and next higher protocols.
- After the second 'Recovery, Higher' is confirmed, press 'Enter'.
- Press 'Input +' to advance to Protocol 2.
- Press 'Next Stage' to go to 'Warm-up, Current' for Protocol 2. (yellow dot will be flashing under this stage).
- Continue this sequence as listed above to confirm speed and grade for all protocols 1 through 8. There are 16 settings to check for each protocol.
- When all settings are confirmed with the settings on the table, press 'Enter' to save the settings.
- Simultaneously press Stop Belt, Cool Down, and Start Course, then release them to exit the programming mode.
- Check that this was completed in the weekly Quality Control Log on the Integrated Survey Information System (ISIS). (Section 1.3 provides a brief explanation of ISIS, the data entry system that is used for this study.)

2.3 Colin STBP-780 Automatic Blood Pressure Monitor

The Colin STBP-780 is an automated electronic heart rate and blood pressure (BP) monitor capable of accurate readings at rest and during exercise. See Figure 2-5. The unit assesses heart rate via wires connected to four electrodes placed on the thorax and abdomen. Blood pressure is assessed during deflation of the cuff via two microphones in the cuff. The front display provides clear, easy to read measurements of heart rate, systolic and diastolic BP, elapsed time, and error messages. When turning the system on, a Self-Check is performed automatically. If an error is detected, see Section 2.9.2.

The field for the systolic blood pressure at the top left corner of the blood pressure monitor displays an estimate of systolic BP during deflation. This number flashes as the measurement is being taken. Once the systolic and diastolic blood pressure and the heart rate are measured, all three values are displayed on the monitor and downloaded to the computer screen and system database during the test.



Figure 2-5. Colin STBP-780 Automatic Blood Pressure Monitor

2.3.1 Maintenance of the Automatic Blood Pressure Monitor

Routine maintenance of the Colin STBP-780 BP monitor involves:

- Checking that hoses and microphone wires are securely connected to the jacks.
- Checking that the computer cable is securely connected to the back of the BP monitor.
- Wiping surfaces clean daily with a clean dry cloth. If necessary, use a mild soap and water solution. DO NOT use alcohol to clean surfaces.
- Clean BP cuffs after each use with disinfectant or alcohol and allow drying.
- Ensure BP cuffs are stored properly. DO NOT store the cuff such that stress is placed on the connection of the microphone to the cuff.

2.3.2 Calibration of the Automatic Blood Pressure Monitor

The calibration of the pressure amplifier inside the STBP monitor will be checked each week with a mercury manometer. The reading on the STBP monitor should be within +/- 5mmHg from the reading on the mercury manometer. If the differences between the readings differ by more than 5 mm Hg, the STBP monitor will need to be returned to the manufacturer for calibration.

- Weekly calibration check procedures are outlined below:
 - Turn off the power switch for the STBP monitor.
 - Get portable mercury manometer, and Y-connector with rubber bulb from the room supply box.
 - Set up the mercury manometer on the computer cart.
 - Connect one end of the Y-connector to the air tubing on the mercury manometer.
 - Connect second end of the Y-connector to the inflation outlet at the back of the STBP monitor.
- Use the following steps to check the calibration of the STBP with the mercury manometer. The calibration will be checked at these five ranges:
 1. 220-250
 2. 180-220
 3. 120-180

4. 70-120

5. 0-70

- While pressing the END switch, turn POWER switch on.
- Press MANUAL switch to set unit in the MANUAL mode. The word “MAN” is displayed on the LED of the ELAPSED TIME indicator.
- While pressing the UP switch, press the START switch.
- Raise the pressure in the mercury manometer to 250 mmHg. Tighten the screw to maintain the mercury at a constant level.
- In about 1 second, the pressure values are displayed on the STBP monitor.
- Compare these values with the value on the mercury manometer.
- The numbers on the STBP LED display should be within +/- 5 of the reading at eye level on the mercury manometer.
- Loosen the screw and slowly deflate the cuff to a reading on the manometer in the next range.
- Press the RESET switch to deflate pressure.
- Press the MANUAL switch to release the MANUAL mode.
- Repeat these steps for each of the ranges listed above.
- If the number is within +/- 5, enter a check beside this item on the Quality Control Screen. See Figures 6-9 to 6-11.

If the BP calibration checks are not within the acceptable range (more than 5):

- Check connections with Y-connector, air hoses, and cuff for leaks.
- If the **BP calibration checks are not within the acceptable range**, enter the difference between the mercury manometer reading and the reading on the STBP monitor in the ‘Result’ field. (Example: Difference = 8mmHg)
- Enter the reading on the manometer and the reading on the STBP monitor in the ‘Comments’ field. (Example: Manometer = 110; STBP=118).
- Notify the MEC manager and Chief Technician. Make a note of this in the ‘Comments’ section. (Example: Reported to the MEC manager.)
- If the differences between the readings differ by more than 5 mm Hg, the STBP monitor will need to be returned to the manufacturer for calibration.
- **NOTE: Notify the Home Office if the difference is greater than 3mmHg.** This is an ‘early warning alert.’ Continue to use the monitor for pressure measurement. This

allows time to have the unit replaced and serviced in a timely manner without compromising data.

2.3.3 Mercury Manometer (Used for Calibration)

The mercury manometer is calibrated when it is manufactured and, once calibrated, re-calibration is unnecessary. However, regular inspection is necessary to eliminate conditions that could cause the blood pressure measurement to be read as erroneously high or low.

- Check to see that the level of mercury in the glass tube is zero. The level of the mercury in the calibrated glass tube should always be at the zero line when the manometer is on a level surface with the inflation system disconnected. Tip the manometer gently to the right and then to the erect position. If the top of the mercury column does not return to zero, replace the equipment.
- Check the shape of the meniscus - it should be a smooth, well-defined curve.
- Check that the mercury rises easily in the tubing and that the mercury column does not bounce noticeably when the valve is closed.
- Check for cracks in the glass tube to make sure it is securely in place.
- Check cap at top of the calibrated glass tube to make sure it is securely in place.
- Check cuff, pressure bulb, and manometer for cracks or tears.
- Check pressure control valve for sticks or leaks.
- Use the following procedure to check for air leaks:
 - Put the male luer plug in the female luer adapter on the Y-connector.
 - Connect the male connector on the manometer to the female connector on the Y-connector.
 - Inflate to 250 mmHg.
 - Open valve and deflate to 200 mmHg and close valve.
 - Wait for 10 seconds. If mercury column falls more than 10 mm Hg, there is an air leak in the system.
- Never attempt to repair the equipment yourself.
- Contact the MEC manager to have equipment sent for repair or replacement. The manometer will be sent to the Home Office and from there will be sent to Baum for repair or replacement.

2.3.4 Mercury Spills and Leaks

A mercury manometer will be used for the calibration check. Mercury is a metallic substance which gives off a toxic vapor when exposed to the atmosphere. Temperature, ventilation, and sunlight affect the level of the vapor's concentration. Mercury vapors will permeate the skin surface and are poisonous when inhaled.

Check the cap at the top of the calibrated glass tubing. If it is not securely closed, the mercury could leak out. Loss of air and mercury will occur if the glass tube is broken. Care should be taken in handling the manometer to prevent this. If the tube appears cracked, check for any spilled mercury and replace the equipment if necessary.

The following procedure is recommended for handling spilled or leaking mercury:

- Contact the MEC manager immediately to report the incident and receive instructions. The Mercury Spill Kit contains all the materials needed for clean up.
- The MEC manager will determine what action should be taken. If necessary, he/she may call the local contact to report such spills. This agency may want to make a followup visit to the facility to check for levels of mercury in the exposed area.
- Do not touch the mercury with your bare hands or attempt to vacuum or clean up the spill yourself.
- The procedures for cleaning up the mercury are different depending on the surface where the spill is detected.
- If the spill is small you may be directed to use the following procedures.
- For leakage or spillage on a hard surface (e.g., hard floor, table, etc.):
 - Get the Mercury Spill Kit from the physician's examination room.
 - Put on the green protective gloves
 - Use the mercury (HG Absorb) sponge.
 - Remove the sponge from the plastic zip closable bag.
 - Dampen the sponge with water.
 - Wipe the area contaminated with mercury. Do this slowly to allow for complete absorption of all free mercury. The chemical layer (Hg Absorb Powder) on the sponge will absorb the mercury droplets.

- After finishing with the Hg Absorb sponge, place it back into its plastic zip closable bag.
- Place the manometer case and mercury sponge into separate plastic bags and secure with tape. Return to the MEC manager.
- Obtain replacement equipment.
- For leakage and/or spillage on a soft or absorbent surface (e.g., rug, etc.):
 - Get the Mercury Spill Kit from the physician's examination room.
 - Put on the green protective gloves.
 - Turn off heat or air conditioning to avoid spreading vapors to other rooms in the MEC.
 - Open a window if possible to ensure adequate ventilation to the outside for the room where the spill occurred.
 - Use the mercury (Hg Absorb) sponge, if possible. See directions above.
 - Wet the MERCSORB™ powder with water. Mercury will react with powder, forming a metal/mercury amalgam.
 - Wipe or sweep wetted powder over all cracks and hard to reach locations for maximum pickup.
 - Pick up amalgam by sponging or by using the small sweep and dustpan provided in the Mercury Spill Kit.
 - Obtain replacement equipment.

2.4 SECA Floor Scale Model 815

The SECA model 815 floor scale is used to measure weight in this component. It has a digital display indicator head fitted at the back of the platform with the connecting cable stored in the compartment underneath the head. The scale is a load cell model with a weighing range of 0-136 kilograms (kg). It operates with a standard 9 V alkaline battery. In order to open the battery compartment, open the battery lid underneath the head. Connect the battery terminals, then insert the battery and close the cover.

The reading of the scale is accurate to 0.25 kg over the entire weighing range. If two values are displayed alternately in the 0.25 kg range, then the exact weight is between the two values and the

intermediate weight should be entered. For example, if the scale is reading between 70 and 70.25 kg, 70.13 kg should be entered in the weight field. Other instructions include:

- If “SUP” appears on the display, it signals that the scale is overloaded and cannot be used for this measurement.
- If no weight display appears when an SP is on the scale, have the SP step off the scale, press the **ON** button and wait for the display 0.0.
- If - - - appears on the display, press the **ON** button and wait for the display 0.0.
- If “Err” appears on the display, have the SP step off the scale, press the **ON** button, and wait for display 0.0.
- If “bAt” appears on the display, change the battery.

2.4.1 Maintenance and Calibration of the Floor Scales

The SECA scale should be calibrated at the beginning of each stand and weekly.

Start of Stand. At the beginning of each stand, check to see that the scale is operating properly with the batteries in place, and then standardize the portable scale to the Toledo scale on the MEC using the calibration weights found in the body measurements room of the MEC. Procedures for standardizing the SECA scale are given below. Any problems with the scale should be reported to the MEC manager immediately.

To standardize each SECA scale:

- Perform the standardization procedures in the body measures room of the MEC.
- Press the ON button of the SECA scale and wait for the display 0.0.
- Place the calibration weights on the scale in increments of 50 pounds, starting with 50 and continuing to 250.
- If the scale is off by 1 pound or more with the 5 weights placed on it, inform the MEC manager and use the other scale.

Weekly. Prior to each examination, place a 50-pound weight on the scale to roughly check the accuracy. Note the check in the QC log. If there is any reason to believe that the scales are not accurate, do a full calibration as described above. Record the calibration in the QC log. If the scale is off by 1 pound or more, inform the MEC manager.

2.5 SECA Portable Stadiometer Model 222

The SECA model 222 portable stadiometer is attached to the wall in the CV Fitness room. The detachable floor piece has been modified to be large enough for the SP to stand on and the plastic headpiece must be lifted and placed on the SP's head. The two vertical measuring bars have both centimeters and inches imprinted on them. Measurements are read in centimeters at the orange plastic line where the two measuring bars connect.

2.5.1 Maintenance and Calibration of the Stadiometer

Calibrate the stadiometer at the beginning of each stand and weekly during the stand. The procedures are as follows:

Start of Stand and weekly: Make sure the measuring column is pushed completely into the floor piece by checking that it is inserted up to the etched line.

- a. Place the calibration rod on the floor of the stadiometer.
- b. Place the horizontal bar of the stadiometer firmly against the top of the calibration rod. The stadiometer should read 100 cm.
- c. Document the performance of the calibration procedure in the QC log.
- d. If the reading is off by > 1 cm, inform the MEC manager. Use the stadiometer until a backup is sent.

2.6 Description of Supplies

2.6.1 ECG Electrodes/Alcohol Prep Pads/Abrasive Pads

Foam stress test V electrodes (1.5 “) from Vermont Medical will be used with the ECG chest leads to monitor the heart rate during the exam. Alcohol prep pads will be used to prepare the site for attachment of electrodes. Marquette abrasive pads will be used to remove old cells from the skin surface to improve adherence of the chest electrodes.

2.6.2 Running Shoes

The criteria that are important to consider when assessing whether or not the SP’s footwear is appropriate for the treadmill test are the presence of arch supports, heel cushions, and traction. If the SP does not bring footwear that is appropriate and safe for this test, acceptable footwear will be provided. A supply of women’s and men’s running shoes are kept in a shoe rack in the cardiovascular fitness room in each mobile examination center (MEC).

Examples of appropriate footwear are running shoes, casual or leisure style shoes, and modern athletic hiking boots. Some examples of inappropriate footwear include pumps, sandals, women’s or men’s dress shoes, hard sole boots or shoes (work boots, traditional hiking boots), any boot or shoe that has a heel for aesthetic purposes, and any shoe with cleats (soccer, golf, baseball, etc.).

2.6.3 Mesh Vests

If a SP has chest hair that covers the sites for electrode placement, the area **will not** be shaved. See Section 3.8.5. In situations where the chest hair makes adherence of the electrode to the skin difficult, a mesh vest will be used to help secure the electrodes and stabilize the leads.

2.7 Contact Information for Equipment Repair or Parts Replacement

Quinton Treadmill and Programmable Controller

Quinton Instrument Company
3303 Monte Villa Parkway, Bothell, WA 98021-8906
Phone: 206-402-2000
Fax: 206-402-2001

STBP Monitor and BP Cuffs

Colin Medical Instruments Corp.
5850 Farinon Drive, University Technology Park
San Antonio, Texas 78249
Phone: 800-829-6427
Fax: 210-696-8808

Mercury Manometer:

W. A. Baum Co. Inc.
Copiague, New York 11726
Phone: 516-226-3940
Fax: 516-226-3969

2.8 Equipment Manuals

The following equipment manuals are kept in the CV Fitness room.

- Medtrack ST 55 & ST 65 Treadmills Operator Manuals (Quinton)
- Medtrack ST 55 & ST 65 Treadmills Service Manuals (Quinton)
- Medtrack ST Programmable Controller Operator and Service Manual (Quinton)
- Stress Test Blood Pressure Monitor Model STBP 780 Operation & Service Manual (Colin Medical Instruments Corporation)

2.9 Troubleshooting

2.9.1 Quinton MedTrack ST 65 Treadmill and Controller

The following table (2-1) may be used to troubleshoot some common problems with the treadmill and programmable controller.

Table 2-1. Troubleshooting problems with the Quinton Med Track ST65 Treadmill and Controller

PROBLEM	POSSIBLE CAUSE	REMEDY
Treadmill will not run.	Treadmill and/or monitor power cord is not plugged in	Plug in power cords.
	Monitor power off	Turn on controller.
	Controller cable disconnected or damaged	Plug in and secure controller cable. Replace if damaged.
	Circuit breaker for treadmill room reads "OFF"	Set the circuit breaker switch on the treadmill to "ON."
	Emergency stop button activated	Turn emergency stop button counter-clockwise to deactivate.
	Internal problem	Call authorized service rep.
Walking belt too far left or right.	Improper belt tracking	Adjust tracking as described in Section 2.1.1.2. If problem persists, call authorized service rep.
Walking belt slips, but front roller turns.	Improper belt tension	Adjust belt tension as described in Section 2.1.1.1. If problem persists, contact authorized service rep.
Walk belt hesitates; adjusting belt tension is ineffective.	Internal drive belt slipping	Service required. Call authorized service rep.
Treadmill will not change grade.	Excess weight on treadmill	Max load is 400 lbs.
	Internal problem	Service required. Call authorized service rep.
Circuit breaker trips during normal operation.	Power fault	Service required. Call authorized service rep.
Treadmill does not run; Exxx appears on display.	Problem with electronic circuitry	Record error number and call service rep.
	Internal problem	Service required. Call authorized service rep.
Other error messages appear on controller.	Varied	Record error message and call service rep.
Treadmill will not reach maxspeed.	Speed limit control is activated	Do not change!

2.9.2 Colin STBP-780 Blood Pressure Monitor

The following table (2-2) may be used to troubleshoot problems with the STBP monitor.

Table 2-2. Troubleshooting problems with the STBP monitor

PROBLEM			PROBABLE CAUSE	REMEDY
DISPLAY	SYMPTOM	VOICE		
Check ECG*	ECG trigger cannot be detected for 5 seconds.	Please check ECG.	A. ECG connectors are disconnected. B. ECG lead-select is not proper.	Check the lead and electrode placement. Check the skin preparation. Call the Home Office if problem

			<p>C. The electrodes are not attached correctly and electrical noise is interfering.</p> <p>D. The ground wire is not attached correctly and electrical noise is interfering.</p> <p>E. Because of the extraordinarily strong ECG signal, ECG is detected as noise.</p>	not resolved.
Check cuff	The inflation exceeds 15 seconds.	Please check cuff.	<p>A. The capability of the pump-up motor has declined.</p> <p>B. Power supply voltage is fluctuating.</p> <p>C. Air is leaking.</p> <p>D. Cuff connector has become loose.</p>	Check the position of the cuff. Check the airhose and microphone connections. Call the Home Office if problem not resolved.
Check SYS	Third reinflation starts during REST MODE.	The safety valve is open.	<p>A. The microphone detects noise.</p> <p>B. There is line noise.</p>	Check the position of the microphone. Call the Home Office if problem not resolved.
<p>Check MIC*</p> <p>Check MIC.*</p> <p>Check MIC!*</p> <p>Check MIC?*</p>	The input signal from the microphones cannot be collected.	Please check microphone.	<p>A. The microphones are not properly located on the artery.</p> <p>B. The signal cannot be detected; cuff may be wrapped upside down.</p> <p>C. The microphone is inferior.</p> <p>D. Microphone cord is not connected or is broken.</p>	Check the position of the microphone. Call the Home Office if the problem is not resolved.

Table 2-2. Troubleshooting problems with the STBP monitor (continued)

PROBLEM			PROBABLE CAUSE	REMEDY
DISPLAY	SYMPTOM	VOICE		
FAIL	In the self-test, an abnormality has been detected.	Machine failure.	Hardware failure.	Service required. Call authorized service representative.
NOISE		None	The increase or decrease between 2 consecutive R wave readings exceeds 25%.	
BATT		None	The voltage of the memory back-up battery goes down below 2 volts.	Service required. Call authorized service representative.

3. PROTOCOL

3.1 Exclusions to the Cardiovascular (CV) Fitness Component

Sampled participants (SPs) are excluded from the Cardiovascular Fitness component based on medical conditions, medications, physical limitations, and edit limits on heart rate, blood pressure, and irregular beats.

3.1.1 Medical Conditions

SPs are asked questions in the Household Questionnaire about their past and current medical conditions. Specific medical conditions are reasons for exclusion from this component for reasons of safety and/or physical limitation. If these conditions are confirmed in the Household Questionnaire, the SP is excluded from the CV component before they enter the mobile examination center (MEC). The system will record the component status as 'Not Done' with the comment 'safety exclusion.' The SP will not need to go to the CV room in the MEC to have this exclusion recorded.

The medical conditions that are a cause for exclusion include history of previous myocardial infarction, coronary heart disease, stroke, emphysema, and certain physical function or development problems that would make it difficult or unsafe to complete this test. See Appendix B for a list of the questions and responses in the Household Interview Questionnaire that would lead to an exclusion to the CV Fitness component.


In addition, SPs who are greater than 12 weeks pregnant will also be excluded from the CV fitness component. A urine test for pregnancy will be completed and a self-report of pregnancy status will be obtained in the MEC. Once the pregnancy condition is reported by the SP or confirmed by the test result, the system will block the SP from the CV fitness exam and set the Component status to 'Not Done' with the comment 'safety exclusion.'

3.1.2 Medications

SPs are excluded from the CV Fitness component if they are taking medications such as calcium channel blockers, anti-arrhythmics, beta blockers, nitrates, nitroglycerin, and digitalis. See Appendix C for a detailed list of exclusion medications. This list is updated every 3 months if necessary to include recently released medications. There are two levels of exclusions for medications.

3.1.2.1 Medications Reported in the Household Interview Questionnaire

Exhibit 3-1. CV exclusions based on medications in the physician's exam

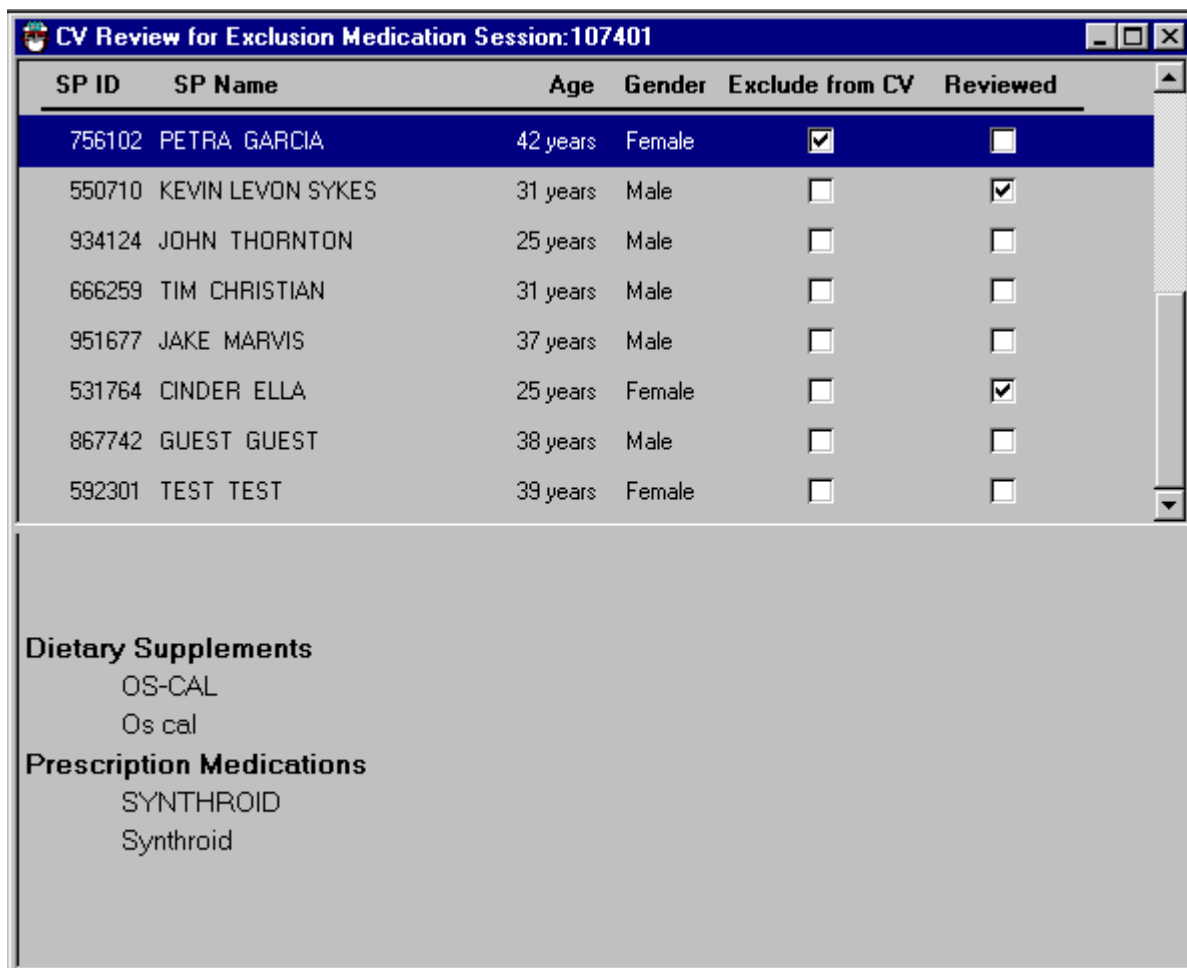

SP ID	SP Name	Age	Gender	Exclude from CV	Reviewed
756102	PETRA GARCIA	42 years	Female	<input checked="" type="checkbox"/>	<input type="checkbox"/>
550710	KEVIN LEVON SYKES	31 years	Male	<input type="checkbox"/>	<input checked="" type="checkbox"/>
934124	JOHN THORNTON	25 years	Male	<input type="checkbox"/>	<input type="checkbox"/>
666259	TIM CHRISTIAN	31 years	Male	<input type="checkbox"/>	<input type="checkbox"/>
951677	JAKE MARVIS	37 years	Male	<input type="checkbox"/>	<input type="checkbox"/>
531764	CINDER ELLA	25 years	Female	<input type="checkbox"/>	<input checked="" type="checkbox"/>
867742	GUEST GUEST	38 years	Male	<input type="checkbox"/>	<input type="checkbox"/>
592301	TEST TEST	39 years	Female	<input type="checkbox"/>	<input type="checkbox"/>

Dietary Supplements
OS-CAL
Os cal

Prescription Medications
SYNTHROID
Synthroid

At the beginning of each session, the MEC physician will review the list of medications the SP reported taking on the Household Interview Questionnaire (See Exhibit 3-1). The physician will review the medication for each SP eligible for the CV Fitness component and will check this with the current list of exclusion medications. If the SP reported taking any of the exclusion medications on the

Household Questionnaire, the SP will be excluded from the CV Fitness component. The physician will check this on the CV Exclusion box. The system will block the SP from CV Fitness and set the Component status for CV Fitness to 'Not Done' with the comment 'safety exclusion.' The SP will not need to answer the CV safety/exclusion questions in the physician exam and will not go to the CV Fitness room.

If the SP did not report taking any exclusion medications on the Household Interview Questionnaire, he/she will not be excluded at this point. An additional safety measure is in place to confirm that the SP has not started on any of the exclusion medications in the time between the Household Interview and the MEC exam visit.

3.1.2.2 Medications Reported in the Physician's Exam in the MEC

The physician will ask the SP about current medications as part of the CV safety/exclusion questions (see the following section). If the SP is taking any of the medications on the exclusion list, the system will exclude him/her from the CV Fitness component and will set the component status to 'Not Done' with the comment 'safety exclusion.' The SP will be blocked from the CV Fitness exam.

3.1.3 Safety/Exclusion Questions

Safety/Exclusion Questions will be asked in the Physician's Exam prior to the CV Fitness exam. This is a second level safety check to confirm that SPs do not have any of the exclusion medical conditions or are not taking exclusion medications. This will avoid including any SPs who may have been recently diagnosed with one of the exclusion medical conditions in the time period since the Household Interview. If the SP responds affirmatively to any of the medical exclusion questions, he/she will be immediately excluded. The remaining exclusion questions will not be answered.

See Appendix D for the CV Exclusion questions asked in the Physician's Exam.

3.1.4 Exclusions Based on Heart Rate and Blood Pressure

The resting heart rate and blood pressure measurements taken in the Physician's Exam will be used as the baseline measurements that are reviewed for exclusion from the treadmill test. SPs will be excluded from the CV Fitness component if one or more of the following criteria are present:

- Heart rate greater than 100 beats per minute;
- Systolic blood pressure greater than 180 mmHg;
- Diastolic blood pressure greater than 100 mmHg; and
- Irregular beats greater than 3 per minute.

If one or more of the above are present, the system will exclude the SP from the treadmill test. The CV Fitness Component Status will be set to 'Not Done' with the comment 'safety exclusion.'

3.2 Introduction to Component

It is important to log on to the computer data entry system as soon as you can after entering the examination room. However, you should provide the SPs with a brief explanation of the treadmill test as you are walking from the reception area to the treadmill test room. Many sampled participants may have never walked on a treadmill and may be a little apprehensive about this component. Make frequent eye contact with the SP when you are explaining the test and make sure he/she is listening and understands what you are saying. Below is an example of the introduction you may use. This is not a standard script so you may use your own words.

Suggested Introduction to the Component (Not a standard script): English Version

"In this exam we want to get some information about your cardiovascular fitness or the fitness of your heart. We will be asking you to walk on the treadmill while we monitor your heart rate and blood pressure. I will explain this in more detail in a few minutes but first I want to ask you a few questions. Please come over here near the computer while I scan your identification number into the system."

Suggested Introduction to the Component (Not a Standard Script): Spanish Version

“En este examen queremos obtener alguna información acerca de su condición cardiovascular o la condición del corazón. Le pediremos que camine sobre la máquina de andar mientras le monitoreamos la frecuencia cardíaca (los latidos del corazón) y la presión de sangre. Explicaré esto con más detalle en unos pocos minutos pero primero quiero hacerle algunas preguntas. Por favor venga aquí cerca de la computadora mientras anoto con el escáner su número de identificación en el sistema.”

3.3 Measurement of Height and Weight to Determine Body Mass Index

Once you are in the CV Fitness room, log on to the Integrated Survey Information System (ISIS) system and confirm the SP's identity. (Section 1.3 provides a brief explanation of ISIS, the data entry system that is used for this study.)

After you have logged on to the system, check the first screen to determine if the height and weight information is displayed on the screen. If the SP was in the anthropometry component prior to this test, the height and weight will be uploaded and displayed on the ISIS screen. If the information is not displayed, you will need to measure the SP's height and weight using the stadiometer and floor scales in the room. The system will use the height and weight measurements to calculate the body mass index and use it in the calculation to determine predicted VO₂ max.

3.3.1 Height

Follow these steps to take the SP's height:

1. Ask the SP to remove his/her shoes if necessary.
2. Place the stadiometer a few inches away from the wall. Check to be sure the measurement column is completely inserted into the floorpiece. Pull the sliding top bar section up and open the head piece to allow the SP to step under the head piece.

3. Ask the SP to stand erect on the floor board with his/her back to the vertical piece of the stadiometer and the wall. The head, shoulder blades, buttocks, and heels should be touching the measurement surface but the SP should not be leaning against the stadiometer.
4. Ask the SP to evenly distribute weight on both feet. The heels are placed together with the feet pointed slightly outward at a 60 degree angle. The arms hang freely by the sides of the trunk with palms facing the thighs. (See Figure 3-1.)
5. Position the head in the Frankfort horizontal plane. The head is in the Frankfort plane when the horizontal line from the ear canal to the lower border of the orbit of the eye is parallel to the floor and perpendicular to the vertical backboard. Many people will assume this position naturally, but for some it may be necessary to make a minor adjustment. If required, gently tilt the head up or down until proper alignment is achieved with the eyes looking straight ahead. Once correctly positioned, ask the SP to inhale deeply and to stand fully erect without altering the position of the heels.
6. Lower the headpiece snugly to the crown of the head with sufficient pressure to compress the hair.
7. Hold the top sliding bar in place at the junction and ask the SP to step out away from the stadiometer.
8. Record the measurement in centimeters (measurements printed on right side of bar) at the orange line on the measuring bar.
9. If an SP has a hairstyle that will interfere with measurement of height or is wearing a hairpiece that he/she declines to remove for the measure, the examiner will measure the hairpiece with a ruler and subtract this from the height measurement.

3.3.2 Weight

Follow these steps to take the SP's weight:

1. Make sure the scale weighs in kilograms by checking the switch on the underside of the digital display.
2. Place the scale on the floor.
3. Switch on the scale by pressing the **ON** button; 888.8 will appear on the digital display.
4. Have the SP remove his/her shoes and any heavy outer clothing such as sweaters, jackets, etc.
5. As soon as the display switches to 0.0, have the SP step on the scale with his/her feet positioned in the center.
6. Ask the SP to stand straight and remain still.

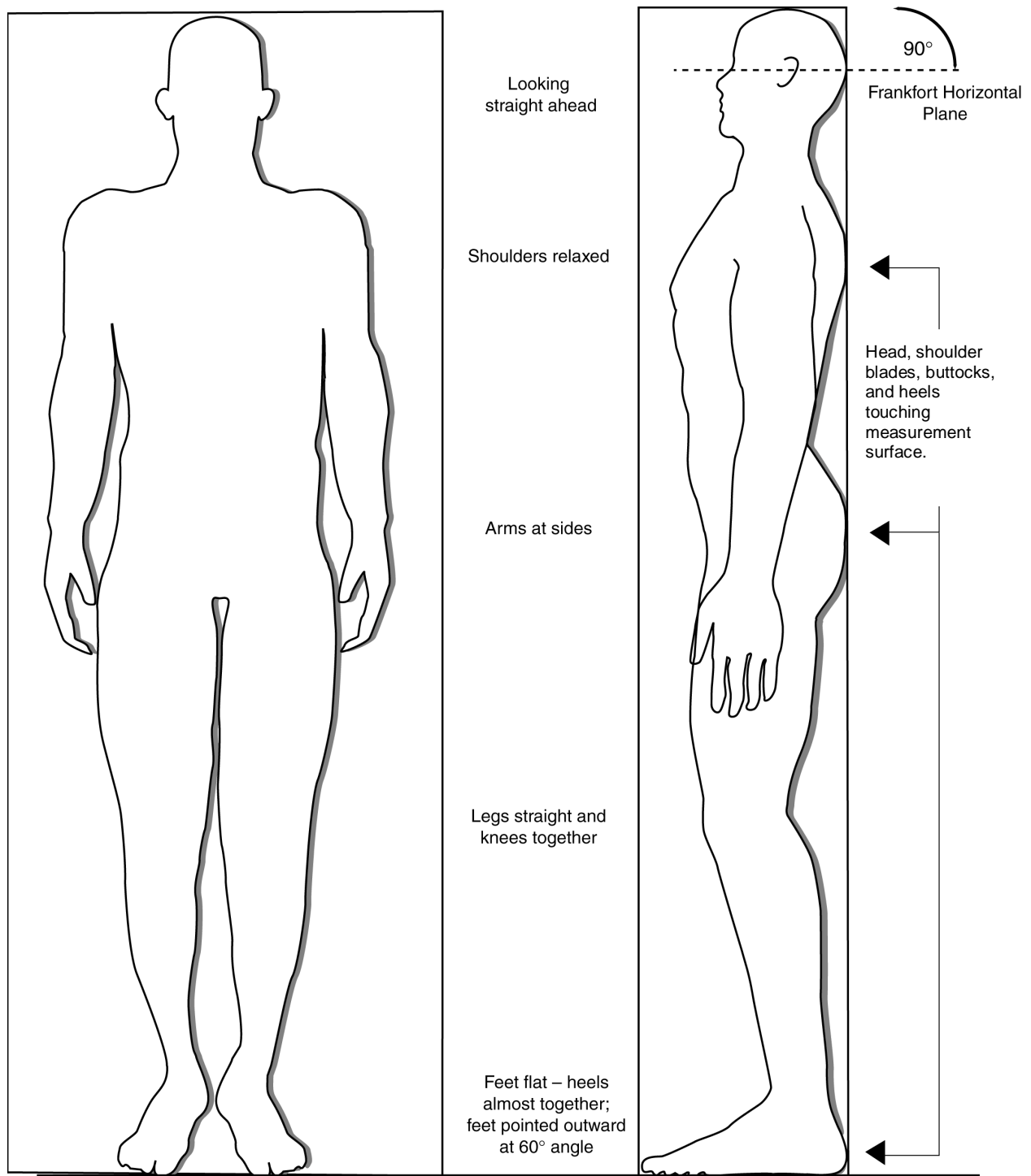


Figure 3-1. Correct body placement for height measurement

7. Wait about 4 seconds for the weight to display on the digital readout.
8. Record the weight in kilograms in the weight field.
9. Ask the SP to step off the scale. The scale switches off automatically after 30 seconds.

If the SP's weight is in excess of 136 kilograms, follow these steps for using 2 scales:

1. Perform steps 1 through 4 (as listed above).
2. Place the 2 scales side by side, but not touching each other. The scales do not need to be parallel to each other. They can be angled to facilitate SP positioning.
3. Instruct the SP to place one foot on the center of each scale platform.
4. Ask the SP to stand straight and remain still as described above in step #6.
5. Note the 2 weights given, following the same procedure as if one scale were used.
6. Add the weights displayed on each scale and enter the sum of the two weights in the weight field.
7. Perform step #9 as listed above.
8. Note: If the SP weighs more than 350 pounds (158.8 kilograms), he/she is excluded from this component due to weight limitation on the equipment.

3.4 Checking SP's Footwear

SPs are asked to bring comfortable walking shoes to the MEC; checking the SP's footwear is an important part of the preparations for the treadmill test. See Section 2.6.2 for a description of appropriate and inappropriate footwear for this test. If the SP does not have appropriate footwear, select a pair of running shoes from the stock supply in the exam room. Do not begin the physical activity readiness questions while the SP is putting on his/her shoes or engaged in any other activity. You should have the SP's full attention when you are asking these questions to ensure that he/she is listening and understands what you are asking.

3.5 Determination of Physical Activity Readiness Code

Each SP will be assigned an exercise protocol based on their sex, age at interview, body mass index, and a Physical Activity Readiness (PAR) Code. The PAR Code is determined by a series of questions to approximate the amount of time spent in physical activity per week over the past year. See Sections 3.5.1 and 3.5.2 for the English and Spanish version of these questions.

- Wait until the SP has his/her shoes on and is ready to listen to your questions. It is important that you read the questions clearly and that this is the only activity you and the SP are doing. Make sure that the SP is listening to you and make frequent eye contact while administering these questions.
- Read the introductory sentence and the first three questions completely before accepting a response.
- If the SP gives a response before you have read all three options, tell them that you are required to read all three responses before they give an answer.
- If the response is refused or 'Don't Know,' the SP is excluded from the exam. The component status would be set to 'Not Done' with the comment 'Unable to calculate protocol.'

3.5.1 Questions to Determine PAR Code: English Version

First Level of Questions

- *I am going to ask you several questions to describe your **typical physical activity level**. Please try to answer as accurately as possible. Please wait until I read all three questions before you give me your answer.*
 1. Would you say that you did **little or no regular** recreation, sport, or physical activity?

OR
 2. Would you say that you participated **regularly** in recreation or work requiring **modest physical activity**?

OR
 3. Would you say that you participated **regularly** in **heavy physical exercise**?

If the SP needs further clarification on the meaning of modest versus heavy, use the following examples:

Modest physical activities include: walking, calisthenics, bowling, weight lifting, yard work.

Heavy physical activities include: jogging, swimming, cycling, rowing, tennis, and basketball.

Enter 1, 2, or 3 in the system.

Second Level of Questions for Little or No Regular Activity

- If the response to the first set of questions is 1 (little or no regular activity), the system will go to a set of questions to further classify the SP's physical activity level. A PAR code will then be assigned based on the response and be used to determine predicted VO₂ max.
- Read both questions before accepting a response.
 1. Would you say that you **avoid walking or exertion**, for example always use the elevator, drive whenever possible instead of walking?

OR
 2. Would you say that you walk for pleasure, routinely use stairs, occasionally exercise sufficiently to cause heavy breathing or perspiration?
- A PAR code will then be assigned as outlined below:
 - If response is 1 (avoid walking), the PAR Code will be assigned as 0.
 - If response is 2 (walk for pleasure), the PAR Code will be assigned as 1.

Second Level of Questions for Regular Participation in Modest Physical Activity

- If the response to the first set of questions is 2 (participated regularly in modest physical activity), the system will go to a set of questions to further classify the SP's physical activity level. A PAR code will then be assigned based on the response and be used to determine predicted VO₂ max.
- Read both options before accepting a response.

Would you say that you participate regularly in recreation or modest physical activity for:

 1. 10 to 60 minutes per week?

OR
 2. More than 60 minutes per week?
- A PAR code will then be assigned as outlined below:
 - If response is 1 (10 - 60 minutes per week), the PAR Code will be assigned as 2.
 - If response is 2 (more than 60 minutes per week), the PAR Code will be assigned as 3.

Second Level of Questions for Participated Regularly in Heavy Physical Activity

- If the response to the first set of questions is 3 (participate regularly in heavy physical exercise), the system will go to a set of questions to further classify the SP's physical activity level. A PAR code will then be assigned based on the response and be used to determine predicted VO₂ max.
- Read all options before accepting a response.

Now I am going to ask for some more detail on the amount of heavy physical exercise you did. I am going to give you four options.

Would you say that you participated regularly in heavy physical activity for:

1. less than 30 minutes per week

OR

2. 30 to 60 minutes per week

OR

3. 1 to 3 hours per week

OR

4. more than 3 hours per week <<do not need “bold”

- The PAR Code will be assigned as outlined below:
 - If response is 1 (less than 30 minutes per week), the PAR Code will be assigned as 4.
 - If response is 2 (30-60 minutes per week), the PAR Code will be assigned as 5.
 - If response is 3 (1-3 hours per week), the PAR Code will be assigned as 6.
 - If response is 4 (more than 3 hours per week), the PAR Code will be assigned as 7.

3.5.2 Questions to Determine PAR Code: Spanish Version

- If Spanish is the first language of the SP, the system can be changed to the Spanish version of the PAR Questions.

Le voy a hacer varias preguntas para describir su nivel típico de actividad física. Por favor trate de responder tan precisamente como sea posible. Por favor espere hasta que yo lea las tres preguntas antes de darme su respuesta.

1. Diría que usted hizo poca o ninguna actividad regular de recreación, deporte, o física?

O

2. ¿Diría que usted participó regularmente en recreación o trabajo que requería actividad física moderada?

O

3. ¿Diría que usted participó regularmente en ejercicios físicos pesados?

Actividades físicas moderadas incluye: caminar, calistenia, bolos, levantamiento de pesas, trabajo en el (patio/jardín).

Actividades físicas pesadas incluye: trotar, nadar, andar en bicicleta, remar, jugar tenis, y baloncesto.

1. ¿Diría que usted evita caminar o hacer esfuerzos (por ejemplo, siempre toma el ascensor, conduce en lugar de caminar cada vez que es posible.)?

O

2. ¿Diría que usted camina por placer, rutinariamente usa las escaleras, ocasionalmente hace ejercicio suficientemente como para que le cause respiración agitada o transpiración?

3. ¿Diría que usted participó regularmente en actividades de recreación o físicas moderadas por:

1. 10 a 60 minutos por semana?

O

2. Más de 60 minutos por semana?

—Ahora le voy a preguntar algunos detalles más sobre la cantidad de ejercicio físico pesado que usted hizo. Le voy a dar cuatro opciones.

—¿Diría que usted participó regularmente en actividades físicas pesadas |
...

1. menos de 30 minutos por semana,

O

2. 30 a 60 minutos por semana,

O

3. 1 a 3 horas por semana,

O

4. más de 3 horas por semana?

3.6 Determination of Predicted VO₂ Max

The system will use the following formula based on the PAR Code, age of SP at interview, body mass index, and sex of SP to calculate the predicted maximal oxygen consumption (VO₂ max).

$$\text{Predicted VO}_2 \text{ max} = 56.363 + [1.921 * (\text{PA-R})] - [0.38 * (\text{Age at interview})] - [0.754 * (\text{BMI})] + [10.987 * (\text{F}=0, \text{M}=1)]$$

- Physical Activity Readiness (PAR) Code = 0-7 (see Section 3.5)
- Body Mass Index (BMI) = body weight(kg) / height squared (m²) (see Section 3.3)
- Age predicted maximal heart rate (HRmax = 220 - Age)

3.7 Determination of Exercise Protocol

The system will use the predicted maximal oxygen consumption (see Section 3.6) to determine the exercise protocol to use for the SP. See Appendix H for the table that is used to determine the assigned exercise protocol for each individual SP. This protocol will be displayed on the Treadmill Test screen. The assigned protocol may be adjusted up or down one level based on the SP's response (see Section 3.12).

3.8 Preparation for Treadmill Test

Prior to beginning the treadmill test, a blood pressure cuff for automatic blood pressure monitoring and chest electrodes for monitoring heart rate will be attached to the SP. The SPs footwear should be checked once again before starting the test to confirm that laces are tied or shoes are on properly.

3.8.1 Preparation of Skin and Placement of Electrodes

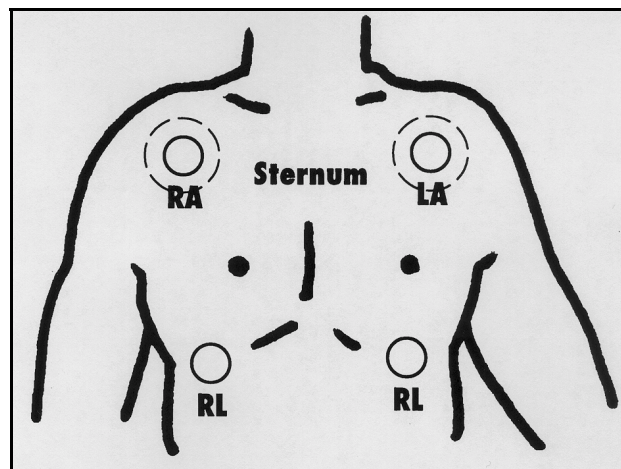


Figure 3-2. Site for placement of electrodes (1)

- The following sites will be used for electrode placement. See Figure 3-2.
 - RA: below the distal end of the right clavicle, on the superior aspect of the pectoralis major.
 - LA: below the distal end of the left clavicle, on the superior aspect of the pectoralis major.
 - RL: right anterolateral surface of the abdomen at the level of the umbilicus.
 - LL: left anterolateral surface of the abdomen at the level of the umbilicus.
 - If the individual has a large amount of body hair, the electrodes can be moved slightly to give better contact with the skin. The site will not be shaved to avoid unintended injury or discomfort to the SP. Move the electrodes approximately 2 inches laterally. See the outer limit circles on RA and LA in Figure 3-2.
- Preparation of the sites for attachment of electrodes:
 - Briskly wipe the skin at placement site of electrode with a skin preparation swab to remove any oils or lotions on the skin.
 - Lightly abrade the skin (no more than 6-8 swipes) with a fine abrasive strip covering an area the size of the electrode.
- Attach electrode to the four sites by pressing around the outer edges of electrode first to maximize adherence and then press gently on the center of the electrode to enhance the contact of the electrode gel to the skin surface.
- Clip the RA, LA, RL, and LL leads to the appropriate electrodes.
- For female SPs, hand them the two top electrodes (RA and LA) and ask them to thread the leads up underneath the gown and hand them to you. After the SP has done this, take the leads and attach them to the appropriate electrodes.

3.8.2 Selection of Arm for Blood Pressure Monitoring

The Colin STBP-780 is an automated electronic heart rate and blood pressure monitor that will be used to measure blood pressure and heart rate during the treadmill test. See Chapter 2 for a description and picture of this monitor.

- The right arm will be used for blood pressure monitoring.
- Do not use the right arm if any of the following are present on that arm: wound, rash, cast, edema, bandages, or any condition that would cause discomfort to the SP or would prevent accurate measurement of blood pressure.
- Do not use the arm if the SP has a dialysis shunt in that arm.

- If the SP has had a radical mastectomy do not use the arm on the side of the mastectomy. Use the opposite arm. If the SP has had a bilateral mastectomy, do not measure the blood pressure in either arm. The SP will be excluded from the exam.
- If the right arm cannot be used, use the left arm.
- Do not use the left arm if any of the conditions listed above for the right arm are present. In this situation, enter “Could Not Obtain” in the system.
- Enter “right, left, or could not obtain” for Arm Selected.
- If “right” or “left” is selected, go to cuff size and enter the cuff size that will be used.
- If “could not obtain” is selected, the SP will be excluded due to inability to monitor blood pressure. The system will display a message to ask you to confirm your selection. If you select yes, another message will be displayed as “SP excluded due to safety considerations/physical limitations” informing you that the SP has been excluded from the test. The system will go to the Component Status screen; the Status will be set to “Partial” with the comment “physical limitation.”

3.8.3 Selection of Cuff Size and Placement of Cuff

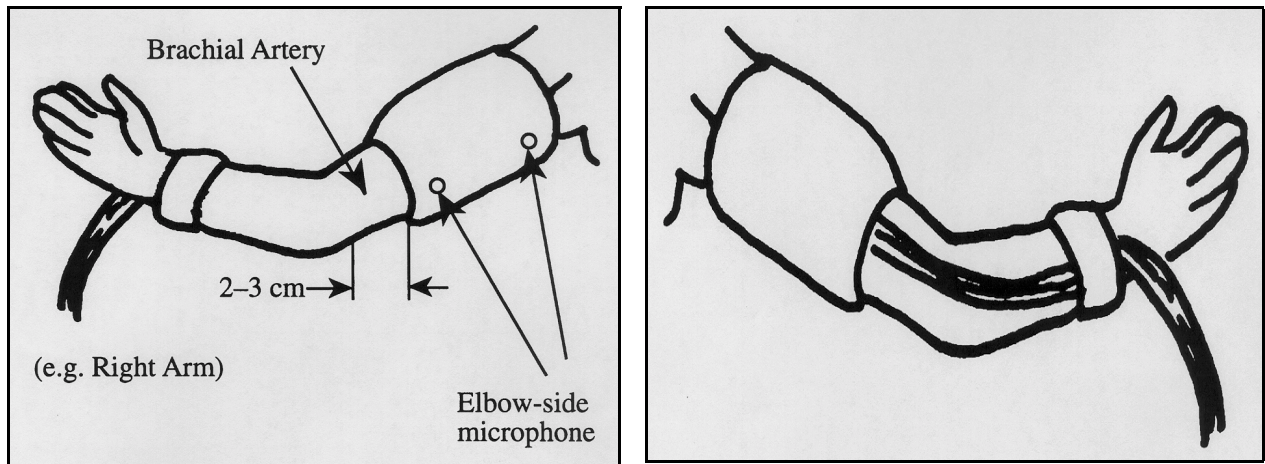


Figure 3-3. Site for placement of electrodes (2)

- See Figure 3-3 for position of cuff on arm.
- The MEC physician will measure the arm to determine the appropriate cuff size and will write the cuff size used on the shoulder or sleeve of the SP's gown. The same size cuff should be used for this exam.
- Attach appropriately sized blood pressure cuff (child, adult, or large adult) on the arm selected. The bladder of the cuff should be about two-thirds of the arm circumference.
- Palpate for the brachial artery and position the cuff so that the two microphone sensors are located directly over the brachial artery.
- Position the cuff snugly around the upper arm, 1 inch above the antecubital space.

3.8.4 Room Temperature

Ideal room temperature for the treadmill test should be 68-72 degrees Fahrenheit. Perspiration on the skin will interfere with the adherence of the electrodes to the skin and higher temperatures may increase the effort needed by the SP to perform the test.

The room temperature will be automatically recorded by integration with a Nimbus digital thermometer. If integration fails, the room temperature will be entered manually by entering the temperature reading on the digital thermometer. The room should be kept within the ideal temperature noted above. If necessary, turn up the air conditioner and turn on the fan. It may be more comfortable for the SP if the fan is turned on after the test has started and he/she is warming up.

3.8.5 Mesh Vest

If male SPs have chest hair that interferes with the attachment of the electrode:

- The electrode can be moved slightly to give better contact with the skin. The electrode may be positioned approximately 2 inches from the optimal site for electrode placement.
- If there is a lot of hair on this spot, try to move the hair aside and place the electrode on the site with as much skin contact as possible.
- In this situation, optimize the situation and use the mesh vest. Record this on the Exam Preparation screen in ISIS. The default for this field is set to 'No' and is disabled for females.

3.9 Treadmill Test

3.9.1 Setup of Colin STBP Automatic BP Monitor

- See Sections 3.8.1 through 3.8.3 for procedure for electrode placement and cuff selection and placement.
- Press Power button to turn on power.
- SELF TEST is displayed on the Elapsed Time/Cuff Pressure Indicator.
- The STBP monitor performs a self-test.
 - If an error is found:
 - a. **FAIL** is displayed on the Elapsed Time/Cuff Pressure Indicator.
 - If **no error** is found, the following are displayed for about 5 seconds:
 - a. The year on the Systolic Indicator,
 - b. The month on the Diastolic Indicator, and
 - c. The date on the Heart Rate Indicator.
- FIRST REST is displayed on the Elapsed Time/Cuff Pressure Indicator if the ECG cable is connected properly to the ECG connector on the rear panel.
- Heart Rate is displayed in the Heart Rate Indicator.
- CHECK ECG is displayed on the Elapsed Time/Cuff Pressure Indicator if electrodes are not connected properly. Check that the leads are in the correct position. Check that the ECG cable is connected tightly to the port.

3.9.2 Setup of Programmable Treadmill Controller

As described in Sections 3.6 and 3.7, the ISIS system will assign and display a set of exercise protocol for each individual SP based on the predicted VO₂ max. This set of predetermined test protocols will be manually programmed into the controller as outlined below.

- Press Power to turn on the power to the controller and the treadmill.
- Press Enter (this enters 150 as weight - do not change this).
- Press Input + twice to bring up the different protocols.
- Press Input + to move forward through the protocols.
- Press Input - to move back through the protocols.
- Two sets of information are indicated on the multifunction display window in the upper middle of the treadmill controller. The protocol number appears to the left and the exercise time appears to the right of the window (see Figure 2-4).
- There are 16 horizontal lights on the lower middle of the treadmill display. The green LED's on the left hand side of the lights indicate the preprogrammed grade. The red LED's on the right hand side indicate the preprogrammed speed. See Figure 3-4 on page 3-32.
 - The first 4 lights represent the current protocol (Warm-up, Stage I, Stage 2, and Recovery).
 - The next 3 lights represent the next lower protocol (Stage 1, Stage 2, and Recovery).
 - The next 3 lights represent the next higher protocol (Stage 1, Stage 2, and Recovery).
 - The next 3 lights represent the next lower protocol again (Stage 1, Stage 2, and Recovery).
 - The next 3 lights represent the next higher protocol again (Stage 1, Stage 2, and Recovery).
 - A blinking light represents the current stage.

As described above, each of the 16 vertical columns of lights in the tricolor display of the treadmill controller represents one stage of the protocol. Because the controller only allows advancing to the next column of protocol, technicians are not able to go back to the previous protocol on the controller without advancing through all 16 levels. Therefore, the next lower and next higher protocols are repeated twice in case the technician overshoots the desired protocol when manually changing the protocol.

3.9.3 Obtaining Pretest Blood Pressure and Heart Rate

- When the Pretest Measurement screen comes up on ISIS, the system will automatically enter the room temperature from an integrated digital thermometer.
- The system will then give an auditory message:

Check footwear, attach chest electrodes and BP cuff. Then explain the treadmill test.

- The SP's footwear is checked when they come in the room and he/she is asked to put on appropriate footwear at this time. This check is just an additional measure to confirm that he/she has proper footwear and that laces are tied.
- Attach electrodes and BP cuff as outlined in Sections 3.8.1 and 3.8.3.
- Explain the treadmill test to the SP following the standard script in Section 3.9.4 when attaching the electrodes and BP cuff.
- When the explanation is completed, have the SP stand quietly with their arms down by the side while the pretest measurements are being obtained.
- Confirm that the BP monitor power button is ON.
- Press CAPTURE on the ISIS screen to tell the system to start the pretest BP measurement. It will take 35 to 45 seconds for the measurements to be taken.
- Pretest blood pressure and heart rate will then be displayed on BP monitor and will be automatically captured in the ISIS system.
- These values are not the official data points for blood pressures for the NHANES. The SP's resting blood pressure and heart rate are measured and recorded in the physician's room prior to the treadmill test. The pretest values obtained in the treadmill testing room simply indicate that the STBP unit is functioning.
- When the blood pressure measurements are captured, the system will display a message:

"Demonstrate and explain how to walk on a treadmill."

- Give a brief demonstration of how to walk on a treadmill. See Section 3.9.5.

3.9.4 Script for Explaining Treadmill Test

While you are attaching the chest electrodes and the blood pressure cuff you can explain the treadmill test to the SP. Below is a standard script that should be read exactly as written to ensure that each SP is given the same explanation for this test.

Explanation of Treadmill Test - Standard Script: English Version

“The purpose of the treadmill test is to estimate your cardiovascular fitness level. You will be walking on the treadmill at various speeds and grades. This is not a maximal test to exhaustion. We will stop when your heart rate reaches a certain level. While this is a very safe test, it is important that you let us know how you feel during the test. Please report any unusual symptoms that you have during the test such as chest pain, dizziness, nausea, or pain in your legs.”

Explanation of Treadmill Test - Standard Script: Spanish Version

“El propósito de la prueba en la máquina de andar es calcular su nivel de capacidad cardiovascular. Usted caminará en la máquina a diferentes velocidades y grados. Esta no es una prueba máxima de esfuerzo. Pararemos cuando la frecuencia cardíaca (los latidos del corazón) alcance cierto nivel. Si bien esta es una prueba muy segura, es importante que nos diga cómo se siente durante la prueba. Por favor informe cualquier síntoma inusual que tenga durante la prueba tal como dolor en el pecho, mareo, náusea, o dolor en las piernas.”

3.9.5 Demonstrating How to Walk on a Treadmill

After the blood pressure measurements have been captured, (see Section 3.9.3) , you should give a brief demonstration of how to walk on a treadmill. The SP should be asked to step off the treadmill and stand to one side. Start the treadmill belt by pressing the Power and Start Course buttons on the treadmill controller. Increase the speed of the treadmill to 2.0 mph. Step on the treadmill and begin a brief demonstration of how to walk on the treadmill. You should hold on to the handrails, straddle the belt, and then gently “paw” the belt with one foot to get a feel for the movement, then step on to the belt and start walking slowly. The demonstration and explanation should include the importance of standing up straight, taking long steps, relaxing, holding the head up and looking forward. You should demonstrate

by removing one hand at a time from the handrail and then walking comfortably on the belt. After this brief demonstration, step off the belt and press Stop Belt on the treadmill controller, then press Clear.

After the data are captured and the demonstration is completed, press the 'Next' button on the ISIS screen to advance to the Treadmill Test screen.

- When the Treadmill Test screen comes up, the ISIS system displays and gives this auditory message:

Choose exercise protocol on treadmill controller. Have SP straddle the belt, then press "Start Belt."

- Select the predetermined protocol (the middle protocol on the ISIS screen) on the treadmill controller as described in section 3.9.2.
- Get the SP to straddle the belt. Press 'START BELT' on the treadmill controller.
- Increase the speed to 2 MPH.
- When you have completed this, click OK to the message on the ISIS screen.
- Have the SP step on the treadmill and begin to walk.
- Tell him/her to use an upright posture and look straight ahead.
- If the SP is holding on to the handrail, tell him/her to drop one arm down by their side, then the other one.
- Press OK on the ISIS screen.
- The SP should now be walking slowly on the treadmill and ready to begin the test.
- If the SP is not able to take his/her hands off the handrails immediately, proceed to the Warm-up stage. (See Section 3.9.6). During the first minute of the Warm-up stage, work with the SP to get him/her to let go of the handrail.

3.9.6 Warm-up (00:00 - 02:00 Minutes)

- When OK was pressed as indicated above, the system will give an auditory message:

“Simultaneously press START on the ISIS screen and press START COURSE on treadmill controller.”

- Click ‘START’ on the ISIS screen to start the timer for the test.

AND

- Press ‘START COURSE’ on the treadmill controller to begin the Warm-up stage for the predetermined protocol selected for the SP.
- The timer will start the time for the 2-minute Warm-up. At the end of each stage, the timer will reset to zero and begin the next stage.
- When START COURSE is pressed on the treadmill controller, the treadmill will increase to the speed and grade for the predetermined Warm-up stage for the selected protocol.
- The system will automatically capture BP and HR at the end of each stage.
- Throughout the Warm-up stage, monitor the %HRmax scale on the ISIS screen. The SP’s HR should be up to the **60 percent range**. The %HR max scale displays the HR in beats per minute for 50 percent, 60 percent, and 100 percent of PMHR for the SP. The SP’s current HR will be displayed on the BP monitor.
- If %HRmax is below 50 percent or above 60 percent, see Section 3.12 for protocol change procedure.
- If %HRmax is within the **50 percent - 60 percent range**, continue with present preprogrammed protocol.

3.9.7 RPE Explanation and Assessment

- Rating of perceived exertion (RPE) will be obtained at the end of Warm-up, Stage 1, and Stage 2. This scale will be explained to the SP during the first minute of Warm-up.
- If the SP is having trouble learning to walk on a treadmill, the explanation can be given after the SP is stable and comfortable. The priority is to get the SP walking on the treadmill with stability as soon as possible.

- At 00:15 minutes into Warm-up, the system will give auditory message:
“Explain RPE”
- The technician will give an explanation of RPE to the SP. See Section 3.9.7.1 for the English version and Section 3.9.7.2 for the Spanish version. This is a standard script and should be read exactly as it is written so that each SP is given the same explanation.
- At 01:30 minutes into Warm-up, the system will give an auditory message:
“Get RPE assessment.”
- The technician will ask the SP to give his/her rating of perceived exertion.
- The technician will enter this number.
- At 1:50 minutes into Warm-up, the system will give an auditory message:
“Press NEXT STAGE on treadmill controller at 1 minute 57 seconds.”
- Technician should press NEXT STAGE on the Treadmill Controller at exactly 01:57 minutes into Warm-up.
- This will change the treadmill to Stage 1.

3.9.7.1 Explanation of Rating of Perceived Exertion (RPE) - Standard Script: English Version

“Please look at the chart on the wall in front of you. Every couple of minutes during the test, I will ask you to rate how hard you feel you are working overall. Please respond with a number from the scale. For example, if you feel that you are working at a very low level of effort, you might reply with a 6 or 7 (very, very light). On the other hand, if you feel you are working at a very high level of effort and cannot continue, you might reply with a 20 (maximal exertion). There is no right or wrong answer, simply rate your perceived level of exertion at the time I ask the question.”

3.9.7.2 Explanation of Rating of Perceived Exertion (RPE) - Standard Script: Spanish Version

“Por favor mire el gráfico en el muro frente a usted. Cada dos minutos durante la prueba, le pediré que califique qué tan fuerte siente que está trabajando en general. Por favor conteste con un número de la escala. Por ejemplo, si usted siente que está trabajando a un nivel de esfuerzo muy bajo, usted puede responder un 6 ó 7 (muy, muy liviano). Por otra parte, si usted siente que está trabajando a un nivel de esfuerzo muy alto y no puede continuar, usted puede contestar con un 20 (máximo esfuerzo). No

hay respuesta correcta o errónea, simplemente califique su nivel de esfuerzo en el momento en que yo le haga la pregunta.”

3.9.8 Stage 1 (00:00 - 03:00 Minutes)

- The system will reset the timer to zero for Stage 1.
- Throughout the Stage 1, monitor the %HRmax scale on ISIS screen. Indicator should be within 60 percent - 80 percent range.
- If the SPs HR goes above the 80 percent range, you must go to the Recovery stage by pressing 'Go to Recovery' on the ISIS screen. See 3.13.2 for detailed procedures.
- At 02:30 minutes into Stage 1, the system will give an auditory message:
“Get RPE assessment”
- The technician will ask the SP to give his/her rating of perceived exertion.
- The technician will enter this number.
- The system will automatically set the BP monitor to start the BP measurement at 02:10 minutes into Stage 1. The end of Stage 1 BP and HR measurement will be captured at approximately 02:40 to 02:50 minutes into Stage 1.
- The system will give the following auditory message:
“Check heart rate on monitor.”
- The ISIS system will then display a visual message:
“Compare heart rate on blood pressure monitor with captured heart rate on ISIS screen. Enter a manual heart rate ONLY if they differ by more than 5 beats.”
- Note the heart rate on the blood pressure monitor and compare this number with the captured heart rate on the ISIS screen. If the captured heart rate on the ISIS screen does not differ by more than 5 beats from the heart rate on the BP monitor, click OK and continue the test.
- If the heart rate is not captured on the ISIS screen or the observed HR on the blood pressure monitor is more than 5 beats from the captured heart rate, enter the observed HR in the manual heart rate field on the ISIS screen.
- Five seconds before the stage ends, all the LED's in that column flash, warning that the speed and grade will change soon.
- When the stage ends, the top and bottom LEDs turn yellow and the bottom LED in the next stage begins to flash.

- The treadmill will automatically change to the speed and grade for Stage 2.

3.9.9 Stage 2 (00:00 - 03:00 Minutes)

After 3 minutes in Stage 1, the speed and grade move to the values preprogrammed into the second stage for the selected protocol.

- Throughout Stage 2, monitor %HRmax scale on ISIS screen. The heart rate should be:
 - below 90 percent for 12- through 18-year-olds, and
 - below 85 percent for 19- through 49-year-olds.
- During the first 2 minutes of Stage 2, if %HR max is above 85 percent (or 90 percent for 12- through 18-year-olds), see Section 3.13 for stopping procedures. (Active cool down and stopping criteria 2).
- During the third minute of Stage 2, if %HRmax is 85 to 100 percent (90 to 100 percent for 12- through 18-year-olds) the technician should observe the SP for signs of exhaustion and or distress. The technician should ask the SP “Are you able to continue for another minute?” If the SP agrees AND the SP does not appear distressed or exhausted, continue with present protocol.
- During the third minute of Stage 2, if %HR max is above 100 percent, then press ‘Go to Recovery’ on the ISIS screen. See Section 3.13 for Stopping Procedures. Select ‘Other’ with comment ‘%HRmax > 100%’ for Priority 2 Stopping Criteria.
- At 02:30 minutes into Stage 2, the system will give an auditory message:

“Get RPE assessment”
- The technician will ask the SP to give his/her rating of perceived exertion.
- The technician will enter this number.
- The system will automatically set the STBP monitor to start the BP measurement at 02:10 minutes into Stage 2. The end of Stage 2 BP and HR measurement will be captured at approximately 02:40 to 02:50 minutes into Stage 2.

- The system will give the following auditory message:
“Check heart rate on monitor.”
- The ISIS system will then display a visual message:
“Compare heart rate on blood pressure monitor with captured heart rate on ISIS screen. Enter a manual heart rate ONLY if they differ by more than 5 beats.”
- Note the heart rate on the blood pressure monitor and compare this number with the captured heart rate on the ISIS screen. If the captured heart rate on the ISIS screen does not differ by more than 5 beats from the heart rate on the BP monitor, click OK and continue the test.
- If the heart rate is not captured on the ISIS screen or the observed HR on the blood pressure monitor is more than 5 beats from the captured heart rate, enter the observed HR in the manual heart rate field on the ISIS screen.
- Five seconds before the stage ends, all the LED’s in that column flash, warning that the speed and grade will change soon.
- When the stage ends, the top and bottom LED’s turn yellow and the bottom LED in the next stage begins to flash.
- After 3 minutes in stage 2, the speed and grade move to the values preprogrammed into Recovery 1 for the selected protocol.

3.9.10 Recovery Stages 1, 2, & 3 (If Necessary)

- The system will automatically set the BP monitor to start the BP measurement at 00:10 minutes into Recovery 1. The end of Recovery 1 BP and HR measurement will be captured at 00:50 minutes into Recovery 1.
- The system will automatically set the BP monitor to start the BP measurement at 00:10 minutes into Recovery 2. The end of Recovery 2 BP and HR measurement will be captured at 00:50 minutes into Recovery 2.
- If $HR \leq 110$, an auditory message should remind technician to stop the timer and stop the treadmill.
“Press COOL DOWN on treadmill controller.”
- Press COOL DOWN as directed above.
- Tell SP to grasp the front handrail and warn that the treadmill is about to stop.
“Press STOP BELT on treadmill controller. Congratulations, you have completed the test.”
- Press STOP BELT as directed.

- If HR > 110, proceed with Recovery 3. The system will give an auditory message:

“Proceed with Recovery 3”

- **Recovery 3:** The system will automatically set the BP monitor to start the BP measurement at 00:10 minutes into Recovery 3. The end of Recovery 3 BP and HR measurement will be captured at 00:50 minutes into Recovery 3. The system will give the following auditory message:

“Press COOL DOWN on treadmill controller.”

- Press COOL DOWN as directed above.
- Tell SP to grasp the front handrail and warn that the treadmill is about to stop.

“Press STOP BELT on treadmill controller. Congratulations, you have completed the test.”

- Press STOP BELT as directed.

3.9.11 End of Test

Assist the SP off the treadmill. It is normal to feel as if you are still moving, especially if you are not accustomed to walking on a treadmill. Remove BP cuff, disconnect the leads, and ask the SP to remove and dispose of the electrodes.

3.10 Procedure if Values Exceed Criteria

The system has the following edit limits set on these measurements:

- HR Stage 1 > 80 percent PMHR
- HR Stage 2 > 85 percent PMHR (>90% for 12 through 18 years). This applies to the first 2 minutes of Stage 2 only. See Section 3.12.3 and Exhibit 4-69.
- SBP > 260 mmHg
- DBP > 115 mmHg
- RPE > 17 (See below “Exception for RPE greater than 17 at end of Stage 2”)
- If one or more of these measurements are above the acceptable criteria, the system will display a message:

“Check message on monitor.”

- The system will then display a message *“(Data item) is not within the acceptable limits. Go to Recovery.”*
- Press the ‘Go to Recovery’ button.
- The system will then give an auditory message:

“Press GO TO RECOVERY on ISIS screen and change treadmill controller to Recovery Stage for current protocol.”

3.10.1 Exception for RPE greater than 17 at end of Stage 2

- If the RPE is greater than 17 in Warm-up, or Stage 1, go to Recovery. There is no exception at these two stages but you should confirm that the SP understands the concept. Confirm that the SP understands what he/she is rating and if he/she still gives a rating of greater than 17, go to Recovery.
- When you get the RPE assessment at 2:30 minutes into Stage 2, there is 15 to 20 seconds left before ISIS captures the heart rate and blood pressure.
- If the SP answers with a number greater than 17 AND their heart rate is within the acceptable limits for the end of stage 2:
 - Determine if he/she is able to continue for another 15 to 20 seconds.
 - Check for signs of distress by looking at their face, color, gait, respiratory rate.
 - Ask if they are able to continue for another 15 to 20 seconds.
 - Tell them the test is almost finished.
 - An example of a script is: “Do you think you can continue for another 15 to 20 seconds? The test will be over then and the treadmill will be slowing down in speed and the incline will become level. Are you sure you can continue?”
 - If the SP says they can continue for another 15-20 seconds, monitor continuously and allow the test to continue until the Stage 2 HR and BP are captured.
- Once the HR and BP are captured, if you think the SP can continue for another 10 seconds, allow the test to finish normally.
- If you think the SP cannot continue for another 10 seconds then:
 - Press ‘Next Stage’ on the treadmill controller to move the treadmill to the Recovery Stage.
 - Do not Press ‘Go to Recovery’ in ISIS.

- NOTE: You may not have to do this because the length of time is so short and the treadmill may have already advanced by the time you are ready to do this.
- Continue Recovery and Cool Down as usual.

3.11 Procedure if SP is Unable to Adhere to the Protocol

- Some people may not be able to walk on the treadmill without holding onto the handrails. The technician should give them further instructions and encouragement. SPs are allowed to hold onto the handrail during Warm-up until he/she gets accustomed to walking on the treadmill. The technician should encourage the SP to remove their hands one at a time. During the treadmill test they are allowed to lightly touch the handrail with a finger but they should not be grasping the handrail with their fingers or hands. If they are unable to do this, the test will not be valid.
- Some of the instructions that you may give to help the SP get accustomed to walking on the treadmill include reminding the SP to stand up straight, to keep his/her head up and look straight ahead, and to take longer steps.
- If the SP is still unable to adhere to the protocol, press GO TO RECOVERY on the ISIS screen.
- The system will record the time in the current stage.
- The system will reset the timer for 00:00 for Recovery.
- The system will give an auditory message:

“Change treadmill controller to Recovery Stage for current protocol”
- Change the setting on the treadmill controller by pressing NEXT STAGE until the Recovery stage is reached.
- The SP will go through the Recovery Stages 1 and 2.
- After the Recovery stages are completed, the system will go to the Priority 2 Stopping Criteria.
- Select “Unable to stop gripping handrails.”
- System will then go to Component Status screen. The component status will be set at “Partial.” Select “Unable to follow protocol” as the comment.

3.12 Procedure to Change Predicted Exercise Protocol

The goal of the predetermined exercise protocol is to elicit a heart rate that is approximately 75 percent of the predicted maximum heart rate (PMHR) by the end of the treadmill test. This goal is

expected to be achieved gradually by stages: to elicit a HR response that is approximately 50 to 60 percent of PMHR during the Warm-up stage, then to elicit a HR response that is approximately 60 to 75 percent of PMHR during Stage 1, and last, to elicit a HR response that is approximately 70 to 80 percent of PMHR during Stage 2. There will be times where the selected protocol results in a HR response that is higher or lower than expected. When this occurs, the treadmill technician will need to change to either an easier or more difficult protocol to help reach the goal. The decision to change the exercise protocol should be made during the beginning of the test – the Warm-up stage. No decision in changing the protocol should be made after the Warm-up stage. If the ISIS “Change Protocol” button is selected during Stage 1 or 2, the system will display this message: “The protocol cannot be changed after Warm-up.”

Throughout the Warm-up stage, monitor the **%HRmax** scale on ISIS screen. The SP’s HR should be within the **50 percent to 60 percent range**. The %HR max scale on the computer screen displays the HR in beats per minute for 50 percent, 60 percent, and 100 percent of PMHR for the SP. The SP’s current HR will be displayed on the STBP monitor. If this reading is above or below the 50 to 60 percent range, follow the procedures in 3.12.1 or 3.12.2 for changing the protocol. The decision to change the protocol must be made by 1 minute 30 seconds into the Warm-up stage to allow time to make all the necessary adjustments. If the protocol needs to be changed, the Change Protocol button on the ISIS screen should be pressed during the last 30 seconds in Warm-up but the setting on the treadmill controller should not be changed until the treadmill automatically moves to the next stage. The treadmill is programmed to move to the next stage at specific times. If it is advanced to a new protocol during one stage, it will advance again to the next stage when the preset time is reached for the next stage.

There are 16 horizontal lights on the treadmill display. See Figure 3-4. The first four lights represent the current protocol (Warm-up, Stage 1, Stage 2, and Recovery). The next three lights represent the next lower protocol (Stage 1, Stage 2, and Recovery) and the next 3 lights represent the next higher protocol (Stage 1, Stage 2, and Recovery). The next lower and next higher protocols are repeated a second time to allow the technician to advance to the correct protocol if they overshoot the protocol the first time. The controller does not allow moving backwards through the protocols. You can only advance to the next protocols. A blinking light represents the stage that the SP is currently performing.



Figure 3-4. Treadmill controller - changing the exercise protocol

3.12.1 Increasing the Exercise Protocol

If %HRmax is below 50 percent PMHR in Warm-up, go to the next higher protocol at the beginning of Stage 1.

- Click on CHANGE PROTOCOL on the ISIS screen after 1 minute 30 seconds into Warm-up.
- Click on INCREASE and APPLY on the ISIS screen to change protocol to the next highest protocol.
- The system will give an auditory message:
“When this stage changes, change stage on treadmill controller to next higher protocol.”
- Do not change the stage on the treadmill controller until the current stage is completed and the treadmill has advanced to the next stage.
- As soon as the treadmill changes to Stage 1, press NEXT STAGE on the treadmill controller to get to Stage 1 for the next higher protocol. Watch the markings on the controller carefully and do not go past the desired setting.
- The system will automatically highlight the new stage and protocol on the ISIS screen.

- Confirm the speed and grade settings on the treadmill controller with the settings on the ISIS screen.

3.12.2 Decreasing the Exercise Protocol

If %HRmax is above 60 percent PMHR in Warm-up, go to next lower protocol at the beginning of Stage 1.

- Click on CHANGE PROTOCOL on the ISIS screen.
- Click on DECREASE and APPLY on the ISIS screen to change protocol to next lowest protocol.
- The system will give an auditory message:
“When this stage changes, change stage on treadmill controller to next lower protocol.”
- Do not change the stage on the treadmill controller until the current stage is completed and the treadmill has advanced to the next stage.
- As soon as the treadmill changes to Stage 1, press NEXT STAGE on the treadmill controller to get to Stage 1 for the next lower protocol. Watch the markings on the controller carefully and do not go past the desired setting.

- The system will automatically highlight the new stage and protocol on the ISIS screen.
- Confirm the speed and grade settings on the treadmill controller with the settings on the ISIS screen.

3.13 Procedure to Stop Treadmill Test Early

Technicians must be aware that not all SPs will be able to complete the treadmill test due to various reasons. From a safety perspective, it is critical that all technicians be very attentive to SPs regarding any signs/symptoms that may occur. There are two sets of stopping criteria and a different procedure will be used for each.

3.13.1 Priority 1 Stopping Criteria

- Pain or pressure in the chest
- Pain in neck, jaw, or teeth
- Pain in the shoulder or between the shoulder blades
- Pain radiating down one or both arms
- Moderate/severe leg pain or cramping
- Lightheadedness, dizziness
- Severe shortness of breath/wheezing
- Change in skin color to bluish/grayish
- Severe headache
- Subject appears confused or agitated
- Nausea
- Unstable gait
- Visual problems

If one or more of the above criteria are present, the technician will terminate the test, stop the treadmill, and initiate the following emergency action plan:

- Press EMERGENCY STOP on the ISIS screen. The system will give the following auditory message:

“Press COOL DOWN on treadmill controller, then press STOP BELT.”
- Assist the SP off the treadmill.
- Place SP in the supine position on the floor if necessary.
- Remain with the SP and call the physician.
- When the SP has stabilized complete the information in the system.
- The system will record the time spent in the current stage at the time the EMERGENCY STOP button was pushed.
- The system will go to the Priority 1 Stopping Criteria (listed above).
- Enter the appropriate response. Check all that apply.
- The remaining stages and Recovery will not be completed.
- After the Stopping Criteria have been entered, press NEXT STAGE to go to the Component Status screen. The Status will be set to “Partial Complete” with the comment “Priority 1 Stopping Criteria.”
- Press the Finish button to exit this test.

3.13.2 Priority 2 Stopping Criteria

- SP HR during Warm-up, Stage 1, or first 2 minutes of Stage 2 exceeds 90 percent for 12- through 18-year-olds or 85 percent PMHR for 19- through 49-year-olds.
- SP exceeds 80 percent PMHR during Stage 1.
- Exercise blood pressure > 260 mmHg systolic and/or >115 mmHg diastolic.
- Significant drop (>20 mmHg) in systolic blood pressure during exercise.
- SP reports RPE > 17.
- Failure of systolic blood pressure to rise with increasing exercise intensity.
- Failure of heart rate to increase by more than 5 beats/minute during a given workload compared to the previous workload.
- Repeated noticed variability during a stage - (**sudden** changes of ± 30 beats/minute) - after checking electrode connections.
- SP requests to stop test, or reports severe fatigue.

- Failure of testing equipment.
- SP is unable to stop gripping handrails (lightly touching handrails for balance purposes is allowed, but only if absolutely necessary). This criterion should cause the system to set the component status to “Partial” with the comment “unable to follow protocol.’
- Technician discretion (technician notices signs of severe exertion).

If one or more of the above criteria are present, terminate test and allow SP to actively cool down:

- Press “GO TO RECOVERY” button on the ISIS screen.
- The system will record the time spent in the current stage.
- The timer will be reset to 00:00 on the ISIS screen for the Recovery stages. The system will give an auditory message:

“ Change treadmill controller to Recovery Stage of current protocol.”
- Press NEXT STAGE on the treadmill controller to get to Recovery for the current protocol. Watch the markings on the controller carefully and do not go past the desired setting. Press NEXT STAGE only once for this change.
- The SP will continue to walk on the treadmill for 1-3 minutes and the system will continue to monitor HR and BP at 1-minute intervals. The technician may decide to terminate the test at any time based on the SP’s condition.
- When the Recovery Stages are completed, the system will give the following auditory message:

“Press STOP on treadmill controller. Congratulations, you have completed the test.”
- The system will go to Priority 2 Stopping Criteria at the end of Recovery.
- Enter the appropriate response. Check all that apply.
- After the Stopping Criteria have been entered, press NEXT STAGE to go to the Component Status screen. The Status will be set to “Partial” with the comment “Priority 2 Stopping Criteria.”
- Press the Finish button to exit this test.

3.14 Procedure to Recover/Restart an Exam

If the treadmill test has not gone past the first 15 seconds in Stage 1, you can restart the test. This might occur in the following situations:

- The SP is not used to walking on a treadmill and jumps off. The SP may be able to continue the test if the instructions are repeated and the test is restarted;
- If the treadmill controller advances or is mistakenly advanced by more than 1 stage when Next Stage is pressed to move from Warm-up to Stage 1; or
- Other problems selecting the protocol on the treadmill controller.

In these situations the treadmill test may be restarted. The restart procedures below are also kept on the quick reference laminated cards in the CV Fitness room.

- Press **Cool Down** and **Stop Belt** to stop the treadmill.
- Confirm that the SP is willing to restart the treadmill test.
- Tell the SP that you are going to repeat the instructions step by step.
- Reset the treadmill controller to the appropriate exercise protocol.
- Repeat the instructions for getting on the treadmill.
- Have SP straddle the belt, press **Start Belt**.
- Have the SP start walking on the belt.
- Press **Start Course** on treadmill controller and '**Restart**' on ISIS screen.
- Continue the test as per protocol.

3.15 Recording Problems with Data Acquisition During the Treadmill Test

Record problems with equipment, protocol, and data measurement or capture on the data acquisition problem screen. Check all that apply. The categories are:

1. No problem.
2. Problem obtaining BP - the BP is not displayed on the BP monitor and not captured in ISIS. Check the stage(s) when this happened.
3. Problem obtaining HR (variability).
4. Problem with data capture - the BP is displayed correctly on the BP monitor but not correctly captured in ISIS.
5. Problem with electrodes (perspiration, chest hair, large chest/abdomen).
6. Problem with BP cuff (large arm, cuff slipped, bending/flexing of arm).

7. Problem with treadmill controller (accidental overshoot when advancing).

3.16 Component Status

The component status is set by the system based on predetermined criteria. If all the data were obtained and captured, the status is set to complete. If the test was stopped early or some of the data were not obtained or captured, the status will be set to partial. If the SP was excluded prior to the exam due to safety/exclusion reasons or refusal, the status will be set to not done. If the status is set to partial or not done, a comment is required. Some comments are set by the system based on predetermined criteria.

4. INTEGRATED SURVEY INFORMATION SYSTEM (ISIS) DATA ENTRY SCREENS

The Integrated Survey Information System (ISIS) is the data entry system that is used in this study. See Section 1.3 for a brief description of ISIS. Detailed procedures for conducting the treadmill test are outlined in Chapter 3. The purpose of this chapter is to provide a step-by-step guideline on the data entry for this component.

4.1 Height and Weight Measurements to Determine BMI

Exhibit 4-1. Weight and height measurements from body composition

Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm

SP ID: 818591 Name: BROWN, MAC Age: 38 years Gender: Male Date: 10/10/2000 Time: 12:46 PM

SP Height and Weight Information

Measured Weight: 77.1 Kgs

Measured Height: 170. Cms

1 of 10

End of Section Close Exam Finish

- As soon as you bring the SP into the CV Fitness room you should log on to the system. Then wand the SP into the system to confirm that you are conducting the test on the correct SP.
- Do not start the height and weight measurements until you have logged on to the system.

- If the SP completed the Body Measures component prior to starting the CV Fitness component, the height and weight measurements from that component will be pulled from the database and displayed on the first screen. The fields for weight and height measurements will be disabled to prevent unnecessary entry. See Exhibit 4-1 above.

Exhibit 4-2. Weight and height measurements in CV Fitness room

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm". Below the title bar, a status bar displays "SP ID: 818591 Name: BROWN, MAC Age: 38 years Gender: Male Date: 10/10/2000 Time: 12:48 PM". The main content area is titled "SP Height and Weight Information" and contains a central box with two input fields: "Measured Weight: 77.1 Kgs" and "Measured Height: 170. Cms". At the bottom of the window, there is a navigation bar with buttons for "End of Section", "Close Exam", and "Finish", along with a "1 of 10" indicator and a small icon in the bottom right corner.

- If the SP didn't complete the body measures component prior to entering the CV Fitness component, the weight and height will have to be measured in the CV Fitness room. The SP should not be wearing his/her shoes for either measurement.
- First, measure the SP's weight using the floor scales. Enter the weight in kilograms in the weight field.
- Next, measure the SP's height using the stadiometer. Enter the height in centimeters in the height field.
- The height and weight measurements will be used to determine Body Mass Index (BMI) which in turn will be used in the calculation of the VO_2 max.

4.1.1 Edits on Data Entry for Weight and Height

Exhibit 4-3. Edit message on incorrect data entry for weight

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm". The header bar contains patient information: "SP ID: 818591 Name: BROWN, MAC Age: 38 years Gender: Male Date: 10/10/2000 Time: 07:57 PM". The main area is titled "SP Height and Weight Information". A modal dialog box titled "Cardiovascular Fitness" is displayed, containing a question mark icon and the text: "The weight is outside the edit range for this age group and gender. Are you sure?". Below the text are "Yes" and "No" buttons. In the background, the "Measured Weight" field contains the value "150." followed by "Kgs", and the "Measured Height" field contains a period "." followed by "Cms". The bottom of the window features navigation controls (back, forward, and a page indicator "1 of 10") and three buttons: "End of Section", "Close Exam", and "Finish". A blue arrow icon is located in the bottom right corner.

The height and weight fields have set edit limits based on the 1999 NHANES data.

- If the weight entered is above the edit limits, the ISIS system will display a message: "The weight is outside the edit range for this age group and gender. Are you sure?"
- If the weight you entered is correct, select 'yes.' If the data was entered incorrectly, select 'no' and re-enter the correct value for weight.

Exhibit 4-4. Edit message on incorrect data entry for weight (below edit range)

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm". The header bar contains the following information: "SP ID: 818591", "Name: BROWN, MAC", "Age: 38 years", "Gender: Male", "Date: 10/10/2000", and "Time: 07:57 PM". Below the header, the section "SP Height and Weight Information" is displayed. A modal dialog box titled "Cardiovascular Fitness" is open, showing a question mark icon and the text: "The weight is outside the edit range for this age group and gender. Are you sure?". Below the text are "Yes" and "No" buttons. In the background, the "Measured Weight" field contains the value "40." followed by "Kgs", and the "Measured Height" field contains a period "." followed by "Cms". At the bottom of the window, there are navigation buttons: "End of Section", "Close Exam", and "Finish", along with a status bar indicating "1 of 10" and a blue arrow icon.

- If the weight entered is below the edit limits, the ISIS system will display a message: "The weight is outside the edit range for this age group and gender. Are you sure?"
- If the weight you entered is correct, select 'yes.' If the data was entered incorrectly, select 'no' and re-enter the correct value for weight.

Exhibit 4-5. Edit message on incorrect data entry for height

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm". The header bar contains the following information: "SP ID: 818591", "Name: BROWN, MAC", "Age: 38 years", "Gender: Male", "Date: 10/10/2000", and "Time: 07:57 PM". Below the header, the section "SP Height and Weight Information" is displayed. A modal dialog box titled "Cardiovascular Fitness" is open, containing a question mark icon and the text: "The height is outside the edit range for this age group and gender. Are you sure?". Below the text are two buttons: "Yes" and "No". In the background, the "Measured Weight" field is set to "75." Kgs and the "Measured Height" field is set to "200." Cms. At the bottom of the window, there are navigation controls including "End of Section", "Close Exam", and "Finish" buttons, along with a status bar indicating "1 of 10" and a blue arrow icon.

- If the height entered is outside the edit limits, the ISIS system will display a message: "The weight is outside the edit range for this age group and gender. Are you sure?"
- If the height you entered is correct, select 'yes.' If the data was entered incorrectly, select 'no' and re-enter the correct value for weight.

Exhibit 4-6. Required data entry for height and weight

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm". Below the title bar, a status bar displays "SP ID: 818591 Name: BROWN, MAC Age: 38 years Gender: Male Date: 10/10/2000 Time: 07:57 PM". The main area is titled "SP Height and Weight Information". A pop-up dialog box titled "Cardiovascular Fitness" is overlaid on the screen. The dialog contains an information icon and the text: "You must enter weight and height information to allow the ISIS system to calculate protocol. If SP refuses to answer, press Cancel to exclude from the exam." Below this text are "OK" and "Cancel" buttons. In the background, the "Measured Weight" field contains "75." and "Kgs", and the "Measured Height" field is empty with "Cms" next to it. At the bottom of the main window, there are navigation buttons: "End of Section", "Close Exam", and "Finish", along with a "1 of 10" indicator and a blue arrow button.

- If the height and/or weight are not entered and the Next button is pressed, the ISIS system will display a message 'You must enter weight and height information to allow the ISIS system to calculate protocol. Press OK to continue to enter the weight and height information.
- If SP refuses to answer or you cannot obtain the height and weight for some other reason such as equipment problems, press the Cancel button to exclude the SP from the exam.

Exhibit 4-7. Cannot calculate VO₂ max

The screenshot shows a dialog box titled "Cardiovascular Fitness". It contains an information icon and the text: "Cannot calculate predicted VO2 max: BMI information not available." Below the text is an "OK" button.

- If the Cancel button is selected, the SP will be excluded. The ISIS system will display a message to indicate that the predicted VO₂ max cannot be calculated if the BMI information is not available. Click OK to this message.

Exhibit 4-8. Partial status with comment 'Other, specify'

Cardiovascular Fitness: Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm

SP ID: 531764 Name: ELLA, CINDER Age: 25 years Gender: Female Date: 07/04/2000 Time: 08:46 AM

Cardiovascular Fitness Component Status

Status

☐ Complete

☒ Partial

☐ Not Done

Comments other, specify

Other text Unable to get ht & wt information - SP refused

9 of 9

End of Section Close Exam Finish

- The ISIS system will go to the Component Status screen and the Status will be set to Partial.
- If the reason the BMI could not be calculated was due to SP refusal to have weight and height measured, select 'Other, specify' for comments and enter 'Unable to get ht and wt information - SP refusal' in the 'Other text' field.
- If the reason the BMI could not be calculated was another reason, enter the appropriate information in the 'Other text' field.

4.2

Questions to Determine Exercise Protocol (PAR Code Questions)

Exhibit 4-9. Questions to determine protocol

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

File View Utilities Reports Window Help

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:04 AM

Questions To Determine Exercise Protocol

I am going to ask you several questions to describe your TYPICAL PHYSICAL ACTIVITY LEVEL. Please try to answer as accurately as possible. Please wait until I read all three questions before you give me your answer.

1. Would you say that you did LITTLE or NO REGULAR recreation, sport, or physical activity?
OR
2. Would you say that you participated REGULARLY in recreation OR WORK requiring MODERST physical activity?
OR
3. Would you say that you participated REGULARLY in HEAVY physical activity?

If necessary use the following prompts to describe modest or heavy physical activity.
 MODEST physical activities include: walking, calisthenics, bowling, weight lifting, yard work.
 HEAVY physical activities include: jogging, swimming, cycling, rowing, tennis, basketball.

If SP's activity level over the last few months is different from their activity level in general, use the activity level which best describes their typical activity level over the last year.

If the SP is still not sure in which category they fall, (little or modest, modest or heavy), select the lower category.

Enter/Select answer:

1 of 8

End of Section Close Exam Finish

Ready

- The questions on the next few screens are designed to determine the exercise protocol that will be assigned to the SP.
- These questions are based on the Physical Activity Readiness (PAR) Questionnaire and will be used to determine a PAR Code that is used to calculate the predicted VO_2 max.
- The predicted VO_2 max is used to determine the exercise protocol for each individual SP. See Appendix F for more details.
- Exhibit 4-9 shows the first level questions used to determine the SP's regular physical activity level. These questions should be read exactly as written.
- If the SP interrupts and gives an answer before you have completed asking the question, tell the SP that you are required to read the entire set of options before you accept a response.
- Read the options clearly and slowly enough for the SP to understand what the options are. Put special emphasis on the words in capital letters.

- If the SP is not sure about the difference between heavy and modest activities, use the examples listed for these activities.
- If the SP's activity level over the last few months is different from his/her activity level in general, use the activity level which best describes the activity level over the last year.
- If the SP is still not sure in which category they fall, (little or modest, modest, or heavy), select the lower category.
- Encourage the SP to answer as honestly as possible to obtain the exercise protocol that will get valid information about their fitness level.
- Select 1, 2, or 3 from the drop-down menu.
- If the SP refuses to answer or is unable to give an answer, select Refused or Don't Know. If 'Refused' or 'Don't Know' are selected, the SP will be excluded from the exam.
- When the Next button is pressed, the ISIS system will go to the Component Status Screen. The Status will be set to 'Not Done' with the comment 'Unable to calculate protocol.'

4.2.1 First Level Questions

Exhibit 4-10. First level questions to determine exercise protocol

The screenshot shows a software window titled "Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm". The window has a menu bar with "File", "View", "Utilities", "Reports", "Window", and "Help". Below the menu bar is a toolbar with various icons. The main content area is titled "Questions To Determine Exercise Protocol". It contains the following text:

I am going to ask you several questions to describe your TYPICAL PHYSICAL ACTIVITY LEVEL. Please try to answer as accurately as possible. Please wait until I read all three questions before you give me your answer.

1. Would you say that you did LITTLE or NO REGULAR recreation, sport, or physical activity?
- OR
2. Would you say that you participated REGULARLY in recreation OR WORK requiring MODEST physical activity?
- OR
3. Would you say that you participated REGULARLY in HEAVY physical activity?

If necessary use the following prompts to describe modest or heavy physical activity.
MODEST physical activities include: walking, calisthenics, bowling, weight lifting, yard work.
HEAVY physical activities include: jogging, swimming, cycling, rowing, tennis, basketball.
If SP's activity level over the last few months is different from their activity level in general, use the activity level which best describes their typical activity level over the last year.
If the SP is still not sure in which category they fall, (little or modest, modest or heavy), select the lower category.

Enter/Select answer: 1

1
2
3
Refused
Don't Know

End of S Finish

Close the exam

Little or no recreation, sport, or physical activity: If the SP is not regularly involved in recreation, sport, or physical activity, select 1 from the drop-down menu. If #1 (little or no regular physical activity) is selected from the first level activity questions in Exhibit 4-10 above, the ISIS system will advance to the screen for the second level activity questions for 'little or no regular physical activity.' See Exhibit 4-11.

Regular modest physical activity: If the SP participates regularly in modest physical activity from recreation or work, select 2 from the drop down menu. The activity should be typical physical activity and should be done on a regular basis rather than sporadically or in an earlier time period. Emphasize the words in capital letters when reading the choices. If #2 (regular, modest physical activity) is selected from the first level activity questions in Exhibit 4-10 above, the ISIS system will advance to the screen for the second level activity questions for 'regular modest physical activity.' See Exhibit 4-12.

Heavy physical activity: If #3 (regular, heavy physical activity) is selected from the first level activity questions in Exhibit 4-10 above, the ISIS system will advance to the screen for the second level activity questions for ‘regular, heavy physical activity.’ See Exhibit 4-13 on page 4-14.

4.2.2 Second Level Questions (Little or No Regular Physical Activity)

Exhibit 4-11. Little or no regular physical activity (second level questions)

The screenshot displays a software window titled "Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm". The menu bar includes File, View, Utilities, Reports, Window, and Help. The toolbar contains various icons for file operations and help. Below the menu bar, a status bar shows "SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:14 AM". The main content area is titled "Questions To Determine Exercise Protocol (Cont...)" and contains two sets of questions. The first set asks: "1 Would you say that you AVOID walking or exertion, (for example, always use the elevator, drive whenever possible instead of walking.)?" followed by "OR" and "2 Would you say that you walk for pleasure, routinely use stairs, occasionally exercise sufficiently to cause heavy breathing or perspiration?". The second set asks: "Would you say that you participated regularly in recreation or modest physical activity for:" followed by "1. 10 TO 60 MINUTES PER WEEK?" and "OR" and "2. MORE THAN 60 MINUTES PER WEEK?". At the bottom, there is a section labeled "Enter/Select answer:" with a dropdown menu showing "1" selected. Below the dropdown are buttons for "End of Section", "Close Exam", and "Finish". A navigation bar at the bottom left shows "2 of 8" and a "Ready" status.

- If 1 is selected from the first level activity questions (see Exhibit 4-10 shown earlier), the ISIS system will advance to the second level activity questions for ‘little or no regular activity.’ The second level questions for this category will be highlighted.
- Read both highlighted options entirely before accepting the SP’s response. As with the previous options, if the SP interrupts with an answer before you have completed reading the options, remind them that you are required to read the options entirely before accepting their response.
- These questions are used to determine if the SP avoids walking or exertion such as always using an elevator or driving whenever possible instead of walking. The other option is if he/she walks for pleasure, routinely uses stairs, or occasionally exercises sufficiently to cause heavy breathing or perspiration.

- If the SP is having difficulty choosing between the options, read the options again to give them time to decide. Do not answer the question for the SP. Ask him/her to choose the option that is closest to their usual activity level. If the SP refuses to select an answer or is unable to give a response, select 'Refused' or 'Don't Know' from the drop-down menu. If 'Refused' or 'Don't Know' is selected, the SP is excluded.
- The ISIS system will go to the Component Status screen. The status will be set to 'Not Done' with the comment 'Unable to calculate protocol.' Most SP's will be able to select one of the options. Select 1 or 2 from the drop-down menu. Press the Next button to advance to the next screen.

4.2.3 Second Level Questions (Regular Modest Physical Activity)

Exhibit 4-12. Regular modest physical activity (second level questions)

The screenshot displays a software window titled "Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm". Below the title bar, a header section contains the following information: "SP ID: 836308", "Name: Anderson, Greg", "Age: 28 years", "Gender: Male", "Date: 01/03/1999", and "Time: 10:14 AM". The main content area is titled "Questions To Determine Exercise Protocol (Cont...)" and contains two sets of questions. The first set includes two options: "1. Would you say that you AVOID walking or exertion, (for example, always use the elevator, drive whenever possible instead of walking.)?" and "2. Would you say that you walk for pleasure, routinely use stairs, occasionally exercise sufficiently to cause heavy breathing or perspiration?". The second set includes two options: "1. 10 TO 60 MINUTES PER WEEK?" and "2. MORE THAN 60 MINUTES PER WEEK?". Below the questions is a label "Enter/Select answer:" followed by a dropdown menu. At the bottom of the window, there is a navigation bar with buttons for "End of Section", "Close Exam", and "Finish", along with a progress indicator showing "2 of 8" and a blue arrow button.

- If 2 is selected from the first level activity questions (see Exhibit 4-10 shown earlier), the ISIS system will advance to the second level activity questions for regular, modest activity. The second level questions for this category will be highlighted.
- Read both highlighted options entirely before accepting the SP's response. As with the previous options, if the SP interrupts with an answer before you have completed reading the options, remind them that you are required to read the options entirely before accepting their response.

- These questions are used to determine the amount of time the SP spends in regular, modest physical activity per week.
- If they spend 10 to 60 minutes per week in regular, modest physical activity, select 1 from the drop-down menu.
- If they spend more than 60 minutes per week in regular, modest physical activity, select 2 from the drop-down menu.
- If the SP is having difficulty choosing between the options, read the options again to give them time to decide. Do not answer the question for them. Ask them to choose the option that is closest to their usual activity level.
- If the SP refuses to select an answer or is unable to give a response, select 'Refused' or 'Don't Know' from the drop-down menu. If 'Refused' or 'Don't Know' is selected, the SP is excluded.
- The ISIS system will go to the Component Status screen. The status will be set to 'Not Done' with the comment 'Unable to calculate protocol.' Most SP's will be able to select one of the options. Select 1 or 2 from the drop-down menu. Press the Next button to advance to the next screen.

4.2.4 Second Level Questions (Regular Heavy Physical Activity)

Exhibit 4-13. Second level questions (regular heavy physical activity)

The screenshot displays the 'Cardiovascular Fitness Subsystem' window. The title bar indicates 'Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm'. The menu bar includes 'File', 'View', 'Utilities', 'Reports', 'Window', and 'Help'. The toolbar contains various icons for file operations and help. The main window area shows the following information:

SP ID: 666259 Name: CHRISTIAN, TIM Age: 31 years Gender: Male Date: 09/26/2000 Time: 03:55 PM

Questions To Determine Exercise Protocol (Cont...)

Now I am going to ask for some more detail on the amount of heavy physical exercise you did. I am going to give you four options.

Would you say that you participated regularly in heavy physical activity for:

1. LESS THAN 30 MINUTES PER WEEK
- OR
2. 30 TO 60 MINUTES PER WEEK
- OR
3. 1 TO 3 HOURS PER WEEK
- OR
4. MORE THAN 3 HOURS PER WEEK

Below the options is a text input field labeled 'Enter/Select answer:' with a dropdown menu showing the selected option '4'. The dropdown menu also lists options 1, 2, 3, and 4, along with 'Refused' and 'Don't Know'. At the bottom of the window, there are navigation buttons: 'End of', 'Refused', 'Don't Know', and 'Finish'. The status bar at the bottom shows 'Ready', 'MEC Layer: 09/21/2000', 'Application: Ver 8.9.26B', 'Not connected to Coordinator', and '04:07 PM'. The taskbar at the very bottom shows the 'Start' button and several open applications: 'Audiocasi', 'Anthropometry Sub...', 'Exploring - cvfitness', and 'Cardiovascular...'.

- If 3 is selected from the first level activity questions (see Exhibit 4-10 shown earlier), the ISIS system will advance to the second level questions for regular, heavy physical activity.
- The second level questions for regular, heavy physical activity have four options to determine the amount of time per week spent in this type of activity.
- These options are used to determine the amount or degree of regular, heavy physical activity the SP does.
- Read all four highlighted options before accepting the SP's response. Read the options carefully and slowly enough for the SP to understand the time period associated with each response.
- If they spend less than 30 minutes/week in regular, heavy physical activity, select 1.
- If they spend 30-60 minutes per week in regular, heavy physical activity, select 2.

- If they spend 1-3 hours per week in regular, heavy physical activity, select 3.
- If they spend more than 3 hours per week in regular, heavy physical activity, select 4.
- If the SP is having difficulty choosing between the options, read the options again to give them time to decide. Do not answer the question for them. Ask them to choose the option that is closest to their usual activity level.
- If the SP refuses to select an answer or is unable to give a response, select 'Refused' or 'Don't Know' from the drop-down menu. If 'Refused' or 'Don't Know' is selected, the SP is excluded. The ISIS system will go to the Component Status screen. The status will be set to 'Not Done' with the comment 'Unable to calculate protocol.' Most SP's will be able to select one of the options. Select 1, 2, 3, or 4 from the drop-down menu.

4.2.5 Required Data Entry for PAR Questions

Exhibit 4-14. Required data entry for first level activity questions

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:33 AM

Questions To Determine Exercise Protocol

I am going to ask you several questions to describe your TYPICAL PHYSICAL ACTIVITY LEVEL. Please try to answer as accurately as possible. Please wait until I read all three questions before you give me your answer.

1. Would you say that you did LITTLE or NO REGULAR recreation, sport, or physical activity?
- OR
2. Would you say that you did MODERATE or HEAVY physical activity?
- OR
3. Would you say that you did MODERATE or HEAVY physical activity?

If necessary use the following definitions:

MODEST physical activities include: walking, canoeing, bowling, weight lifting, yard work.

HEAVY physical activities include: jogging, swimming, cycling, rowing, tennis, basketball.

If SP's activity level over the last few months is different from their activity level in general, use the activity level which best describes their typical activity level over the last year.

If the SP is still not sure in which category they fall, (little or modest, modest or heavy), select the lower category.

Enter/Select answer:

Cardiovascular Fitness

You must enter selection to proceed to the next slide.

OK

1 of 8

End of Section Close Exam Finish

Ready

- If a response is not selected for the physical activity questions and the Next button is pressed, the ISIS system will display a reminder message 'You must enter selection to proceed to the next slide.'
- Click OK and enter the response.
- Press the Next button to advance to the next screen.

4.3 Pretest Measurements

Exhibit 4-15. Pretest measurements

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm". The header bar contains the following information: "SP ID: 618864", "Name: PARKER, JACK", "Age: 21 years", "Gender: Male", "Date: 08/14/1999", and "Time: 07:15 AM". The main area is titled "Pretest Measurements". It features a "Room Temperature:" label followed by a text box containing ".00" and the word "Degrees". Below this is a section titled "Exam Preparation" containing three dropdown menus: "Arm:" (empty), "Cuff Size:" (empty), and "Mesh Vest Used:" (set to "No"). Another section titled "Pretest Measurements" contains three text boxes for "Heart Rate:", "Systolic BP:", and "Diastolic BP:". To the right of these boxes is a "Capture" button. At the bottom of the window, there is a navigation bar with left and right arrow buttons, a "4 of 9" indicator, and three buttons: "End of Section", "Close Exam", and "Finish". A blue arrow icon is located in the bottom right corner.


- A digital thermometer is directly integrated with the ISIS system. The system will capture the room temperature and display it on the screen. The system will then give an auditory message:

“Check footwear, attach chest electrodes and BP cuff. Then explain treadmill test.”

- The initial footwear check is done as soon as the SP comes in the room. This is another reminder to check the footwear again for proper type and to make sure laces are tied if appropriate. If the footwear is not appropriate, get the shoe size and give the SP a pair of running shoes from the stock supply. (See Section 2.6.2 for discussion of appropriate footwear and Section 3.4 for discussion of checking footwear prior to the test.)
- Attach the chest electrodes and blood pressure cuff. (See Section 3.8.1 to 3.8.3 for discussion of proper placement of electrodes and cuff.) If the SP's chest hair interferes with the attachment of the electrodes, do not shave the site to avoid unintended injury or discomfort to the SP. If the SP has chest hair that makes the adherence of the electrodes on the chest difficult, ask the SP to put on a mesh vest that will assist in keeping the electrodes attached. If the mesh vest is used, check 'yes' to Mesh Vest. This field is defaulted to 'No' and is disabled for females. See Section 3.8.5.

4.3.1 Digital Thermometer for Measuring Room Temperature

Exhibit 4-16. Unable to obtain room temperature

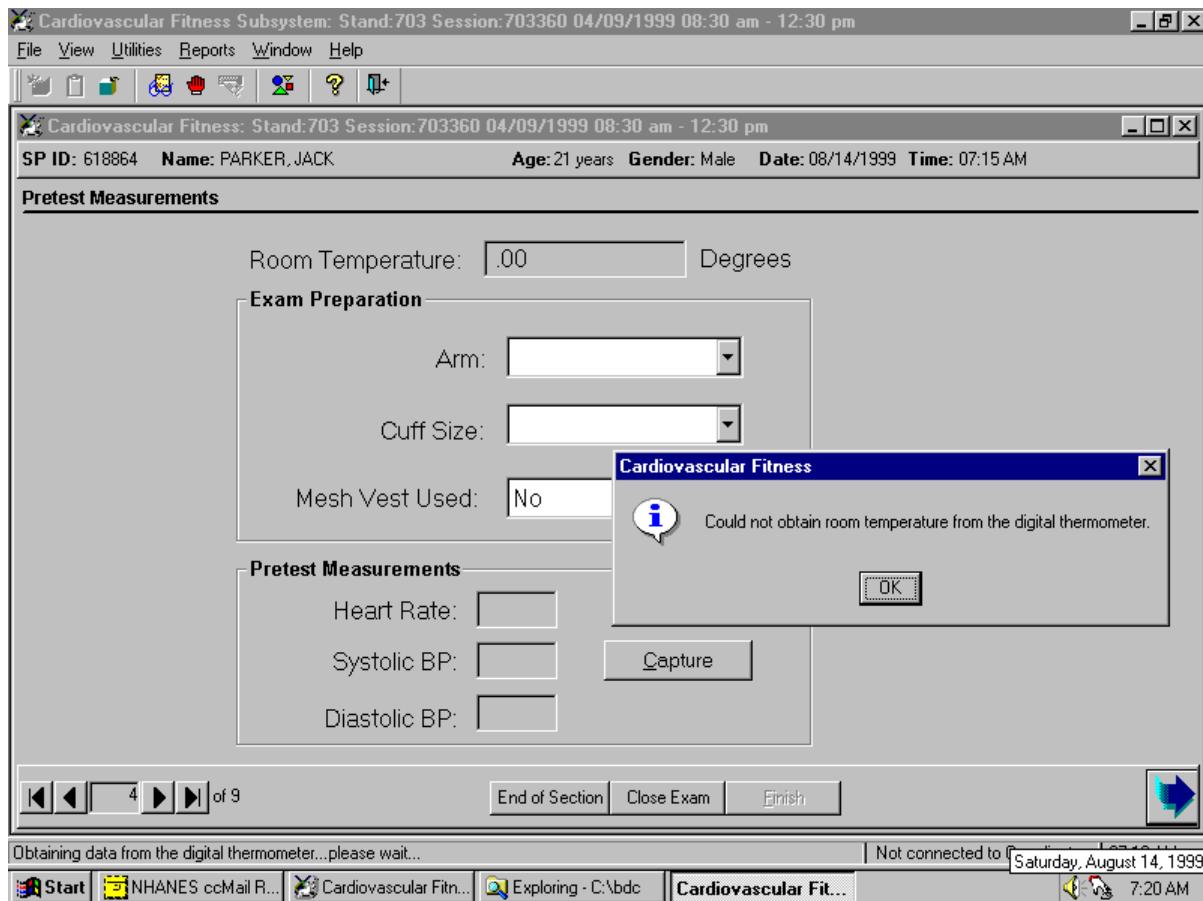

- If the digital room thermometer does not capture the room temperature, the ISIS system will display a message: 'Could not obtain room temperature from the digital thermometer.'
- Click OK.
- Check the room temperature on the digital display and type in the temperature in the appropriate field on the screen.

4.3.2 Exam Preparation –Arm Used, Cuff Size, Mesh Vest Used

Exhibit 4-17. Exam preparation – arm selected

Cardiovascular Fitness Subsystem: Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm

File View Utilities Reports Window Help

SP ID: 934124 Name: THORNTON, JOHN Age: 25 years Gender: Male Date: 07/08/2000 Time: 01:03 PM

Pretest Measurements

Room Temperature: 68.00 Degrees

Exam Preparation

Arm: Right
Cuff Size: Left
Mesh Vest Used: No

Pretest Measurements

Heart Rate:
Systolic BP:
Diastolic BP:

Capture

4 of 9

End of Section Close Exam Finish

Ready Not connected to Coordinator 01:25 PM

Start Cardiovascular Fitness... Exploring - Beryl 13:25


- Select the arm used for blood pressure measurements from the drop-down menu. The options are right, left, and could not obtain.
- The right arm will be used for all BP measurements unless there is a reason for not using this arm. See Section 3.8.2 for discussion of selection of arm to use for blood pressure measurement.

Exhibit 4-18. SP excluded if neither arm can be used for BP measurement

The screenshot shows the 'Cardiovascular Fitness' software window. At the top, it displays 'Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm'. Below this, a header bar contains 'SP ID: 531764', 'Name: ELLA, CINDER', 'Age: 25 years', 'Gender: Female', 'Date: 07/04/2000', and 'Time: 09:27 AM'. The main area is titled 'Pretest Measurements'. Under 'Exam Preparation', there are three dropdown menus: 'Arm:' set to 'Could Not Obtain', 'Cuff Size:', and 'Mesh Vest Used:' set to 'No'. A 'Room Temperature:' field shows '68.00' Degrees. A 'Diastolic BP:' field is empty. A warning dialog box is overlaid on the screen with the title 'Cardiovascular Fitness'. It contains a yellow warning icon and the text: 'SP will be excluded from CV Fitness if BP cannot be measured. Are you sure?'. There are 'Yes' and 'No' buttons. At the bottom of the main window, there are navigation buttons: 'End of Section', 'Close Exam', and 'Finish'. A status bar at the very bottom shows '4 of 9' and a blue arrow icon.

- If you are unable to put a cuff on either arm, select 'could not obtain.' If the blood pressure could not be monitored, the SP will be excluded due to physical limitation.
- The ISIS system will display a message 'SP will be excluded from CV Fitness if BP cannot be measured. Are you sure?' Select Yes if this is true.

Exhibit 4-19. SP excluded due to physical limitation

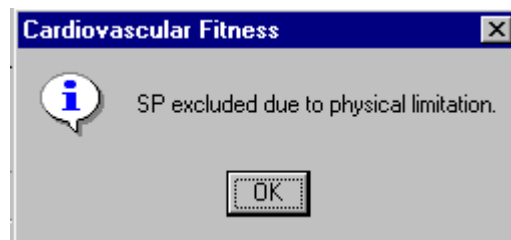

- If Yes is selected, the ISIS system will display a message: "The SP will be excluded due to physical limitation." Select OK.
- The Component Status screen will be displayed.

Exhibit 4-20. Partial status with physical limitation if CNO selected for arm

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm". Below the title bar, a header section contains the following information: "SP ID: 427581", "Name: DONOHUE, MAURIE", "Age: 23 years", "Gender: Female", "Date: 07/12/2000", and "Time: 03:33 PM". The main content area is titled "Cardiovascular Fitness Component Status". Inside this area, there is a "Status" section with three radio buttons: "Complete", "Partial" (which is selected), and "Not Done". Below the status section, there is a "Comments" section with a dropdown menu showing "physical limitation" and an "Other text" input field. At the bottom of the window, there is a navigation bar with buttons for "End of Section", "Close Exam", and "Finish", along with a "9 of 9" indicator and navigation arrows.

- If you selected 'Could Not Obtain' for arm used, the Status will be set to 'Partial' with the comment 'physical limitation.'
- Press the Finish button to end the exam.

Exhibit 4-21. Selection of cuff size

Cardiovascular Fitness Subsystem: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

SP ID: 427581 Name: DONOHUE, MAURIE Age: 23 years Gender: Female Date: 07/11/2000 Time: 05:43 PM

Pretest Measurements

Room Temperature: 68.00 Degrees

Exam Preparation

Arm: Right

Cuff Size: child
adult
large arm

Mesh Vest Used:

Pretest Measurements

Heart Rate:

Systolic BP: Capture

Diastolic BP:

4 of 9

End of Section Close Exam Finish

Ready Not connected to Coordinator 05:52 PM

- Select the cuff size used from the 'Cuff size' drop-down menu. The sizes are 'child,' 'adult,' and 'large arm.'
- Most problems with cuff placement will be a problem with the arm rather than the cuff and 'Could not obtain' should be selected for arm as described above. There may be some situations where the arm can be used but for some reason the cuff cannot be positioned as described in Section 3.8.3. Make every effort to attempt to position the cuff correctly before deciding to exclude the SP for this reason.
- If you are unable to position a cuff correctly on the SP's arm, the SP will be excluded.
- Press the 'Close Exam' button at the bottom of the screen. The status screen will be displayed.

Exhibit 4-22. Problem with cuff size

The screenshot displays the 'Cardiovascular Fitness' software interface. At the top, a status bar shows 'Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm'. Below this, a header bar contains 'SP ID: 531764 Name: ELLA, CINDER', 'Age: 25 years', 'Gender: Female', 'Date: 07/04/2000', and 'Time: 09:27 AM'. The main window is titled 'Pretest Measurements'. It features a 'Room Temperature' field set to '68.00' Degrees. Below this is an 'Exam Preparation' section with 'Arm' set to 'Right' and 'Cuff Size' set to an empty dropdown. A 'Mesh Vest Used' checkbox is also present. To the right, a 'Pretest Measurements' section includes fields for 'Heart Rate:', 'Systolic BP:', and 'Diastolic BP:'. A modal dialog box titled 'Cardiovascular Fitness' is open, showing 'SP ID:531764 SP Name:ELLA, CINDER'. It has an 'Exam Status' field set to 'partial' and a 'Comments' dropdown menu with 'Problem with BP cuff' selected. Below the dropdown is an 'Other text' field. At the bottom of the dialog, a blue message reads: 'This will close the exam, saving all results and the exam status. Do you wish to close the exam?'. There are 'OK' and 'Cancel' buttons. The main window's bottom bar shows navigation controls, a page indicator '4 of 9', and buttons for 'End of Section', 'Close Exam', and 'Finish'.

- Select 'Close Exam' and the Status screen will be displayed.
- Select 'Problems with BP cuff' or "Other, specify" and enter the reason in the 'Other text' field.
- Select OK to close the exam.

Exhibit 4-23. Mesh vest used

Cardiovascular Fitness Subsystem: Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness: Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm

SP ID: 934124 Name: THORNTON, JOHN Age: 25 years Gender: Male Date: 07/08/2000 Time: 01:03 PM

Pretest Measurements

Room Temperature: 68.00 Degrees

Exam Preparation

Arm: Right

Cuff Size: adult

Mesh Vest Used: No

Yes
No
Refused

Pretest Measurements

Heart Rate:

Systolic BP: Capture

Diastolic BP:

4 of 9

End of Section Close Exam Finish

Ready Not connected to Coordinator 01:21 PM

Start Cardiovascular Fitness... Exploring - Beryl 13:23

- If the SP is male, the field for Mesh Vest Used will be enabled. A mesh vest is used to optimize adherence of the chest electrodes in men with a large amount of chest hair. See Section 3.8.5.
- Select Yes or No from the drop-down menu.
- If the SP is female, this field will be disabled.

4.3.3 Capturing Blood Pressure Measurements

Exhibit 4-24. Blood pressure data capture

The screenshot shows a software window titled "Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm". The header bar contains patient information: "SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM". The main area is titled "Pretest Measurements" and contains two sections. The "Exam Preparation" section has three dropdown menus: "Arm:" set to "Right", "Cuff Size:" set to "adult", and "Mesh Vest Used:" set to "No". The "Pretest Measurements" section has three input fields: "Heart Rate:" with the value "60", "Systolic BP:" with the value "118", and "Diastolic BP:" with the value "78". A "Capture" button is located to the right of the Systolic BP field. At the bottom, there is a navigation bar with left and right arrow buttons, a "4 of 9" indicator, and three buttons: "End of Section", "Close Exam", and "Finish". A blue arrow icon is in the bottom right corner.

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Pretest Measurements

Room Temperature: 68.00 Degrees

Exam Preparation

Arm: Right

Cuff Size: adult

Mesh Vest Used: No

Pretest Measurements

Heart Rate: 60

Systolic BP: 118

Diastolic BP: 78

Capture


4 of 9

End of Section Close Exam Finish

- When the cuff is attached to the arm, ask the SP to relax his/her arm by the side.
- Press the 'Capture' button to initiate the blood pressure readings.
- The ISIS system will initiate the blood pressure measurement and capture the reading.
- The readings will be displayed on the screen.

4.3.4 Demonstrating and Explaining How to Walk on a Treadmill

Exhibit 4-25. Demonstrate and explain how to walk on a treadmill

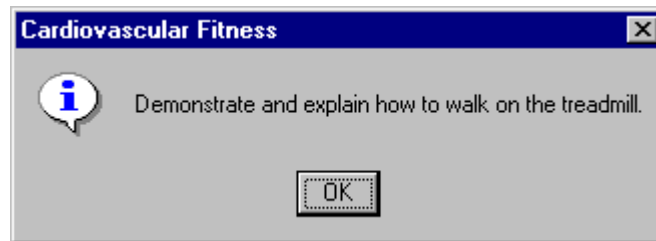

- When the pretest heart rate and blood data are captured, the ISIS system displays a second visual message:

‘Demonstrate and explain how to walk on a treadmill.’

- Give a brief demonstration of how to walk on a treadmill. See Section 3.9.5.
- Press “Start Belt” and increase the speed of the treadmill to 2.0 mph.
- Explain the procedures step-by-step as you demonstrate.
- Straddle the belt and hold onto the handrail with both hands.
- Tell them to ‘paw’ the belt with one foot to get a feel for the movement of the belt and then step on the treadmill and begin to walk.
- Place emphasis on the importance of an upright posture, looking straight ahead nor down at the feet, and taking long steps with a relaxed gait. Drop one arm by the side and then the other one.
- Step off the belt and press 'Stop Belt' and then 'Clear' on treadmill controller.
- Press OK to the message on the ISIS screen.
- Press the Next button to advance to the treadmill test screen.

4.3.5 Required Data Entry for Arm and Cuff Size

Exhibit 4-26. Required data entry for arm and cuff used

The screenshot displays the 'Cardiovascular Fitness Subsystem' window. At the top, it shows 'Stand: 703 Session: 703360 04/09/1999 08:30 am - 12:30 pm'. Below this is a menu bar (File, View, Utilities, Reports, Window, Help) and a toolbar. A patient information bar contains: 'SP ID: 618864', 'Name: PARKER, JACK', 'Age: 21 years', 'Gender: Male', 'Date: 08/14/1999', and 'Time: 07:15 AM'. The main area is titled 'Pretest Measurements'. It includes a 'Room Temperature' field with '68.00' and 'Degrees'. An 'Exam Preparation' section contains 'Arm:' and 'Cuff Size:' dropdown menus, and a 'Mesh Vest Used:' field with 'No'. Below this is another 'Pretest Measurements' section with 'Heart Rate: 60', 'Systolic BP: 118', and 'Diastolic BP: 78'. A 'Capture' button is next to the BP fields. A pop-up dialog box titled 'Cardiovascular Fitness' is overlaid, displaying an information icon and the message: 'You must select arm and cuff size to proceed to the next slide.' with an 'OK' button. At the bottom of the main window are navigation buttons: 'End of Section', 'Close Exam', and 'Finish'. The status bar at the very bottom shows 'Ready', 'Not connected to Coordinator', and '07:24 AM'. The Windows taskbar at the bottom includes the Start button, 'NHANES ccMail Router', 'Cardiovascular Fitness...', and 'Exploring - C:\bdc'.

- You cannot advance to the next screen if all data fields are not complete. If you do not select the arm or cuff size used, the ISIS system will display a message:

'You must select arm and cuff size to proceed to the next slide.'

- Click OK.
- Select arm and cuff size.

4.3.6 Required Data Entry for BP Measurements

Exhibit 4-27. Required data entry for blood pressure

The screenshot displays the 'Cardiovascular Fitness' software window. The title bar reads 'Cardiovascular Fitness: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm'. The main window has a header with patient information: 'SP ID: 427581 Name: DONOHUE, MAURIE Age: 23 years Gender: Female Date: 07/11/2000 Time: 05:43 PM'. Below this is a section titled 'Pretest Measurements'. It contains a 'Room Temperature' field with the value '68.00' and a 'Degrees' label. An 'Exam Preparation' box includes three dropdown menus: 'Arm' set to 'Right', 'Cuff Size' set to 'adult', and 'Mesh Vest Used' set to 'No'. Below this is another 'Pretest Measurements' box with three input fields for 'Heart Rate', 'Systolic BP', and 'Diastolic BP'. A 'Capture' button is positioned to the right of these fields. A small message dialog box titled 'Cardiovascular Fitness' is overlaid on the right, containing an information icon, the text 'Obtain Pretest BP Measurements.', and an 'OK' button. The bottom of the window features a navigation bar with buttons for 'End of Section', 'Close Exam', and 'Finish', along with a status indicator showing '4 of 9'.

- See Section 3.9.1 for setup of blood pressure monitor. If the Next button is pressed before the pretest data are captured, the ISIS system will display a message: 'Obtain Pretest BP measurements.'
- Click OK.
- Click on the 'Capture' button to obtain the heart rate and blood pressure readings.

4.4 Treadmill Test

Exhibit 4-28. Treadmill test data entry screen

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:04 PM

Treadmill Test

00:00
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE:

Captured HR:

Manual HR:

PMHR: 184.00 beats/min

%HR max: 92 119.6 184

50% 60% 100%

End of Section Close Exam Finish

4.4.1 Function of Treadmill Test Screen

The treadmill test is completed using the screen in Exhibit 4-28 above. The following sections will describe the buttons and functionality of this data entry screen.

4.4.1.1 Timer

The timer in the upper left of the screen will be used to time the exercise stages. When the 'Start' button is pressed, the timer begins for the Warm-up Stage. The timer counts up to 2 minutes for Warm-up, then resets to zero for the beginning of Stage 1. It then counts up to 3 minutes for Stage 1, resets to zero for the beginning of Stage 2, counts up to 3 minutes for Stage 2, then resets to zero for the beginning of Recovery 1. The timer counts up to 1 minute for Recovery 1 and one minute for Recovery 2

and stops. If the heart rate is greater than 110 after Recovery 2, the timer will count up to 1 minute for Recovery 3 and then stop.

4.4.1.2 Start Button

The 'Start' button is used to start the timer for the exercise stages.

4.4.1.3 Restart Button

The Restart' button is used to restart the timer at zero if the test needs to be restarted. This would be used in a nonemergency situation. For example, if the test has started and the SP's shoe laces become untied, the test can be interrupted and the timer restarted. This function is only used in the Warm-up stage and during the first 15 seconds in Stage 1. The test should not be restarted for any reason after the first 15 seconds in Stage 1. See Section 3.14.

4.4.1.4 Emergency Stop Button

The Emergency Stop is used only to stop the test immediately if the condition of the SP warrants this action. This will be described in greater detail in subsequent screens.

4.4.1.5 Predetermined Protocols

The box to the right of the 'Start' button displays the predetermined protocol with the speed and grade to be used for the Warm-up, Stages 1 and 2, and for the Recovery Stages. The middle protocol is the predetermined protocol that will be used. When the Start button is pressed, the middle protocol will be highlighted in red for the current stage. The next higher and next lower protocols are also displayed in this box.

4.4.1.6 Change Protocol Button

The 'Change Protocol' button is used to increase or decrease the predetermined protocol. This button can only be used during the Warm-up stage.

4.4.1.7 Go To Recovery Button

The 'Go to Recovery' button is used to stop the exercise test at any point and go directly to Recovery. The SP will continue through the Recovery Stages but will not complete the remainder of time in the exercise stages.

4.4.1.8 Show Data Button

The 'Show Data' button is used to look at the heart rate and blood pressure from the previous stages. See Exhibit 4-33 for an example of the 'Show Data' window on the treadmill test screen.

4.4.1.9 Estimated VO₂ Max Field

The 'Estimated VO₂ max' field is used to view this variable after the treadmill test is completed. The VO₂ max is calculated and displayed after Recovery 2 has been completed.

4.4.1.10 Continuum Bar


The continuum bar has the predicted maximal heart rate (PMHR) in beats per minute on the top and the percent of the maximum heart rate (%HR) along the bottom.

- During Warm-up, the SP's heart rate should increase to approximately 60 percent of the maximum heart rate. This range is displayed on the continuum bar. If the heart rate goes higher than this, the protocol should be decreased for Stages 1 and 2.
- During Stage 1, the SP's heart rate should be about 60-80 percent of the maximum heart rate. If it goes higher than this, you should go to Recovery. When the ISIS system changes to Stage 1, the continuum changes to reflect this change.

- During Stage 2, the SP's heart rate should be below 85 percent of the predicted maximal heart rate. The continuum changes to reflect these numbers. EXCEPTION: If the SP is 12 through 18 years old the heart rate should be below 90 percent of the predicted heart rate. The ISIS system will display 85 percent or 90 percent based on the age of the SP.
- During Recovery, the SP's heart rate is expected to decrease back to approximately 110 beats per minute or less. The continuum changes to reflect the 50 percent PMHR. This is used as a guideline.

4.4.2 Preparing the Treadmill and Controller

Exhibit 4-29. Preparing the treadmill and controller

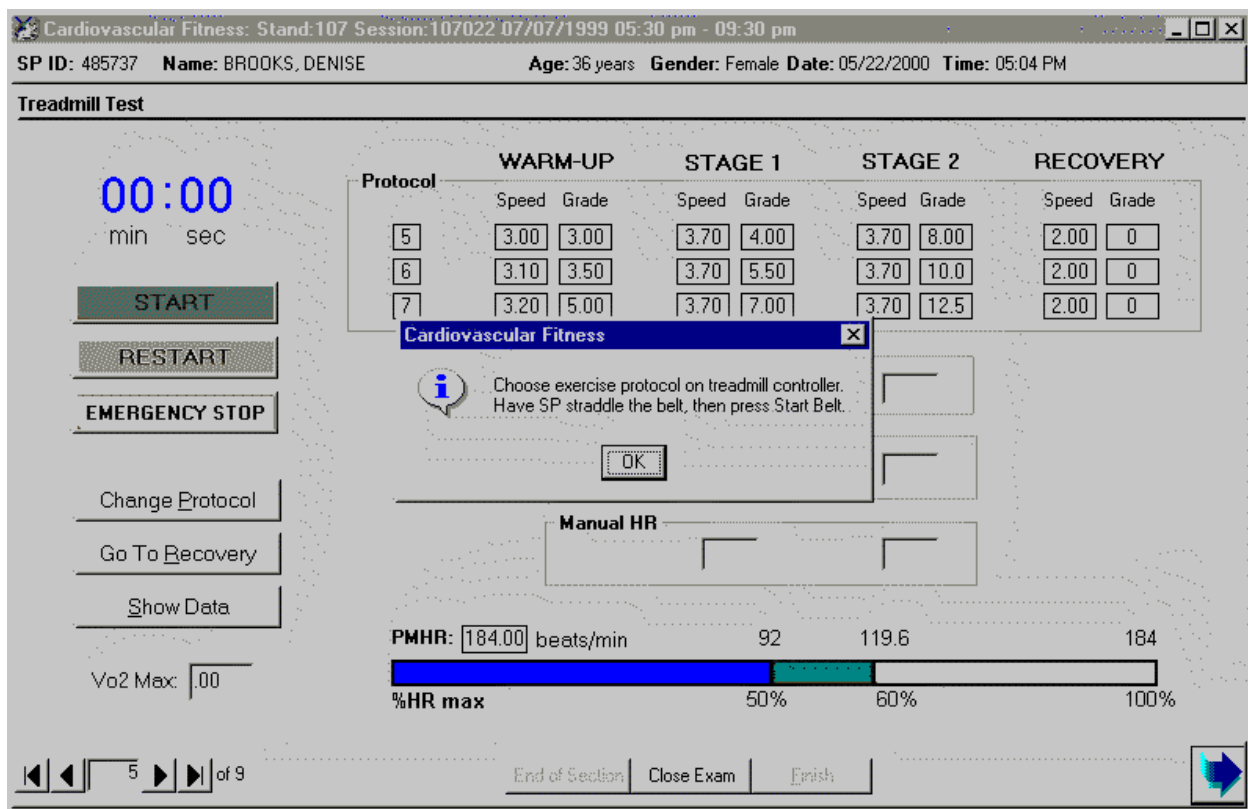

When this screen initially comes up, the ISIS system displays a visual message and gives the same auditory message:

“Choose exercise protocol on treadmill controller. Have SP straddle the belt, then press ‘Start Belt’.”

- The middle protocol is the predetermined protocol for this SP. (This would be protocol # 6 in this example.)
- Press the Clear button on the treadmill controller to reset the treadmill setting.

- Press Input plus (+) or minus (–) to move through protocols 1 to 8. Select the predetermined protocol (the middle protocol on the ISIS screen) on the treadmill controller. (See Section 3.9.2 for details)

4.4.3 Getting the SP on the Treadmill

- You just demonstrated and explained the treadmill test after the BP and HR data were captured. Now you want to get the SP on the treadmill to start the Warm-up phase. You will need to repeat the explanation step-by-step. Don't assume that the SP will know how to walk on a treadmill.
- Have the SP straddle the treadmill belt on treadmill controller.
- Press "Start Belt" and increase the speed of the treadmill to 2.0 mph.
- Repeat the explanation as outlined below.
- Get the SP to straddle the belt and hold onto the handrail with both hands.
- Tell them to 'paw' the belt with one foot to get a feel for the movement of the belt and then instruct them to step on the treadmill and begin to walk.
- Tell them to use an upright posture and look straight ahead. Once they have started, get them to drop one arm by their side and then the other one. Click OK.
- They should now be walking slowly on the treadmill and ready to begin the test.
- At this point they are walking very slowly and you should go to Warm-up as soon as he/she seems to be getting somewhat comfortable walking on the belt.

4.4.4 Warm-up

Exhibit 4-30. Starting the timer

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:04 PM

Treadmill Test

00:01
min sec

START
RESTART
EMERGENCY STOP
Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 184.00 beats/min

%HR max

92 119.6 184

50% 60% 100%

5 of 9

End of Section Close Exam Finish

- The ISIS system will give the following auditory message:
"Simultaneously press 'Start' on the ISIS screen and press 'Start Course' on the treadmill controller."
- Click 'Start' on the ISIS screen. The timer will start timing the 2-minute Warm-up Stage.
- Press 'Start Course' on the treadmill controller.
- The treadmill will increase to the speed and grade for the predetermined Warm-up Stage for the selected protocol.

Exhibit 4-31. Explaining the rating of perceived exertion (RPE)

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:04 PM

Treadmill Test

00:15
min sec

START
RESTART
EMERGENCY STOP
Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 184.00 beats/min

%HR max

92 119.6 184

50% 60% 100%

5 of 9

End of Section Close Exam Finish

- When the timer starts, the Warm-up Stage for the predetermined protocol (the middle protocol on the screen) will be highlighted in red.
- At 15 seconds into Warm-up, the ISIS system will give the auditory message:
"Explain RPE."
- Explain the Rating of Perceived Exertion (RPE) using the standard script. See Section 3.9.7 for procedures for explaining RPE and Section 3.10 for procedures if the RPE exceeds the criteria.
- Throughout this stage, monitor the % HR scale on the ISIS screen and the SP's current heart rate, which is displayed digitally on the BP monitor.
- The SP's heart rate should be below 60 percent of the PMHR.

Exhibit 4-32. Getting RPE for warm-up

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:04 PM

Treadmill Test

01:30
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.70	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 184.00 beats/min 92 119.6 184

%HR max 50% 60% 100%


5 of 9

End of Section Close Exam Finish

- At 1 minute 30 seconds into Warm-up, the ISIS system will give an auditory message:
"Get RPE Assessment."
- Enter RPE score in the box under Warm-up.
- Note the heart rate on the monitor. The ISIS system is capturing the heart rate at this time for the Warm-up Stage. When the stage changes you can check the heart rate that was captured. At 1 minute 50 seconds into Warm-up, the ISIS system will give an auditory message:
"Press 'Next Stage' on treadmill controller at 1 minute 57 seconds."
- Watch the timer and at 1:57 in Warm-up, press 'Next Stage' on treadmill controller. This will advance the treadmill to Stage 1 of the current protocol.

4.4.5 Show Data


Exhibit 4-33. Show data window

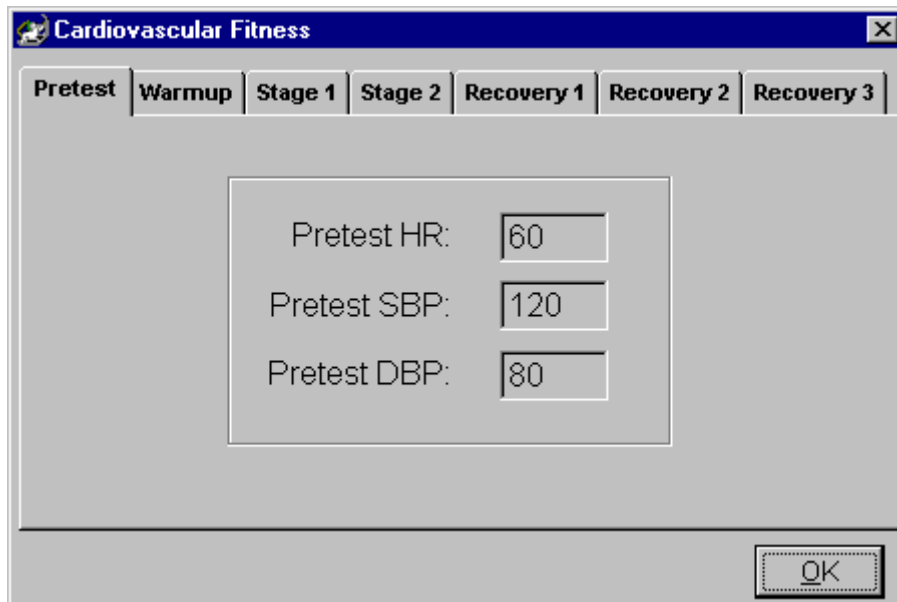

The screenshot shows a window titled "Cardiovascular Fitness" with a tabbed interface. The tabs are "Pretest", "Warmup", "Stage 1", "Stage 2", "Recovery 1", "Recovery 2", and "Recovery 3". The "Pretest" tab is selected. In the center of the window, there is a box containing three labels and their corresponding values in input fields: "Pretest HR:" with the value "60", "Pretest SBP:" with the value "120", and "Pretest DBP:" with the value "80". An "OK" button is located at the bottom right of the window.

- The captured data can be displayed by clicking on the 'Show Data' button. The data for Pretest, Warm-up, Stage 1, Stage 2, and Recovery stages can be viewed by clicking on the tab of the stage you want to see.

Exhibit 4-34. Show data message

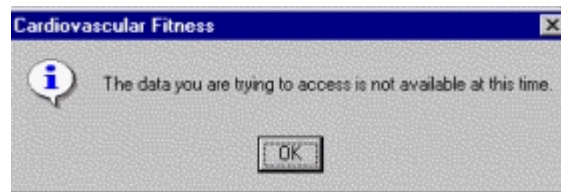

- Only the data from the previous stages can be accessed. If you click on the current stage the ISIS system will display a message: 'The data you are trying to access is not available at this time.' Click OK to this message and select a tab from the previous stage. Once the current stage is completed, the data from that stage can be viewed.

4.4.6 Stage 1

Exhibit 4-35. Treadmill test stage 1 – routine test

Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm

SP ID: 818591 Name: BROWN, MAC Age: 38 years Gender: Male Date: 10/10/2000 Time: 07:57 PM

Treadmill Test

00:20
min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0
8	3.60	5.00	3.70	8.50	3.70	14.5	2.00	0

RPE

9

Captured HR

Manual HR

PMHR: 182.00 beats/min

109.2 145.6 182

%HR max

60% 80% 100%

6 of 10

End of Section Close Exam Finish

- The ISIS system will advance to Stage 1. The timer will be reset to zero and will start timing the 3-minute Stage 1. Throughout this stage, monitor the BP monitor for the SP's heart rate. Watch the continuum bar to confirm that the SP's heart rate is within the 60 – 80 percent range.

Exhibit 4-36. Get rating of perceived exertion (RPE)

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

02:30
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9

Captured HR

Manual HR

PMHR: 184.00 beats/min

%HR max: 60% 80% 100%

110.4 147.2 184


5 of 9

End of Section Close Exam Finish

- At 2 minutes 30 seconds into Stage 1 the ISIS system will give the auditory message:
"Get RPE Assessment."
- Get the SP's assessment of RPE and enter the number in the appropriate field below Stage 1.

4.4.6.1 Manually Entered Heart Rate

Exhibit 4-37. Compare heart rate on monitor and screen

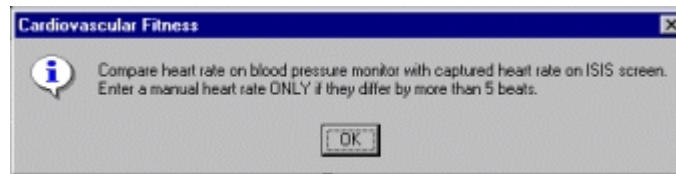

- The BP monitor will automatically measure BP and HR. These values will be captured by the ISIS system. At approximately 2 minutes 40 seconds into Stage 1 the ISIS system will give the auditory message:

“Check heart rate on monitor.”

- The ISIS system will then display a visual message:

“Compare heart rate on blood pressure monitor with captured heart rate on ISIS screen. Enter a manual heart rate ONLY if they differ by more than 5 beats.”

- Note the heart rate on the blood pressure monitor and compare this number with the captured heart rate on the ISIS screen. If the captured heart rate on the ISIS screen does not differ by more than 5 beats from the heart rate on the BP monitor, click OK and continue the test.

Exhibit 4-38. Example 1 – manually entered heart rate - not necessary

Cardiovascular Fitness: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm

SP ID: 819357 Name: JOHNSON, SUE Age: 39 years Gender: Female Date: 07/12/2000 Time: 02:20 PM

Treadmill Test

01:27 min sec

START RESTART EMERGENCY STOP

Change Protocol Go To Recovery Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE 10 11

Captured HR 110

Manual HR

PMHR: 181.00 beats/min 153.8 181

%HR max 85% 100%

End of Section Close Exam Finish

If the heart rate on the blood pressure monitor is within 5 beats of the captured heart rate, it should not be entered manually. The ISIS system will accept the captured heart rate.

- Example 1: In the example above, the ISIS system captured the heart rate on the ISIS screen as 110.

Exhibit 4-39. Show data screen for Example 1

Cardiovascular Fitness

Pretest Warmup Stage 1 Stage 2 Recovery 1 Recovery 2 Recovery 3

End Stage Heart Rate: 110 beats/min

Stage 1 Systolic BP: 140 mmHg

Stage 1 Diastolic BP: 80 mmHg

Stage 1 Time Completed: 3.00 minutes

OK

- If the heart rate on the BP monitor was 112 and the captured heart rate on the ISIS screen is 110 (as in the above example), it is not necessary to enter 112 in the manual heart rate field. The ISIS system will store 110 to the "show data" field and use it as the heart rate for calculation of VO₂max.

Exhibit 4-40. Example 2 – manually entered heart rate - necessary

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

00:06
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9 11

Captured HR: 110

Manual HR: 120

PMHR: 184.00 beats/min 156.4 184

%HR max: 85% 100%

5 of 9 End of Section Close Exam Finish

- Example 2: In the above example, the ISIS system captured the heart rate as 110. The technician noticed that the heart rate was actually 120 at the time of data capture. This is not within 5 beats of the captured value, so the technician manually entered 120 in the heart rate field.

Exhibit 4-41. Show data window for Example 2.

Cardiovascular Fitness

Pretest Warmup Stage 1 Stage 2 Recovery 1 Recovery 2 Recovery 3

End Stage Heart Rate: 120 beats/min

Stage 1 Systolic BP: 140 mmHg

Stage 1 Diastolic BP: 80 mmHg

Stage 1 Time Completed: 3.00 minutes

OK

- The ISIS system will save both values but will store 120 to the "show data" field and use it as the heart rate for calculation of VO_2 max. This is an example of a correct data entry for manual heart rate.

Exhibit 4-42. Example 3 - manually entered heart rate - incorrect data entry

Cardiovascular Fitness: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm

SP ID: 819357 Name: JOHNSON, SUE Age: 39 years Gender: Female Date: 07/12/2000 Time: 02:20 PM

Treadmill Test

00:38 min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

10 11

Captured HR

110 130

Manual HR

128

PMHR: 181.00 beats/min 90.5 181

%HR max 50% 100%

Vo2 Max: 00

End of Section Close Exam Finish

- Example 3: In the above example for Stage 2, the ISIS system captured the heart rate as 130. The technician noticed that the heart rate was actually 128 at the time of data capture. This is within 5 beats of the captured value, but the technician manually entered 128 in the heart rate field. This is an example of an incorrect data entry for this field. The heart rate did not have to be entered manually because the captured heart rate was within 5 beats of the heart rate noted on the blood pressure monitor.

Exhibit 4-43. Show data window for Example 3

Cardiovascular Fitness

Warmup Stage 1 Stage 2 Recovery 1 Recovery 2 Recovery 3

End Stage Heart Rate: 130 beats/min

Stage 2 Systolic BP: 160 mmHg

Stage 2 Diastolic BP: 80 mmHg

Stage 2 Time Completed: 3:00 minutes

OK

- In the example above, the manually entered HR (128 for Stage 2) is not necessary and should not have been entered. In this situation, the ISIS system will save both values but will store 130 to the 'show data' field, and will use this for calculation of the VO_{2max} .

Exhibit 4-44. Manual entry fields disabled

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 06:41 PM

Treadmill Test

00:09
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9 13 15

Captured HR: 110 130

Manual HR: 120 128

PMHR: 184.00 beats/min 92 184

Vo2 Max: 60.00

%HR max 50% 100%

5 of 9

End of Section Close Exam Finish

- The fields for manual entry of heart rate are disabled after Recovery 2 is completed. At this point, the VO₂ max is calculated and the ISIS system will not allow a change to be made to the heart rate.

4.4.7 Stage 2

- After 3 minutes in Stage 1 the ISIS system will advance to Stage 2. The treadmill controller will automatically increase the speed and grade of the treadmill. The timer will be reset to zero and will start timing the 3-minute Stage 2.

4.4.7.1 First Two Minutes in Stage 2 – For SPs 19 Through 49 Years Old

Exhibit 4-45. Treadmill test stage 2 – (19 through 49 years)

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

00:30 min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
3	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9 12

Captured HR: 100

Manual HR:

PMHR: 184.00 beats/min

%HR max: 156.4 85% 184 100%

5 of 9

End of Section Close Exam Finish

- Throughout this stage, monitor the BP monitor for the SP's heart rate. Watch the continuum bar to confirm that the SP's heart rate is below 85 percent of PMHR.
- At 2 minutes 30 seconds into Stage 2 the ISIS system will give the auditory message:
"Get RPE Assessment."
- Get the SP's assessment of RPE and enter the number in the appropriate field below Stage 2.
- The BP monitor will automatically measure BP and HR. These values will be captured by the ISIS system. At approximately 2 minutes 40 seconds into Stage 2 the ISIS system will give the auditory message:
"Check heart rate on monitor."
- The ISIS system will then display a visual message:
"Compare heart rate on blood pressure monitor with captured heart rate on ISIS screen. Enter a manual heart rate ONLY if they differ by more than 5 beats."
- Note the heart rate on the blood pressure monitor and compare this number with the captured heart rate on the ISIS screen. Follow the same rules in Section 4.4.6.1 for manual entry of heart rate.

4.4.7.2 First Two Minutes in Stage 2 – For SPs 12 Through 18 Years Old

Exhibit 4-46. Treadmill test stage 2 – (12 through 18 years)

Cardiovascular Fitness: Stand:107 Session:107261 07/31/1999 01:30 pm - 05:30 pm

SP ID: 567817 Name: KLEIN, CYNTHIA Age: 12 years Gender: Female Date: 05/22/2000 Time: 07:11 PM

Treadmill Test

00:17 min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
4	2.70	2.00	3.10	4.00	3.10	8.50	2.00	0
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0

RPE

9 13

Captured HR

130

Manual HR

PMHR: 208.00 beats/min

%HR max

187.2 208

90% 100%

5 of 9

End of Section Close Exam Finish

- Throughout this stage, monitor the BP monitor for the SP's heart rate. Watch the continuum bar to confirm that the SP's heart rate is below 90 percent of PMHR.
- At 2 minutes 30 seconds into Stage 2 the ISIS system will give the auditory message:
"Get RPE Assessment."
- Get the SP's assessment of RPE and enter the number in the appropriate field below Stage 2.
- The BP monitor will automatically measure BP and HR. These values will be captured by the ISIS system. At approximately 2 minutes 40 seconds into Stage 2 the ISIS system will give the auditory message:
"Check heart rate on monitor."
- The ISIS system will then display this visual message:
"Compare heart rate on blood pressure monitor with captured heart rate on ISIS screen. Enter a manual heart rate ONLY if they differ by more than 5 beats."

Note the heart rate on the blood pressure monitor and compare this number with the captured heart rate on the ISIS screen. Follow the same rules in Section 4.4.6.1 for manual entry of heart rate.

4.4.7.3 The Last Minute of Stage 2

Exhibit 4-47. Treadmill test – last minute of stage 2

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

02:02
min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.00	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9 12

Captured HR: 100

Manual HR: 1

PMHR: 184.00 beats/min

%HR max: 85% 100%

5 of 9

End of Section Close Exam Finish

- In the last minute of stage 2, the red area of the continuum bar changes from red to yellow. This is the 85 – 100 percent mark for 19-through 49-year-olds and the 90-100 percent for the 12-through 18-year-olds.
- If the SP's HR is in the yellow zone, ask the SP if s/he feels OK. If the SP feels and looks OK, allow the test to continue. Monitor the SP very closely for the remaining seconds and stop the test if you think s/he is exhausted and cannot continue.
- After 3 minutes in Stage 2 the ISIS system will advance to Recovery 1. The treadmill controller will automatically decrease the speed and grade of the treadmill.

4.4.8 Recovery Stages

Exhibit 4-48. Treadmill test, recovery – routine test

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

00:30 min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: 0.00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9 12

Captured HR: 100

Manual HR:

PMHR: 184.00 beats/min 156.4 184

%HR max 85% 100%

5 of 9

End of Section Close Exam Finish

- The ISIS system will advance to Recovery 1. The timer will be reset to zero and will start timing the 1-minute Recovery Stage. The BP monitor will automatically measure BP and HR. These values will be captured by the ISIS system.
- After 1 minute in Recovery 1, the ISIS system will advance to Recovery 2 for an additional minute. If the SP's HR is above 110 at the end of Recovery 2, the ISIS system will give an auditory message:

"Proceed with Recovery 3."

- Recovery will continue for 1 minute, then stop as described below. If the SP's heart rate is 110 or below, the ISIS system will give an auditory message:

"Press 'Cool Down' on treadmill controller."

- Tell SP to grasp the front of the handrail and warn him/her that the treadmill is about to stop. The ISIS system will give the following auditory message:

"Press 'Stop' on treadmill controller. Congratulations, you have completed the test."

- The timer will stop and return to zero.

4.4.9 Checking Estimated VO₂ Max

Exhibit 4-49. Checking estimated VO₂ max


Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9, 12, 14
Captured HR: 100, 130
Manual HR: [] []

PMHR: 184.00 beats/min
%HR max: 50% (92 bpm) to 100% (184 bpm)

- The estimated VO₂ max is calculated after Recovery 2 and displayed in the VO₂Max field.
- The SP will be given the results of this exam in the Report of Findings. This estimate of VO₂ max will be used to determine how the SP's fitness level compares to other people of similar age and gender.


Exhibit 4-50. Checking estimated VO₂ max

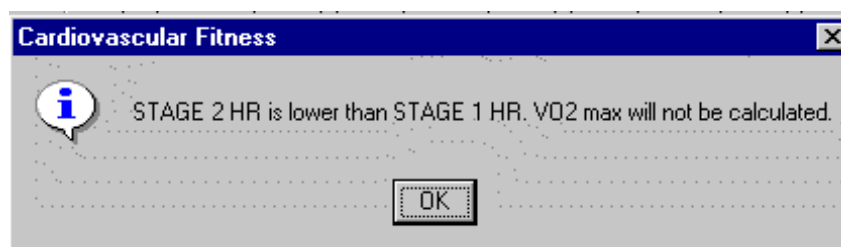

- If the heart rate at Stage 2 is lower than the heart rate at Stage 1, the ISIS system will display the above message and the VO₂ will not be calculated.

4.4.10 Component Status for Routine Test

Exhibit 4-51. Component status - complete

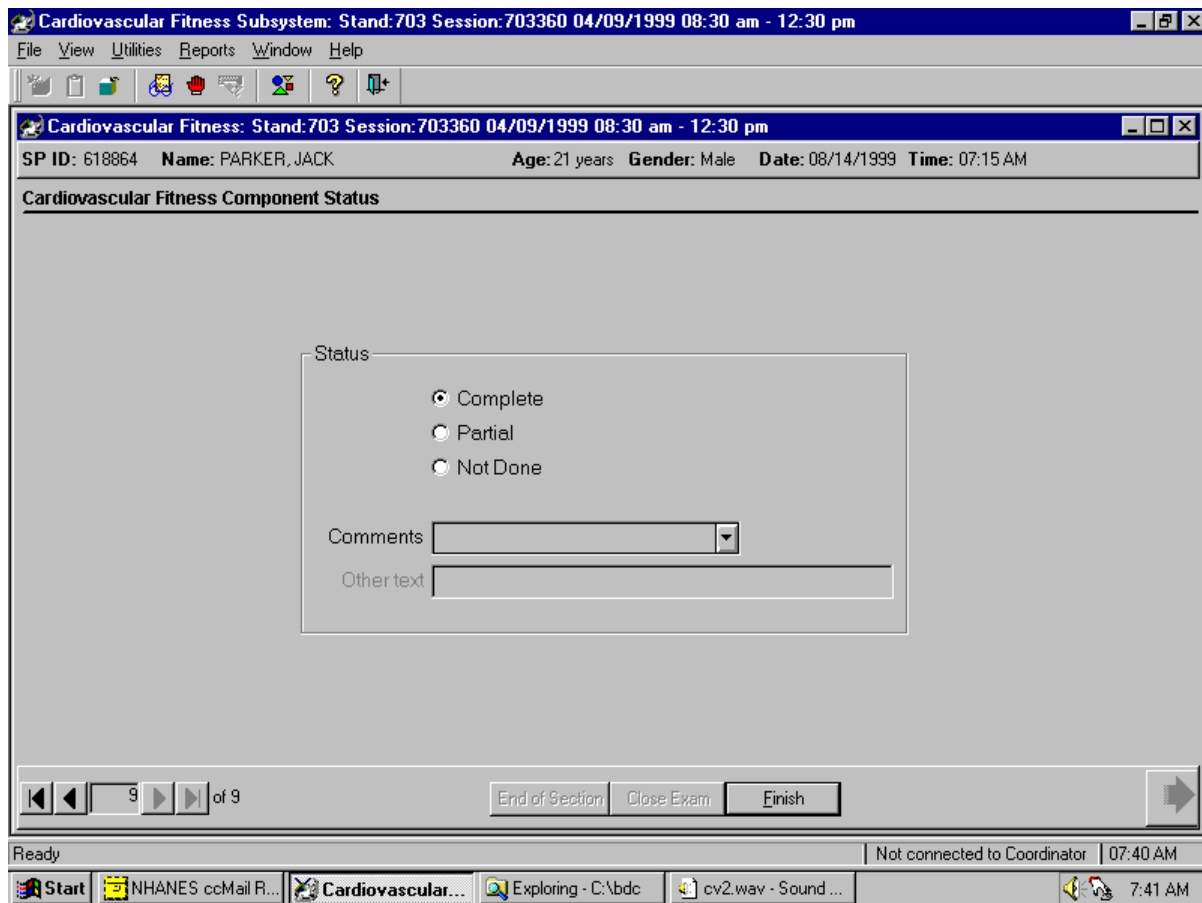


- At the end of the test, the ISIS system will advance to the Component Status screen.
- The Component Status will be set to 'Complete' if all the data are successfully obtained and captured in the ISIS system.
- Press the 'Finish' button to close the exam.

4.5 Changing the Exercise Protocol

4.5.1 Changing the Treadmill Controller

Exhibit 4-52. Treadmill controller

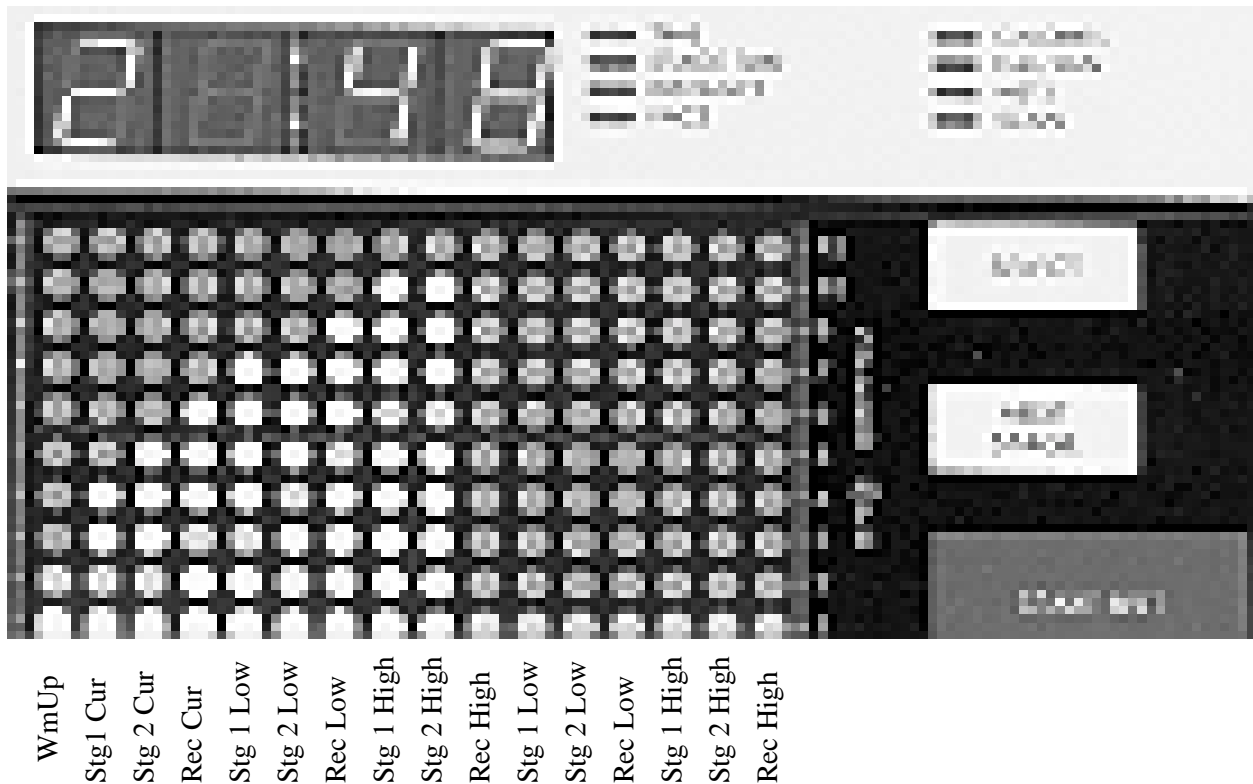

- There are 16 horizontal lights on the treadmill controller display (See Exhibit 4-52) to represent the predetermined protocols, including the currently selected protocols and the ones in the next lower and next higher levels.
- ‘WmUp’ is the Warmup Stage of the current protocol. ‘Stg 1 Cur’ is Stage 1 of the current protocol, ‘Stg 2 Cur’ is Stage 2 of the current protocol and ‘Rec Cur’ is the Recovery Stage of the current protocol. If the protocol does not have to be changed during the test, the controller will advance through these four stages only.
- ‘Stg 1 Low’ is Stage 1 for the next lower protocol, ‘Stg 2 Low’ is Stage 2 for the next lower protocol, and ‘Rec Low’ is the Recovery Stage for the next lower protocol. If the protocol has to be decreased, when the stage changes from Warmup to Stage 1, Press ‘Next Stage’ **slowly** until you get to ‘Stg 1 Low.’
- ‘Stg 1 High’ is Stage 1 for the next higher protocol, ‘Stg 2 High’ is Stage 2 for the next higher protocol, and ‘Rec High’ is the Recovery Stage for the next higher protocol. If the protocol has to be increased, when the stage changes from Warmup to Stage 1, Press ‘Next Stage’ **slowly** until you get to ‘Stg 1 High.’

4.5.2 Decreasing the Exercise Protocol

Exhibit 4-53. Decreasing the exercise protocol

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 06:41 PM

Treadmill Test

00:39
min sec

START
RESTART
EMERGENCY STOP
Change Protocol
Go To Recovery
Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
3	3.70	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE 9

Change Protocol: 5
Apply Cancel

%HR max 50% 60% 100%

5 of 9
End of Section Close Exam Finish

- If the heart rate is greater than 60 percent of PMHR during Warm-up, the protocol should be decreased for the remainder of the exam. Watch the heart during Warm-up. You need to make a decision to change the protocol by 1 minute 30 seconds into Warm-up to allow time to initiate the changes
- Click the 'Change Protocol' button. A box with change protocol options will be displayed. Click the down arrow to decrease the protocol. Click 'Apply.' The ISIS system will give an auditory message:

"When this stage changes, change treadmill controller to the next lower protocol."

Exhibit 4-54. Decreased protocol

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 06:41 PM

Treadmill Test

00:01
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9

Captured HR

Manual HR

PMHR: 184.00 beats/min

%HR max: 60% 80% 100%

110.4 147.2 184


5 of 9

End of Section Close Exam Finish

- When the protocol is changed, the ISIS system will advance to Stage 1 and will automatically select the next lower protocol.
- The protocol will be highlighted in red below Stage 1. The same protocol will be used for Stage 2.
- In contrast, the treadmill controller does not move to the next lower protocol automatically, even though ISIS was changed. It is necessary to manually change the setting on the treadmill controller to advance to Stage 1 on the next lower protocol.
- When the timer in the ISIS screen is reset to zero and starts timing the 3-minute Stage 1, press 'Next Stage' on the treadmill controller **slowly** until you get to 'Stg 1 Low' for the next lower protocol. It is important to wait until the stage changes before changing the controller.
- If you press too quickly and go past the desired protocol, keep pressing 'Next Stage' to advance to the second set of protocols. The next lower and next higher protocols are programmed twice in sequence. Look for Stage 1 of the next lower protocol and press 'Next Stage' **slowly** until you get to the desired protocol.
- When the protocol changes to Stage 2, the treadmill controller will now go automatically to Stage 2 of the current protocol (the lower protocol).


4.5.3 Increasing the Exercise Protocol

Exhibit 4-55. Increasing the exercise protocol

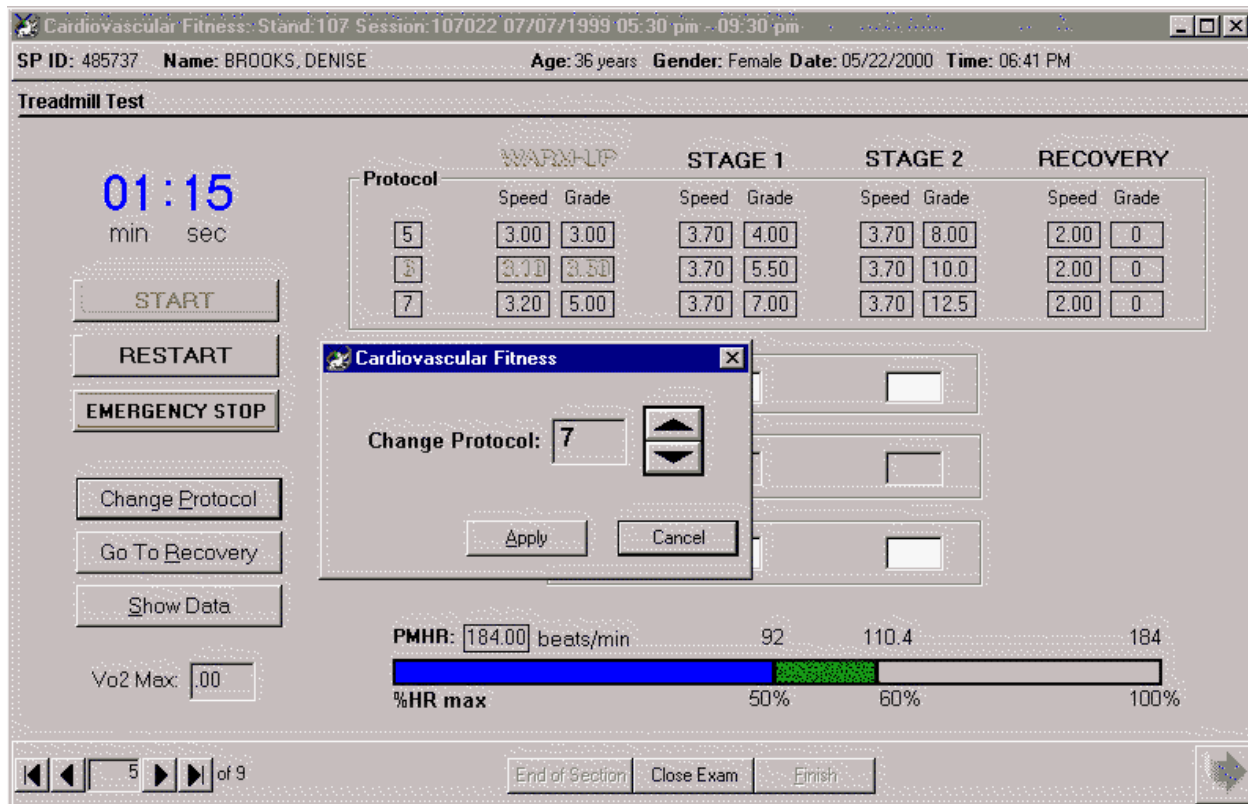

- The method for increasing the protocol is similar to decreasing the protocol. If the heart rate is less than 50 percent during Warm-up, you should increase the protocol for the remainder of the exam.
- Click 'Change Protocol' button. Click the up arrow to increase the protocol. Click Apply.
- The ISIS system will give an auditory message:
"At the beginning of the next stage, change treadmill controller to the next higher protocol."
- When the Warm-up is finished, the ISIS system will advance to Stage 1 and will automatically select the next higher protocol. This protocol will be highlighted in red below Stage 1. The same protocol will be used for Stage 2.
- In contrast, the treadmill controller does not move to the next higher protocol automatically even though ISIS was changed. Press 'Next Stage' on the treadmill controller to advance to Stage 1 on the next higher protocol.
- When the protocol changes to Stage 2, the treadmill controller will now go automatically to Stage 2 of the current protocol (the higher protocol).

Exhibit 4-56. Changing the protocol in stage 1 or 2 - message

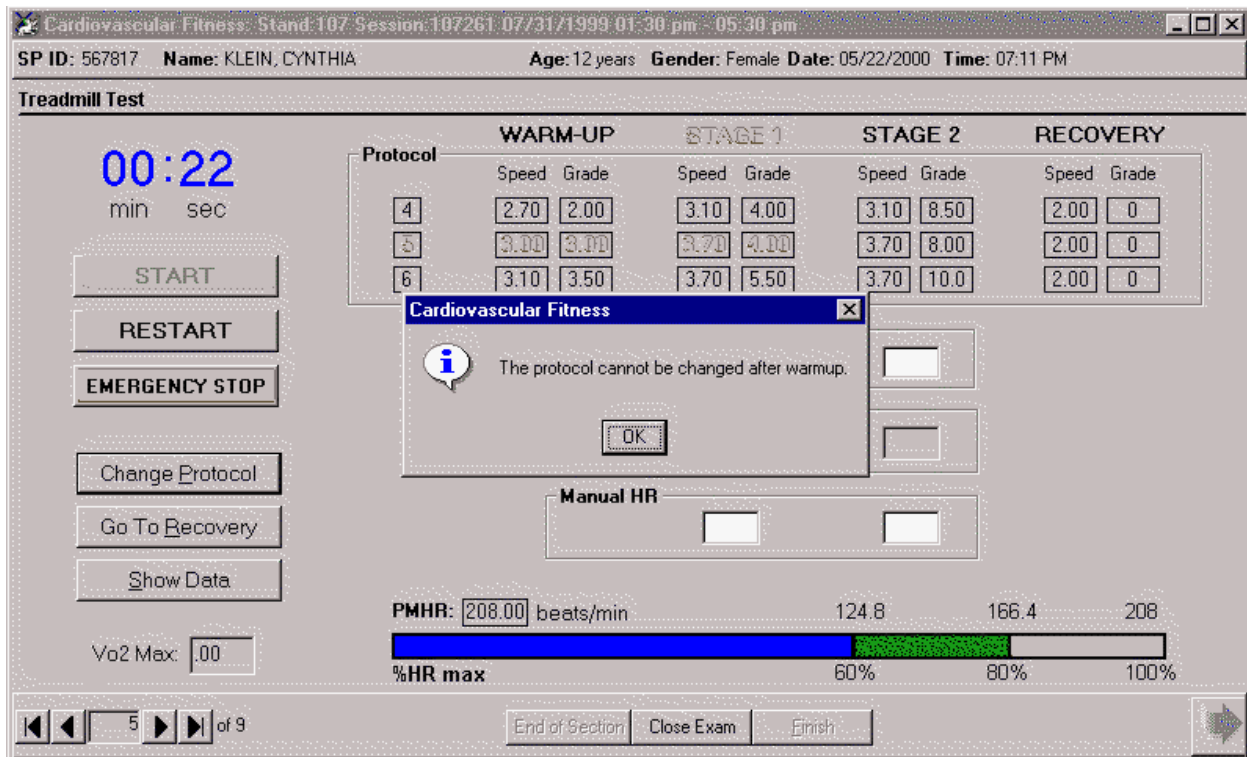

- The exercise protocol may only be changed during Warm-up.
- If you click on the 'Change Protocol' button during Stage 1 or Stage 2, the ISIS system will display a message:
 - *'The protocol cannot be changed after warm-up.'*
- Click OK.
- If there is an indication to change the protocol during Stage 1 or Stage 2, use 'Go to Recovery' or 'Emergency Stop.' (See Sections 4-6 to 4-8 for detailed procedures).

4.5.4 Lowest and Highest Protocol

Exhibit 4-57. Protocol 0

Cardiovascular Fitness Stand 107 Session 107280 06/02/1999 08:30 am - 12:30 pm

SP ID: 988340 Name: GARCIA, ANA MARIA Age: 44 years Gender: Female Date: 05/22/2000 Time: 07:06 PM

Treadmill Test

00:00 min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
1	1.70	0	2.10	.50	2.10	4.50	2.00	0
2	2.00	1.00	2.30	2.00	2.30	6.50	2.00	0

RPE

Captured HR

Manual HR

PMHR: 176.00 beats/min

%HR max


88 114.4 176

50% 60% 100%

End of Section Close Exam Finish

- If Protocol 1 is the predetermined protocol, the next lower protocol (Protocol zero) is grayed out on the ISIS screen.

Exhibit 4-58. Protocol 0 message

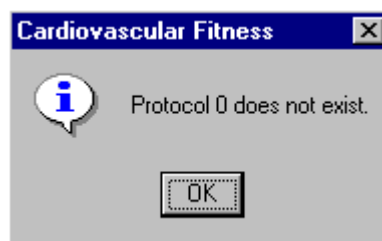

- If the predetermined protocol is 1, and the request for changing protocol is made by pressing the "Change Protocol" button, the ISIS system will display a message: 'Protocol 0 does not exist.' This is a reminder that the protocol can only be increased to Protocol 2.
- If the predetermined protocol is 1 (the lowest protocol in the program) and the heart rate is greater than 60% of PMHR during Warm-up, the protocol will keep unchanged due to the unavailability of the next lower protocol.

Exhibit 4-59. Protocol 9

Cardiovascular Fitness: Stand:107 Session:107020 07/07/1999 08:30 am - 12:30 pm

SP ID: 998919 Name: ROCHE, BRIAN Age: 30 years Gender: Male Date: 05/22/2000 Time: 07:08 PM

Treadmill Test

00:00
min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0
8	3.60	5.00	3.70	8.50	3.70	14.5	2.00	0
9								

RPE

Captured HR

Manual HR

PMHR: 190.00 beats/min

%HR max

95 123.5 190


50% 60% 100%

5 of 9

End of Section Close Exam Finish

- If Protocol 9 is the predetermined protocol, the next higher protocol (Protocol 9) is grayed out.

Exhibit 4-60. Protocol 9 message

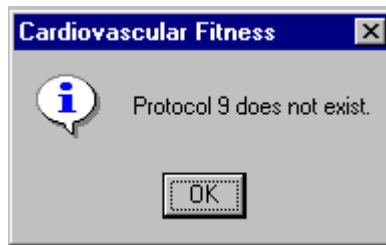

- If the predetermined protocol is 8, and the request for changing protocol is made by pressing the "Change Protocol" button the ISIS system will display a message: 'Protocol 9 does not exist.' This is a reminder that the protocol can only be decreased to Protocol 7.
- If the predetermined protocol is 8 (the highest protocol in the program) and the heart rate is lower than 50% of PMHR during Warm-up, the protocol will remain unchanged due to the unavailability of the next higher protocol.

4.6 Emergency Stop

Exhibit 4-61. Emergency stop

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

02:52
min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

9 13

Captured HR

170

Manual HR

1 1

PMHR: 184.00 beats/min

%HR max

110.4 147.2 184

60% 80% 100%

5 of 9

End of Section Close Exam Finish

- During the treadmill test you should monitor the SP constantly for any signs of discomfort or distress. If any of the conditions listed in the Priority 1 Stopping Criteria are present (or if for any other reason you think the test should be stopped immediately), follow the procedures for 'Emergency Stop.'
- Press 'Emergency Stop' on the ISIS screen. The ISIS system will go to Priority 1 Stopping Criteria and will give this auditory message:
"Press Cool Down on Treadmill Controller. Press Stop Belt."
- Press 'Cool Down' first. This will lower the grade on the treadmill. Then press Stop Belt. Assist the SP off the treadmill.
- Depending on the situation, you may call for assistance immediately when you notice signs of distress. If the situation is not critical but you want to terminate the test, follow the procedures and get the SP to the chair or floor and then get help. Try not to alarm the SP unnecessarily and maintain a calm, controlled atmosphere.
- Call the MEC physician and have the SP evaluated.
- When the situation is under control and the SP is stabilized, enter the appropriate information on the ISIS screens.

4.7 Priority 1 Stopping Criteria

Exhibit 4-62. Priority 1 stopping criteria - other

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:57 AM

Priority 1 Stopping Criteria

You have stopped the exam early. You must select one or more of the following stopping criteria:

- ☐ Pain or pressure in the chest
- ☐ Pain in neck, jaw, or teeth
- ☐ Pain in the shoulder or between the shoulder blades
- ☐ Pain radiating down one or both arms
- ☐ Moderate/severe leg pain or cramping
- ☐ Lightheadedness, dizziness
- ☐ Severe shortness of breath/wheezing
- ☐ Change in skin color to bluish/grayish
- ☐ Severe headache
- ☐ Subject appears confused or agitated
- ☐ Nausea
- ☐ Unstable gait
- ☐ Visual problems
- ☒ Other Priority 1

Navigation: 6 of 8 | End of Section | Close Exam | Finish

- If the 'Emergency Stop' button is pressed, the ISIS system will go immediately to Priority 1 Stopping Criteria. When the SP is stabilized and the situation is under control, check the appropriate reason for terminating the test for the criteria on the screen.

- The Priority 1 Stopping Criteria are:
 - Pain or pressure in the chest
 - Pain in the neck, jaw, or teeth
 - Pain in the shoulder or between the shoulder blades
 - Pain radiating down one or both arms
 - Moderate/severe leg pain or cramping
 - Lightheadedness, dizziness
 - Severe shortness of breath/wheezing
 - Change in skin color to bluish/grayish
 - Severe headache
 - Subject appears confused or agitated
 - Nausea
 - Unstable gait
 - Visual problems
 - Other
- If you notice any or the above criteria or if the SP complains of any of the above, follow the procedures for 'Emergency Stop' and terminate the test immediately.

Exhibit 4-63. Priority 1 stopping criteria - required entry

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:57 AM

Priority 1 Stopping Criteria

You have stopped the exam early. You must select one or more of the following stopping criteria:

- ☐ Pain or pressure in the chest
- ☐ Pain in neck, jaw, or teeth
- ☐ Pain in the shoulder or between the shoulder blades
- ☐ Pain radiating down one or both arms
- ☐ Moderate/severe leg pain or cramping
- ☐ Lightheadedness, dizziness
- ☐ Severe shortness of breath/wheezing
- ☐ Change in skin color to bluish/grayish
- ☐ Severe headache
- ☐ Subject appears confused or agitated
- ☐ Nausea
- ☐ Unstable gait
- ☐ Visual problems
- ☐ Other Priority 1

Cardiovascular Fitness

You must select a stopping criteria.

OK

6 of 8

End of Section Close Exam Finish

- If the 'Emergency Stop' button is pressed, the Priority 1 Stopping Criteria screen appears. You are required to select one or more criteria.
- If you do not select a criterion and the Next button is pressed, the ISIS system will display a message: 'You must select a stopping criteria.'
- Click OK to this message. Select on or more criteria as appropriate and press the Next button to advance to the Component Status screen.
- The Component status will be set to 'Partial' with the comment 'Priority 1 Stopping Criteria.'
- Press the Finish button to close the exam.

4.8 Early Termination of Test – Go to Recovery

Exhibit 4-64. Go to recovery

Cardiovascular Fitness: Stand:107 Session:107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

01:02
min sec

START
RESTART
EMERGENCY STOP
Change Protocol
Go To Recovery
Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE: 9

Captured HR

Manual HR

PMHR: 184.00 beats/min

%HR max: 60% 80% 100%

End of Section Close Exam Finish

- If any of the situations listed in 3.13.2 Priority 2 Stopping Criteria occur, the test should be terminated and moved to the Recovery stages. Some of the criteria are data related and the ISIS system will display a message that you should go to the Recovery Stage (See exhibits in 4.8.1 for examples). Other criteria will be noted by observation. You should be monitoring the SP closely for any signs of these criteria.
- If a Priority 2 Stop is warranted, Press 'Go to Recovery' on the ISIS screen.
- The ISIS system will then give an auditory message:
"Change treadmill controller to 'Recovery Stage' for current protocol."
- Change the setting on the treadmill controller by pressing 'Next Stage' until the Recovery Stage is reached. Recovery stages 1 and 2 (and 3 if necessary) will be completed as usual. When Recovery is complete, press the 'Next' button to advance to the next screen.

4.8.1 Data Not Within Acceptable Limits

Exhibit 4-65. Data not within acceptable limits



One of the reasons for early termination is some of the data items exceed set limits.

- There are edit limits set on the following measurements:
 - Heart rate greater than 85 percent PMHR (>90 percent PMHR for 12 through 18 years)
 - SBP greater than 260 mmHg
 - DBP greater than 115 mmHg
 - RPE greater than 17
- If these limits are exceeded, the ISIS system will display a message:
"(Data item) is not within acceptable limits. Go to Recovery."
- Click OK.
- The ISIS system will then give an auditory message:
"Press 'Go to Recovery' on ISIS screen."
- Press 'Go to Recovery' on the ISIS screen. The ISIS system will then give an auditory message:
"Change treadmill controller to Recovery Stage for current protocol."
- Change the setting on the treadmill controller by pressing Next Stage until the Recovery Stage is reached.
- Recovery Stages 1 and 2 (and 3 if necessary) will be completed as usual. When Recovery is complete, press the Next button to advance to the next screen.

Exhibit 4-66. SBP not within acceptable limits

Cardiovascular Fitness: Stand: 107 Session: 107022 07/07/1999 05:30 pm - 09:30 pm

SP ID: 485737 Name: BROOKS, DENISE Age: 36 years Gender: Female Date: 05/22/2000 Time: 05:32 PM

Treadmill Test

02:57
min sec

START
RESTART
EMERGENCY STOP

Change Protocol
Go To Recovery
Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.00	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

Cardiovascular Fitness

SBP is not within acceptable limits. Go to RECOVERY.

OK

Manual HR

PMHR: 184.00 beats/min

%HR max

110.4 147.2 184


60% 80% 100%

5 of 9

End of Section Close Exam Finish

- If the systolic blood pressure is greater than 260 mmHG at any time during Warm-up, Stage 1 or Stage 2, the ISIS system will display a message: 'SBP is not within acceptable limits. Go to Recovery.'
- Click OK. Follow the steps outlined in Section 3.13.2.

Exhibit 4-67. DBP not within acceptable limits

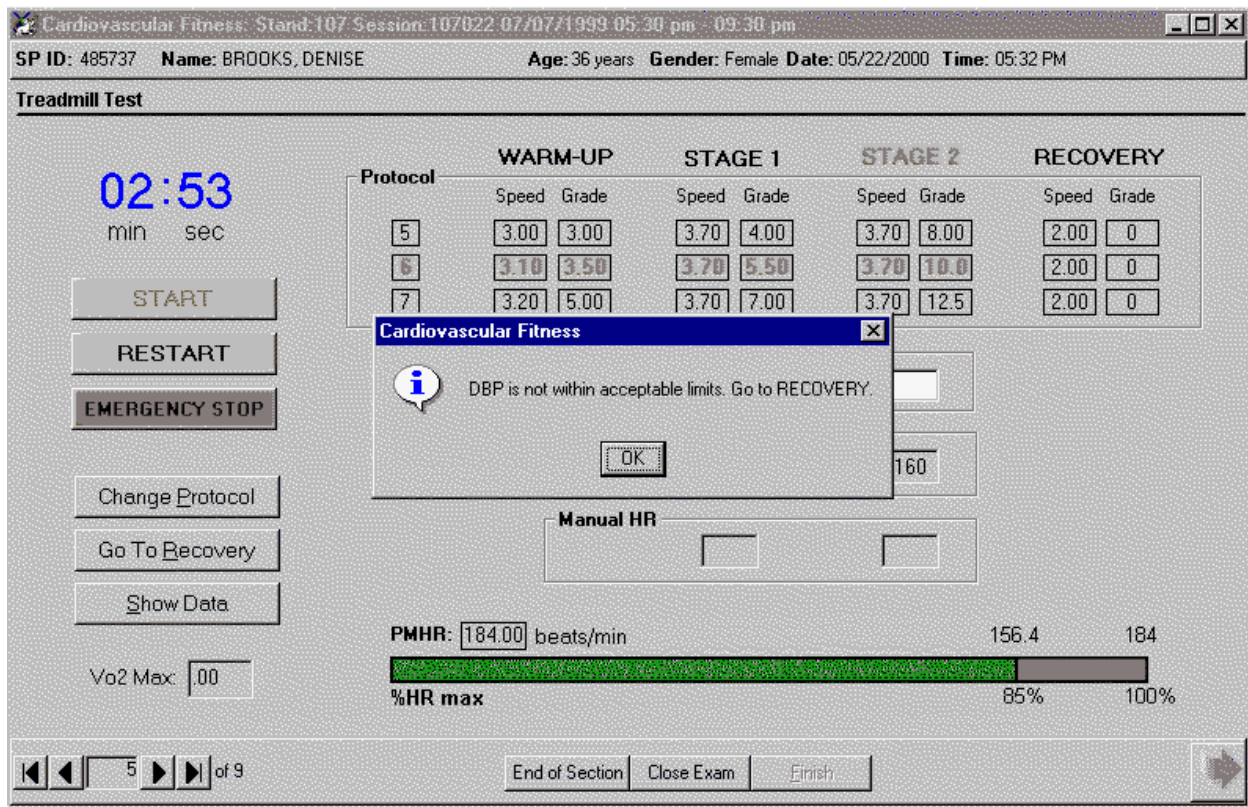

- If the diastolic blood pressure is greater than 115 mmHg at any time during Warm-up, Stage 1 or Stage 2, the ISIS system will display a message: 'DBP is not within acceptable limits. Go to Recovery.'
- Click OK. Follow the steps outlined in Section 3.13.2.

Exhibit 4-68. Rating of perceived exertion not within acceptable limits

Cardiovascular Fitness: Stand:107 Session:107020 07/07/1999 08:30 am - 12:30 pm

SP ID: 884264 Name: DUCK, DAFFY Age: 35 years Gender: Female Date: 05/31/2000 Time: 03:34 PM

Treadmill Test

01:09
min sec

START

RES Cardiovascular Fitness

EMERGE

Change

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
4	2.70	2.00	3.10	4.00	3.10	8.50	2.00	0
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0

RPE is above acceptable criteria, go to recovery.

OK

18

Manual run

PMHR: 185.00 beats/min


%HR max

50% 60% 100%

End of Section Close Exam Finish

- If the Rating of Perceived Exertion (RPE) is greater than 17 at any time during Warm-up, Stage 1 or Stage 2, the ISIS system will display a message: 'RPE is not within acceptable limits. Go to Recovery.'
- Click OK. Follow the steps outlined in Section 3.13.2 if this happened during Warm-up or Stage 1. Follow the steps outlined in Section 3.10.1 if this happened during Stage 2.

Exhibit 4-69. Heart rate greater than 85 percent PMHR

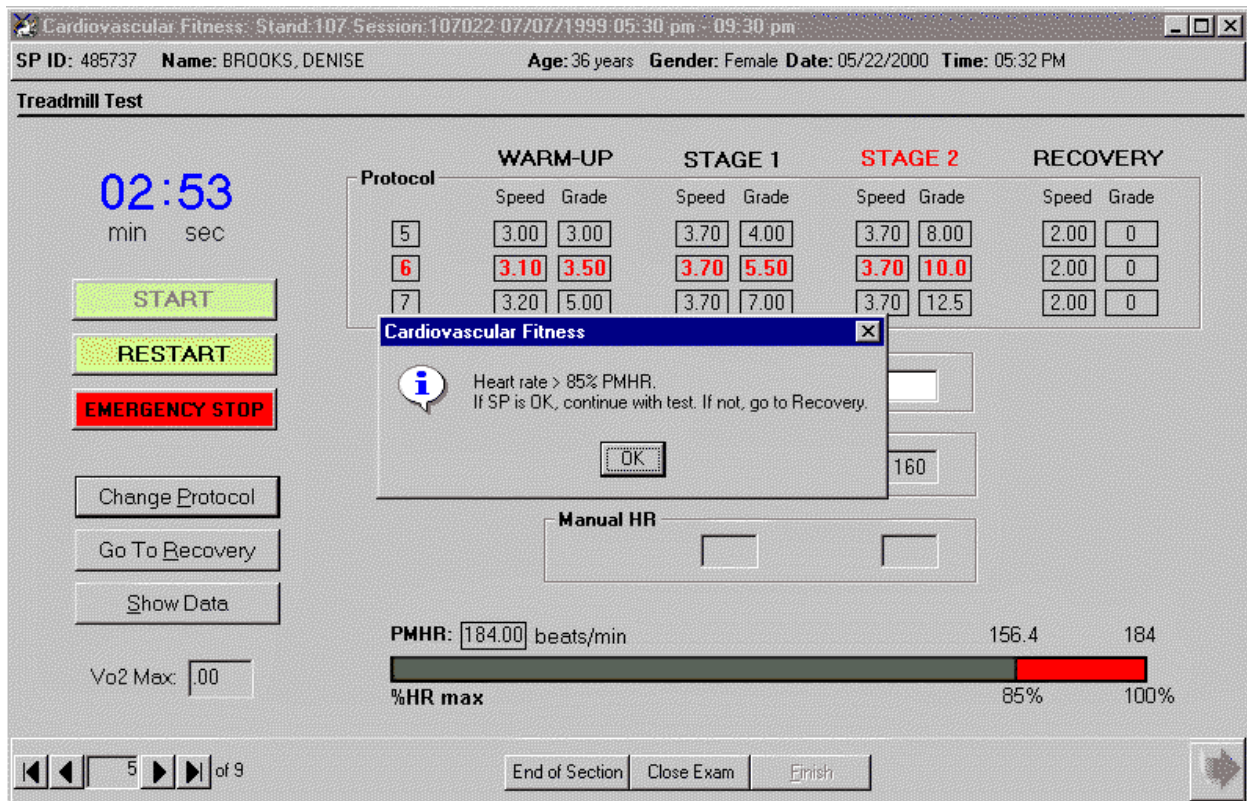

- In Stage 2, when the heart rate is captured by the ISIS system, if the heart rate is greater than 85 percent (or 90 percent for 12- through 18-year-olds), the following message will be displayed.

Heart Rate > 85 percent PMHR. If SP is OK, continue with test, if not, go to Recovery.

- In the last minute of Stage 2, the SP can continue with the exam if the HR is greater than 85 percent PMHR as long as the SP feels OK and looks OK. The SP should be monitored closely in these last few seconds and the exam should be stopped if s/he looks exhausted or does not want to continue.

4.9 Priority 2 Stopping Criteria

Exhibit 4-70. Priority 2 stopping criteria

The screenshot shows a software window titled "Cardiovascular Fitness Subsystem: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm". The menu bar includes File, View, Utilities, Reports, Window, and Help. Below the menu is a toolbar with various icons. The main display area has a header with patient information: "SP ID: 427581 Name: DONOHUE, MAURIE Age: 23 years Gender: Female Date: 07/11/2000 Time: 05:43 PM". The title of the current screen is "Priority 2 Stopping Criteria". The main text reads: "You have stopped the exam early. You must select one or more of the following stopping criteria:". Below this is a list of 13 criteria, each preceded by an unchecked checkbox. The criteria are: 1. SP HR during WARM-UP, STAGE 1, or first 2 minutes of STAGE 2 exceeds 90% PMHR for 12 through 18 yr olds or 85% PMHR for 19 through 49 yr olds; 2. SP exceeds 80% of PMHR during STAGE 1; 3. Exercise blood pressure > 260 mmHg systolic and/or 115 mmHg diastolic; 4. Failure of systolic blood pressure to rise with increasing exercise intensity; 5. Significant drop (>20 mmHg) in systolic blood pressure during exercise; 6. Failure of heart rate to increase with increasing exercise intensity; 7. Sudden, noticeable variability in heart rate (sudden change of + 30 beats/minute during a STAGE); 8. SP requests to stop test, or reports severe fatigue; 9. Failure of testing equipment; 10. SP is unable to stop gripping hand rails (lightly touching hand rails for balance purposes is allowed, but only if absolutely necessary); 11. SP reports RPE > 17; 12. Technician discretion (tech notices signs of severe exertion); 13. Other. Below the "Other" checkbox is a text field labeled "Enter Comment:". At the bottom of the screen is a navigation bar with buttons for "End of Section", "Close Exam", and "Finish", along with a "Next" button represented by a blue arrow. The status bar at the very bottom shows "Ready", "Not connected to Coordinator", and "06:01 PM".

Cardiovascular Fitness Subsystem: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

SP ID: 427581 Name: DONOHUE, MAURIE Age: 23 years Gender: Female Date: 07/11/2000 Time: 05:43 PM

Priority 2 Stopping Criteria

You have stopped the exam early. You must select one or more of the following stopping criteria:

- ☐ SP HR during WARM-UP, STAGE 1, or first 2 minutes of STAGE 2 exceeds 90% PMHR for 12 through 18 yr olds or 85% PMHR for 19 through 49 yr olds
- ☐ SP exceeds 80% of PMHR during STAGE 1.
- ☐ Exercise blood pressure > 260 mmHg systolic and/or 115 mmHg diastolic
- ☐ Failure of systolic blood pressure to rise with increasing exercise intensity
- ☐ Significant drop (>20 mmHg) in systolic blood pressure during exercise
- ☐ Failure of heart rate to increase with increasing exercise intensity
- ☐ Sudden, noticeable variability in heart rate (sudden change of + 30 beats/minute during a STAGE)
- ☐ SP requests to stop test, or reports severe fatigue
- ☐ Failure of testing equipment
- ☐ SP is unable to stop gripping hand rails (lightly touching hand rails for balance purposes is allowed, but only if absolutely necessary)
- ☐ SP reports RPE > 17
- ☐ Technician discretion (tech notices signs of severe exertion).
- ☐ Other Enter Comment:

7 of 9 End of Section Close Exam Finish

Ready Not connected to Coordinator 06:01 PM

- If the treadmill test is terminated early due to one or more of the Priority 2 Stopping Criteria, the SP will go to the Recovery Stage of the current protocol.
- After completing the normal cool-down, the ISIS system will advance to the Priority 2 Stopping Criteria screen.
- Select one or more of the appropriate stopping criteria.
- You are required to select one or more of these criteria if 'Go to Recovery' was selected on the previous screen.
- If you do not select a criterion, the ISIS system will display a message when the Next button is selected: 'You must select a stopping criterion.'
- Click OK. Select one or more of the criteria and press the Next button to advance to the next screen.

Exhibit 4-71. Component status for priority 2

Cardiovascular Fitness Subsystem: Stand:806 Session:808571 12/11/1998 05:30 pm - 09:30 pm

SP ID: 553247 Name: Adams, Henry Age: 44 years Gender: Male Date: 01/24/1999 Time: 05:58 PM

Cardiovascular Fitness Component Status

Status

☐ Complete

☒ Partial

☐ Not Done

Comments Priority 2 stopping criteria

Other text

8 of 8

End of Section Close Exam Finish

- The ISIS system will advance to the Component Status Screen.
- The Component Status will be set to 'Partial.'
- The comment will be set to 'Priority 2 Stopping Criteria'.
- Press the 'Finish' button to close the exam.
- The ISIS system will not calculate the VO2 max.

Exhibit 4-72. SP unable to stop gripping handrails

Cardiovascular Fitness Subsystem: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness: Stand:990 Session:990300 10/05/1999 08:30 am - 12:30 pm

SP ID: 427581 Name: DONOHUE, MAURIE Age: 23 years Gender: Female Date: 07/11/2000 Time: 05:43 PM

Priority 2 Stopping Criteria

You have stopped the exam early. You must select one or more of the following stopping criteria:

- ☐ SP HR during WARM-UP, STAGE 1, or first 2 minutes of STAGE 2 exceeds 90% PMHR for 12 through 18 yr olds or 85% PMHR for 19 through 49 yr olds
- ☐ SP exceeds 80% of PMHR during STAGE 1.
- ☐ Exercise blood pressure > 260 mmHg systolic and/or 115 mmHg diastolic
- ☐ Failure of systolic blood pressure to rise with increasing exercise intensity
- ☐ Significant drop (>20 mmHg) in systolic blood pressure during exercise
- ☐ Failure of heart rate to increase with increasing exercise intensity
- ☐ Sudden, noticeable variability in heart rate (sudden change of + 30 beats/minute during a STAGE)
- ☐ SP requests to stop test, or reports severe fatigue
- ☐ Failure of testing equipment
- ☒ SP is unable to stop gripping hand rails (lightly touching hand rails for balance purposes is allowed, but only if absolutely necessary)
- ☐ SP reports RPE > 17
- ☐ Technician discretion (tech notices signs of severe exertion).
- ☐ Other Enter Comment:

7 of 9

End of Section Close Exam Finish

Ready Not connected to Coordinator 06:02 PM

- SP's are required to walk on the treadmill without holding on to the handrail. Grasping the handrail throughout the test will tend to overestimate the SP's VO₂ max.
- Some SP's may hold the handrail for support at the beginning of the exam for support. If this occurs, repeat the instructions for walking on a treadmill and try to coach him/her to walk on the treadmill without holding. Remind the SP to stand up straight, look ahead not down at their feet and take long steps as the belt moves.
- Lightly touching the handrails for balance purposes is allowed but only if absolutely necessary.
- Give them lots of encouragement and try to keep the SP walking on the belt. If they continue to have problems and are holding on or briefly grabbing the handrails after you have repeated these instructions, use the method for stopping the test early and go to Recovery.
- Do not allow the SP to complete the exam if he/she is holding on to the handrails for a large part of the test as this will result in an overestimation of the VO₂ max.
- After the Recovery stage is completed, the ISIS system will advance to the Priority 2 Stopping Criteria.

- Select 'SP is unable to stop gripping handrails.' See Exhibit 4-72.

4.10 Problems With Data Acquisition During Treadmill Test

There are several potential sources for problems with acquisition of data during the treadmill test. After the test is completed a screen will be displayed to allow entry of problems that may occur during the test with equipment, integration, and technician error. You are able to check all that apply. It is extremely important to check off the correct problem.

Exhibit 4-73 - Problems with data acquisition - No problems

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Problems with data acquisition during treadmill test

No Problem ☒

Problem obtaining BP - BP not displayed on the BP monitor, not captured in ISIS:

☐ BP not obtained in Baseline ☐ BP not obtained in Warmup ☐ BP not obtained in Stage 1

☐ BP not obtained in Stage 2 ☐ BP not obtained in Recovery 1, 2, and /or 3

☐ Other

Problem with data capture - BP displayed correctly on STBP monitor but not correctly captured in ISIS

☐ BP on STBP OK, not captured in ISIS ☐ BP displayed normally on BP monitor, different data captured in ISIS

☐ Other

Problem obtaining HR

☐ Too much variability on HR on monitor to determine accurate HR

☐ Other

Problem with electrodes

☐ adherence problem due to excessive perspiration ☐ adherence problem due to excessive chest hair

☐ large chest/abdomen ☐ Other

Problem with BP cuff

☐ large arm caused difficulty with microphone placement ☐ cuff slipped off during test

☐ SP bending/flexing arm during test ☐ Other

Problem with treadmill controller

☐ Controller advanced more than one stage - Restarted test ☐ Other

8 of 9

End of Section Close Exam Finish

- **No Problem** - If there have been no problems with data acquisition during the test, check 'No Problem.' The remaining fields on the screen will be disabled. If you try to advance to the next screen without checking a problem or 'No problem,' the ISIS system will display a message to remind you to check this.

Exhibit 4-74. Problems with data acquisition - BP, data capture, and HR

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Problems with data acquisition during treadmill test

No Problem ☐

Problem obtaining BP - BP not displayed on the BP monitor, not captured in ISIS:

☐ BP not obtained in Baseline ☐ BP not obtained in Warmup ☐ BP not obtained in Stage 1

☐ BP not obtained in Stage 2 ☐ BP not obtained in Recovery 1, 2, and /or 3

☐ Other

Problem with data capture - BP displayed correctly on STBP monitor but not correctly captured in ISIS

☐ BP on STBP OK, not captured in ISIS ☐ BP displayed normally on BP monitor, different data captured in ISIS

☐ Other

Problem obtaining HR

☐ Too much variability on HR on monitor to determine accurate HR

☐ Other

Problem with electrodes

☐ adherence problem due to excessive perspiration ☐ adherence problem due to excessive chest hair

☐ large chest/abdomen ☐ Other

Problem with BP cuff

☐ large arm caused difficulty with microphone placement ☐ cuff slipped off during test

☐ SP bending/flexing arm during test ☐ Other

Problem with treadmill controller

☐ Controller advanced more than one stage - Restarted test ☐ Other

8 of 9

End of Section Close Exam Finish

- **Problem obtaining BP - BP not displayed on the BP monitor, not displayed in ISIS** - This problem can occur for several reasons. If the BP cuff is not positioned correctly on the arm, if the SP has a large arm and the signal to the microphone is not reaching the microphone, or if the SP is bending or flexing the arm at the time of data acquisition, select this option. When this problem occurs, the technician should check the position of the microphone. This problem can occur in one or all stages of the treadmill test. If this occurs at the pretest measurement, you should reposition the cuff and recapture the BP. You should always get a pretest BP before starting the test.
- **Problem with data capture - BP not displayed correctly on BP monitor** - This is not a problem with cuff position or the BP monitor. This is a problem with the transfer of data from the BP monitor to the ISIS system. This may be caused by interference with the data transfer from an outside source. If this happens, check the connections of the BP monitor to the computer and to the power source at every connection. The BP monitor should be on a different circuit than the treadmill. This could also be a problem with the ISIS system in receiving the data transfer.
- **Problems obtaining HR** - This can occur if there is something interfering with the electrodes picking up the signal. This can be caused by such things as inadequate preparation of the skin, poor placement of the electrodes, perspiration, chest hair interfering with the adherence of the electrode to the skin, an irregular heart rate, faulty electrodes, leads, or cables. If this occurs, troubleshoot to determine the cause.

Exhibit 4-75. Problems with data acquisition - electrodes, BP cuff, and treadmill controller

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Problems with data acquisition during treadmill test

No Problem ☐

Problem obtaining BP - BP not displayed on the BP monitor, not captured in ISIS:

☐ BP not obtained in Baseline ☐ BP not obtained in Warmup ☒ BP not obtained in Stage 1

☐ BP not obtained in Stage 2 ☒ BP not obtained in Recovery 1, 2, and /or 3

☐ Other

Problem with data capture - BP displayed correctly on STBP monitor but not correctly captured in ISIS

☐ BP on STBP OK, not captured in ISIS ☐ BP displayed normally on BP monitor, different data captured in ISIS

☐ Other

Problem obtaining HR

☐ Too much variability on HR on monitor to determine accurate HR

☐ Other

Problem with electrodes

☐ adherence problem due to excessive perspiration ☐ adherence problem due to excessive chest hair

☐ large chest/abdomen ☐ Other

Problem with BP cuff

☐ large arm caused difficulty with microphone placement ☐ cuff slipped off during test

☐ SP bending/flexing arm during test ☒ Other

Problem with treadmill controller

☐ Controller advanced more than one stage - Restarted test ☐ Other

8 of 9

End of Section Close Exam Finish

- **Problem with electrodes** - Causes of problems with electrodes include inadequate skin preparation, chest hair, faulty electrodes, and a large chest or abdomen.
- **Problem with BP cuff** - The SP may have a large arm that makes the positioning of the cuff difficult. If the cuff is not attached snugly to the arm, it may slip during the exam. In addition, if the SP is flexing when the cuff is applied, the cuff may slip out of position when the arm is relaxed. Bending or flexing the arm may also cause some problems. The SP should be encouraged to let his/her arm hang in a relaxed manner and should avoid extreme flexing although they do not have to hold their arms in a stiff uncomfortable position.
- **Problem with treadmill controller** - If the buttons on the treadmill are inadvertently pushed in error, the test may be restarted if this occurs before the first 15 seconds in Stage 1. If this happens after warm-up when changing protocols or going to recovery, you have the option of going to the next set of protocols on the controller.
- Check all data acquisition problems that apply. Be very careful that the selection you make accurately describes the problem you had.

4.11 Component Status Categories

4.11.1 Complete

Exhibit 4-76. Component status - complete

Cardiovascular Fitness Subsystem: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Cardiovascular Fitness Component Status

Status

☒ Complete
☐ Partial
☐ Not Done

Comments

Other text

9 of 9

End of Section Close Exam Finish

Ready Not connected to Coordinator 07:40 AM

Start NHANES ccMail R... Cardiovascular... Exploring - C:\bdc cv2.wav - Sound ... 7:41 AM

- Component status is set by the ISIS system as complete, partial, or not done. If all data were collected and entered, the system will set the status to complete.

4.11.2 Partial

Exhibit 4-77. Component status - partial (1)

Cardiovascular Fitness Subsystem: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Cardiovascular Fitness Component Status

Status

☐ Complete

☒ Partial

☐ Not Done

Comments

Other text

safety exclusion

SP refusal

no time

physical limitation

communication problem

equipment failure

9 of 9

End of Section Close Exam Finish

Finish slide show

Not connected to Coordinator 07:45 AM

Start NHANES ccMail Router Cardiovascular Fitness... Exploring - C:\bdc 7:45 AM

- If one or more pieces of data are missing (RPE measurement, blood pressure, PAR questions), the status will be set to partial.
- If the component status is partial, a comment is required. In some situations, the comment will be set to default. The comments that may be selected are listed below with a description of the situation where this comment is applicable. See Exhibits 4-77 to 4-80.
 - **Safety exclusion.** The SP is excluded from the exam due to a medical condition, current medications, or other situation that might cause the SP harm. The technician should not have to select this comment after the SP is in the CV Room because all the safety exclusion checks are completed prior to being assigned to this component. This would be set by the coordinator or physician applications.
 - **SP refusal.** The SP refused to do the exam.
 - **No time.** There was not sufficient time in the session to complete this exam.

- **Physical limitation.** The SP has a physical limitation that would prevent him/her from doing the exam. The technician should not have to select this comment after the SP is in the CV Room because all the physical limitation checks are completed prior to being assigned to this component. This would be set as default by the coordinator or physician applications.
- **Equipment failure.** The equipment malfunctioned making it impossible to complete the exam. This will not be selected because more specific information will be collected from the data acquisition screen.
- **Interrupted.** The exam or session had to be interrupted due to an emergency in this room or other area in the MEC making it impossible to complete the exam.

Exhibit 4-78. Component status - partial (2)

Cardiovascular Fitness Subsystem: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Cardiovascular Fitness Component Status

Status

☐ Complete

☒ Partial

☐ Not Done

Comments

Other text

SP ill/emergency

interrupted

Priority 1 stopping criteria

Priority 2 stopping criteria

SP pregnant

Unable to calculate protocol

9 of 9

End of Section Close Exam Finish

Finish slide show

Not connected to Coordinator 07:46 AM

Start NHANES ccMail Router Cardiovascular Fitness... Exploring - C:\bdc 7:46 AM

Other comments for 'partial' include:

- **Priority 1 Stopping Criteria.** The treadmill test was stopped early due to one or more of these criteria. See Section 3.13.1.
- **Priority 2 Stopping Criteria.** The treadmill test was stopped early due to one or more of these criteria. (Recovery Stages completed). See Section 3.13.2.

- **Protocol change.** The protocol was increased or decreased.
- **SP pregnant.** The SP was more than 12 weeks pregnant. See Section 3.1.1.
- **Unable to calculate protocol.** The ISIS system was unable to calculate the protocol because the data for height and weight for PAR Code data were not entered.

Exhibit 4-79. Component status - partial (3)

Cardiovascular Fitness Subsystem: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

File View Utilities Reports Window Help

Cardiovascular Fitness: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm

SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM

Cardiovascular Fitness Component Status

Status

☐ Complete

☒ Partial

☐ Not Done

Comments

Other text

Unable to follow protocol

Weight limitation on equipment

Problem obtaining BP

Problem obtaining HR

Problem with data capture

Problem with electrodes

9 of 9

End of Section Close Exam Finish

Finish slide show Not connected to Coordinator 07:46 AM

Start NHANES ccMail Router Cardiovascular Fitness... Exploring - C:\bdc 7:47 AM

Other comments for 'partial' include:

- **Unable to follow protocol.** The SP is unable to follow the directions for walking on the treadmill and are unable to get accustomed to the exercise.
- **Weight limitation on equipment.** The weight limitation on the treadmill is 400 pounds. If the weight is greater than 350 pounds, SPs are excluded from this exam.

The following comments will be selected on the data acquisition screen and are described in Section 4-10. If one or more of these problems are selected in the data acquisition screen, the comment 'data acquisition' will be set to default for this comment field.

- **Problem obtaining BP** – data acquisition problem set to default
- **Problem obtaining HR** – data acquisition problem set to default
- **Problem with data capture** – data acquisition problem set to default
- **Problem with electrodes** – data acquisition problem set to default
- **Problem with BP cuff** – data acquisition problem set to default
- **Problem with BP controller** – data acquisition problem set to default

Exhibit 4-80. Component status – partial (4)

The screenshot shows a Windows-style application window titled "Cardiovascular Fitness Subsystem: Stand:703 Session:703360 04/09/1999 08:30 am - 12:30 pm". The menu bar includes File, View, Utilities, Reports, Window, and Help. Below the menu is a toolbar with various icons. The main window area displays patient information: "SP ID: 618864 Name: PARKER, JACK Age: 21 years Gender: Male Date: 08/14/1999 Time: 07:15 AM". Below this is the "Cardiovascular Fitness Component Status" section. A dialog box is open with the title "Status". It contains three radio buttons: "Complete", "Partial" (which is selected), and "Not Done". Below the radio buttons is a "Comments" section with a text area labeled "Other text" and a list box. The list box contains the following items: "Problem with data capture", "Problem with electrodes", "Problem with BP cuff", "Problem with treadmill controlle", "Data acquisition problems", and "other, specify". The list box is currently empty. At the bottom of the dialog box, there are navigation buttons: "End of Section", "Close Exam", and "Finish". The bottom of the main window shows a status bar with "Finish slide show", "Not connected to Coordinator", and "07:46 AM". The Windows taskbar at the very bottom shows the Start button, several open applications including "NHANES ccMail Router", "Cardiovascular Fitnes...", and "Exploring - C:\bdc", and the system clock showing "7:47 AM".

- **Data acquisition problems.** If one or more of the problems selected on the data acquisition screen are selected, 'data acquisition problems' will be selected as the default comment. The specific information about the problem can be obtained from the data acquisition screen, which allows multiple entries as opposed to the comments for component status which allow only one comment to be selected.
- **Other.** Other situations that cause the SP to be excluded from this exam at some point during the exam should be noted.

Exhibit 4-81. Component status - partial - required comments

Cardiovascular Fitness Subsystem: Stand:806 Session:808571 12/11/1998 05:30 pm - 09:30 pm

SP ID: 342322 Name: Dee, Mike Age: 31 years Gender: Male Date: 01/24/1999 Time: 07:14 PM

Cardiovascular Fitness Component Status

Status

☐ Complete

☒ Partial

☐ Not Done

Comments

Other text

Cardiovascular Fitness

Please, select comments

OK

- If a comment is not selected from the drop-down menu, a message will be displayed: 'Please select comments.'

4.11.3 Not Done

Exhibit 4-82. Component status - not done

The screenshot displays the 'Cardiovascular Fitness Subsystem' window. The title bar indicates 'Stand: 806 Session: 808571 12/11/1998 05:30 pm - 09:30 pm'. The menu bar includes 'File', 'View', 'Utilities', 'Reports', 'Window', and 'Help'. The toolbar contains various icons for file operations and help. The main window area shows patient information: 'SP ID: 977952', 'Name: Ball, Britney', 'Age: 37 years', 'Gender: Female', 'Date: 01/24/1999', and 'Time: 07:26 PM'. Below this, the 'Cardiovascular Fitness Component Status' section is visible. It features a 'Status' group box with three radio buttons: 'Complete', 'Partial', and 'Not Done' (which is selected). To the right of the 'Status' group box is a 'Comments' text area with a dropdown menu. The dropdown menu is open, showing a list of options: 'safety exclusion', 'SP refusal', 'no time', 'physical limitation', 'communication problem', 'equipment failure', and 'SP ill/emergency'. At the bottom of the window, there is a navigation bar with buttons for 'End of Section', 'Close Exam', and 'Finish'. A status bar at the very bottom indicates 'Finish slide show'.

- The comments for 'Not Done' are the same as for 'partial.'


5. REFERRALS AND REPORT OF FINDINGS

5.1 Observation Referral

Observation referrals are nonemergency situations that may arise in any of the examination rooms in the MEC. Technicians may send an observation referral to the MEC physician if they notice any condition that may be abnormal or that may warrant further assessment. This type of referral may be sent at any time from any of the exam rooms. The referral may or may not have anything to do with the current exam being performed.


Once a technician sends a referral to the physician, the Integrated Survey Information System (ISIS) will flag the referral for the SP in the Physician Referral Review Box. (ISIS is the data entry system that will be used in this study; see Section 1.3 for a brief description.) The SP will not be checked out of the MEC until the physician has reviewed this referral. The physician will make a decision whether further action is warranted and a physician referral to the SP's health care provider may be given to the SP.

Exhibit 5-1. Menu to select observation referral

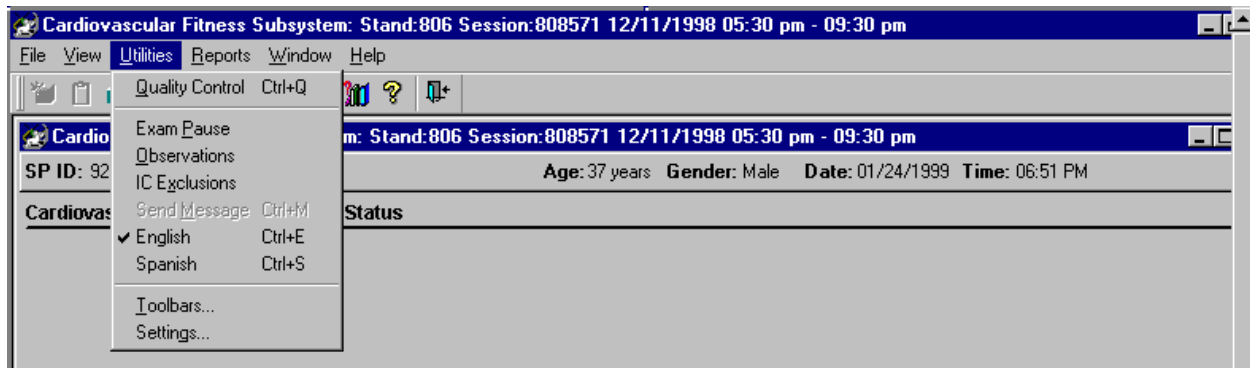

- Observation referrals can be sent during the exam or after the exam has been closed.
- Under 'Utilities,' select 'Observations.' If an exam is already opened, the observation referral box for that SP will be displayed. See Exhibit 5-1.
- If the exam has been closed, select 'Observation' from the 'Utilities' menu. A pick list with the names of the SPs in the current session will be displayed. This list is only displayed if 'Utilities/Observations' is selected when an exam is not open. See Exhibit 5-2.

Exhibit 5-2. Pick list of SPs in current session

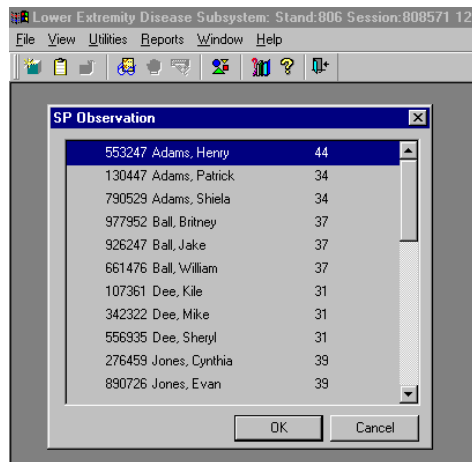

- Select the name of the SP for whom an observation referral should be sent. Click OK.

Exhibit 5-3. Observation referral in CV Fitness

Cardiovascular Fitness Subsystem: Stand:806 Session:808571 12/11/1998 05:30 pm - 09:30 pm

SP ID: 926247 Name: Ball, Jake Age: 37 years Gender: Male Date: 01/24/1999 Time: 06:51 PM

Cardiovascular Fitness Component Status

Observations

SP ID: 926247 Name: Ball, Jake Age: 37 Gender: M

This is the place to write an observation referral to the MEC physician about this SP. Type your message here. The system will send the message to the physician. The SP will not be checked out of the MEC until the physician has reviewed this message.

OK Cancel

8 of 8 End of Section Close Exam Finish

- Type in the message you would like to send to the physician. When you are finished, click OK.

Exhibit 5-4. Observation referral from CV Fitness in physician's referral review box

The screenshot shows a software window titled "Referral Review Stand:806 Session: 808571". It contains a table with columns for "SP", "Pregnancy", "CBC", "BP", "LED", and "Observation". Each column has sub-columns for "Name", "ID", "Age", "Gender", and "Reviewed Status". Below the table, a patient's details are displayed: "J. Ball", "926247", "37 Male". A "Review in Box Details" pop-up window is open, showing "SP ID: 926247", "Name: Jake, Ball", "Age: 37", and "Gender: M". This pop-up has tabs for "Pregnancy", "CBC", "BP", "LED", and "Observations". The "Observations" tab is selected, showing a section titled "OBSERVATIONS Cardiovascular Fitness" with a text area containing a referral message: "This is the place to write an observation referral to the MEC physician about this SP. Type your message here. The system will send the message to the physician. The SP will not be checked out of the MEC until the physician has reviewed this message."

- This is the referral as it appears in the physician's referral review box. The message typed in the observation referral box in CV Fitness appears in this box in the physician referral review.
- The physician will review this referral and make a decision about further action if warranted.
- The SP cannot be checked out of the MEC until the physician has reviewed this referral.

5.2 Report of Findings for CV Fitness

Each SP will be given a report of the results or findings for the exams performed. The Report of Findings for the Cardiovascular Fitness Exam will include a report on the SP's level of cardiovascular fitness. The level of cardiovascular fitness will be based on the estimated maximal oxygen consumption achieved by the SP in the treadmill test. SPs will be categorized as having a low, moderate, or a high level of cardiovascular fitness compared to other people of their age and sex.

Reference ranges for estimated maximal oxygen consumption (VO_2 max) are listed in Appendix H. The reference ranges used for this report are based on studies completed at the Cooper Institute for Aerobic Research and published in the Fifth Edition of the American College of Sports Medicine (ACSM) *Guidelines for Exercise Testing and Prescription*. There are no national norms based on random samples of the population (until data is obtained from NHANES).

General Statement

"Your fitness test was done on a treadmill. The test consisted of a warm-up, two exercise periods (Stage 1 and Stage 2), and a recovery period. The table below shows your test results.

<u>Stage</u>	<u>Required time</u>	<u>Your time</u>	<u>Your heart rate</u>
Warm-up	2 min.		
Stage 1	3 min.		
Stage 2	3 min.		
Recovery	2-3 min.		

During this exercise your maximum incline on the treadmill was ____ percent and your maximum speed was _____ mph.

Compared with other people of your age and sex, your cardiovascular fitness level is *<low, moderate, high>*."

In addition to the general statement, the report will include the speed and grade of the treadmill at each of the stages in the treadmill test and the length of time spent in each of the stages. SPs will be given this report just before they leave the MEC. The report includes a phone number to call if they have any questions about this or any part of the exam in the MEC.

6. QUALITY CONTROL

6.1 Overview

Quality assurance/quality control procedures for the Cardiovascular Fitness component includes development of standard procedures for the treadmill test, and intensive training and certification of technicians prior to the start of the study. After study start-up, technician adherence to procedures is monitored regularly through observational visits and specific data analysis for evaluation of technician performance. The validity and reliability of specific data is also monitored routinely in addition to regular equipment calibration and maintenance procedures.

6.2 Training

A 5-day training period for the Cardiovascular Fitness component is conducted at the beginning of the study with a 2-day annual retraining. Consultants from the Cooper Institute and Westat development staff supervise training for the treadmill test, data entry, and integration issues with the biomedical equipment and the computer system.

6.2.1 Initial Training

The initial training period includes a discussion of the importance of collecting data on the cardiovascular fitness level of the Nation, followed by a brief overview of exercise physiology focusing on normal and abnormal responses to a graded exercise test. Methods for excluding participants from this component (medical history, medications, heart rate and blood pressure criteria, and safety exclusion questions) are reviewed. Technicians are taught how to explain and demonstrate the treadmill test to a participant and how to monitor a participant's comfort level and to observe signs of distress that would require stopping the test early. Operation, maintenance, and calibration procedures for the treadmill, programmable controller, and automatic blood pressure monitor are reviewed and practiced. Procedures for a normal treadmill test are taught in addition to special procedures for changing the exercise protocol, interrupting the test to go to the recovery stage early, or to stop the test immediately if the condition of the participant warrants it. Data entry and integration of the biomedical equipment is demonstrated and practice sessions are conducted.

A modular approach to training is used where the technicians are given classroom lectures using slide presentations and familiarization with each piece of equipment in a classroom session. The technician learns how to use the blood pressure monitor and to obtain heart rate and blood pressure readings in one session. Another session is used to familiarize the technicians with the operation of the treadmill and programmable controller. Other classroom sessions are used to train technicians in procedures for asking questions, giving explanations, and reading scripts. Once the technicians are exposed to modular aspects of the treadmill test procedures in the classroom setting, practice sessions are moved to the mobile examination center (MEC). Here the technicians are given an opportunity to complete a practice treadmill on a classmate. This provides an opportunity for the technician to experience a treadmill test from the participant's perspective and to practice the test in a low stress situation. Volunteer participants are used for the latter part of training to expose technicians to a variety of individuals from diverse backgrounds and age groups and who are not familiar with the protocol.

At the end of the training session, each technician is given a written exam and a technician competency form is completed. The technicians must pass the written exam and complete all the competencies on the form before being cleared to conduct the treadmill test. Technicians who do not pass the exam or have not completed all competencies will be given additional training and testing on aspects of the procedures that need more attention. A debriefing session is held at the end of training and technicians are given individual and group feedback on performance.

6.2.2 Annual Retraining

Annual retraining is conducted over a 2-day period on site in the MECs. Topics for these sessions included protocol changes that have been started within the last year, staff requests, and procedures that have caused problems or have been noted on quality control visits. A CIAR consultant and the Westat project staff person conduct these sessions. The retraining is given on both an individual and group level. Each technician is observed for a minimum of two treadmill tests and is individually questioned on specific protocol issues and is questioned on responses or actions that should be taken to specific scenarios.

A debriefing session is held at the end of the retraining session to provide an opportunity for group discussion and mutual feedback on this component including problem areas, unusual situations, and any modifications in procedures. A technician competency form is completed for each technician and a Site Visit Report form is completed for each team.

Additional training is conducted throughout the year if necessary, depending on protocol changes, technician performance, and quality control reports.

6.3 Equipment and Room Set-up Checks

The equipment, room supplies, and room set-up need to be checked on a regular basis. Some checks are completed daily and others need only be completed on a weekly basis or at the beginning of each stand. These checks include calibration checks, maintenance inspection of equipment and supplies, and preparation of the room and equipment for the session exams.

Each time you log on to the application, the system will remind you to do Quality Control (QC) checks if the checks have not been completed for that time period. The checks are to be completed daily, weekly, and/or every stand. If you do not have time to do the checks when you log on, you can bypass this message and complete the checks at a later time. However, this message will be displayed each time you log on until you have completed the checks for that time period. Once you have completed the checks and entered this in the system, the message box with the reminder will not be displayed again until the appropriate time period has passed.

The daily, weekly, and once-a-stand checks are listed in the following sections.

6.3.1 Daily

- Confirm that the daily equipment checks are complete.
 - Turn on the power on the treadmill controller.
 - Turn on the power on the blood pressure (BP) monitor. (The Colin STBP-780 is the automated electronic heart rate and blood pressure monitor that will be used to measure blood pressure and heart rate during the treadmill test. See Chapter 2 for a description and picture of this monitor.)
 - Confirm that the BP monitor power cord is plugged in at the power source and on the back of the BP monitor.
 - Confirm that the computer cable for the BP monitor is securely connected at the back of the BP monitor.
 - Confirm that the BP microphone and airhose are securely connected to the back of the BP monitor.

- Confirm that the cable from the treadmill is securely connected to the back of the treadmill controller.
- Confirm that the treadmill power cord is securely plugged in to the power source.
- Check that the digital thermometer has registered the room temperature.
- Check that the treadmill belt is centered in the roller.
- Visually inspect walking belt for hazards.
- Clean surface of belt with disinfectant wipes.
- Clean both sides of all BP cuffs with disinfectant wipes. (Cuffs are also wiped after each exam.)
- Visually inspect BP cuffs and hoses for cracks and tears.
- Check that each BP cuff is hanging properly from the clip on the grid. The microphone wires (taped areas on air hose below cuff) should be hanging free - not twisted or bent.
- Check that all shoes are in the shoe rack.
- Confirm that supplies (electrodes, alcohol, and abrasive pads) are stocked for today's session.

6.3.2 Weekly

- Complete all daily checks.
- Check speed of treadmill. See Section 2.1.2.1.
- Check grade of treadmill. See Section 2.1.2.2.
- Check belt tension on treadmill. See Section 2.1.1.1.
- Check belt tracking on treadmill. See Section 2.1.1.2
- Check rear roller guards on treadmill. See Section 2.1.1.3
- Check calibration on BP monitor with mercury manometer (0-70 mmHg within +/- 5 mmHg). See Section 2.3.2.
- Check calibration on BP monitor with mercury manometer (70-120 mmHg within +/- 5 mmHg).
- Check calibration on BP monitor with mercury manometer (120-180 mmHg within +/- 5 mmHg).

- Check calibration on BP monitor with mercury manometer (180-220 mmHg within +/- 5 mmHg).
- Check calibration on BP monitor with mercury manometer (220-255 mmHg within +/- 5 mmHg).
- Check that calibration mercury manometer is intact. See Section 2.3.3.
- Calibrate the SECA scale - weekly weight check. If the scale is off, do a full calibration and if the scale is off by a pound or more with the 5 weights on it, inform the MEC manager and a replacement scale will be sent to the MEC. See Section 2.4.1 for calibration procedures.
- Calibrate the SECA stadiometer. If the reading is off by greater than 1 cm, inform the MEC manager. Use the stadiometer until a replacement is sent. See Section 2.5.1 for calibration procedures.


6.3.3 Stand

- Complete all daily checks.
- Complete all weekly checks.
- Remove all BP cuffs from bladder and wash cuffs with soap and water. Allow drying.
- Wipe down BP monitor with cloth dampened with water. Do not use alcohol.
- Wipe down treadmill controller with cloth dampened with water. Do not use alcohol.
- Confirm treadmill controller settings for speed and grade for all stages are accurate (for each of the eight protocols) (Chief Technician).
- Confirm with the data manager that the BP monitor is plugged in to an uninterruptible power supply (UPS). (Enter 'Confirmed with {initials of data manager}' in Comments field).
- Confirm with the data manager that the treadmill and the UPS for the BP monitor are on different power circuits. (Enter 'Confirmed with {initials of data manager}' in Comments field).
- Calibrate SECA scale - Full. If the scale is off by 1 pound or more with the 5 weights placed on it, inform the MEC manager. A replacement scale will be sent to the MEC. See Section 2.4.1 for calibration procedures.
- Calibrate SECA stadiometer. If the reading is off by greater than 1 cm, inform the MEC manager. Use the stadiometer until a replacement is sent. See Section 2.4.2 for calibration procedures.

6.4 Data Entry Screens for QC on Equipment

6.4.1 Quality Control Reminder Message Box


Exhibit 6-1. Quality control reminder message box

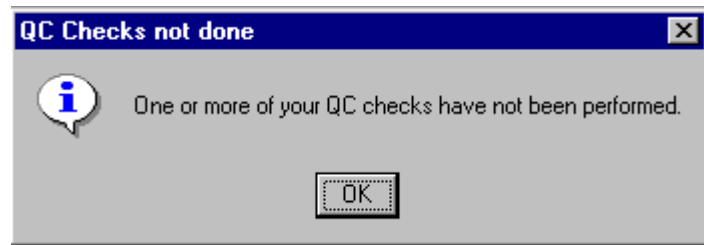

- When you log on to the application before the quality control checks are performed, the system displays a message: 'One or more of your QC checks have not been performed.' See Exhibit 6-1.
- Click OK to this message.

6.4.2 Utilities Menu to Select Quality Control


Exhibit 6-2. Utilities menu to select quality control

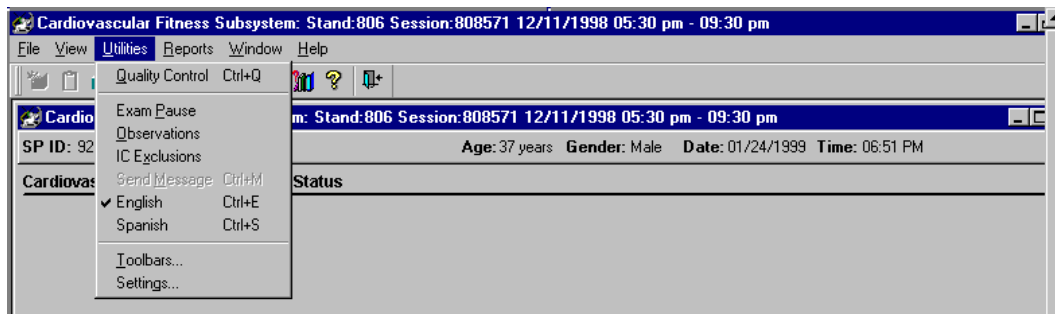

- When you want to complete the QC checks, select 'Utilities.' Then select 'Quality Control' from the menu. See Exhibit 6-2.
- Clicking on the QC icon from the Toolbar can also access the QC screens.

6.4.3 Quality Control Log-on Box

Exhibit 6-3. Quality control log-on

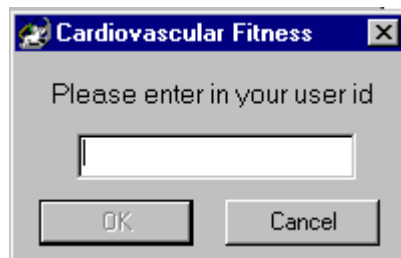

When QC is selected from the Utilities menu, the User ID entry box will be displayed. See Exhibit 6-3. Each technician will have a personal ID. This ID will be used to identify the person who completed the QC checks for this time period.

- Enter your User ID and click OK.
- If you do not want to do the QC checks at this time, click Cancel.

6.4.4 Data Entry Screen for QC Daily Checks

Exhibit 6-4. Quality control daily checks (1)

QC Check	Done	Result	Comment
Turn power on treadmill controller.	<input type="checkbox"/>		
Turn power on BP monitor.	<input type="checkbox"/>		
Confirm BP monitor power cord is plugged in at UPS and at back of BP monitor.	<input type="checkbox"/>		
Confirm computer cable for BP monitor is securely connected at back of BP monitor.	<input type="checkbox"/>		

On the QC screens, check 'Done' for the listed items when that item has been completed. You are not required to enter anything in the 'Result' or 'Comment' fields unless there is a problem. The 'Result' field is used to enter values for selected QC items if required. The 'Comment' field is used to enter information about problems encountered with the QC item check. This field is also used to indicate if something has been reported or confirmed with another staff member (Chief Technician, MEC manager, data manager).

- Turn on the power to the treadmill controller and the BP monitor.
- Confirm that the power cord on the BP monitor is plugged in securely at the back of the monitor and at the UPS.
- Confirm that the computer cable connecting the BP blood pressure monitor to the ISIS monitor is connected securely to the back of the blood pressure monitor.

Exhibit 6-5. Quality control daily checks (2)

QC Check	Done	Result	Comment
Confirm BP microphone and airhose are securely connected at back of BP monitor.	<input type="checkbox"/>		
Confirm cable from treadmill is securely connected to back of treadmill controller.	<input type="checkbox"/>		
Confirm treadmill power cord is securely plugged in to power source.	<input type="checkbox"/>		
Check digital thermometer has registered room temperature.	<input type="checkbox"/>		

Use the scroll bar to move to the next items on the list. See Exhibit 6-5.

- Confirm that the microphone and air hose are securely connected at the back of the BP monitor.
- Confirm that the cable from the treadmill is securely connected to the back of the treadmill controller.
- Confirm that the treadmill power cord is securely plugged into the power source.
- Check that the digital thermometer has registered the room temperature.

Exhibit 6-6. Quality control daily checks (3)

QC Check	Done	Result	Comment
Check treadmill belt is centered in roller.	<input type="checkbox"/>		
Visually inspect walking belt for hazards.	<input type="checkbox"/>		
Clean surface of belt with disinfectant wipes.	<input type="checkbox"/>		
Clean both sides of all BP cuffs with disinfectant wipes. (Cuffs are also wiped after each exam.)	<input type="checkbox"/>		

See Exhibit 6-6 for the additional items.

- Check that the treadmill belt is centered in the rollers.
- Visually inspect the walking belt for objects that may have been left on the belt. Remove all objects and keep the belt clear for the entire session.
- Wipe down the belt with disinfectant wipes.
- Clean both sides of the BP cuffs with disinfectant wipes. The BP cuffs should be wiped down with disinfectant wipes after each use.

Exhibit 6-7. Quality control daily checks (4)

QC Check	Done	Result	Comment
Visually inspect BP cuffs and hoses for cracks and tears.	<input type="checkbox"/>		
Check BP cuffs are hanging from clip on grid with no twisting on microphone wires.	<input type="checkbox"/>		
Check that all shoes are in the shoe rack.	<input type="checkbox"/>		
Confirm electrodes, alcohol and abrasive pads are stocked for this session.	<input type="checkbox"/>		

- See Exhibit 6-7 for the remaining items.
- Check all BP cuffs and hoses for cracks and tears.
- Check that the cuffs are hanging properly from the grid. The cuffs should be clipped to the grid with the weight of the microphone tubing and air hose supported in the basket beneath. There should be no twisting on the microphone wires.
- Check that all shoes are in the shoe rack.
- Confirm that there are adequate supplies for the entire session and restock as necessary.
- When you are finished with the daily item checks, click OK to close the QC box.

6.4.5 Data Entry Screen for QC Weekly Checks

Exhibit 6-8. Quality control weekly checks (1)

QC Check	Done	Result	Comment
Complete all daily checks.	<input type="checkbox"/>		
Check speed of treadmill.	<input checked="" type="checkbox"/>	Difference=0.8	Controller speed=2.5.Belt speed=3.3: Reported to MEC Manager
Check grade of treadmill.	<input type="checkbox"/>		
Check belt tracking on treadmill.	<input type="checkbox"/>		

- Complete all daily checks. Check 'Done' for each item on the weekly checks when complete. See Exhibit 6-8.
- Check speed of the treadmill. See Section 2.1.2.1 for calibration procedures for speed of treadmill. If the speed readout on the treadmill controller exceeds +/- 0.2 mph of actual speed, repeat the procedure. If the treadmill still exceeds the range, enter the difference between the belt speed on the controller and the calculated belt speed in the 'Result' field. (Example: Difference = 0.8 mph). Enter the speed of the treadmill controller and the calculated belt speed in the 'Comment' field. (Example: Controller speed = 2.5 mph; Belt speed = 3.3 mph). Report this to the Chief Technician. Note this report in the 'Comment' field. (Example: Reported to Chief Technician). The Chief Technician or MEC manager will call a service representative.
- Check grade of the treadmill. See Section 2.1.2.2 for calibration procedures for grade of treadmill. If the elevation error exceeds +/- 1% grade, repeat the calibration check. If the elevation error still exceeds +/- 1% grade, enter the elevation error in the 'Result' field. (Example: Elevation error = 1.5%). Enter the treadmill grade setting and the measured grade in the 'Comment' section. (Example: Treadmill grade setting=20%; Measured grade=21.5%). Report this to the Chief Technician. Note the report in the 'Comment' section. (Example: Reported to Chief Technician.). The Chief Technician will call a service representative.
- Check belt tracking on treadmill. See Section 2.1.1.2.

Exhibit 6-9. Quality control weekly checks (2)

QC Check	Done	Result	Comment
Check belt tension on treadmill.	<input type="checkbox"/>		
Check rear roller guards on treadmill.	<input type="checkbox"/>		
Check calibration on STBP monitor with mercury manometer (0-70 mmHg within +/- 5 mmHg).	<input type="checkbox"/>		
Check calibration on STBP monitor with mercury manometer (70-120 mmHg within +/- 5 mmHg).	<input type="checkbox"/>		

- Check belt tension on treadmill. See Section 2.1.1.1 for a description of how to check treadmill belt tension.
- Check rear roller guards. See Section 2.1.1.3 for a description of how to check rear roller guards.
- Check calibration on the BP monitor with the mercury manometer (0-70 mmHg within +/- 5 mmHg). See Section 2.3.2 for BP calibration procedures. If the BP calibration checks are not within the acceptable range, enter the difference between the mercury manometer reading and the reading on the BP monitor in the 'Result' field. (Example: Difference = 8mmHg). See Figure 6-10. Enter the reading on the manometer and the reading on the BP monitor in the 'Comment' field. (Example: Manometer = 110; BP=118). Notify the Chief Technician. Make a note of this in the 'Comment' section. (Example: Reported to the Chief Technician.)
- Check calibration on BP monitor with mercury manometer (70-120 mmHg within +/- 5 mmHg). If the BP calibration checks are not within the acceptable range, follow the same procedures as in the preceding bullet.

Exhibit 6-10. Quality control weekly checks (3)

QC Check	Done	Result	Comment
Check calibration on STBP monitor with mercury manometer (120-180 mmHg within +/- 5 mmHg).	<input checked="" type="checkbox"/>	Difference=8	Manometer=110; STBP=118. Reported to MEC Manager.
Check calibration on STBP monitor with mercury manometer (180-220 mmHg within +/- 5 mmHg).	<input checked="" type="checkbox"/>	Difference=8	Manometer=110; STBP=118. Reported to MEC Manager.
Check calibration on STBP monitor with mercury manometer (220-255 mmHg within +/- 5 mmHg).	<input checked="" type="checkbox"/>	Difference=8	Manometer=110; STBP=118. Reported to MEC Manager.
Check calibration mercury manometer is intact.	<input checked="" type="checkbox"/>		

- Check calibration on the BP monitor with the mercury manometer (120-180 mmHg within +/- 5 mmHg). If the BP calibration checks are not within the acceptable range, follow the same procedures as in the first bullet for calibration of the BP monitor.
- Check calibration on the BP monitor with mercury manometer (180-220 mmHg within +/- 5 mmHg). If the BP calibration checks are not within the acceptable range follow the same procedures as in the first bullet for calibration of the BP monitor.
- Check calibration on the BP monitor with mercury manometer (220-255 mmHg within +/- 5 mmHg). If the BP calibration checks are not within the acceptable range, follow the same procedures as in the first bullet for calibration of the BP monitor.
- Check that the calibration of the mercury manometer is intact. The procedures for this are in Section 2.3.3 and are also on the laminated Quick Reference cards in the CV Fitness room.

Exhibit 6-11. Quality control weekly checks (4)

QC Check	Done	Result	Comment
Check calibration on STBP monitor with mercury manometer (220-255 mmHg within +/- 5 mmHg).	<input type="checkbox"/>		
Check calibration mercury manometer is intact.	<input type="checkbox"/>		
Calibrate SECA Scale - weekly weight check	<input type="checkbox"/>		
Calibrate SECA stadiometer	<input type="checkbox"/>		

- Calibrate the SECA scale - weekly weight check. If the scale is off, do a full calibration and if the scale is off by a pound or more with the 5 weights on it, inform the MEC manager and a replacement scale will be sent to the MEC. See Section 2.4.1 for calibration procedures.
- Calibrate the SECA stadiometer. If the reading is off by greater than 1 cm, inform the MEC manager. Use the stadiometer until a replacement is sent. See Section 2.5.1 for calibration procedures.

6.4.6 Data Entry Screen for QC Start of Stand Checks

Exhibit 6-12. Quality control 'start of stand' checks (1)

The screenshot shows a software window titled "Cardiovascular Fitness Quality Control Checks". It has three tabs: "Start of Stand", "Daily", and "Weekly". The "Start of Stand" tab is selected. Below the tabs is a table with four columns: "QC Check", "Done", "Result", and "Comment". There are four rows of checks, each with a checkbox in the "Done" column, a text input field in the "Result" column, and a larger text input field in the "Comment" column. The checks are: "Complete all daily checks.", "Complete all weekly checks.", "Remove all BP cuffs from bladder and wash cuffs with soap and water.", and "Wipe down STBP monitor with cloth dampened with water. Do not use alcohol." At the bottom right of the window are "OK" and "Cancel" buttons.

QC Check	Done	Result	Comment
Complete all daily checks.	<input type="checkbox"/>		
Complete all weekly checks.	<input type="checkbox"/>		
Remove all BP cuffs from bladder and wash cuffs with soap and water.	<input type="checkbox"/>		
Wipe down STBP monitor with cloth dampened with water. Do not use alcohol.	<input type="checkbox"/>		

- Complete all daily checks.
- Complete all weekly checks.
- Remove all BP cuffs from bladder and wash cuffs with soap and water. Allow time to dry before using again.
- Wipe down the BP monitor with a damp cloth. Do not use alcohol.

Exhibit 6-13. Quality Control 'start of stand' checks (2)

The screenshot shows the same software window as Exhibit 6-12, but with different checks listed under the "Start of Stand" tab. The table structure is identical, with four columns: "QC Check", "Done", "Result", and "Comment". There are four rows of checks, each with a checkbox in the "Done" column, a text input field in the "Result" column, and a larger text input field in the "Comment" column. The checks are: "Wipe down treadmill controller with cloth dampened with water. Do not use alcohol.", "Confirm treadmill controller settings for speed and grade for all protocols.", "Confirm with Data Manager that BP monitor is plugged in to UPS. (Enter 'Confirmed with DM' in Comments field)", and "Confirm with DM that UPS for BP monitor and treadmill are on different power circuits. (Enter 'Confirmed with DM')". At the bottom right of the window are "OK" and "Cancel" buttons.

QC Check	Done	Result	Comment
Wipe down treadmill controller with cloth dampened with water. Do not use alcohol.	<input type="checkbox"/>		
Confirm treadmill controller settings for speed and grade for all protocols.	<input type="checkbox"/>		
Confirm with Data Manager that BP monitor is plugged in to UPS. (Enter 'Confirmed with DM' in Comments field)	<input type="checkbox"/>		
Confirm with DM that UPS for BP monitor and treadmill are on different power circuits. (Enter 'Confirmed with DM')	<input type="checkbox"/>		

- Wipe down treadmill controller with cloth dampened with water. Do not use alcohol.

- Confirm treadmill controller settings for speed and grade for all protocols. The Chief Technician is responsible for doing this check. See Section 2.2.3 for details and Appendix J for the Speed and Grade Table by Protocol and Stage. These procedures and the table are also on the laminated Quick Reference cards in the CV Fitness room.
- Confirm with the data manager that the BP monitor is plugged into the UPS. Enter 'Confirmed with the data manager' in the Comment field.
- Confirm with the data manager that the UPS for BP monitor and treadmill are on different power circuits. Enter 'Confirmed with the data manager' in the Comment field. If not, ask him/her to do this and then enter 'Confirmed with the data manager' in the Comment field. The BP monitor should be plugged into the UPS in the room to minimize the 'noise' to the monitor from the treadmill and other sources.

Exhibit 6-14. Quality Control 'start of stand' checks (3)

The screenshot shows a software window titled "Cardiovascular Fitness Quality Control Checks". It has three tabs: "Start of Stand", "Daily", and "Weekly". The "Start of Stand" tab is selected. Below the tabs is a table with four columns: "QC Check", "Done", "Result", and "Comment". There are four rows of checks, each with a checkbox in the "Done" column, a text input field in the "Result" column, and a larger text input field in the "Comment" column.

QC Check	Done	Result	Comment
Confirm with Data Manager that BP monitor is plugged in to UPS. (Enter 'Confirmed with DM' in Comments field)	<input type="checkbox"/>		
Confirm with DM that UPS for BP monitor and treadmill are on different power circuits. (Enter 'Confirmed with DM')	<input type="checkbox"/>		
Calibrate SECA Scale - Full	<input type="checkbox"/>		
Calibrate SECA stadiometer	<input type="checkbox"/>		

At the bottom right of the window are "OK" and "Cancel" buttons.

- Calibrate SECA scale - Full. If the scale is off by 1 pound or more with the 5 weights placed on it, inform the MEC manager. A replacement scale will be sent to the MEC. See Section 2.4.1 for calibration procedures.
- Calibrate SECA stadiometer. If the reading is off by greater than 1 cm, inform the MEC manager. Use the stadiometer until a replacement is sent. See Section 2.5.1 for calibration procedures.
- When the 'Start of Stand' checks are complete, click OK to close the QC checks.

6.4.7 Incomplete Entry of Quality Control Data

Exhibit 6-15. Quality control incomplete entry

The screenshot shows a window titled "Cardiovascular Fitness Quality Control Checks". It has three tabs: "Start of Stand", "Daily", and "Weekly". The "Daily" tab is selected. Below the tabs is a table with four columns: "QC Check", "Done", "Result", and "Comment".

QC Check	Done	Result	Comment
Remove all BP cuffs from bladder and wash cuffs with soap and water.	<input checked="" type="checkbox"/>		
Wipe down STBP monitor with cloth dampened with water. Do not use alcohol.	<input checked="" type="checkbox"/>		
Wipe down treadmill controller with cloth dampened with water. Do not use alcohol.	<input checked="" type="checkbox"/>		
Confirm treadmill controller settings for speed and grade for all protocols.	<input type="checkbox"/>		

Below the table, there is a message box with a yellow warning icon and the text: "Not all of the QC items were done. Do you wish to exit?". The message box has "Yes" and "No" buttons. The "QC Check" table has "OK" and "Cancel" buttons at the bottom right.

- When you have completed all checks, click OK to close the QC box.
- If you do not check that all items are complete, the system will display this message: 'Not all the QC items were done. Do you wish to exit?' See Exhibit 6-15.
- If you want to complete the items before exiting, click 'No' to this message and complete the items.
- If you do not wish to exit without completing all the QC checks, click 'Yes' to this message.
- If all QC items were not complete, the system will remind you each time you log on that the QC checks are not complete.

6.5 Observational Visits

The Cardiovascular Fitness component has two types of quality control visits, expert visits conducted by consultants from the Cooper Institute and nonexpert visits conducted by Westat project development staff. Consultants from the CIAR will make two field visits per team per year. Each visit will last approximately 2 days. Whenever possible, visits will be made during peak testing days in order to observe as many CV tests as possible. The consultants will complete the Technician Competency Evaluation form (See Section 6.5.1) during each test and review the completed form with the technician. Any necessary improvements that need to be made by the technician will be documented on the form, and the technician will sign the form after reading it. Technicians who are observed having CV testing problems during field visits will be observed more frequently than other technicians. A secondary activity

during each site visit will be to monitor the physician's examination to ensure that the CV exclusion questions are being asked according to the script.

Nonexpert quality control visits are conducted twice a year for each team. Westat development staff for this component will conduct these visits to monitor procedures and to conduct re-training for changes in procedures if necessary.

A Site Visit Report form is completed by the consultants and development staff following each field visit. The Site Visit Report form is used to identify, document, and follow up on issues or problems noted during observational field visits. In contrast to the Technician Competency Evaluation form, the Site Visit Report form is used to give an overall assessment of the status of the component and team performance. The intent is to report general problems and issues rather than problems and issues regarding the specific performance of specific technicians. The latter is reported on the Technician Competency Evaluation form.

6.5.1 Technician Competency Evaluation Form

This form is used by consultants and trainers to evaluate health technicians conducting the exam. Items are checked as Pass/Fail. If the technician does not pass an item, the observer will discuss this with the technician, offer suggestions for improvement, and re-evaluate performance in a subsequent exam later in that visit or at the next visit. See template form in Exhibit 6-16.

Exhibit 6-16. Technician Competency Evaluation Form

NHANES CARDIOVASCULAR FITNESS COMPONENT
Quality Control
Technician Competency Evaluation Form

Technician Name: _____ Stand Number: _____

Cooper Institute Examiner Name: _____ Date: _____

	Satisfactory	Unsatisfactory
Introduction to Component		
<ul style="list-style-type: none"> ■ Greets SP and brings into treadmill room ■ Provides brief overview of treadmill component ■ Scans SP identification number into system 		
Obtaining Weight and Height Measurement		
<ul style="list-style-type: none"> ■ Position SP on the weight scale and instruct SP to stand still prior to taking weight measurement ■ Position SP on the stadiometer and instruct SP to inhale prior to taking height measurement ■ Enters answers correctly into system 		
Determination of PAR Code		
<ul style="list-style-type: none"> ■ Reads questions clearly and quickly ■ Encourages SP to answer questions quickly ■ Enters answers correctly into system 		
Check Appropriate Footwear		
<ul style="list-style-type: none"> ■ No heels, cleats, or dress shoes ■ Shoes tied securely 		
Preparation of Skin and Placement of Electrodes		
<ul style="list-style-type: none"> ■ Correctly identifies electrode sites ■ Thoroughly wipes each site with alcohol swab ■ Lightly abrades skin at each site (6-10 swipes) ■ Applies electrode at each site and presses firmly ■ Clips RA, LA, RL, and LL wires to appropriate electrodes 		
Blood Pressure Cuff Placement		
<ul style="list-style-type: none"> ■ Selects appropriately sized cuff ■ Palpates for brachial artery ■ Places cuff 1 inch (2-3 cm) above antecubital space ■ Places cuff so that 2 microphones are directly over brachial artery ■ Cuff placed snugly, no slipping 		

Exhibit 6-16. Technician Competency Evaluation Form (continued)

	Satisfactory	Unsatisfactory
■ Explanation of Treadmill Test		
■ Follows script, speaks clearly		
Obtaining Pretest Blood Pressure and Heart Rate		
■ Instructs SP to stand quietly on treadmill with arms by sides		
■ Selects arm and cuff size on ISIS screen		
■ Selects “Capture” on ISIS screen		
Selecting Protocol On Treadmill Controller		
■ Presses ENTER, then INPUT + twice on controller		
■ Presses INPUT + or INPUT – until the correct protocol number is displayed on the controller		
Tech Demonstrates How to Walk on Treadmill		
■ Straddles belt near front, holding onto front rail		
■ Presses START BELT on treadmill controller		
■ Begins walking, stands up straight, takes long steps, looks ahead, drops right hand first, then left hand		
■ Presses STOP BELT on treadmill controller		
■ Presses CLEAR on treadmill controller		
Getting SP on to the Treadmill		
■ Instructs SP to straddle belt near front, holding onto front rail		
■ Presses START BELT on treadmill controller		
■ Instructs SP to step on belt and begin walking (big steps, eyes ahead, stand up straight)		
Starting the Treadmill Test		
■ Presses START on ISIS screen and press START COURSE on treadmill controller at the same time		
■ Encourages SP to let go of front rail one hand at a time		
■ A light touch with fingertips is OK, but no gripping rail		
■ If SP continues to grip front rail by end of warm-up, press GO TO RECOVERY on ISIS screen, press NEXT STAGE on treadmill controller so that recovery stage is displayed		
Explanation of RPE		
■ Reads script quickly and clearly when prompted by ISIS		
Obtaining RPE – Warm Up		
■ Asks SP for RPE at 1:30 of warm-up following auditory prompt from ISIS		
■ Enters RPE correctly into system		

Exhibit 6-16. Technician Competency Evaluation Form (continued)

	Satisfactory	Unsatisfactory
Normal Transition From Warm-Up to Stage I		
<ul style="list-style-type: none"> ■ Presses NEXT STAGE on treadmill controller at 1:59 following auditory prompt from ISIS 		
Obtaining RPE – Stage I		
<ul style="list-style-type: none"> ■ Asks SP for RPE at 2:30 of Stage I following auditory prompt from ISIS ■ Enters RPE correctly into system 		
Check HR – Stage I		
<ul style="list-style-type: none"> ■ Check heart rate on BP monitor and captured data in ISIS ■ Manually enter correct HR into system if the difference is > 5 beats. ■ HR checked _____ HR captured _____ HR entered _____ 		
Obtaining RPE – Stage II		
<ul style="list-style-type: none"> ■ Asks SP for RPE at 2:30 of Stage II following auditory prompt ■ Enters RPE correctly into system 		
Check HR – Stage II		
<ul style="list-style-type: none"> ■ Check heart rate on BP monitor and captured data in ISIS ■ Manually enter correct HR into system if the difference is > 5 beats. ■ HR checked _____ HR captured _____ HR entered _____ 		
Normal Recovery		
<ul style="list-style-type: none"> ■ Informs SP that treadmill will be slowing down and grade will level out ■ Informs SP that it is OK to lightly touch front railing on treadmill during recovery 		
Normal Test Termination		
<ul style="list-style-type: none"> ■ Presses COOL DOWN on treadmill controller after auditory prompt from ISIS ■ Tells SP to grab front rail because treadmill is about to stop ■ Presses STOP BELT on treadmill, be attentive to SP ■ Detaches lead wires from electrodes ■ Selects FINISH on ISIS screen 		
Changing Protocols (example scenario for decreasing protocol)		
<ul style="list-style-type: none"> ■ At 1:30 of warm-up, examiner presents scenario that warm up heart rate is >60% of predicted max ■ Selects CHANGE PROTOCOL on ISIS screen ■ Selects DECREASE and APPLY on ISIS screen ■ At 1:59 presses NEXT STAGE 4 times on treadmill controller so that next easiest protocol (Stage I) is selected 		

Exhibit 6-16. Technician Competency Evaluation Form (continued)

	Satisfactory	Unsatisfactory
Terminating a Test Early (example scenario for priority 2 stop)		
<ul style="list-style-type: none"> ■ At 2:00 into Stage I, examiner presents scenario that heart rate is 100% of predicted max ■ Selects GO TO RECOVERY on ISIS screen ■ Presses NEXT STAGE 2 times on treadmill controller so that recovery stage for current protocol is selected ■ When recovery stage is completed, alerts SP that treadmill is about to stop and tells SP to grab front rail ■ Presses STOP BELT on treadmill controller ■ Selects appropriate reason for Priority 2 stop on ISIS screen 		
Terminating a Test Early (priority 1)		
<ul style="list-style-type: none"> ■ Examiner presents a scenario that SP is having chest pain ■ Selects EMERGENCY STOP on ISIS screen ■ Alerts SP that treadmill is about to stop and tells SP to grab front rail ■ Presses COOL DOWN and STOP BELT on treadmill controller ■ Places SP in supine position on treadmill belt ■ Remains with SP and calls physician ■ Selects appropriate reason for Priority 1 stop in ISIS screen 		
Examiner Comments:		
My signature below confirms that Cooper Institute staff has reviewed this completed document with me, and made me aware of any and all cardiovascular component procedures that I need to improve upon.		
Signature _____	Date _____	

6.5.2 Quality Control Review of CV Fitness Exclusions in Physician's Exam

Two activities in the physician's exam are related to the cardiovascular fitness component: reviewing medications recorded at the SP household questionnaire for safety exclusions and administering CV fitness exclusion questions. Monitoring these activities is included as part of CV fitness component quality control visits.

6.5.2.1 Medication Review

Physicians are responsible for reviewing the medications reported by the SP at the household interview. Physicians will then determine if the SP should be excluded from the CV fitness exam based on a list of exclusion medications that are updated quarterly by the consultants from the Cooper Institute. This review process is expected to take place at the beginning of the session. The following form (Exhibit 6-17) lists the observation points for this part of the activity.

Exhibit 6-17. Checklist for medication review

Observation points	Yes	No	Specific comments and notes
Reported medications are reviewed with the exclusion list			
Result of the review is recorded correctly			

6.5.2.2 CV Fitness Exclusion Questions

A set of exclusion questions will be asked in the Physician Exam for the CV fitness component. This is a second level of safety check to prevent including SPs who may have been recently diagnosed with one of the exclusion medical conditions or have started taking one of the exclusion medications in the time period since the household interview. If the SP responds affirmatively to any of the exclusion questions, he/she will be immediately excluded and will not need to answer the remaining questions. Physicians are expected to read these questions exactly as written. Below is a checklist (Exhibit 6-18) used by the trainers in their quality control visits. The purposes of this list are to monitor the administration of the CV fitness exclusion questions and to aggregate specific comments for these questions, especially the necessity and appropriateness for extra prompt use.

Exhibit 6-18. Checklist for CV Fitness exclusion questions

CV Fitness exclusion questions	Follow script		SP receptive		Prompt use		Specific comments and notes
	Yes	No	Yes	No	Yes	No	
Have you been hospitalized in the past 3 months?							
Has a doctor ever said you should not participate in sports or other activities because of a health condition? (12-19 years only)							
Has a doctor ever said you have a heart condition and that you should only do physical activity recommended by a doctor?							
Do you feel pain in your chest when you do physical activity? (20-49 years only)							
In the past month, have you had chest pain when you were not doing physical activity? (20-49 years only) <i>Probe: Have you seen a medical doctor about your chest pain? Did the doctor tell you that the chest pain was related to your heart?</i>							
Do you lose your balance because of dizziness? <i>Probe: Is this an isolated incident or does it occur on a regular basis?</i>							
Do you ever lose consciousness? <i>Probes: Did this occur as a result of illness or was it unexplained? Is this an isolated incident or does it occur on a regular basis?</i>							
Do you have a bone or joint problem that could be made worse by walking? <i>Probe: Do you think you can do the test without injuring yourself?</i>							
Are you currently taking any prescription medications?							
Are you currently taking any prescription medications for your blood pressure?							
What is the name of this medication?							
Are you currently taking prescription medications for the following conditions: heart condition, or prescription eye drops for glaucoma?							
Do you know of any other reason why you should not do a treadmill test?							

6.5.3 Site Visit Report Form

The Site Visit Report form (Exhibit 6-19) is an online method to identify, document, and followup on issues or problems noted during observational field visits. This form is completed by consultants and development staff following each field visit and is posted to the project intranet site.

In contrast to the Technician Competency Evaluation form, the Site Visit Report form is used to give an overall assessment of the status of the component and team performance. The intent is to report general problems and issues rather than problems and issues regarding the specific performance of specific technicians. The latter is reported on the Technician Competency Evaluation form.

Exhibit 6-19. Site Visit Report form

Date(s) of Visit:		Stand #	
Name of Observer:		Visit # (this year to this team)	
Component:	Cardiovascular Fitness	# of SP exams observed:	
ID number of examiners/recorders observed			
<div>1. Outstanding issues from last visit or QC meeting (if any):</div> <div>2. Observations/issues:</div> <div>3. Issues that need to be addressed at the next visit:</div> <div>4. Comments:</div>			

APPENDIX A

SCRIPTS

APPENDIX A. SCRIPTS

English Version

Introduction to the component: “In this exam we want to get some information about your cardiovascular fitness or the fitness of your heart. We will be asking you to walk on the treadmill while we monitor your heart rate and blood pressure. I will explain this in more detail in a few minutes but first I want to ask you a few questions about your health in general. Please come over here near the computer while I scan your identification number into the system.”

Explanation of the test: “The purpose of the treadmill test is to estimate your cardiovascular fitness level. You will be walking on the treadmill at various speeds and grades. This is not a maximal test to exhaustion. We will stop when your heart rate reaches a certain level. While this is a very safe test, it is important that you let us know how you feel during the test. Please report any unusual symptoms that you have during the test such as chest pain, dizziness, nausea, or pain in your legs.”

Explanation of Rating of Perceived Exertion (RPE): “Please look at the chart on the wall in front of you. Every couple of minutes during the test, I will ask you to rate how hard you feel you are working overall. Please respond with a number from the scale. For example, if you feel that you are working at a very low level of effort, you might reply with a 6 or 7 (very, very light). On the other hand, if you feel you are working at a very high level of effort and cannot continue, you might reply with a 20 (maximal exertion). There is no right or wrong answer, simply rate your level of effort at the time I ask the question.”

Spanish Version

Introduction to the component: “En este examen queremos obtener alguna información acerca de su condición cardiovascular o la condición del corazón. Le pediremos que camine sobre la máquina de andar mientras le monitoreamos la frecuencia cardíaca (los latidos del corazón) y la presión de sangre. Explicaré esto con más detalle en unos pocos minutos pero primero quiero hacerle algunas preguntas acerca de su salud en general. Por favor venga aquí cerca de la computadora mientras anoto con el escáner su número de identificación en el sistema.”

Explanation of the test: “El propósito de la prueba en la máquina de andar es calcular su nivel de capacidad cardiovascular. Usted caminará en la máquina a diferentes velocidades y grados. Esta no es una prueba máxima de esfuerzo. Pararemos cuando la frecuencia cardíaca (los latidos del corazón) alcance cierto nivel. Si bien esta es una prueba muy segura, es importante que nos diga cómo se siente durante la prueba. Por favor informe cualquier síntoma inusual que tenga durante la prueba tal como dolor en el pecho, mareo, náusea, o dolor en las piernas.”

Explanation of Rating of Perceived Exertion (RPE): “Por favor mire el gráfico en el muro frente a usted. Cada dos minutos durante la prueba, le pediré que califique qué tan fuerte siente que está trabajando en general. Por favor conteste con un número de la escala. Por ejemplo, si usted siente que está trabajando a un nivel de esfuerzo muy bajo, usted puede responder un 6 ó 7 (muy, muy liviano). Por otra parte, si usted siente que está trabajando a un nivel de esfuerzo muy alto y no puede continuar, usted puede contestar con un 20 (máximo esfuerzo). No hay respuesta correcta o errónea, simplemente califique su nivel de esfuerzo en el momento en que yo le haga la pregunta.”

APPENDIX B
EXCLUSIONS BASED ON MEDICAL CONDITIONS

APPENDIX B. EXCLUSIONS BASED ON MEDICAL CONDITIONS

From Responses to Household Questionnaire

Medical Conditions and Health Status (MCQ)

- | | | |
|------------|-------------------------|--------------------------|
| ■ MCQ.160b | (if 1, 7, or 9 exclude) | Congestive heart failure |
| ■ MCQ.160c | (if 1, 7, or 9 exclude) | Coronary heart disease |
| ■ MCQ.160d | (if 1, 7, or 9 exclude) | Angina pectoris |
| ■ MCQ.160e | (if 1, 7, or 9 exclude) | Myocardial infarction |
| ■ MCQ.160f | (if 1, 7, or 9 exclude) | Stroke |
| ■ MCQ.160g | (if 1, 7, or 9 exclude) | Emphysema |

Physical Functioning - (PFQ) -

- | | | |
|---------------|-----------------------------|--|
| ■ PFQ.060b | (if 3, 4, 7, or 9, exclude) | Difficulty walking for a quarter mile (2-3 blocks) |
| ■ PFQ.060c | (if 3, 4, 7, or 9, exclude) | Difficulty walking up 10 steps without resting |
| ■ PFQ.060h | (if 3, 4, 7, or 9, exclude) | Walking from one room to another on the same level |
| ■ PFQ.060i | (if 3, 4, 7, or 9, exclude) | Standing up from an armless straight chair |
| ■ PFQ.067 =11 | (if 2,3, or 4, exclude) | Back or Neck Problem |
| =15 | (if 2,3, or 4, exclude) | Developmental Problems (Cerebral Palsy) |
| =17 | (if 2,3, or 4, exclude) | Fractures, Bone/Joint Injury |
| =19 | (if 2,3, or 4, exclude) | Heart Problem |
| =21 | (if 2,3, or 4, exclude) | Lung/Breathing Problem |
| =25 | (if 2,3, or 4, exclude) | Stroke Problem |
| =77 | (if 2,3, or 4, exclude) | Refused |
| =99 | (if 2,3, or 4, exclude) | Don't Know |
| ■ PFQ.090 | (if 1,7, or 9, exclude) | Use of a device such as a cane or wheelchair |

Diabetes (DIQ)

- | | | |
|-----------|-------------------------|-------------|
| ■ DIQ.080 | (if 1,7, or 9, exclude) | Retinopathy |
|-----------|-------------------------|-------------|

Cardiovascular (CDQ)-

- | | | |
|-----------|--------------------------|---|
| ■ CDQ.030 | (if 1, 7, or 9, exclude) | Stop when walking at own pace on the level |
| ■ CDQ.040 | (if 1, 7, or 9, exclude) | SOB after walking 100 yards or few minutes on the level |
| ■ CDQ.050 | (if 1, 7, or 9, exclude) | PND |
| ■ CDQ.060 | (if 1, 7, or 9, exclude) | PND relieved by sitting on side of bed |
| ■ CDQ.070 | (if 1, 7, or 9, exclude) | Orthopnea |

Respiratory Health (RDQ) –

- | | | |
|-----------|--|--|
| ■ RDQ.080 | (if 12 or more attacks, 77, or 99 exclude) | Wheezing in past 12 months |
| ■ RDQ.110 | (if 1,7,or 9, exclude) | Wheezing that limits speech (last 12 months) |

Vision (VIQ) -

- | | | |
|-----------|-------------------------|--------------------|
| ■ VIQ.020 | (if 1,7, or 9, exclude) | Blind |
| ■ VIQ.030 | (if 5,7, or 9, exclude) | Very poor eyesight |

DSQ.240

Medications will be available in SP History in Physician's Exam. Physician will check this list and exclude based on medication on the Exclusion list.

See Exclusionary Medication List:

Interviewer will have participant get medication bottles and will enter the complete product name of the medication in CAPI. This list of medications will be compared to the medications on the exclusion list - see Appendix C. Participants taking a medication that is on the exclusion medication list will be excluded from the CV Fitness component.

APPENDIX C
EXCLUSIONS BASED ON MEDICATIONS

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS

SP medications are reviewed by the Physician and the SP is excluded from CV Fitness if he/she is taking any of the following medications.

Updated 06-01-2000

Calcium Channel-Blockers

Amlodipine (Norvasc)
Bepridil (Vascor)
Diltiazem (Cardizem, Dilacor, Tiazac)
Felodipine (Plendil)
Isradipine (Dyna Circ)
Norvasc (Amlodipine)
Nicardipine (Cardene)
Nifedipine (Procardia, Adalat)
Nimodipine (Nimotop)
Nisoldipine (Sular)
Posicor (Mibefradil)
Tiazak (Diltiazem HCl)
Verapamil (Covera, Verelan, Calan, Isoptin)

Anti Arrhythmics

Amiodarone (Cordarone)
Bretylium (Bretylol)
Disopyramide (Norpace)
Encainide (Enkaid)
Ethmozine (Moricizine)
Flecainide (Tambocor)
Lidocaine (Xylocaine, Xylocard)
Metoprolol Succinate (Toprol-XL)
Mexiletine (Mexitil)
Moricizine (Ethmozine)
Procainamide (Pronestyl, Procan SR)
Propafenone (Rhythmol)
Sotalol (Betapace)
Tocainide (Tonocard)
Quinidine (Quinora, Quinalan, Cardioquin, Quinidex, Quinaglute)

Beta Blockers

Acebutolol (Sectral)
Atenolol (Tenormin)
Betaxolol (Kerlone)
Bisoprolol (Zebeta)
Carteolol (Cartrol)
Corzide *
Esmolol (Brevibloc)
Inderide *
Labetalol (Normodyne)
Lopressor Hydrochlorothiazide *

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Metoprolol tartrate (Lopressor)
Nadolol (Corgard)
Penbutolol (Levatol)
Pindolol (Visken)
Propranolol (Inderal)
Sotolol (Betapace)
Tenorectic *
Timolide *
Timolol (Blocadren)
Ziac *

* Beta blocker/diuretic combinations

Eye Drops

Betoptic Eyedrops
Timoptic Eyedrops

Nitrates and Nitroglycerin

Isosorbide dinitrate (Isordil, Diltrate)
Isosorbide mononitrate (Ismo, Monoket)
Nitroglycerin, translingual (Nitrostat, Nitrolingual spray)
Nitroglycerin, transmucosal (Nitrogard)
Nitroglycerin, topical (Nitrol, Nitro-Bid, Transderm Nitro, Nitro-Dur II, Nitrodisc, Minitran, Deponit, Nitroderm)
Nitroglycerin, sustained release (Nitrong, Nitrocine, Nitroglyn)
Pentaerythritol tetranitrate (Cardilate)

Digitalis

Digoxin (Lanoxin)

CNS Stimulant

Ma Huang (ephedrine)

APPENDIX D
CARDIOVASCULAR SAFETY AND
EXCLUSION QUESTIONS

APPENDIX D. CARDIOVASCULAR SAFETY AND EXCLUSION QUESTIONS

Exclusions based on questions: These questions will be asked as part of shared questions in the MEC exam or during the physician's examination that always occurs before the fitness procedure.

Shared questions asked in MEC Exam (relevant to CV Fitness):

1. Do you have any amputations of your legs and feet other than toes?
Yes/No (If Yes, exclude from CV Fitness)
2. How much do you weigh without shoes?
_____ pounds (If > 350, exclude from CV Fitness)
3. Do you have a pacemaker or automatic defibrillator?
Yes/No/Don't Know (If Yes or Don't Know, exclude from CV Fitness)
4. Are you currently pregnant? (Self-Report question for females 12-49 years)
Yes/ No (If Yes, exclude from Body Composition, ask followup question)
- 4b. How many weeks? _____ weeks (If >12 weeks, exclude from CV Fitness)

CV Safety/Exclusion questions asked in the Physician's Exam:

5. Have you been hospitalized in the past 3 months?
Yes/No/Don't Know (If Yes, and condition on exclusion list or Don't Know, exclude)
6. (12-19 years only) Has a doctor ever said you should not participate in sports or other activities because of a health condition?
Yes/No/Don't Know (If Yes or Don't Know, exclude)
7. Has a doctor ever said you have a heart condition and that you should only do physical activity recommended by a doctor?
Yes/No/Don't Know (If Yes or Don't Know, exclude)
8. (20-49 years only) Do you feel pain in your chest when you do physical activity?
Yes/No/Don't Know (If Yes or Don't Know, exclude)
9. (20-49 years only) In the past month, have you had chest pain when you were not doing physical activity? (*Probe: Have you seen a medical doctor about your chest pain? Did the doctor tell you that the chest pain was related to your heart?*)
Yes/No/Don't Know (If Yes or Don't Know, exclude)
10. Do you lose your balance because of dizziness? (*Probe: Is this an isolated incident or does it occur on a regular basis?*)
Yes/No/Don't Know (If Yes or Don't Know, exclude)
11. Do you ever lose consciousness? (*Probes: Did this occur as a result of illness or was it unexplained? Is this an isolated incident or does it occur on a regular basis?*)
Yes/No/Don't Know (If Yes or Don't Know, exclude)
12. Do you have a bone or joint problem that could be made worse by walking? (*Probe: Do you think you can do the test without injuring yourself?*)
Yes/No/Don't Know (If Yes or Don't Know, exclude)
13. Are you currently taking any prescription medications?
Yes/No/Don't Know (If Yes, go to 13b. If No, go to 14. If Don't Know, exclude.)
- 13b. Are you currently taking any prescription medications for your blood pressure?
Yes/No/Don't Know (If Yes, go to 13c. If No go to 13f. If Don't Know, exclude.)
- 13c. What is the name of this medication? _____.
(If medication is on list, exclude. If not on list, go to 13d.)

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

- 13d. Are you taking any other medication for your blood pressure?
Yes/No/Don't Know (Yes, go to 13e/ No go to 13f./Don't Know, exclude.)
- 13e. What is the name of this medication? _____.
(If medication is on list, exclude. If not on list, go to 13f.)
- 13f. Are you currently taking prescription medications for the following conditions:
- heart condition Yes/No/Don't Know (If Yes or Don't Know, exclude)
 - prescription eye drops for glaucoma? (If Yes or Don't Know, exclude)
14. Do you know of any other reason why you should not do a treadmill test?
Yes/No/Don't Know (If Yes or Don't Know, exclude)

See Safety/exclusion Question 5: List of reasons for exclusion based on hospitalization from *ACSM Guidelines, 5th edition, page 42.*

- A recent significant change in the resting ECG suggesting infarction or other acute cardiac event
- Recent complicated myocardial infarction
- Unstable angina
- Uncontrolled ventricular arrhythmia
- Uncontrolled atrial arrhythmia that compromises cardiac function
- Third degree AV heart block
- Acute congestive heart failure
- Severe aortic stenosis
- Suspected or known dissecting aneurysm
- Active or suspected myocarditis or pericarditis
- Thrombophlebitis or intracardiac thrombi
- Recent systemic or pulmonary embolus
- Significant emotional distress (psychosis)
- Moderate valvular heart disease
- Known electrolyte abnormalities
- Fixed rate pacemaker
- Frequent or complex ventricular ectopy
- Ventricular aneurysm
- Uncontrolled metabolic disease (diabetes, thyrotoxicosis, myxedema, etc)
- Chronic infections disease (mononucleosis, hepatitis, AIDS)
- Neuromuscular, musculoskeletal, or rheumatoid disorders that are exacerbated by exercise
- Complicated or late stage pregnancy

Exclusions based on Data:

Resting BP and HR taken in Physician's Exam.

These measurements are used for exclusion from CV Fitness (not Pretest measurements taken in the CV Fitness Room)

Excluded based on any of the following criteria:

SBP > 180 mmHg

DBP > 100 mmHg

HR > 100 beats per minute

Irregular Beats > 3 beats per minute

APPENDIX E
DETERMINATION OF PHYSICAL ACTIVITY
READINESS CODE

APPENDIX E. DETERMINATION OF PHYSICAL ACTIVITY READINESS CODE

English Version Questions

Physical Activity Questions

I am going to ask you several questions to describe your typical physical activity level. Please try to answer as accurately as possible. Please wait until I read all three questions before you give me your answer.

1. Would you say that you did **little or no regular** recreation, sport, or physical activity?
OR
2. Would you say that you participated **regularly** in recreation or work requiring **modest physical activity**?
OR
3. Would you say that you participated **regularly** in **heavy physical exercise**?

Modest physical activities include: walking, calisthenics, bowling, weight lifting, yard work.

Heavy physical activities include jogging, swimming, cycling, rowing, tennis, and basketball.

1. Would you say that you avoid walking or exertion (for example, always use the elevator, drive whenever possible instead of walking.)?
OR
2. Would you say that you walk for pleasure, routinely use stairs, occasionally exercise sufficiently to cause heavy breathing or perspiration?

Would you say that you participated regularly in recreation or modest physical activity for:

1. **10 to 60 minutes per week**
OR
2. **More than 60 minutes per week?**

Now I am going to ask for some more detail on the amount of heavy physical exercise you did. I am going to give you four options.

Would you say that you participated regularly in heavy physical activity for

1. **less than 30 minutes per week**
OR
2. **30 to 60 minutes per week**
OR
3. **1 to 3 hours per week**
OR
4. **more than 3 hours per week**

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Spanish Version Questions

Le voy a hacer varias preguntas para describir su **nivel típico de actividad física**. Por favor trate de responder tan precisamente como sea posible. Por favor espere hasta que yo lea las tres preguntas antes de darme su respuesta.

4. Diría que usted hizo **poca o ninguna** actividad **regular** de recreación, deporte, o física?
☐
5. ¿Diría que usted participó **regularmente** en recreación o trabajo que requería **actividad física moderada**?
☐
6. ¿Diría que usted participó **regularmente** en **ejercicios físicos pesados**?
☐

Actividades físicas moderadas incluye: caminar, calistenia, bolos, levantamiento de pesas, trabajo en el (patio/jardín).

Actividades físicas pesadas incluye: trotar, nadar, andar en bicicleta, remar, jugar tenis, y baloncesto.

2. ¿Diría que usted **evita** caminar o hacer esfuerzos (por ejemplo, siempre toma el ascensor, conduce en lugar de caminar cada vez que es posible.)?
☐
2. ¿Diría que usted camina por placer, rutinariamente usa las escaleras, ocasionalmente hace ejercicio suficientemente como para que le cause respiración agitada o transpiración?
☐

¿Diría que usted participó regularmente en actividades de recreación o físicas moderadas por:

1. **10 a 60 minutos por semana**?
☐
2. **Más de 60 minutos por semana**?
☐

Ahora le voy a preguntar algunos detalles más sobre la cantidad de ejercicio físico pesado que usted hizo. Le voy a dar cuatro opciones.

¿Diría que usted participó regularmente en actividades físicas pesadas ...

1. **menos de 30 minutos por semana**,
☐
2. **30 a 60 minutos por semana**,
☐
3. **1 a 3 horas por semana**,
☐
4. **más de 3 horas por semana**?
☐

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Determination of PAR code from responses

The PAR code is one of the determinants of which treadmill protocol will be assigned to the SP. PAR codes range from 0 to 7. The code is assigned as follows:

If the SP reports little or no regular recreation, sport, or physical activity and avoids walking or exertion, PAR code = 0

If the SP reports little or no regular recreation, sport, or physical activity and walks for pleasure, routinely uses stairs, occasionally exercises sufficiently to cause heavy breathing or perspiration, PAR code = 1

If the SP reports participating regularly in recreation or work requiring modest physical activity for 10-60 minutes per week, PAR code = 2

If the SP reports participating regularly in recreation or work requiring modest physical activity for more than 60 minutes per week, PAR code = 3

If the SP reports participating regularly in heavy physical exercise for less than 30 minutes per week, PAR code = 4

If the SP reports participating regularly in heavy physical exercise for 30-60 minutes per week, PAR code = 5

If the SP reports participating regularly in heavy physical exercise for 1-3 hours per week, PAR code = 6

If the SP reports participating regularly in heavy physical exercise for more than 3 hours per week, PAR code = 7

Sources:

Jackson, Blair, Mahar, Wier, Ross, and Stuteville. (1990). Prediction of functional aerobic capacity without exercise testing. *Medicine and Science in Sports and Exercise*, 22 (6), 863-870.

Ross and Jackson. (1990). *Exercise Concepts, Calculations and Computer Applications*. Benchmark Press.

APPENDIX F
CALCULATION OF PREDICTED VO_2 MAX

APPENDIX F. CALCULATION OF PREDICTED VO₂ MAX

Predicted VO₂ max = $56.363 + [1.921 * (\text{PA-R})] - [0.381 * (\text{Age at interview})] - [0.754 * (\text{BMI})] + [10.987 * (\text{F}=0, \text{M}=1)]$

PAR Code (as determined from Appendix E)	Code 0-7
Use Age at Interview	years
Measured body weight	Kg.
Measured body height	cm.
BMI (body mass index) = body weight (kg)/height squared (meters)	(kg/m ²)

Use Table in Appendix H to determine protocol based on predicted VO₂ max calculated above.

APPENDIX G
CALCULATION OF ESTIMATED VO_2 MAX

APPENDIX G. CALCULATION OF ESTIMATED VO₂ MAX

Because the relation between heart rate and oxygen consumption is linear during exercise, we can estimate maximal oxygen consumption (VO₂ max) by measuring the heart rate response to known levels of submaximal work. This relation can be summed up by the equation for a line: $y = mx + b$, where y is heart rate, m is slope, x is VO₂, and b is the intercept. By rearranging these terms, we can derive the following equation:

$$\text{Estimated VO}_2\text{max} = \frac{\text{PMHR} - \text{Intercept}}{\text{Slope}}$$

$$\begin{aligned} \text{Slope} &= n(\sum x_i y_i) - (\sum x_i)(\sum y_i) / n \sum x_i^2 - (\sum x_i)^2 \\ \text{PMHR} &= 220 - \text{Age at interview} \\ \text{Intercept} &= \bar{Y} - b \bar{X} \end{aligned}$$

\bar{Y} = mean of End of Stage 1 Heart Rate and End of Stage 2 Heart Rate

\bar{X} = mean of End of Stage 1 submax VO₂ and End of Stage 2 submax VO₂

b = slope = $n(\sum x_i y_i) - (\sum x_i)(\sum y_i) / n \sum x_i^2 - (\sum x_i)^2$

$$\text{Slope} = n(\sum x_i y_i) - (\sum x_i)(\sum y_i) / n \sum x_i^2 - (\sum x_i)^2$$

$$n = 2$$

$$x_i = x_1, x_2$$

$$x_1 = \text{submax VO}_2 \text{ at end of Stage 1}$$

$$x_2 = \text{submax VO}_2 \text{ at end of Stage 2}$$

$$y_i = y_1, y_2$$

$$y_1 = \text{end of Stage 1 HR}$$

$$y_2 = \text{end of Stage 2 HR}$$

Get end of stage VO₂ from Appendix H.

Example: A 40 year old man was assigned Protocol 6. Looking in Appendix H, we see that the corresponding submax VO₂ values for these stages are 23.3 and 31.4 ml/kg/min respectively. His heart rate at the end of Stage 1 is 120 and at the end of Stage 2 is 140 beats/min respectively. Plugging these numbers into the formulas above, we then solve our equation.

Can you double check the formula – I plot the number into the formula but didn't get the same slope as you put here....

$$\text{Estimated VO}_2 \text{ max} = \frac{180 - 62.7}{2.46} = 47.7 \text{ ml/kg/min}$$

APPENDIX H

DETERMINATION OF EXERCISE PROTOCOL

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

APPENDIX H. DETERMINATION OF EXERCISE PROTOCOL

	WARMUP				STAGE 1			STAGE 2			RECOVERY		
Exercise Protocol Number	Predicted VO ₂ max (ml/kg/min)	Speed (mph)	Grade (%)	% predicted VO ₂ Max	Speed (mph)	Grade (%)	% predicted VO ₂ Max (ml/kg/min)	Speed (mph)	Grade (%)	% predicted VO ₂ Max (ml/kg/min)	Speed (mph)	Grade (%)	% predicted VO ₂ Max (ml/kg/min)
1	<20	1.7	0	45(8.1)	2.1	0.5	55(9.6)	2.1	4.5	76(13.6)	2.0	0	49(8.9)
2	20-24	2.0	1	45(10)	2.3	2.0	55(11.9)	2.3	6.5	75(16.9)	2.0	0	41(8.9)
3	25-29	2.2	2	45(11.5)	2.7	3.0	55(14.6)	2.7	7.5	75(20.6)	2.0	0	33(8.9)
4	30-34	2.7	2	45(13.3)	3.1	4.0	55(17.8)	3.1	8.5	75(24.4)	2.0	0	28(8.9)
5	35-39	3.0	3	45(16)	3.7	4.0	55(20.6)	3.7	8.0	75(27.8)	2.0	0	24(8.9)
6	40-44	3.1	3.5	45(17)	3.7	5.5	55 (23.3)	3.7	10	75(31.4)	2.0	0	21(8.9)
7	45-49	3.2	5	45(19.8)	3.7	7.0	55(26.0)	3.7	12.5	75(35.9)	2.0	0	19(8.9)
8	>50	3.6	5	45(21.8)	3.7	8.5	55(28.7)	3.7	14.5	75(39.5)	2.0	0	17(8.9)

Protocol Number

Speed

Grade

% predicted VO2 max

number to identify protocol being used (see table)

miles per hour

percent incline

predicted VO₂ max as % (submaximal VO₂ in ml/kg/min in parenthesis)

APPENDIX I
REFERENCE RANGES FOR VO₂ MAX

APPENDIX I. REFERENCE RANGES FOR VO₂ MAX

Males ages 12-19

Low	≤ 42.13
Moderate	42.14-53.18
High	≥ 53.19

Females ages 12-19

Low	≤ 35.63
Moderate	35.64 - 45.96
High	≥ 45.97

Males ages 20-29

Low	≤ 37.13
Moderate	37.14-44.22
High	≥ 44.23

Females ages 20-29

Low	≤ 30.63
Moderate	30.64-36.64
High	≥ 36.64

Males ages 30-39

Low	≤ 35.35
Moderate	35.36-42.41
High	≥ 42.42

Females ages 30-39

Low	≤ 28.70
Moderate	28.71-34.59
High	≥ 34.60

Males ages 40-49

Low	≤ 33.04
Moderate	33.05-39.88
High	≥ 39.89

Females ages 40-49

Low	≤ 26.54
Moderate	26.55-32.30
High	≥ 32.31

These norms are based on studies completed at the Cooper Institute for Aerobic Research and published in the fifth edition of the American College of Sports Medicine (ACSM) Guidelines for Exercise Testing and Prescription. There are no national norms based on random samples of the population (until NHANES is completed). The CIAR norms are well-established norms used by military, law enforcement, and ASCM.

Report of Findings statement:

Compared with other people of your age and sex, your cardiovascular fitness level is <low, moderate, high>.

APPENDIX J

SPEED & GRADE TABLE BY PROTOCOL & STAGE

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

APPENDIX J. SPEED AND GRADE TABLE BY PROTOCOL AND STAGE

Protocol #	Stage	Current/ Lower/ Higher	Grade (%)	Speed (mph)	Protocol #	Stage	Current/ Lower/ Higher	Grade (%)	Speed (mph)
1	Warm-up	Current (1)	0%	1.7	2	Warm-up	Current (2)	1.0%	2.0
1	Stage 1	Current (1)	0.5%	2.1	2	Stage 1	Current (2)	2.0%	2.3
1	Stage 2	Current (1)	4.5%	2.1	2	Stage 2	Current (2)	6.5%	2.3
1	Recovery	Current (1)	0%	2.0	2	Recovery	Current (2)	0%	2.0
1	Stage 1	Lower (0)	0%	0.8	2	Stage 1	Lower (1)	0.5%	2.1
1	Stage 2	Lower (0)	0%	0.8	2	Stage 2	Lower (1)	4.5%	2.1
1	Recovery	Lower (0)	0%	0.8	2	Recovery	Lower (1)	0%	2.0
1	Stage 1	Higher (2)	2.0%	2.3	2	Stage 1	Higher (3)	3.0%	2.7
1	Stage 2	Higher (2)	6.5%	2.3	2	Stage 2	Higher (3)	7.5%	2.7
1	Recovery	Higher (2)	0%	2.0	2	Recovery	Higher (3)	0%	2.0
1	Stage 1	Lower (0)	0%	0.8	2	Stage 1	Lower (1)	0.5%	2.1
1	Stage 2	Lower (0)	0%	0.8	2	Stage 2	Lower (1)	4.5%	2.1
1	Recovery	Lower (0)	0%	0.8	2	Recovery	Lower (1)	0%	2.0
1	Stage 1	Higher (2)	2.0%	2.3	2	Stage 1	Higher (3)	3.0%	2.7
1	Stage 2	Higher (2)	6.5%	2.3	2	Stage 2	Higher (3)	7.5%	2.7
1	Recovery	Higher (2)	0%	2.0	2	Recovery	Higher (3)	0%	2.0
3	Warm-up	Current (3)	2.0%	2.2	4	Warm-up	Current (4)	2.0%	2.7
3	Stage 1	Current (3)	3.0%	2.7	4	Stage 1	Current (4)	4.0%	3.1
3	Stage 2	Current (3)	7.5%	2.7	4	Stage 2	Current (4)	8.5%	3.1
3	Recovery	Current (3)	0%	2.0	4	Recovery	Current (4)	0%	2.0
3	Stage 1	Lower (2)	2.0%	2.3	4	Stage 1	Lower (3)	3.0%	2.7
3	Stage 2	Lower (2)	6.5%	2.3	4	Stage 2	Lower (3)	7.5%	2.7
3	Recovery	Lower (2)	0%	2.0	4	Recovery	Lower (3)	0%	2.0
3	Stage 1	Higher (4)	4.0%	3.1	4	Stage 1	Higher (5)	4.0%	3.7
3	Stage 2	Higher (4)	8.5%	3.1	4	Stage 2	Higher (5)	8.0%	3.7
3	Recovery	Higher (4)	0%	2.0	4	Recovery	Higher (5)	0%	2.0
3	Stage 1	Lower (2)	2.0%	2.3	4	Stage 1	Lower (3)	3.0%	2.7
3	Stage 2	Lower (2)	6.5%	2.3	4	Stage 2	Lower (3)	7.5%	2.7
3	Recovery	Lower (2)	0%	2.0	4	Recovery	Lower (3)	0%	2.0

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS
(continued)

3	Stage 1	Higher (4)	4.0%	3.1		4	Stage 1	Higher (5)	4.0%	3.7
3	Stage 2	Higher (4)	8.5%	3.1		4	Stage 2	Higher (5)	8.0%	3.7
3	Recovery	Higher (4)	0%	2.0		4	Recovery	Higher (5)	0%	2.0

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

APPENDIX J. SPEED AND GRADE TABLE BY PROTOCOL AND STAGE (continued)

Protocol #	Stage	Current/ Lower/ Higher	Grade (%)	Speed (mph)	Protocol #	Stage	Current/ Lower/ Higher	Grade (%)	Speed (mph)
5	Warm-up	Current (5)	3.0%	3.0	6	Warm-up	Current (6)	3.5%	3.1
5	Stage 1	Current (5)	4.0%	3.7	6	Stage 1	Current (6)	5.5%	3.7
5	Stage 2	Current (5)	8.0%	3.7	6	Stage 2	Current (6)	10.0%	3.7
5	Recovery	Current (5)	0%	2.0	6	Recovery	Current (6)	0%	2.0
5	Stage 1	Lower (4)	4.0%	3.1	6	Stage 1	Lower (5)	4.0%	3.7
5	Stage 2	Lower (4)	8.5%	3.1	6	Stage 2	Lower (5)	8.0%	3.7
5	Recovery	Lower (4)	0%	2.0	6	Recovery	Lower (5)	0%	2.0
5	Stage 1	Higher (6)	5.5%	3.7	6	Stage 1	Higher (7)	7.0%	3.7
5	Stage 2	Higher (6)	10.0%	3.7	6	Stage 2	Higher (7)	12.5%	3.7
5	Recovery	Higher (6)	0%	2.0	6	Recovery	Higher (7)	0%	2.0
5	Stage 1	Lower (4)	4.0%	3.1	6	Stage 1	Lower (5)	4.0%	3.7
5	Stage 2	Lower (4)	8.5%	3.1	6	Stage 2	Lower (5)	8.0%	3.7
5	Recovery	Lower (4)	0%	2.0	6	Recovery	Lower (5)	0%	2.0
5	Stage 1	Higher (6)	5.5%	3.7	6	Stage 1	Higher (7)	7.0%	3.7
5	Stage 2	Higher (6)	10.0%	3.7	6	Stage 2	Higher (7)	12.5%	3.7
5	Recovery	Higher (6)	0%	2.0	6	Recovery	Higher (7)	0%	2.0
7	Warm-up	Current (7)	5.0%	3.2	8	Warm-up	Current (8)	5.0%	3.6
7	Stage 1	Current (7)	7.0%	3.7	8	Stage 1	Current (8)	8.5%	3.7
7	Stage 2	Current (7)	12.5%	3.7	8	Stage 2	Current (8)	14.5%	3.7
7	Recovery	Current (7)	0%	2.0	8	Recovery	Current (8)	0%	2.0
7	Stage 1	Lower (6)	5.5%	3.7	8	Stage 1	Lower (7)	7.0%	3.7
7	Stage 2	Lower (6)	10.0%	3.7	8	Stage 2	Lower (7)	12.5%	3.7
7	Recovery	Lower (6)	0%	2.0	8	Recovery	Lower (7)	0%	2.0
7	Stage 1	Higher (8)	8.5%	3.7	8	Stage 1	Higher (9)	0%	0.8
7	Stage 2	Higher (8)	14.5%	3.7	8	Stage 2	Higher (9)	0%	0.8
7	Recovery	Higher (8)	0%	2.0	8	Recovery	Higher (9)	0%	0.8
7	Stage 1	Lower (6)	5.5%	3.7	8	Stage 1	Lower (7)	7.0%	3.7
7	Stage 2	Lower (6)	10.0%	3.7	8	Stage 2	Lower (7)	12.5%	3.7
7	Recovery	Lower (6)	0%	2.0	8	Recovery	Lower (7)	0%	2.0

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS
(continued)

7	Stage 1	Higher (8)	8.5%	3.7		8	Stage 1	Higher (9)	0%	0.8
7	Stage 2	Higher (8)	14.5%	3.7		8	Stage 2	Higher (9)	0%	0.8
7	Recovery	Higher (8)	0%	2.0		8	Recovery	Higher (9)	0%	0.8

**APPENDIX K
REVIEW OF EDIT CRITERIA**

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

APPENDIX K. REVIEW OF EDIT CRITERIA

Entire Exam (these criteria apply throughout the entire exam):

- | | | |
|---|--|---------------------------------------|
| ■ %PMHR greater than 100% | ■ G
o

t
o

r
e
c
o
v
e
r
y | ■ (no
exceptions
) |
| ■ %PMHR greater than 85%
(19 through 49 years) | ■ G
o

t
o

r
e
c
o
v
e
r
y | ■ (exception
last min.
Stage 2) |
| ■ %PMHR greater than 90%
(12 through 18 years) | ■ G
o

t
o

r
e
c
o
v
e | ■ (exception
last min.
Stage 2) |

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

	r y	
▪ SBP greater than 260 mmHg	▪ G o t o r e c o v e r y	▪ (no exceptions)
▪ DBP greater than 115 mm Hg	▪ G o t o r e c o v e r y	▪ (no exceptions)
▪ RPE greater than 17	▪ G o t o r e c o v e r y	▪ (exception at end of Stage 2)

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Warmup

- | | |
|---|------------------------------|
| ■ %PMHR less than 50% | ■ Increase exercise protocol |
| ■ %PMHR greater than 60% | ■ Decrease exercise protocol |
| ■ Criteria for 'Entire Exam' also apply | ■ Go to recovery |

Stage 1

- | | |
|---|------------------|
| ■ %PMHR greater than 80% | ■ Go to recovery |
| ■ Criteria for 'Entire Exam' also apply | ■ Go to recovery |
| ■ RPE >17 | ■ Go to recovery |

Stage 2

- | | |
|---|--|
| ■ %PMHR greater than 85%
(19 through 49 years) | ■ Go to recovery (exception
last minute) |
| ■ %PMHR greater than 90%
(12 through 18 years) | ■ Go to recovery (exception
last minute) |
| ■ Criteria for 'Entire Exam' also apply | ■ Go to recovery |
| ■ RPE >17 | ■ Go to recovery (exception
last 15-20 sec) |

RPE Exception:

- | | |
|-----------------------|--|
| ■ Warm-up and Stage 1 | ■ No Exception, if RPE is greater than 17 in Warm-up or Stage 1, go to Recovery <u>BUT</u> confirm that SP understands RPE. If RPE is still greater than 17, go to recovery. |
| ■ Stage 2 | ■ EXCEPTION, if RPE >17 in Stage 2, ask SP if he/she can continue for another 15-20 seconds until the exam is completed. |

%PMHR Exception in last minute of Stage 2:

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

- % PMHR 85% – 100% for 12-18 years
- If SP looks OK, says she/he is OK, and agrees to continue, complete the treadmill test.
- %PMHR 90% - 100% for 19-49 years
- If SP looks OK, says she/he is OK, and agrees to continue, complete the treadmill test.

APPENDIX L
SAMPLE DATA ENTRY SCREENS

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

APPENDIX L. SAMPLE DATA ENTRY SCREENS

Cardiovascular Fitness: Stand:107 Session:107012 07/06/1999 05:30 pm - 09:30 pm

SP ID: 818591 Name: BROWN, MAC Age: 38 years Gender: Male Date: 10/10/2000 Time: 12:48 PM

SP Height and Weight Information

Measured Weight: Kgs

Measured Height: Cms

1 of 10

End of Section Close Exam Finish

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

File View Utilities Reports Window Help

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:04 AM

Questions To Determine Exercise Protocol

I am going to ask you several questions to describe your TYPICAL PHYSICAL ACTIVITY LEVEL. Please try to answer as accurately as possible. Please wait until I read all three questions before you give me your answer.

1. Would you say that you did LITTLE or NO REGULAR recreation, sport, or physical activity?
OR
2. Would you say that you participated REGULARLY in recreation OR WORK requiring MODEST physical activity?
OR
3. Would you say that you participated REGULARLY in HEAVY physical activity?

If necessary use the following prompts to describe modest or heavy physical activity.
MODEST physical activities include: walking, calisthenics, bowling, weight lifting, yard work.
HEAVY physical activities include: jogging, swimming, cycling, rowing, tennis, basketball.

If SP's activity level over the last few months is different from their activity level in general, use the activity level which best describes their typical activity level over the last year.

If the SP is still not sure in which category they fall, (little or modest, modest or heavy), select the lower category.

Enter/Select answer:

1 of 8

End of Section Close Exam Finish

Ready

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

File View Utilities Reports Window Help

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:14 AM

Questions To Determine Exercise Protocol (Cont...)

1 Would you say that you AVOID walking or exertion, (for example, always use the elevator, drive whenever possible instead of walking.)?
OR
2 Would you say that you walk for pleasure, routinely use stairs, occasionally exercise sufficiently to cause heavy breathing or perspiration?

Would you say that you participated regularly in recreation or modest physical activity for:
1. 10 TO 60 MINUTES PER WEEK?
OR
2. MORE THAN 60 MINUTES PER WEEK?

Enter/Select answer: 1

1
2
Refused
Don't Know

2 of 8

End of Section Close Exam Finish

Ready

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:14 AM

Questions To Determine Exercise Protocol (Cont...)

1 Would you say that you AVOID walking or exertion, (for example, always use the elevator, drive whenever possible instead of walking.)?
OR
2 Would you say that you walk for pleasure, routinely use stairs, occasionally exercise sufficiently to cause heavy breathing or perspiration?

Would you say that you participated regularly in recreation or modest physical activity for:
1. 10 TO 60 MINUTES PER WEEK?
OR
2. MORE THAN 60 MINUTES PER WEEK?

Enter/Select answer:

2 of 8

End of Section Close Exam Finish

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm

SP ID: 666259 Name: CHRISTIAN, TIM Age: 31 years Gender: Male Date: 09/26/2000 Time: 03:55 PM

Questions To Determine Exercise Protocol (Cont...)

Now I am going to ask for some more detail on the amount of heavy physical exercise you did. I am going to give you four options.

Would you say that you participated regularly in heavy physical activity for:

1. LESS THAN 30 MINUTES PER WEEK
OR
2. 30 TO 60 MINUTES PER WEEK
OR
3. 1 TO 3 HOURS PER WEEK
OR
4. MORE THAN 3 HOURS PER WEEK

Enter/Select answer:

4 of 10

End of Section Close Exam Finish

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:33 AM

Pretest Measurements

Room Temperature: Degrees

Exam Preparation

Arm:

Cuff Size:

Pretest Measurements

Heart Rate:

Systolic BP:

Diastolic BP:

4 of 8

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm
_ □ ×

SP ID: 132135 **Name:** MCCULLOUGH, KRISTINA MARIE **Age:** 16 years **Gender:** Female **Date:** 08/04/2000 **Time:** 01:31 PM

Treadmill Test

00:00

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 204.00 beats/min

102 132.6 204

%HR max

50% 60% 100%

⏮ ⏪ 5 ⏩ ⏭

End of Section Close Exam Finish

➡

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm

SP ID: 132135 **Name:** MCCULLOUGH, KRISTINA MARIE **Age:** 16 years **Gender:** Female **Date:** 08/04/2000 **Time:** 01:31 PM

Treadmill Test

00:30

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 204.00 beats/min

102
132.6
204

%HR max
50%
60%
100%

⏮ ⏪ 5 ⏩ ⏭

End of Section Close Exam Finish

➡

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm

SP ID: 132135 **Name:** MCCULLOUGH, KRISTINA MARIE **Age:** 16 years **Gender:** Female **Date:** 08/04/2000 **Time:** 01:31 PM

Treadmill Test

00:30

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
Protocol	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 204.00 beats/min

122.4

163.2

204

%HR max
60%
80%
100%

⏮ ⏪ 5 ⏩ ⏭ of 9

End of Section
Close Exam
Finish

➡

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm

SP ID: 142069 **Name:** MC CARRON, MARY ANN **Age:** 34 years **Gender:** Female **Date:** 08/04/2000 **Time:** 01:53 PM

Treadmill Test

00:10

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
3	2.20	2.00	2.70	3.00	2.70	7.50	2.00	0
4	2.70	2.00	3.10	4.00	3.10	8.50	2.00	0
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0

RPE

Captured HR

Manual HR

PMHR: 186.00 beats/min

158.1 186

%HR max
85%
100%

◀

◀

5

▶

▶

of 9

End of Section

Close Exam

Finish

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm

SP ID: 132135 **Name:** MCCULLOUGH, KRISTINA MARIE **Age:** 16 years **Gender:** Female **Date:** 08/04/2000 **Time:** 01:31 PM

Treadmill Test

00:06

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 204.00 beats/min

183.6 204

%HR max
90%
100%

⏮ ⏪ 5 ⏩ ⏭ of 9

End of Section Close Exam Finish

➡

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm

SP ID: 132135 **Name:** MCCULLOUGH, KRISTINA MARIE **Age:** 16 years **Gender:** Female **Date:** 08/04/2000 **Time:** 01:31 PM

Treadmill Test

02:04

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: 00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0
7	3.20	5.00	3.70	7.00	3.70	12.5	2.00	0

RPE

Captured HR

Manual HR

PMHR: 204.00 beats/min

183.6 204

%HR max
90%
100%

⏮ ⏪ 5 ⏩ ⏭

of 9

End of Section Close Exam Finish

➡

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107030 07/08/1999 08:30 am - 12:30 pm

SP ID: 855853 **Name:** HALL, JEFF D. **Age:** 12 years **Gender:** Female **Date:** 08/21/2000 **Time:** 03:04 PM

Treadmill Test

00:19

min sec

START

RESTART

EMERGENCY STOP

Change Protocol

Go To Recovery

Show Data

Vo2 Max: .00

Protocol	WARM-UP		STAGE 1		STAGE 2		RECOVERY	
	Speed	Grade	Speed	Grade	Speed	Grade	Speed	Grade
4	2.70	2.00	3.10	4.00	3.10	8.50	2.00	0
5	3.00	3.00	3.70	4.00	3.70	8.00	2.00	0
6	3.10	3.50	3.70	5.50	3.70	10.0	2.00	0

RPE

Captured HR

Manual HR

PMHR: 208.00 beats/min

104

208

%HR max
50%
100%

◀◀

5

▶▶

of 9

End of Section

Close Exam

Finish

➡

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness Subsystem: Stand:806 Session:808559 12/10/1998 08:30 am - 12:30 pm

SP ID: 836308 Name: Anderson, Greg Age: 28 years Gender: Male Date: 01/03/1999 Time: 10:57 AM

Priority 1 Stopping Criteria






You have stopped the exam early. You must select one or more of the following stopping criteria:

- ☐ Pain or pressure in the chest
- ☐ Pain in neck, jaw, or teeth
- ☐ Pain in the shoulder or between the shoulder blades
- ☐ Pain radiating down one or both arms
- ☐ Moderate/severe leg pain or cramping
- ☐ Lightheadedness, dizziness
- ☐ Severe shortness of breath/wheezing
- ☐ Change in skin color to bluish/grayish
- ☐ Severe headache
- ☐ Subject appears confused or agitated
- ☐ Nausea
- ☐ Unstable gait
- ☐ Visual problems
- ☐ Other Priority 1

6 of 8

End of Section Close Exam Finish

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107341 08/08/1999 01:30 pm - 05:30 pm				
SP ID: 132135	Name: MCCULLOUGH, KRISTINA MARIE	Age: 16 years	Gender: Female	Date: 08/04/2000 Time: 01:31 PM
Priority 2 Stopping Criteria				
You have stopped the exam early. You must select one or more of the following stopping criteria:				
<input type="checkbox"/> SP HR during WARM-UP, STAGE 1, or first 2 minutes of STAGE 2 exceeds 90% PMHR for 12 through 18 yr olds or 85% PMHR for 19 through 49 yr olds				
<input type="checkbox"/> SP exceeds 80% of PMHR during STAGE 1.				
<input type="checkbox"/> Exercise blood pressure > 260 mmHg systolic and/or 115 mmHg diastolic				
<input type="checkbox"/> Failure of systolic blood pressure to rise with increasing exercise intensity				
<input type="checkbox"/> Significant drop (>20 mmHg) in systolic blood pressure during exercise				
<input type="checkbox"/> Failure of heart rate to increase with increasing exercise intensity				
<input type="checkbox"/> Sudden, noticeable variability in heart rate (sudden change of + 30 beats/minute during a STAGE)				
<input type="checkbox"/> SP requests to stop test, or reports severe fatigue				
<input type="checkbox"/> Failure of testing equipment				
<input type="checkbox"/> SP is unable to stop gripping hand rails (lightly touching hand rails for balance purposes is allowed, but only if absolutely necessary)				
<input type="checkbox"/> SP reports RPE > 17				
<input type="checkbox"/> Technician discretion (tech notices signs of severe exertion).				
<input type="checkbox"/> Other Enter Comment: <input type="text"/>				
<div><div>  <input type="text" value="7"/>  </div> of 9</div> <div><div>End of Section</div><div>Close Exam</div><div>Finish</div></div> <div></div>				

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness: Stand:107 Session:107401 08/14/1999 01:30 pm - 05:30 pm

SP ID: 531764 **Name:** ELLA, CINDER **Age:** 25 years **Gender:** Female **Date:** 07/04/2000 **Time:** 08:46 AM

Data acquisition problems

No Problem ☒

Problem obtaining BP - BP not displayed on the BP monitor, not captured in ISIS

☐ BP not obtained in Baseline ☐ BP not obtained in Warmup ☐ BP not obtained in Stage 1

☐ BP not obtained in Stage 2 ☐ BP not obtained in Recovery 1, 2, and /or 3

☐ Other

Problem with data capture - BP displayed correctly on STBP monitor but not correctly captured in ISIS

☐ BP on STBP OK, not captured in ISIS ☐ BP displayed normally on BP monitor, different data captured in ISIS

☐ Other

Problem obtaining HR

☐ Too much variability on HR on monitor to determine accurate HR

☐ Other

Problem with electrodes

☐ adherence problem due to excessive perspiration ☐ adherence problem due to excessive chest hair

☐ large chest/abdomen ☐ Other

Problem with BP cuff

☐ large arm caused difficulty with microphone placement ☐ cuff slipped off during test

☐ SP bending/flexing arm during test ☐ Other

Problem with treadmill controller

☐ Controller advanced more than one stage - Restarted test ☐ Other

◀ ◀ 8 ▶ ▶ of 9

End of Section Close Exam Finish

APPENDIX C. EXCLUSIONS BASED ON MEDICATIONS (continued)

Cardiovascular Fitness Subsystem: Stand:806 Session:808571 12/11/1998 05:30 pm - 09:30 pm

SP ID: 926247 Name: Ball, Jake Age: 37 years Gender: Male Date: 01/24/1999 Time: 06:51 PM

Cardiovascular Fitness Component Status

Status

☒ Complete
☐ Partial
☐ Not Done

Comments

Other text

8 of 8

End of Section Close Exam Finish