


How to Export Data from NHSN

Description

There are three ways to export your data from NHSN.

1. To export ALL data entered from your facility, click on the “Import/Export” option on the left navigation bar, then select “Export Facility Data” from the drop-down menu. You can then export the data into Excel or other available formats (SAS, .txt, etc.).

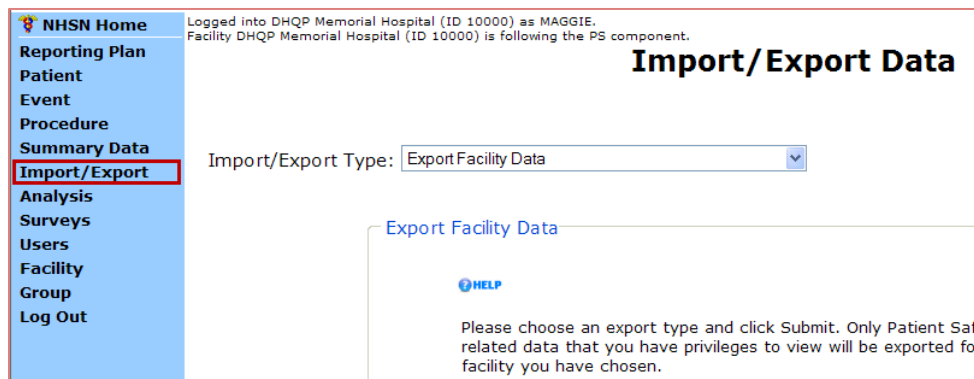

NHSN Home

Reporting Plan

Patient

Event

Procedure

Summary Data

Import/Export

Analysis

Surveys

Users

Facility

Group

Log Out

Logged into DHQP Memorial Hospital (ID 10000) as MAGGIE.
Facility DHQP Memorial Hospital (ID 10000) is following the PS component.

Import/Export Data

Import/Export Type: Export Facility Data

[Export Facility Data](#)


[HELP](#)

Please choose an export type and click Submit. Only Patient Safe related data that you have privileges to view will be exported for facility you have chosen.

Note: This option is only available to facility users with administrative rights. The export will include multiple data tables that mirror the structure in NHSN.

2. To export a specific analysis data set, you must first generate datasets by selecting Analysis > Generate Data Sets on the left navigation bar (*this may take a few minutes*). Once this is done, go to Analysis > Output Options where you can find the specific output option of the data set you would like to export (e.g., Line Listing-All SSI Events). Click “Modify” next to the desired option; at the top of the modification page, click the button to “Export Analysis Data Set”. You can then select the data export format, such as Excel or SAS.


Note: This option does not allow you to export calculated rates or SIRs. For rates or SIRs, export the output dataset, described below.

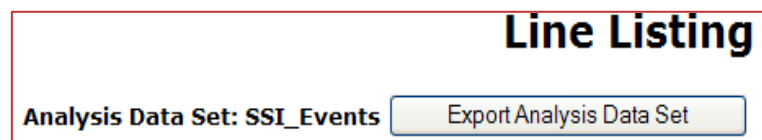

Line Listing

Analysis Data Set: SSI_Events Export Analysis Data Set

3. You can also export an output data set which would include any filtering options you have made, such as time period or location. On the bottom of the modification page, click “Export Output Data Set”. You can then select the data export format, such as Excel or SAS.



Run Save As Reset Back **Export Output Data Set**

