

# PPE CASE



## *Personal Protective Equipment Conformity Assessment Studies and Evaluations*

## **Performance of Stockpiled Air-Purifying Respirators, Facility Five of Ten: Inhalation and Exhalation Resistance and Filtration Efficiency Performance**

*National Institute for Occupational Safety and Health (NIOSH)  
National Personal Protective Technology Laboratory (NPPTL)  
Lee A. Greenawald, Susan M. Moore, Patrick L. Yorio  
March 1, 2020*

In the event of a national emergency, eighteen million U.S. healthcare workers may face high-consequence infectious disease exposures [NIOSH 2017]. Personal protective equipment (PPE), such as gowns, gloves, goggles, and respirators, is an important measure within the infection prevention hierarchy of controls. During public health emergencies, the sudden increase in PPE demand may exceed supplies for upwards of three months while manufacturers increase production [ASTHO 2013; Carias et al. 2015]; [Patel et al. 2017]. For example, during the 2009 H1N1 pandemic, local respirator shortages were reported and, during the 2016 Ebola outbreak and the first U.S. fatality, there was a 10-200 fold increase in PPE orders [DHHS 2012; NIOSH 2018]. To prepare for these shortages, large quantities of PPE are strategically stockpiled at hospital, local, state, and federal facilities [NIOSH 1997].

*NIOSH found that 602 of 602 N95 filtering facepiece respirators and P95 filters stockpiled at Facility Five that were 8-12 years old maintained their inhalation and exhalation resistance and filtration performance in accordance with NIOSH performance standards.*

Due to the decision to stockpile PPE, stockpile personnel and decision makers have sought to understand if stockpiled PPE is still viable following long-term storage. NIOSH does not require approval holders (i.e. those granted the approval from NIOSH) to designate a shelf life for particulate-only air-purifying respirators (APR), although some choose to do so and may provide this information on product packaging or online. There is limited published data to understand the viability of respirators that have undergone long-term storage with or without a designated shelf life. Over the past decade, the Strategic National Stockpile (SNS) and state and local stockpile personnel asked NIOSH to evaluate the performance of stockpiled PPE as well as better understand storage conditions in U.S. stockpile facilities that store PPE.

In 2017, NIOSH established a PPE Stockpile Partnership consisting of 1) federal entities and stockpiles; 2) state, county, and city stockpiles; 3) hospital-related stockpile entities; and 4) a manufacturer trade association to inform the design and execution of an empirical study to evaluate stockpiled APRs. NIOSH obtained samples of PPE from geographically dispersed stockpiles with varying storage conditions.

***This report details the inhalation/exhalation resistance and filtration performance of N95 filtering facepiece respirators and P95 filters collected from Facility Five of Ten. This facility is a state stockpile facility.***

---

# How NIOSH Evaluated Respirators and Storage Conditions

## Description of Facility Five

- NIOSH researchers visited Facility Five in March 2018 (**Figure 1**). This facility was located within the U.S. Department of Health and Human Services Region 9, representing Arizona, California, Hawaii and Nevada.

## Assessment of Storage Conditions

- NIOSH, in conjunction with the PPE Partnership members, developed checklists to document site and packaging (i.e. pallet, case, and box) conditions that may impact respirator performance.
- NIOSH documented the following storage conditions: 1) the PPE packaging presence of dust, shrink-wrapping, chemicals, and moisture, 2) exposure to sunlight and direct light; 3) proximity to fans, windows, doors, and ventilation systems; 4) damage to pallet and product packaging; and 5) location of pallet on storage rack (e.g., top, bottom) and location of PPE product on pallet (e.g., top/not load-bearing, bottom/load-bearing).



**Figure 1:** NIOSH researchers documented storage practices at Facility Five such as location and type of lighting, pallet stacking practices, and conditions of the flooring, roofing, and exterior walls.

- NIOSH collected facility temperature and percent relative humidity (%RH) data by placing two data loggers in the facility. This data was collected in 60-minute intervals from October 2017 to October 2018.

---

## Collection of Respirator Samples

- Facility Five’s inventory included APRs that are classified as N95 filtering facepiece respirators (FFRs), as well as one P95 filter cartridge model. Samples were collected from six different manufacturing models<sup>1</sup>: 1) 3M 1860 (two different manufacturing years); 2) 3M 1870; 3) 3M 8210; 4) 3M 9010; 5) 3M 2071; and 6) Kimberly Clark (KC) 46827 (**Table 1**).
- Upon reviewing the detailed APR inventories and storage location by lot within Facility Five, two different manufacturing lots for each model were identified and sampled within Facility Five. Two lots were sampled to evaluate and attempt to account for inter-lot variation. More than two lots were sampled from a model if conditions presented “worst-case scenario” conditions within the facility. Products were sampled and shipped to the NIOSH facility overnight to reduce exposure to non-climate-controlled conditions.
- Forty-three respirators were tested from each manufacturing lot for inhalation and exhalation resistance (n=3) and filtration performance testing (n=40)<sup>2</sup>.

## Selection of Control Respirators

- Control respirators of the same model as those sampled from the facility were purchased from the open market to be used as a comparison between stockpiled and new respirators.

## Characteristics of Sampled Respirators

- **Table 1** provides a summary of the respirator models sampled from Facility Five.

---

<sup>1</sup> Based on the other nine collaborating stockpiles’ inventories, these six models were sampled in order to compare performance within common respirator models when stored under disparate conditions.

<sup>2</sup> NIOSH testing requirements state that a minimum of three respirator units must be tested for inhalation and exhalation resistance and a minimum of 20 must be tested for filtration efficiency [NIOSH 2018].

**Table 1.** APRs Sampled from Stockpile Facility Five

Model	Lot #	Year of Manufacture	Shelf Life on Packaging?	Respirator Age at Time of Testing <sup>3</sup>	Shelf Life Status at Time of Testing
3M 1860	Lot A	2007	No	11 years	Past 5-year shelf life <sup>4</sup>
3M 1860	Lot A	2007	No		
3M 1860	Lot A	2010	No	8 years	Past 5-year shelf life <sup>4</sup>
3M 1860	Lot B	2010	No		
3M 1870	Lot A	2011	No	7 years	Past 5-year shelf life <sup>4</sup>
3M 1870	Lot B	2011	No		
3M 8210	Lot A	2006	No	12 years	Past 5-year shelf life <sup>4</sup>
3M 8210	Lot B	2006	No		
3M 9010	Lot A	2006	No	12 years	Past 5-year shelf life <sup>4</sup>
3M 9010	Lot B	2006	No		
3M 2071	Lot A	2006	No	12 years	Past 5-year shelf life <sup>4</sup>
3M 2071	Lot B	2006	No		
KC 46827	Lot A	2006	No	12 years	Past 5-year shelf life <sup>5</sup>
KC 46827	Lot B	2006	No		

### Evaluation of Inhalation and Exhalation Resistance and Filtration Performance

- Twenty-three control respirators were tested for inhalation and exhalation resistance and filtration performance. The 3M 1860 and 3M 2071 controls were manufactured in 2018, the 3M 8210 controls were manufactured in 2015, the 3M 1870 controls were manufactured in 2014, the 3M 9010 controls were purchased in 2018 but have an unknown manufacturing date, and the KC 46827 controls were manufactured in 2017. NIOSH testing requirements state that a minimum of three respirator units must be tested for inhalation and exhalation resistance. The same three respirators can be used for both inhalation and exhalation resistance testing [NIOSH 2018].
- Inhalation and exhalation resistance and filtration performance of the stockpiled and control respirators were evaluated using the same Standard Test Procedures (STPs) NIOSH uses for approving respirators under 42 Code of Federal Regulations Part 84, “Approval of Respiratory Protective Devices” [NIOSH 2018] (**Table 2**).
- **Table 2** describes the method for evaluating the inhalation and exhalation resistance and filtration performance of sampled respirators and control respirators.

<sup>3</sup> Testing was completed in 2018.

<sup>4</sup> 3M designated a five-year shelf life for these models [3M 2018]. As of February 2020, these models still have a five-year shelf life.

<sup>5</sup> KC designated a five-year shelf life for this model [KC 2018]. As of February 2020, this model still has a five-year shelf life.

**Table 2.** NIOSH Tests Conducted to Evaluate Inhalation and Exhalation Resistance and Filtration Performance.

NIOSH Standard Test Procedures (STPs)	Pass/Fail Criteria for APRs	Stockpiled Respirators Tested Per Manufacturing Lot	Control Respirators Tested
STP 3: Exhalation Resistance	<25 mm H <sub>2</sub> O column @ 85 liters per minute (LPM)	3 <sup>6</sup>	3 <sup>6</sup>
STP 7: Inhalation Resistance	<35 mm H <sub>2</sub> O column @ 85 LPM	3 <sup>6</sup>	3 <sup>6</sup>
STP 53: Liquid Particulate Filter Efficiency for P95	≤5.0% particulate penetration (≥95.0% filter efficiency)	40 <sup>7</sup>	20
STP 59: Particulate Filter Efficiency for N95	≤5.0% particulate penetration (≥95.0% filter efficiency)	40 <sup>7</sup>	20

## What NIOSH Found Through Inspection, Testing, and Evaluation

### Storage Conditions

- Visual Inspections—Dust and damage to product packaging was limited or not observed at Facility Five; examples of the most amount of dust and damage to product packaging are shown in **Figures 2 and 3**. No product boxes showed damage. Of the 602 respirators visually inspected, no concerns were noted.



**Figure 2:** Most amount of dust observed on product cases from Facility Five.

<sup>6</sup> NIOSH testing requirements state that a minimum of three respirator units must be tested for inhalation and exhalation resistance. The same three respirators can be used for both inhalation and exhalation resistance testing.


<sup>7</sup> An increased sample size was used for the stockpiled respirators as opposed to the control respirators to increase the precision of the performance estimates investigated.




**Figure 3:** Most amount of product case damage observed from Facility Five.

- Temperature nor %RH were controlled. Facility lights were off when not in use. A skylight allowed sunlight to enter the facility in specific places, including indirectly onto pallets storing respirators. No evidence of excess moisture or chemical spills that persisted beyond immediate mitigation were observed. Pallets were generally shrink wrapped around the four pallet sides but not across the top or bottom. Generally, all pallets were stacked three (or more) high causing some weight/load to be applied to the bottom pallet.
- Percent RH (**Figure 4**) and Temperature (**Figure 5**)
  - **At the time of publication**, the recommended storage requirements for %RH and temperature are
    - 3M 1860, 3M 1870, 3M 8210, 3M 9010, and 3M 2071: remain under 80 %RH; remain within -4°F to 86°F [3M 2017]
    - KC 46827: remain under 60 %RH; remain within 68°F to 77°F [KC 2020]
  - Data was aggregated across the two data loggers. The average temperature between the 2017-2018 time period was 70.2°F. The average %RH between 2017-2018 was 42.3%; these averages are within the 3M and KC recommended temperature and %RH storage conditions.
  - For the 3M models, 7.4% of the total temperature data points deviated from the recommended storage conditions. For the KC 46827 model, 44.5% of data points were below 68°F and 30.9% were above 77°F. For %RH, no data points deviated from 3M's recommended storage conditions, and 2.2% of data points deviated from KC's recommended conditions.

## Data Logger 1

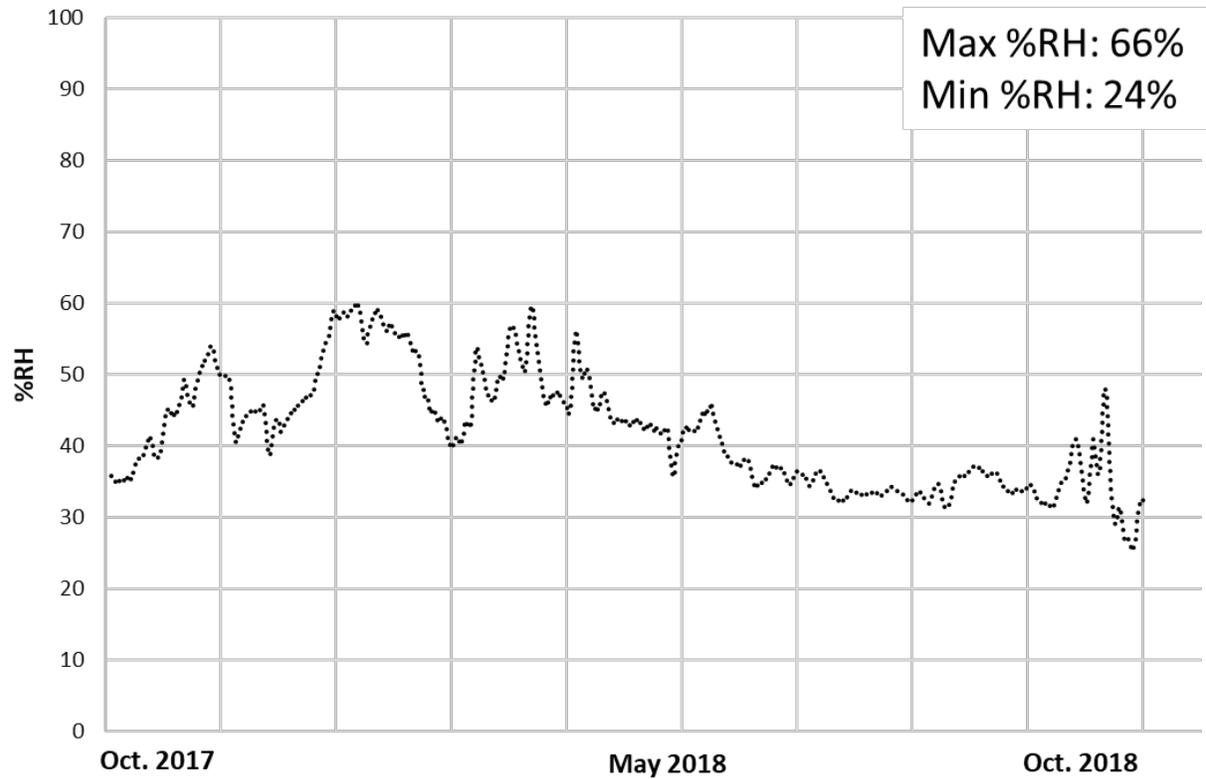


## Data Logger 2

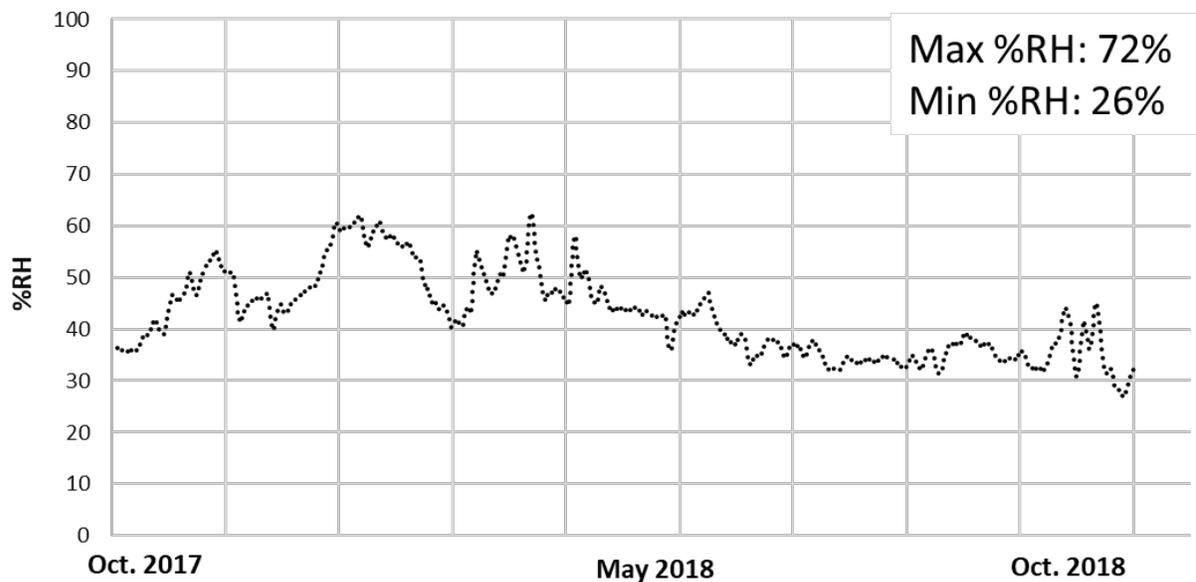


**Figure 4:** Percent Relative Humidity (% RH) from October 2017 – October 2018 for two data loggers stored at Facility Five. Data is plotted as a 50-point moving average for visualization purposes. Maximum and minimum temperatures reported are noted for each data logger.

### Data Logger 1

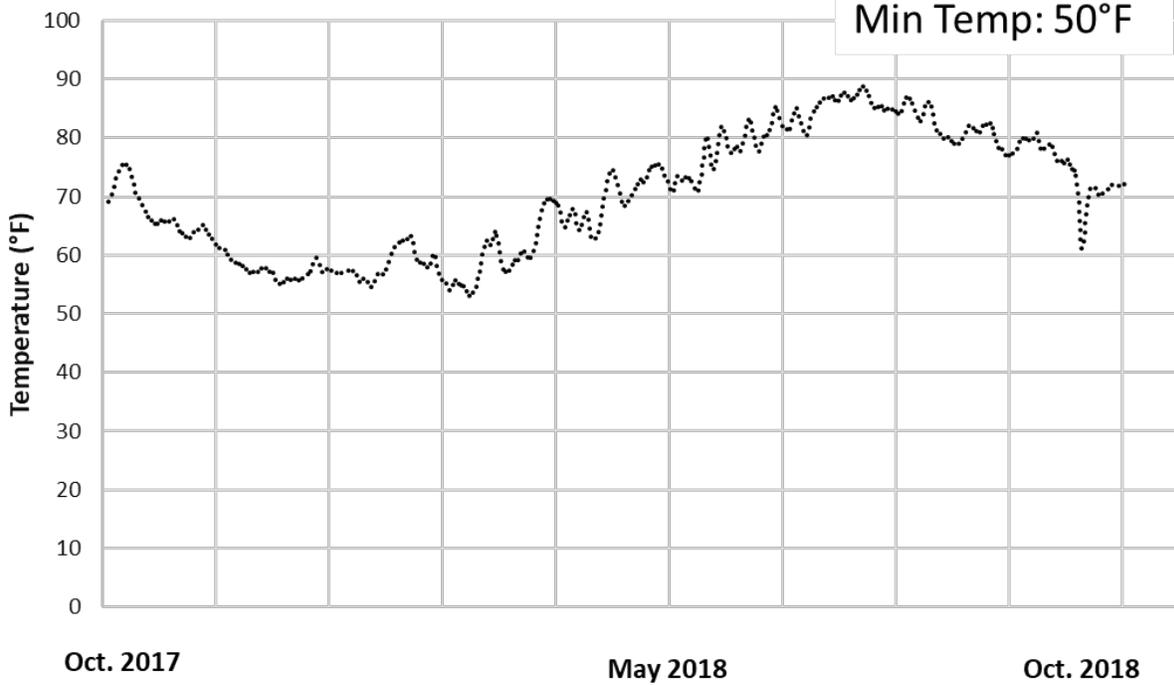

### Data Logger 2

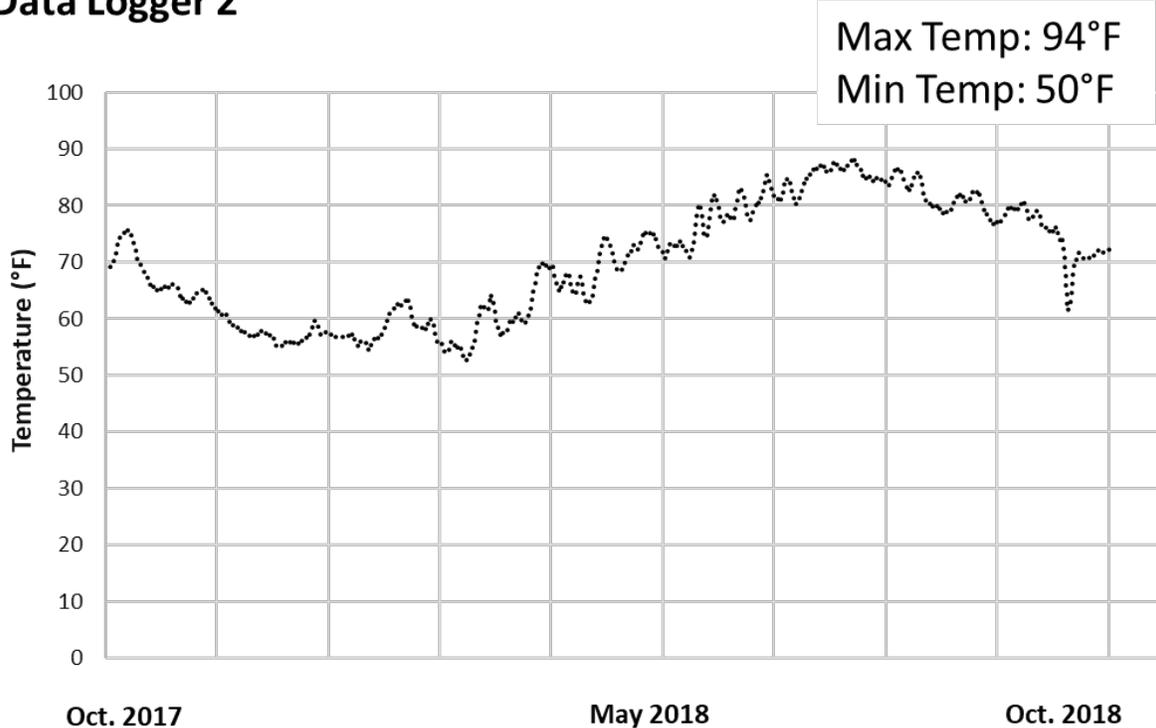


**Figure 5:** Temperatures from October 2017 – October 2018 for two data loggers stored at Facility Five. Data is plotted as a 50-point moving average for visualization purposes. Maximum and minimum temperatures reported are noted for each data logger.


---

## Inhalation and Exhalation Resistance


- NIOSH evaluated the inhalation and exhalation resistance for a total of 42 stockpiled and 18 control respirators. **All stockpiled and control respirators from each model passed these tests.** The 3M 1860, 3M 1870, 3M 9210, and 3M 9010 data are shown in **Figure 6**. The 3M 2071 (P95 filter cartridge) data is shown in **Figure 7**. The KC 46827 data is shown in **Figure 8**.
- Using an analysis of variance (ANOVA), there were no statistically significant differences (defined as  $\alpha \leq 0.05$ ) between the controls and stockpiled respirators for inhalation and exhalation resistance when averaging across models.
- When comparing the individual respirator models to their respective controls through an ANOVA with adjusted, post-hoc multiple comparisons, the following differences were observed for exhalation resistance: 1) 3M 1870 stockpiled 2006 Lots A and B were lower; 2) the 3M 8210 stockpiled 2006 Lot B was higher; 3) and the KC 46827 stockpiled 2006 Lot B was lower. For inhalation resistance the following differences were observed for inhalation resistance: 1) 3M 1870 stockpiled 2006 Lots A and B were lower; and 2) the 3M 8210 stockpiled 2006 Lot B was higher.
- For inhalation resistance, the individual stockpiled respirator with the highest resistance (16.00 mm H<sub>2</sub>O) was below the NIOSH maximum limit for product approval (35 mm H<sub>2</sub>O allowable maximum). For exhalation resistance, the individual stockpiled respirator with the highest resistance (13.46 mm H<sub>2</sub>O) was below the NIOSH maximum limit for product approval (25 mm H<sub>2</sub>O allowable maximum).

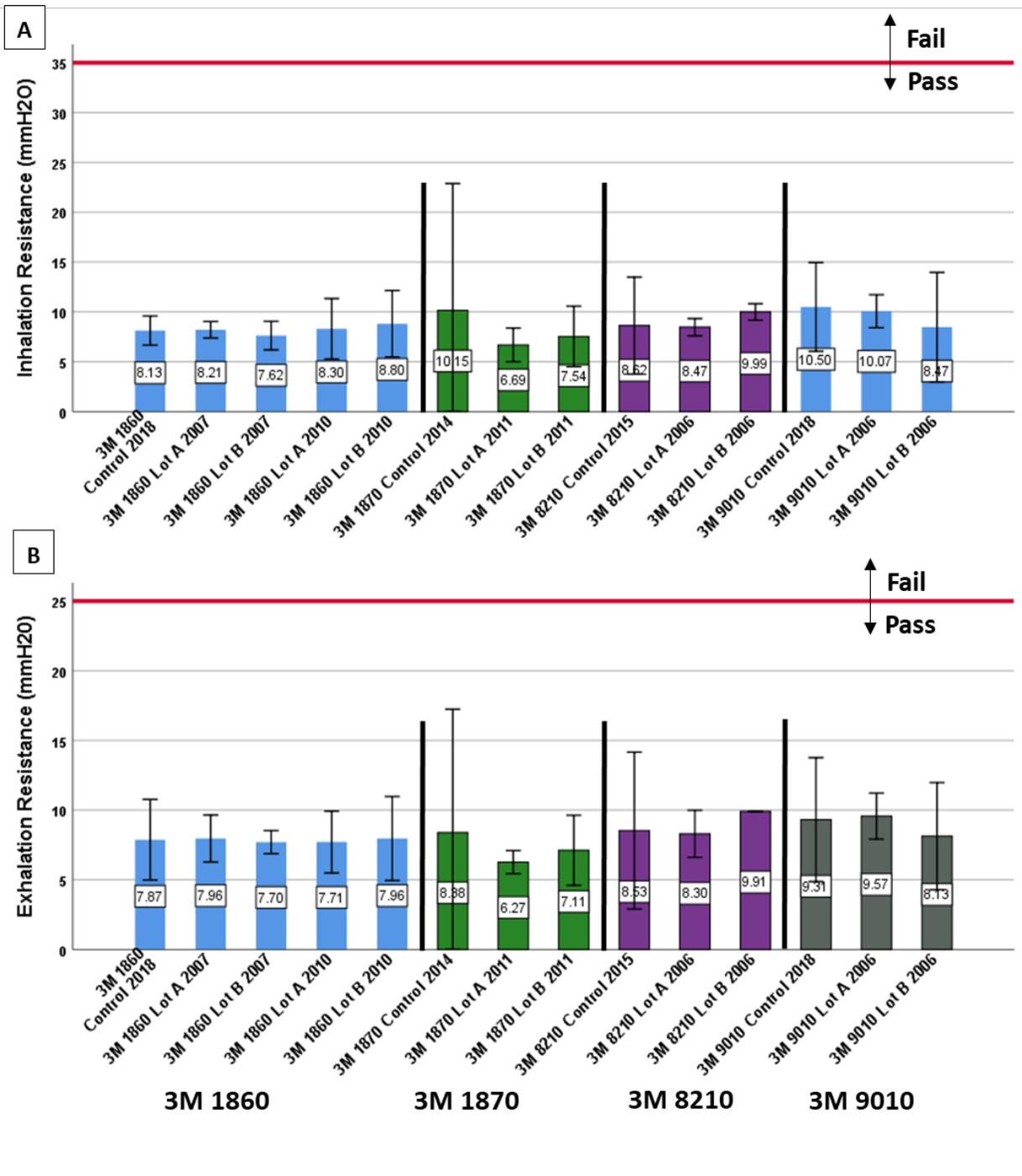

**Figure 6:** Control and stockpiled respirator inhalation (A) and exhalation (B) resistance data for the 3M 1860, 3M 1870, 3M 8210, and 3M 9010 models. N95 FFRs must have an inhalation resistance less than 35 mmH<sub>2</sub>O and an exhalation resistance less than 25 mmH<sub>2</sub>O. The pass/fail threshold for inhalation (A) and exhalation (B) resistance is shown by the red line. Error bars represent the 99% confidence interval and estimate the population parameters. This confidence interval suggests that 99% of any repeated samples tested and evaluated from this lot will have a mean between the upper and lower bounds.

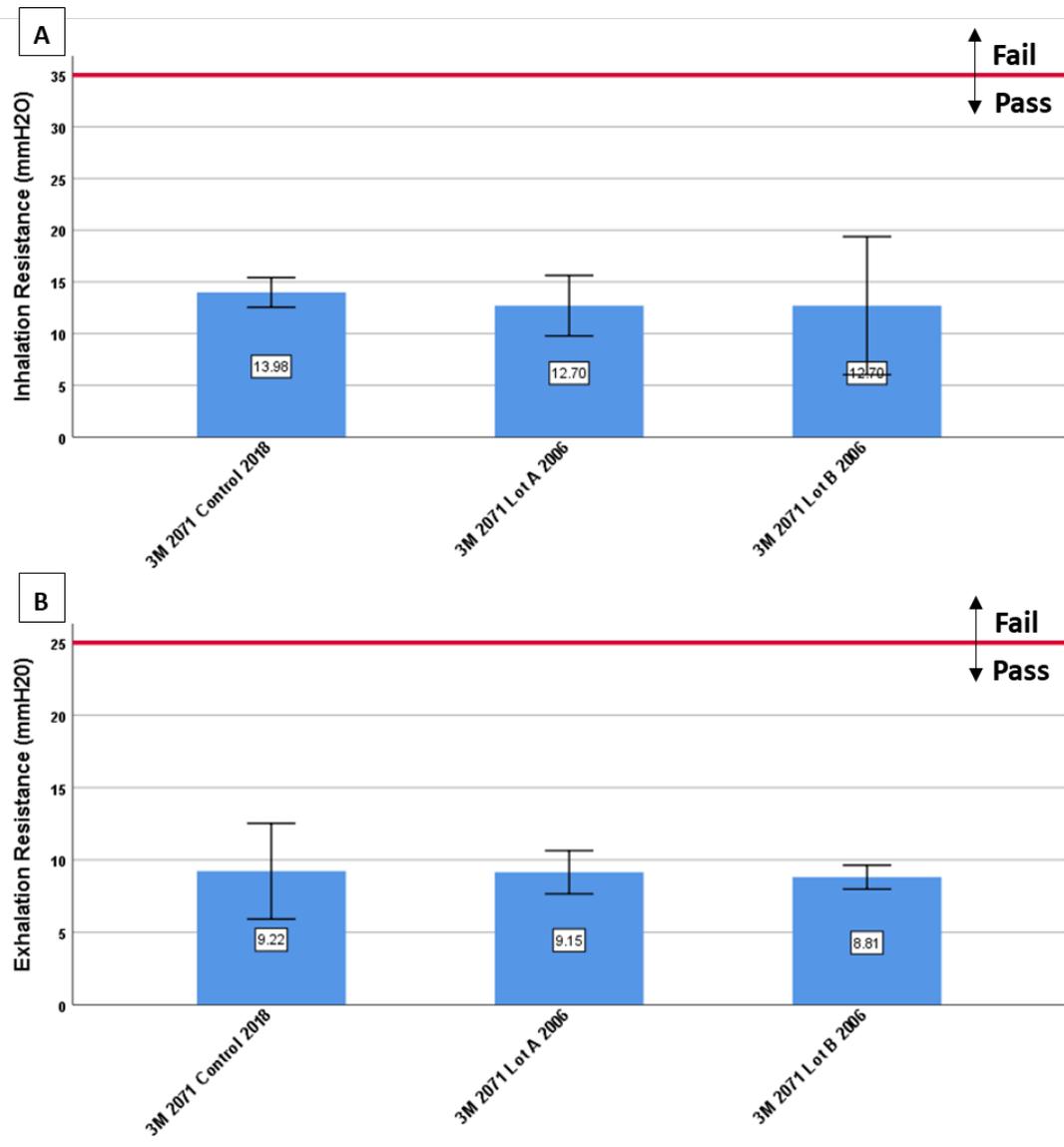


**Figure 7:** Control and stockpiled respirator inhalation (A) and exhalation (B) resistance data for the 3M 2071 model. P95 filters must have an inhalation resistance less than 35 mmH<sub>2</sub>O and an exhalation resistance less than 25 mmH<sub>2</sub>O. The pass/fail threshold for inhalation (A) and exhalation (B) resistance is shown by the red line. Error bars represent the 99% confidence interval and estimate the population parameters. This confidence interval suggests that 99% of any repeated samples tested and evaluated from this lot will have a mean between the upper and lower bounds.

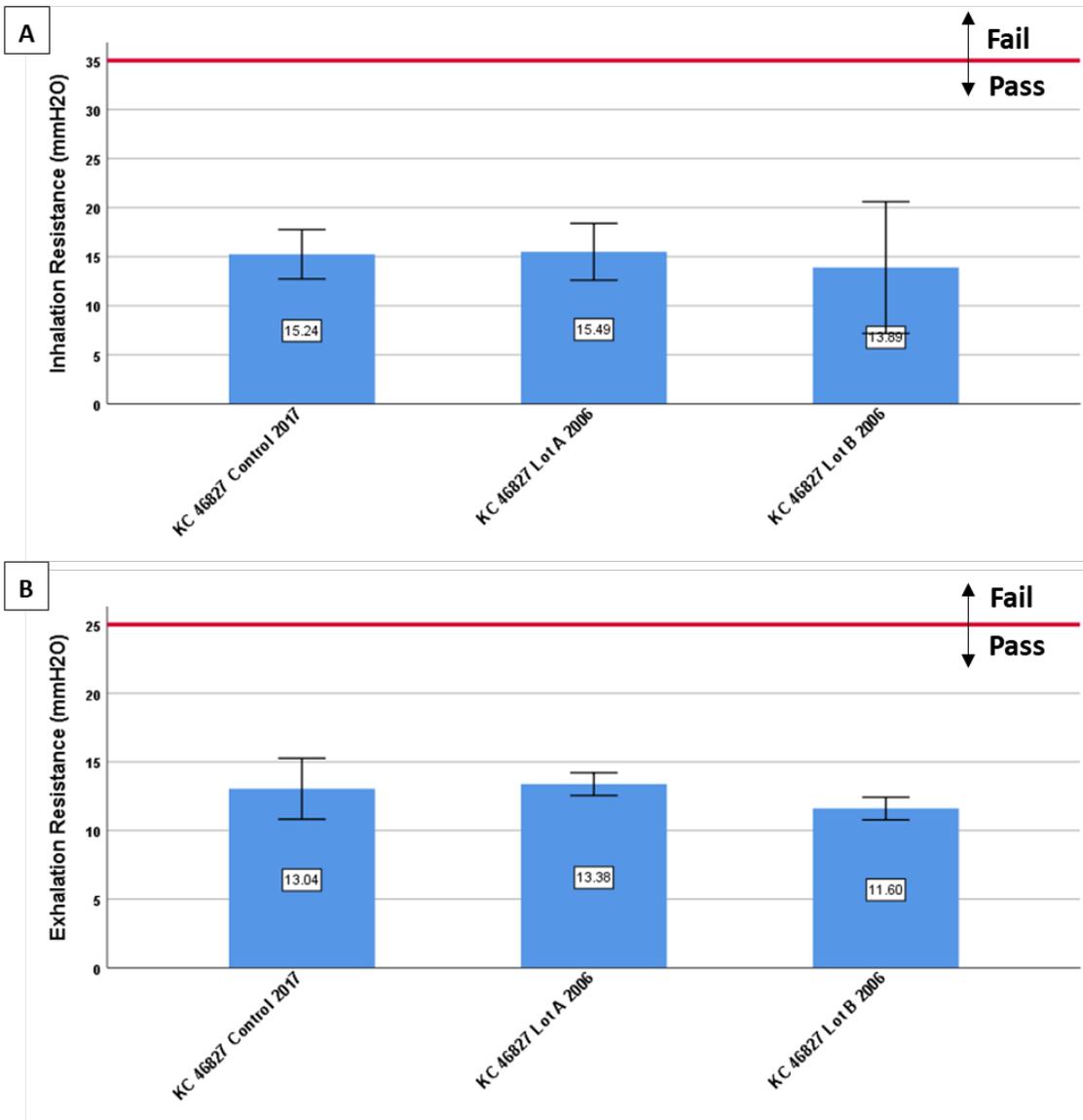


**Figure 8:** Control and stockpiled respirator inhalation (A) and exhalation (B) resistance data for the KC 46827 model. P95 filters must have an inhalation resistance less than 35 mmH<sub>2</sub>O and an exhalation resistance less than 25 mmH<sub>2</sub>O. The pass/fail threshold for inhalation (A) and exhalation (B) resistance is shown by the red line. Error bars represent the 99% confidence interval and estimate the population parameters. This confidence interval me suggests that 99% of any repeated samples tested and evaluated from this lot will have a mean between the upper and lower bounds.

## Filtration Performance


- NIOSH evaluated the particulate penetration efficiency for 560 stockpiled respirators and 120 controls. The 3M 1860, 3M 1870, 3M 9210, and 3M 9010 data are shown in **Figure 9**. The 3M 2071 (P95 filter cartridge) data is shown in **Figure 10**. The KC 46827 data is shown in **Figure 11**.
- None of the individual respirators tested exceeded the 5.0% maximum. The highest penetration for an individual stockpiled respirator was 4.53% (KC 46827) and the highest penetration for an individual control respirator was 2.84% (KC 46827).
- Using an analysis of variance (ANOVA), **there were no statistically significant differences (defined as  $\alpha \leq 0.05$ )** between the 3M controls and 3M stockpiled respirators for filtration when averaging across models. There was a statistically significant difference between the KC controls (mean percent penetration=1.77%, SD=0.45) and KC stockpiled respirators (mean percent penetration=2.73%, SD=0.68),  $p < 0.001$ .
- When comparing the individual respirator models to their respective controls through an ANOVA with adjusted, post-hoc multiple comparisons, the following statistically significant differences were detected: 1) 3M 1860 stockpiled 2007 Lot A had a higher penetration; 2) 3M 1860 stockpiled 2010 Lot A had a higher penetration; 3) 3M 1870 stockpiled 2011 Lot A had a higher penetration; 4) Both 3M 8210 stockpiled 2006 Lot A and Lot B had a lower penetration; and 5) Both KC 46827 stockpiled 2006 Lots A and B had a higher penetration.

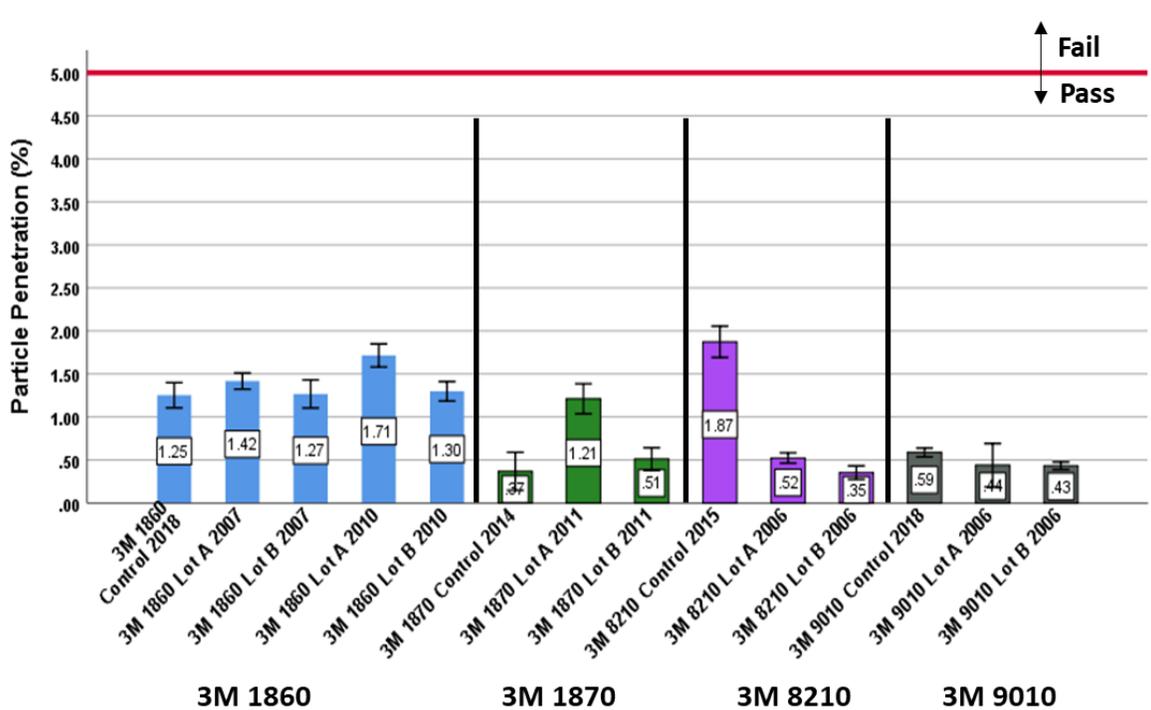

**Figure 9:** Control and stockpiled respirator particle filtration performance data for the 3M 1860, 3M 1870, 3M 8210, and 3M 9010 models. N95 FFRs must have a particle penetration of less than 5.0%. Error bars represent the 99% confidence interval and estimate the population parameters. This confidence interval suggests that 99% of any repeated samples tested and evaluated from this lot will have a mean between the upper and lower bounds.

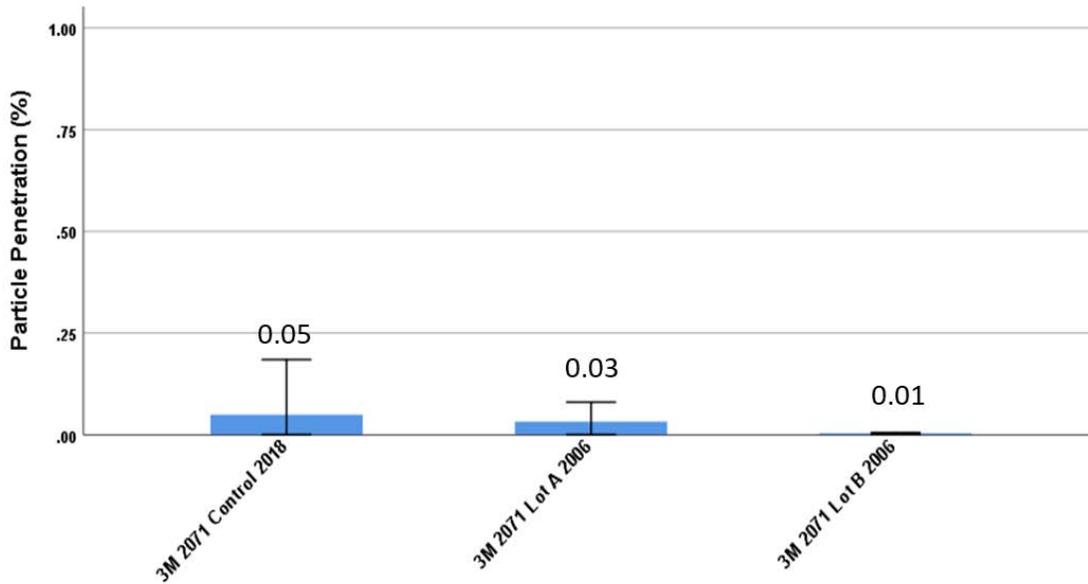

**Figure 10:** Control and stockpiled respirator particle filtration performance data for the 3M 2071 model. P95 FFRs must have a liquid particle penetration of less than 5.0%. Error bars represent the 99% confidence interval and estimate the population parameters. This confidence interval means that 99% of any repeated samples tested and evaluated from this lot will have a mean between the upper and lower bounds.

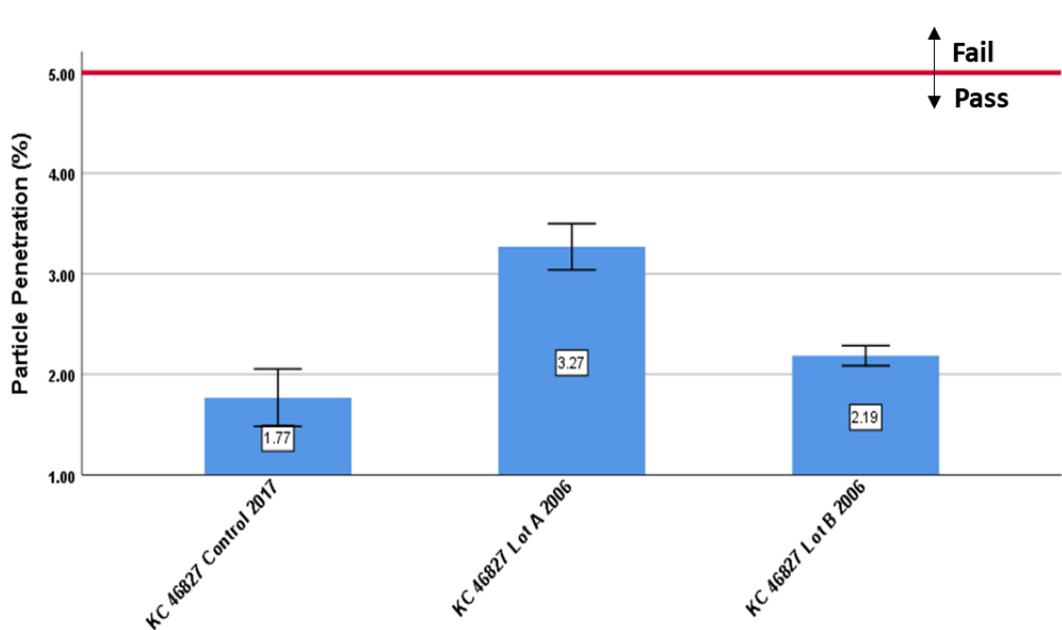

**Figure 11:** Control and stockpiled respirator particle filtration performance data for the KC 46827 model. N95 FFRs must have a particle penetration of less than 5.0%. Error bars represent the 99% confidence interval and estimate the population parameters. This confidence interval suggests that 99% of any repeated samples tested and evaluated from this lot will have a mean between the upper and lower bounds.

---

## CASE Findings

### [Findings for the KC 46827 Model:](#)

**No failures for inhalation resistance, exhalation resistance, or filtration performance were observed—i.e., the performance data suggests that these units would be protective so long as a proper fit is achieved.** This model currently has a five-year recommended shelf life; Appendix 1 shows a KC letter to end users with shelf life information, which states respirators past their shelf life should be discarded [KC 2018]. Thus, these respirators tested are past their recommended shelf life. These findings pertain to KC units from Facility Five and may not be applicable to other stockpile facilities and/or under different environmental storage conditions.

### [Findings for the 3M 1860, 3M 1870, 3M 8210, 3M 9010, and 3M 2071 Models:](#)

**No failures for inhalation resistance, exhalation resistance, or filtration performance were observed—i.e., the performance data suggests that these units would be protective so long as a proper fit is achieved.** These five models currently have a five-year recommended shelf life; Appendix 2 shows two 3M letters to end users with shelf life and recommended storage condition information [3M 2018, 3M 2020]. Thus, all respirators tested are past their recommended shelf life. These findings pertain to 3M units from Facility Five and may not be applicable to other stockpile facilities and/or under different environmental storage conditions.

### [Stockpile Storage Conditions:](#)

For the 3M models, 7.4% of the total temperature data points deviated from the previously described recommended storage conditions. For the KC 46827 model, 75.4% of the total temperature data points deviated from the previously described recommended storage conditions. For %RH, no data points deviated from 3M's recommended storage conditions, and 2.2% of data points deviated from KC's recommended conditions. Stored under these conditions, NIOSH found that 602 APRs evaluated in this study, which were 8-12 years old, maintained their inhalation and exhalation resistance and filtration performance (i.e. all sampled respirators were below the NIOSH maximum limit as defined by 42 CFR Part 84).

NIOSH regulation sets the minimum quality and performance requirements for the approval of respirators [NIOSH 1997]. NIOSH does not have requirements for shelf life or storage conditions for particulate-only APRs. The approval holder<sup>8</sup> (i.e. the entity that is granted the approval from NIOSH) is responsible for understanding how their products' design or performance may be affected by various use or storage conditions and must provide instruction for establishing the proper use, storage, and maintenance procedures for their approved products, which may include designating a shelf life [NIOSH 2019]. FFR or particulate filter packaging (such as the box) often includes NIOSH-approved user instructions, label information, and recommendations on shelf life. Additionally, some approval holders also disseminate recommendations related to storage and shelf life through resources such as user and web notices. The respirators tested in this study were generally not designed for long-term storage.

At this time, we do not have enough information to definitively know the level of protection that may be provided by respirators that 1) are stored for prolonged periods of times; 2) are stored under various storage conditions; or 3) have exceeded the approval holder's designated shelf life. Users of respirators that have

---

<sup>8</sup> An approval may be granted to a non-manufacturing entity.

---

exceeded the designated shelf life should be forewarned to avoid a false sense of confidence; these devices may not provide the same level of protection as those that have not exceeded the designated shelf life. We recommend contacting the approval holder(s) of the respirators in the stockpile with specific questions regarding the use of product beyond the designated shelf life.

---

## What Can Stockpile Personnel Do to Learn More about the Respirators in their Stockpile?

- Stockpile personnel should check the product information from the approval holder as well as the NIOSH Certified Equipment List to remain up-to-date on product storage conditions, shelf-life information, and NIOSH approval status. Check NIOSH's Certified Equipment List to verify the respirator model currently maintains its NIOSH approval at <https://www.cdc.gov/niosh/npptl/topics/respirators/cel/default.html>
- Stockpile personnel should work with the approval holder(s) of the stockpiled products with specific questions regarding the use of expired product.
- Sign up for NPPTL's Listserv at <https://www.cdc.gov/niosh/npptl/sub-NPPTL.html> to receive email notifications relevant to PPE.

For more information related to personal protective equipment, visit the [NIOSH NPPTL website](https://www.cdc.gov/niosh/npptl/) <https://www.cdc.gov/niosh/npptl/>

### Get More Information

Find NIOSH products and get answers to workplace safety and health questions:

1-800-CDC-INFO (1-800-232-4636) | TTY: 1-888-232-6348

CDC/NIOSH INFO: [cdc.gov/info](https://www.cdc.gov/info) | [cdc.gov/niosh](https://www.cdc.gov/niosh)

Monthly NIOSH eNews: [cdc.gov/niosh/eNews](https://www.cdc.gov/niosh/eNews)

---

All photos courtesy of NIOSH NPPTL.

### Disclaimer

The recommendations in this report are made based on the findings at the stockpile evaluated and may not be applicable to other stockpile facilities.

Mention of any company or product does not constitute endorsement by the National Institute for Occupational Safety and Health (NIOSH). In addition, citations to websites external to NIOSH do not constitute NIOSH endorsement of the sponsoring organizations or their programs or products. Furthermore, NIOSH is not responsible for the content of these websites. All web addresses referenced in this document were accessible as of the publication date.



Centers for Disease Control  
and Prevention  
National Institute for Occupational  
Safety and Health

---

## Suggested Citation

NIOSH [2020]. PPE CASE: Inhalation and Exhalation Resistance and Filtration Performance of Stockpiled Air-Purifying Respirators: Facility Five of Ten. By Greenawald, L., Moore, S., and Yorio, P. Pittsburgh, PA U.S. Department of Health and Human Services, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health, NPPTL Report Number P2020-0105.

## References

- 3M [2017]. Frequently Asked Questions: 3M Filtering Facepiece/Disposable Respirator Storage Conditions and Shelf Life. St. Paul, MN. <https://multimedia.3m.com/mws/media/10158530/faq-3m-filtering-facepiece-disposable-respirator-storage-conditions-and-shelf-life.pdf> [accessed February 27, 2020].
- 3M [2018]. 3M Personal Safety Division Letter. St. Paul, MN: The 3M Company. <https://multimedia.3m.com/mws/media/14303830/3m-filtering-facepiece-shelf-life-letter.pdf>.
- 3M [2020]. Frequently Asked Questions: 3M Health Care Particulate Respirator and Surgical Masks Storage Conditions and Shelf Life. St. Paul, MN. <http://multimedia.3m.com/mws/media/8692380/3m-health-care-particulate-respirator-and-surgical-masks-storage-conditions-and-shelf-life-faq.pdf>
- ASTHO [2013]. Federal Emergency Preparedness Directives and Systems- Emergency Authority and Immunity Toolkit. Arlington, VA: Association of State and Territorial Health Officials. <http://www.astho.org/Programs/Preparedness/Public-Health-Emergency-Law/Emergency-Authority-and-Immunity-Toolkit/Federal-Emergency-Preparedness-Directives-and-Systems-Fact-Sheet/>.
- Carias C, Rainisch G, Shankar M, Adhikari BB, Swerdlow DL, Bower WA, Piliai SK, Meltzer MI, Koonin LM [2015]. Potential demand for respirators and surgical masks during a hypothetical influenza pandemic in the United States. *Clinical Infectious Diseases* 60 (1): S42-S51 <https://doi.org/10.1093/cid/civ141>.
- DHHS [2012]. An HHS retrospective on the 2009 H1N1 influenza pandemic to advance all hazards preparedness. Washington D.C: U.S. Department of Health and Human Services, <https://www.phe.gov/Preparedness/mcm/h1n1-retrospective/Documents/h1n1-retrospective.pdf>.
- KC [2018]. Kimberly Clark Letter to Customers (Appendix 1)
- KC [2020]. Kimberly Clark representatives, email correspondence, February 25, 2020.
- NIOSH [1997]. 42 CFR Part 84 respiratory protective devices. Cincinnati, OH: U.S. Department of Health and Human Services, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health, <https://www.cdc.gov/niosh/npptl/topics/respirators/pt84abs2.html>.
- NIOSH [2017]. Workplace safety and health topics: healthcare workers. Cincinnati, OH: U.S. Department of Health and Human Services, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health, <https://www.cdc.gov/niosh/topics/healthcare/default.html>.
- NIOSH [2018]. Standard respirator testing procedures. Cincinnati, OH: U.S. Department of Health and Human Services, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health, <https://www.cdc.gov/niosh/npptl/stps/apresp.html>.
- NIOSH [2019]. NIOSH conformity assessment notice: NIOSH CA 2019-1012, NIOSH respirator approval contents and meaning. Pittsburgh, PA: U.S. Department of Health and Human Services, Centers for Disease Control and Prevention, National Institute for Occupational Safety and Health, <https://www.cdc.gov/niosh/npptl/resources/pressrel/letters/conformitynotice/CA-2019-1012.html>.
- Patel A, D'Alessandro M, Ireland KJ, Burel WG, Wencil EB, Rasmussen SA [2017]. Personal protective equipment supply chain: lessons learned from recent public health emergency responses. *Health Security*. 15 (3): 244-252 <https://doi.org/10.1089/hs.2016.0129>.

---

Rottach DR, Lei Z [2017] Stockpiled N95 filtering facepiece respirator polyisoprene strap performance. J Int Soc Respir Prot. 34 (2): 69-80 <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC6198819/>.

## Appendix 1 [KC 2018]



Date: June 7, 2018

Subject: Kimberly-Clark\* N95 Particulate Filter Respirator and Surgical Mask Shelf Life (Codes: 62355, 62126, 46827, 46727, 46867, and 46767)

Dear Valued KCP Customer,

Since 2014, all Kimberly-Clark\* N95 Particulate Filter Respirator and Surgical Mask† packaging has included the storage conditions within the user instructions and the expiration date printed on each dispenser. If you have product in inventory produced prior to 2014 and without a printed expiration date, confirm that the product is within the recommended five year shelf life prior to using. Verify either through your purchase records or by contacting us with the printed lot number to determine the date of manufacture. We recommend disposing of any product that is beyond the established shelf life, has not been stored according to the user instructions, is damaged, does not provide a proper fit, or has missing parts.

For further information regarding the shelf life or interpreting the lot number to determine the expiration date, please contact us via the Kimberly-Clark Professional\* Technical or Quality hotline at 888-346-4652, email [kcpinfo@kcc.com](mailto:kcpinfo@kcc.com)

Thank you for your continued business and support of Kimberly-Clark Professional\*.

†Kimberly-Clark\* N95 Particulate Filter Respirator and Surgical Mask codes: 62355, 62126, 46827, 46727, 46867, and 46767

## 3M Personal Safety Division

St. Paul, MN 55144-1000

June, 2018

Dear Valued Customer:

Thank you for your inquiry regarding the shelf life of 3M filtering facepiece respirators. 3M is currently in the process of establishing a shelf life for various filtering facepiece respirators. The table below provides a list of models that currently have or will shortly have storage conditions and shelf life information communicated in either the User Instructions and/or packaging in the form of symbols, printed use by dates, etc. As some models were not introduced with storage conditions/shelf life markings, the year listed indicates when these packaging updates were implemented for that particular model.

**Table A.** 3M Filtering Facepiece Respirators

Model	Years of Implementation on Package	Shelf Life from Date of Manufacture*
8110S	2016	5 years
8200 (AAD#07023)	2009	5 years
8210	2016	5 years
8210Plus (AAD#07048)	2016	5 years
8210V	2011	5 years
8211	2014	5 years
8212	2012	3 years
8214 (AAD#07187)	2012	3 years
8233	2018	5 years
8240	Future	5 years
8246	2017	3 years
8247 (AAD#07186)	2017	3 years
8271	Future	5 years
8293	2018	5 years
8510**	2010	5 years
8511 (AAD#07185)	2014	5 years
8512	2012	3 years
8514	2012	3 years
8515 (AAD#07189)	2012	2 years
8516	Future	3 years
8576	2017	3 years
8577	2017	3 years
9105, 9105S	2010	5 years
9210** (AAD#37021)**	2010	5 years
9211** (AAD#37022)**	2010	5 years
9210+ (AAD#37192)	2013	5 years
9211+ (AAD#37193)	2013	5 years
<b>Medical</b>		
1804, 1804S	2018	5 years
1860, 1860S	2013	5 years
1870**	2013	5 years
1870+	2013	5 years

\* Please refer to respirator user instructions and packaging for specific storage conditions and use by date information.

\*\*Discontinued.

---

## 3M Personal Safety Division

**i IMPORTANT NOTE**

OSHA requires that all respirators be stored to protect them from damage, contamination, dust, sunlight, extreme temperatures, excessive moisture, and damaging chemicals, and that they be stored to prevent deformation of the facepiece and exhalation valve. Always follow the product's User Instructions, including that respirators should always be inspected prior to use and discarded if damage to any component is observed.

Additional information, including updates regarding shelf life and storage conditions of 3M Filtering Facepiece Respirators can be found at [www.3M.com/workersafety](http://www.3M.com/workersafety). The following resource documents are offered for your reference:

- [Shelf Life FAQ Industrial Filtering Facepiece - Disposable Respirators](#)
- [Shelf Life FAQ Health Care Particulate Respirators and Surgical Mask](#)

Please call 3M Personal Safety Division's Technical Service at 1-800-243-4630 if you have further questions. Thank you for using 3M products.

Sincerely,

3M Personal Safety Division (PSD)

## Frequently Asked Questions: 3M Health Care Particulate Respirator and Surgical Masks Storage Conditions and Shelf Life

### Why is 3M adding shelf life information for the 3M™ Health Care Particulate Respirator and Surgical Masks\* 1804/1804S, 1860/1860S, 1870, 1870+?

The addition of shelf life information to our 3M NIOSH-approved respirators is a way to communicate to our customers the storage conditions and potential longevity of our respirators. Traditionally the life cycle of these respirators commonly used in health care workplace applications, from date of manufacture to use by the customer, has been short in duration as they are disposable. However, with the increased attention to respirator stockpiling, many customers have requested information on storage conditions and shelf life. We hope that by adding this information to the respirator packaging it will encourage our customers to employ good practices such as appropriate long term storage, rotation of stock and inventory management.

In the United States, per 29 CFR 1910.134, OSHA has required that respirators be stored in the original packaging and away from contaminated areas, dust, sunlight, extreme temperatures, excessive moisture and damaging chemicals. Canada's CSA Standard Z94.4 has a similar requirement.

\*Models that are both NIOSH approved N95 filtering facepiece respirators and FDA cleared as a surgical mask.

### What 3M Health Care Particulate Respirator and Surgical Masks have a shelf life?

The 3M Health Care Particulate Respirator and Surgical Mask models 1804/1804S, 1860/1860S, 1870, 1870+ have an established 5 year shelf life when respirators are stored in their original packaging within climatic conditions ranging from -4 °F (-20 °C) to +86 °F (+30 °C) and not exceeding 80% RH.

### Why does the packaging for some 3M Health Care respirators have shelf life information and other respirator packaging does not?

The transition to updated packaging/labeling in relation to the storage conditions and shelf life has been initiated. However, for a period of time, you may see product packaging in the market place with and without storage condition and shelf life information included/incorporated.

### How is the respirator's shelf life communicated?

The shelf life information is usually found on the side or bottom of the primary box. Storage conditions are included in the instructions for use (IFU). The shelf life for the health care NIOSH-approved respirators is in the form of a "use by" date such as "YYYY-MM-DD" (year-month-day) and should be located near the hourglass icon. This information is also located on the label of the shipper case or corrugated box. An explanation of the icons and additional information regarding shelf life and storage conditions can be found in the IFU provided with the respirator. Please refer to the respirator packaging as shelf life is specific to each model.

Here is an example of how storage conditions and shelf life will be depicted in the IFU and primary box respectively (this is an example only):

*When stored in original packaging between temperatures from -4 °F (-20 °C) to +86 °F (+30 °C) and not exceeding 80% RH, the respirator may be used until the date specified on packaging located next to the "Use by Date" symbol.*

## 3M Personal Safety Division



Use by Date

Here are some additional symbols that you will see in the updated instructions for use.



Date of Manufacture



Manufacturer's Lot Number relevant to the device bearing the symbol



Manufacturer

### What happens if storage conditions are not met?

3M's goal is to help our customers ensure that filtering facepiece respirators stored for extended periods of time will meet the performance requirements to which they were approved and function as intended. When establishing a shelf life, 3M takes into account the filter media as well as the component parts of the respirator such as the strap and any staples. Therefore, we are confident that the respirators will meet performance requirements when the identified conditions are met.

However, when respirators are maintained outside of the established storage conditions, 3M cannot ensure that the respirators will meet performance requirements. In this event, many different kinds of changes can occur to the respirator including cosmetic changes and degradation of components such as headbands, nose foam and noseclips. Examples of cosmetic changes include discoloration of materials. Examples of degradation include crumbling of nose foam or breaking of headbands.

It is always critical that the respirator be inspected and a user seal check be conducted by the wearer per the IFU. If the person wearing the respirator cannot achieve a proper seal the respirator should not be used.

### How do we know when not to use the respirator?

First refer to the packaging for a "use by" date. 3M's recommendation is that respirators be disposed of after the stated use by date. Always inspect the respirator and conduct a user seal check before use per the IFU. If the person wearing the respirator cannot achieve a proper seal, then the respirator should not be used. Even for respirators within the stated shelf life, the respirator should be disposed of immediately upon observation of damaged or missing parts. For those respirators that have established shelf life but which packaging is not yet marked with a "use by" date, 3M recommends they no longer be used if 5 years has passed since the date of manufacture.

### If the respirator is not marked with shelf life information, how can I determine the age of the respirator?

For respirators that are not currently labeled with shelf life information, the date of manufacture can be determined from the label or printed information located on the primary box as well as the shipper case or corrugated box. For assistance in interpreting the date of manufacture, please call 3M Health Care Helpline at 1-800-228-3957 if in the U.S. In Canada call 1-800-267-4414. Release 5, February 2020. Other countries please contact your local 3M office.

### Is it okay to exceed storage conditions and, if so, for how long?

It is recognized that recommended storage conditions may be exceeded for short periods of time during transportation. This has been accounted for in the shelf life determination. However, storage outside the recommended conditions should be avoided when possible.

## 3M Personal Safety Division

### Should the respirator be disposed of after the shelf life has expired?

3M's recommendation is that the respirator be disposed of after the stated use by date has expired.

### Will 3M take back respirators that have reached the end of their stated shelf life?

No, 3M will not accept returns of respirators on the basis of shelf life.

### Will all 3M respirators have the same shelf-life?

No, not all 3M respirators will have the same shelf life. In making shelf life determinations, 3M takes into account the filter media as well as the components of the respirator. Components vary from model to model. See the [3M Filtering Facepiece Shelf Life](#) document for model specific information.

**Personal Safety Division**  
3M Center, Building 235-2W-70  
St. Paul, MN 55144-1000

3M PSD products are  
occupational use only.

**In United States of America**  
Technical Service: 1-800-243-4630  
Customer Service: 1-800-328-1667  
[3M.com/workersafety](http://3M.com/workersafety)

**In Canada**  
Technical Service: 1-800-267-4414  
Customer Service: 1-800-364-3577  
[3M.ca/Safety](http://3M.ca/Safety)

© 3M 2020. All rights reserved.  
3M is a trademark of 3M Company and  
its affiliates.  
Used under license in Canada.  
All other trademarks are property  
of their respective owners.  
Please recycle.  
Release 5, February 2020

